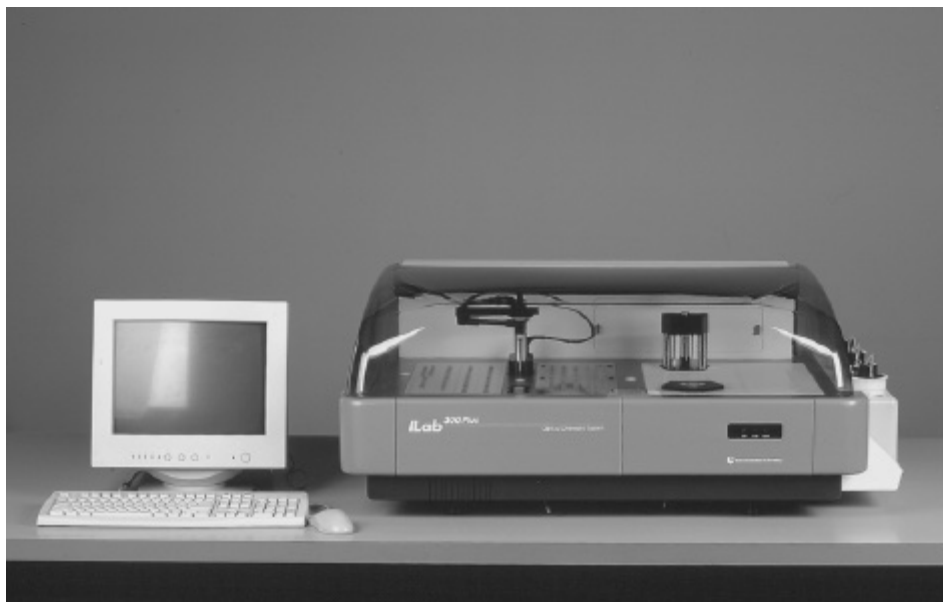


# ILab 300 Plus Service Manual

---

**P/N xxxxxx**

**Revision 02      -      January 2004**



This publication and all materials (including software) concerning the products of ILab 300 Plus System are of proprietary nature and are communicated on a strictly confidential basis.; they may not be reproduced, recorded, stored in a retrieval system, transmitted or disclosed in any way and by any means whatsoever, whether electronic, mechanical through photocopies or otherwise, without IL's prior written consent. Information contained herein is believed by IL to be accurate: in any event, no responsibility, whether express or implied, is assumed hereby by IL for or in connection with the use thereof, or for infringement of any third party right which arise therefrom, or from any representation or omission contained herein. Information is subject to change and/or updating without notice.

## **ILab 300 Plus SERVICE MANUAL**

---

### **GENERAL INDEX**

◆ CHAPTER 1	INTRODUCTION
◆ CHAPTER 2	SYSTEM DESCRIPTION
◆ CHAPTER 3	INSTALLATION
◆ CHAPTER 4	ELECTRICAL SCHEMES AND DRAWINGS
◆ CHAPTER 5	DIAGNOSTIC PROGRAM
◆ CHAPTER 6	SETTINGS AND ADJUSTMENTS
◆ CHAPTER 7	MAINTENANCE
◆ CHAPTER 8	HOST INTERFACING
◆ CHAPTER 9	TROUBLESHOOTING (to be done)
◆ CHAPTER 10	OPTIONS (to be done)
◆ CHAPTER 11	BAR CODE READER (to be done)
◆ CHAPTER 12	ISE MODULE
◆ CHAPTER 13	SPARE PARTS LIST

# CHAPTER 01

---

## - INTRODUCTION -

### INDEX

1	INTRODUCTION .....	1
	1.1 THE AIM OF THE TECHNICAL MANUAL .....	2
1.2	SYSTEM INTRODUCTION .....	2
1.3	PRECAUTIONARY MEASURES.....	2
	1.3.1 CHEMICAL RISKS .....	2
	1.3.2 ELECTRICAL RISKS.....	2
	1.3.3 MECHANICAL RISKS.....	3
1.4	TECHNICAL OPERATING FEATURES.....	3
	1.5.1 COMPUTER & SOFTWARE FEATURES .....	4
	1.5.2 OPTIONS.....	5
	1.5.3 DIMENSIOIS, WEIGHT & ENVIRONMENT .....	5
	1.5.4 INSTALLATION REQUIREMENTS.....	6

# **1 INTRODUCTION**

## **1.1 THE AIM OF THE TECHNICAL MANUAL**

This manual has been written in order to supply the technical staff, the persons who are responsible for the maintenance and for resolving instrument failures, a complete and detailed guide of the ILab 300 Plus analyzer, in accordance with the standard UNI EN591 (which requires a manual to be supplied with vitro diagnostic instruments for professional use).

## **1.2 SYSTEM INTRODUCTION**

The “ILab 300 Plus” system is a continuous-load, random access, bench top instrument for performing chemical and immuno-turbo-dymetric clinical analysis. It is totally automatic and computer controlled.

## **1.3 PRECAUTIONARY MEASURES**

### **1.3.1 Chemical risks**

The individual operator is responsible for assuring that all possible precautionary measures are taken against eventual risks associated with the use of the “ILab 300 Plus” instrument in clinical laboratory settings. The manufacturer will provide the reagents kit and specific written information on the use of each of the reagents.

It is important that the samples be well coagulated and then carefully centrifuged.

Samples which contain fibrinogen clots can obstruct the probe and lead to inexact sampling.

If blood samples containing gel are used, it is suggested that the manufacturer’s recommendations be followed.

Immediately clean and remove any accidental leakage of reagent or other liquid.

### **1.3.2 Electrical risks**

As with any electrical device, the risk of electric shock exist.

Is therefore necessary to take every precautionary measure possible when working with this, or any other, electrical instrument to avoid contact with power supply wires, electrical components or electronic boards.



### 1.3.3 Mechanical risks

Several precautionary measures should be taken when operating the analyzer:

**avoid** wearing very loose clothing or jewelry that could become tangled in the instrument's moving parts (e.g. the sample probe); whenever possible, operate the instrument with the main cover panel lowered.

**WARNING: Never attempt to service or substitute any part(s) of the analyzer when the instrument is turned on.**

Any and all technical repairs or servicing must be performed by specialized personnel only.

## 1.4 TECHNICAL OPERATING FEATURES

<b>DESCRIPTION</b>	<ul style="list-style-type: none"> <li>◆ Fully automatic, random access, continuous loading, benchtop analyzer for clinical chemistry and immunoturbidimetric assays</li> </ul>
<b>ASSAY TYPE</b>	<ul style="list-style-type: none"> <li>◆ End Point, Initial Rate, Kinetic, Bichromatic, Differential,</li> </ul>
<b>TEST ENTRY MODE</b>	<ul style="list-style-type: none"> <li>◆ Selective, Batch, Profiles, STAT</li> </ul>
<b>THROUGHPUT</b>	<ul style="list-style-type: none"> <li>◆ 200 test per hour without ISE module.</li> </ul>
<b>WORKING TEMPERATURE</b>	<ul style="list-style-type: none"> <li>◆ 37° C</li> </ul>
<b>ON LINE REAGENTS</b>	<ul style="list-style-type: none"> <li>◆ 4 removable racks (two racks cooled by Peltier cells and two racks at room temperature).               <ul style="list-style-type: none"> <li>• 33+2 reagent positions of 40 ml or 7 ml each (15 positions at room temperature and 20 cooled)</li> <li>• 14 positions for Standards and Controls</li> </ul> </li> </ul>
<b>SAMPLE LOADING</b>	<ul style="list-style-type: none"> <li>◆ 5 racks for continuous loading of samples split in:               <ul style="list-style-type: none"> <li>• 4 racks for 16 positions each.</li> <li>• 1 for 14 positions (STAT)</li> </ul> </li> <li>◆ Continuous loading</li> <li>◆ Positive barcode reader</li> </ul>

<b>REACTION VOLUME</b>	<ul style="list-style-type: none"> <li>◆ Minimum 300 µl</li> <li>◆ Maximum 670 µl</li> </ul>
<b>SAMPLING ARM</b>	<ul style="list-style-type: none"> <li>◆ A single mechanical arm provides all sampling operations with:               <ul style="list-style-type: none"> <li>• Capacitive liquid level sensing</li> <li>• Reagent pre-warming at 37° C</li> <li>• Automatic probe washing</li> </ul> </li> </ul>
<b>DILUTER</b>	<ul style="list-style-type: none"> <li>◆ Integrated module without syringe with the following specifications:               <ul style="list-style-type: none"> <li>• Sample volume: 3.0 µl ÷ 99 µl (1 µl incr.)</li> <li>• Reagent 1 volume: 3.0 µl ÷ 500 µl (1 µl incr.)</li> <li>• Reagent 2 volume: 3.0 µl ÷ 330 µl (1 µl incr.)</li> <li>• Reagent 3 volume: 3.0 µl ÷ 330 µl (1 µl incr.)</li> </ul> </li> </ul>
<b>PRECISION</b>	<ul style="list-style-type: none"> <li>◆ CV &lt; 1 % at 3 µl</li> </ul>

<b>READING SYSTEM</b>	<ul style="list-style-type: none"> <li>◆ Direct Photometry</li> </ul>
<b>OPTIC SYSTEM</b>	<ul style="list-style-type: none"> <li>◆ Photometer: interferential double channel filter wheel</li> <li>◆ Wavelength: 8 narrow band interferential filters from 340 nm to 620 nm</li> <li>◆ Light source: halogen lamp 6V/10 W</li> <li>◆ Linear range: 0.001 ÷ 2.500 Abs</li> <li>◆ Resolution: 0.0001 Abs</li> </ul>
<b>OPTICAL PATH</b>	<ul style="list-style-type: none"> <li>◆ 10 mm.</li> </ul>
<b>WASHING STATION</b>	<ul style="list-style-type: none"> <li>◆ Composed of five probes that empty, wash and dry the reaction cuvettes.</li> </ul>
<b>REACTION PLATE</b>	<ul style="list-style-type: none"> <li>◆ 60 semi-disposable cuvettes</li> <li>◆ Cuvettes Q.C. continuously computer controlled</li> <li>◆ Incubation temperature: 37°C</li> </ul>

#### 1.4.1. Computer & Software feature

<b>TYPE</b>	<ul style="list-style-type: none"> <li>◆ IBM Compatible</li> </ul>
<b>CPU</b>	<ul style="list-style-type: none"> <li>◆ Pentium III 500 MHz, 512 Kb Cache or plus</li> </ul>
<b>MEMORY</b>	<ul style="list-style-type: none"> <li>◆ RAM 128 Mb or plus</li> <li>◆ Hard Disk 10 Gb or plus</li> </ul>

<b>MONITOR</b>	<ul style="list-style-type: none"> <li>◆ Floppy Disk 3 1/2" 1.44 Mb</li> <li>◆ Colour SVGA 15" low radiation</li> <li>Resolution 800 x 600 pixels;</li> <li>max number of colors 65536 (16 bit)</li> </ul>
<b>PRINTER</b>	<ul style="list-style-type: none"> <li>◆ 80 Columns impact graphic (EPSON LX 300)</li> <li>◆ PS2</li> </ul>
<b>INTERFACE</b>	<ul style="list-style-type: none"> <li>◆ Two Bi-directional RS 232C serial ports and one parallel</li> </ul>
<b>SOFTWARE AVAILABLE LANGUAGES</b>	<ul style="list-style-type: none"> <li>◆ Multitasking WINDOWS 98 II E</li> <li>◆ Italian, English, Chinese, Czech. Software to be released soon in these languages: Russian, Portuguese, French, Polish. Upon request it is possible to release the software in other languages.</li> </ul>
<b>SETTINGS</b>	<ul style="list-style-type: none"> <li>◆ Disable all the energy saving options</li> <li>◆ Disable the screen saver</li> <li>◆ Select English "USA" as language</li> <li>◆ Select date and time in Regional setting</li> </ul>

**NOTE:** Even though the computers demonstrate the same technical and operative characteristics, some of these could have different hardware installed.

This could cause problems for the ILab 300 Plus software when running tests (A message appears indicating "Random" error or blocks the program).

Therefore, if the PC is bought separately/locally, it is highly recommended to test the system at your offices before proceeding with the installation at the final client's.

Consequently, AMS denies any responsibility for software problems that are due to buying the computer separately from the instrument.

### 1.4.2 Options

- ◆ **ISE MODULE**
- ◆ **POSITIVE BARCODE READER**

### 1.4.3 Dimensions, Weight & Environment

<b>DIMENSION</b>	<ul style="list-style-type: none"> <li>◆ Height: 42 cm</li> <li>◆ Depth: 65 cm</li> <li>◆ Length: 100 cm</li> </ul>
<b>WEIGHT</b>	<ul style="list-style-type: none"> <li>◆ 65 Kg</li> </ul>
<b>OPERATING ENVIRONMENT</b>	<ul style="list-style-type: none"> <li>◆ Temperature: 18°C ÷ 32°C.</li> <li>◆ Relative humidity: 20% ÷ 85%</li> </ul>

<b>1.4.4 Installation requirements</b>	
<b>POWER REQUIREMENTS</b>	<ul style="list-style-type: none"> <li>◆ Input Voltage 90 ÷ 250 Vac</li> <li>◆ Input Frequency: 47 ÷ 63 Hz</li> <li>◆ Power consumption: <ul style="list-style-type: none"> <li>◆ 300 W for the analytical unit</li> <li>◆ 400 W for the work station</li> </ul> </li> </ul>
<b>SAFETY REGULATIONS</b>	<ul style="list-style-type: none"> <li>◆ EN 61010-1:1993 +A2:1995 (in compliance with the main European safety directives 73/23/CEE e 3/68/CEE)</li> </ul>

<b>ELECTROMAGNETIC COMPATIBILITY</b>	<ul style="list-style-type: none"> <li>◆ EMC 89/336/CEE – 92/31/CEE Directives</li> <li>◆ EN 55011, Class B, Group 1</li> <li>◆ EN 50081-1:1992 EMC</li> <li>◆ EN 55022</li> <li>◆ ENV 50140 – ENV 50141</li> <li>◆ EN 60601-1-2</li> <li>◆ EN 61000-4</li> </ul>
--------------------------------------	---

**Warning:** A steady power supply (+ 10%) must be provided for the instrument.

If it is not, the manufacturer highly recommends the use of:

- ◆ **UPS** Uninterruptible Power Supply ( No-break module)
- ◆ **ELECTRONIC STABILISER**

## CHAPTER 2

---

### - GENERAL DESCRIPTION -

#### INDEX

2	GENERAL DESCRIPTION OF THE SYSTEM .....	3
2.1	ANALYTICAL CYCLE.....	4
2.1.1	REACTION PLATE.....	4
2.1.2	REACTION PLATE CYCLE.....	4
2.2	SAMPLING ARM - OPERATIONAL SEQUENCE.....	5
2.3	SAMPLING SYSTEM .....	6
2.4	WASHING STATION.....	7
2.4.1	WASHING STATION CYCLE .....	7
	APPENDIX.....WASHING STATION CYCLE	
2.5	ELECTRICAL SCHEMES DESCRIPTION.....	8
2.5.1	ANALYTICAL CONTROL BOARD (CPU) .....	8
2.5.2	STEPPER MOTORS DRIVER BOARD .....	18
2.5.2.1	TA 8435 H.....	21
2.5.3	SAMPLING INTERFACE ASSY .....	23
2.5.3.1	SAMPLING INTERFACE ASSY .....	24
2.5.4	LEVEL SENSOR BOARD 1-2 .....	25
2.5.5	REACTION TRAY INTERFACE .....	26
2.5.5.1	REACTION TRAY INTERFACE .....	27

2.5.6	REACTION CHAMBER MOTOR INTERFACE .....	28
2.5.7	PRE-AMPL/ADC .....	29
2.5.8	PHOTOMETER LAMP PWS .....	30
2.5.9	SAMPLE RACKS IDENTIFICATION BOARD .....	31
2.5.10	CONTROL PANEL .....	32
2.5.11	HYDRAULICS INTERFACE BOARD .....	33

## CHAPTER 2

---

### GENERAL DESCRIPTION

#### 2 GENERAL DESCRIPTION OF THE SYSTEM

*"ILab 300 Plus"* is a random access, counter-top, computer controlled clinical analysis instrument.

The system can perform 200 tests per hour and has a machine cycle of 18 sec. Its execution time ranges from a minimum of 48 sec to a maximum of 756 sec, depending on the analysis method chosen.

The first time the system is used for normal laboratory analyses, the operator must configure the system based on the specific needs of that laboratory; i.e.: the chemical parameters; the reagents racks and the normal, calibrated and control values must all be defined.

The daily routine analyses will be carried out according to patient sample arrival in a sequential and continuous, non-stop manner.

The work list is organized using a loading rack holding 16 patient samples plus a STAT rack for 14 patient samples. Rack loading is non-stop.

The racks can accommodate both test tubes and micro caps. The bar code for primary tubes is optional.

When the system is first turned on, the analytical unit, the computer and the color-meter lamp are supplied with low voltage power (1.2 volts). The sampling arm pre-heater is turned off while the reaction plate heater, the reagents refrigerating unit and the electronic components are turned on. The "Stand-by" light, placed on the front panel will flash until the reaction plate reaches a temperature of 36° C. At this point, the "Stand-by" light will stop flashing and will remain constantly lit.

In the case of system failure or malfunction, the "Ready" light, situated on the front panel of the instrument, will light up red.

In order to enter the main program, double click on the *"ILab 300 Plus"* icon on the computer desktop.

The main menu - "System Monitor" - will appear.

Whenever any system function is launched, the color-meter lamp and the sampling arm pre-heater will receive regular power.

## **2.1 ANALYTICAL CYCLE**

### **2.1.1 REACTION PLATE**

The reaction plate of the *"ILab 300 Plus"* system contains 60 washable and reusable, plastic cuvettes. The cuvettes can be removed individually.

Cuvette washing takes place under the reaction plate washing station. There are five positions which alternate in the washing and drying of the cuvettes.

The basic operating cycle of the reaction plate takes 18 seconds. Said cycle includes: optic reading of the cuvettes in incubation, aspiration and dispensation of the reagents and samples by the arm, repositioning of the plate and washing of the cuvettes.

The reactions take place at 37° C. this temperature is constantly maintained by the heater placed under the reaction plate.

### **2.1.2 REACTION PLATE CYCLE**

After reagents and samples have been placed in cuvette #1, the reaction plate will rotate counter-clockwise to position #31 so as to bring the first cuvette to be analyzed in front of the color-meter for reading with one or two wavelengths as required.

The plate will then, moving counter-clockwise, carry out all the readings of any other prepared cuvettes. After having effected all the readings, the plate will move clockwise to its initial position minus one cuvette, ready for a new dispensing.

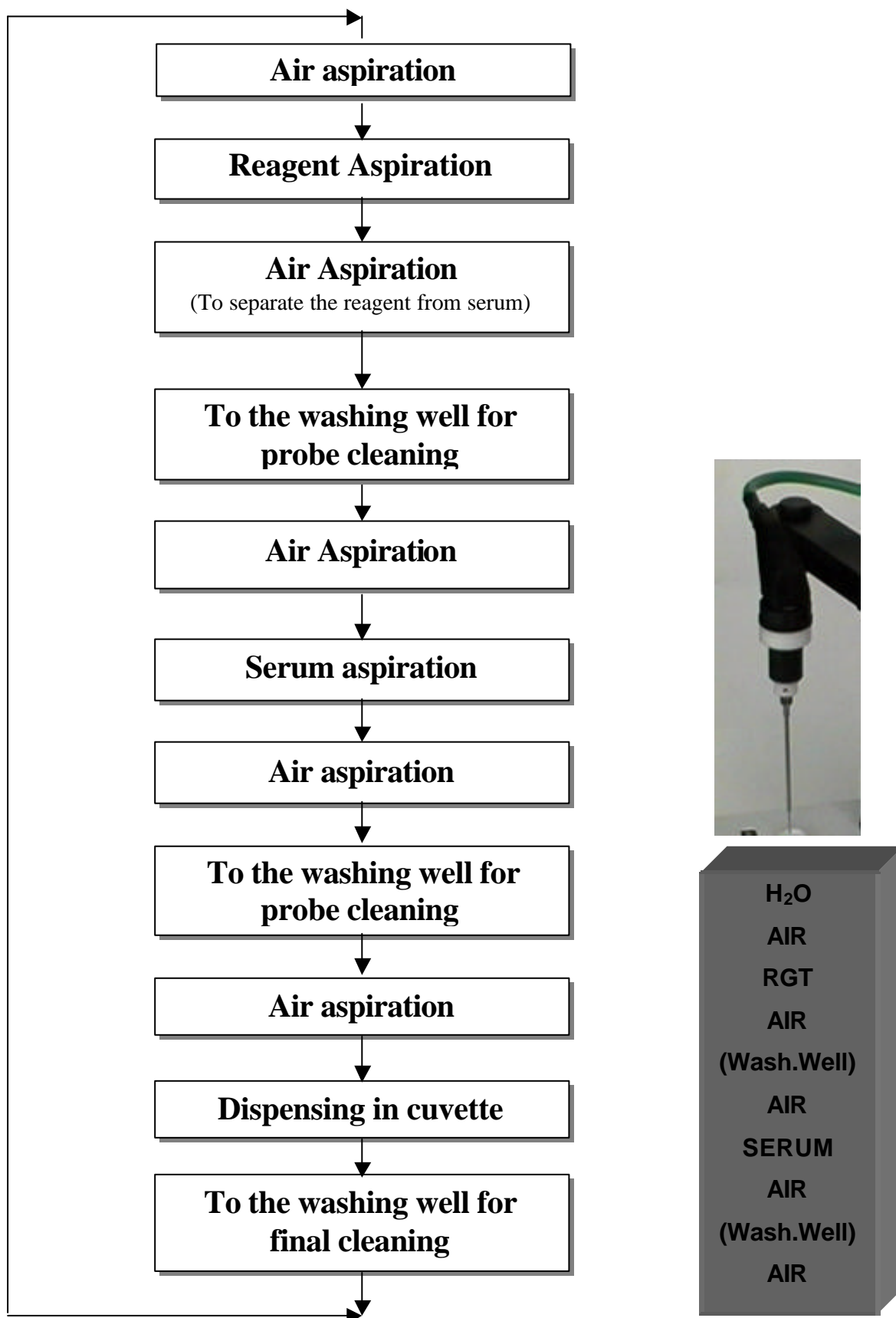
Therefore, the reaction cuvettes will move clockwise 1 - 2 - 3 - 4 for their dispensation and washing, and counter-clockwise 4 - 3 - 2 - 1 for their reading.



## 2.2 SAMPLING ARM - OPERATIONAL SEQUENCE

1. The sampling arm lifts up from the washing well and carries out a washing cycle.
2. The arm moves toward the specific reagent container, while the diluter aspirates an air bubble to separate the rinse column from the reagent.
3. The arm descends into the reagent to the level indicated by the sensor and aspirates the required volume of reagent.  
If the method requires a line prime volume, an extra amount of reagent (not used in the analysis) will be aspirated before the quantity of reagent necessary for the analysis along with a smaller air bubble for separation.
4. While the diluter aspirates a second air bubble, the arm rises and then lowers into the wash well so that it can be washed externally, to minimize cross contamination.
5. The arm moves to the specified sample and aspirates a third air bubble.
6. When the level sensor senses the liquid, the arm stops and aspirates the sample.
7. The arm rises once again, while the diluter aspirates a fourth air bubble to prevent sample loss.
8. At this point, the arm returns to the wash well in order to clean the outside of the probe and it aspirates a fifth air bubble.
9. The arm moves to the reaction plate and dispenses reagent and sample in the reaction cuvette for incubation and reading.
10. The arm returns to the wash well and carries out a probe wash cycle.

### 2.3 SAMPLING SYSTEM



## 2.4 WASHING STATION

The reaction plate washing station is made up of a series of five probes situated to one side of the reaction plate. The probes are connected to the valve and pump system for emptying, washing and drying operations (please see the hydraulic diagram).

### 2.4.1 WASHING STATION CYCLE

The washing station carries out its operation by alternating upward and downward movement. In its downward movement phase the probes are suitably guided to carry out the following operations.

- The first probe, using its central cannula, removes the read sample while the external cannula dispenses, shower fashion, the wash solution (this operation is repeated twice).
- The second probe operates exactly like the first but uses distilled water to wash and repeats the operation three times.
- The third probe dispenses distilled water into the cuvette so that an optics check can be performed (if the results are negative the cuvette is discarded).
- The fourth probe aspirates and removes the previously dispensed distilled water.
- The fifth probe is a drying pad which removes any residual moisture from the sides of the cuvette.

All these operations are part of the routine operation of the instrument. Every reaction cuvette is washed at the end of each round of analysis.

The reusability of each reaction cuvette is always tested before the next round of analysis.

## **2.5 ELECTRICAL SCHEMES DESCRIPTION**

---

### **2.5.1 ANALYTICAL CONTROL BOARD (CPU)**

#### **INTRODUCTION**

The ANALYTICAL CONTROL BOARD (C. P. U.), integrated in the ILab 300 Plus system is the heart of the low level, real time, processing management.

#### **GENERAL CHARACTERISTICS**

The board has been designed to be a generic control of input/output, governed by an external C.P.U. or by an internal PC board based platform, in light of possible different future instrument configurations.

#### **CIRCUIT DESCRIPTION**

The following description refers to the [SE]00195-01 Rev. C. electrical wiring diagram

#### **POWER SUPPLY**

- The main power supply VM (24V), filtered by C10 and C12, is taken up via connector J6 and used to generate the +5 volts; moreover, it is brought back on connectors J9 and J10 to 25 pins
- The 5V is obtained through U4 (LM2575), the Schottky D5 diode, the L1 inductance and condensers C14 and C15: the entire circuit makes up a step down converter - from 24V to 5V - with 1 Ampere of current.
- Said power supply is brought back on connector J6 forecasting the possibility that it could come directly from an external power source (in this case the step-down components are not mounted, with the exclusion of the C14 and C15 condensers).
- A second power supply, 12.5V, filtered via C9 and C11 is taken up via connector J6, in addition to being brought back on connectors J9 and J10 to 25 pins and on connectors J13 and J14 to 10 pins and utilized to generate the 12V for the two flash memories inside the two U1 and U2 microprocessors and the AVCC.

AVCC (5V) used for the operational U23 and U24 and also for the two microprocessors' internal A/D converters.

## MICRO-CONTROLLERS AND BUFFER

In order to quickly and most efficiently manage all the devices, there are two independent micro-controllers, U1 and U2. Each sends input and output signals to the relative connectors.

All signals, which begin with the letter 'P' refer to the Micro 'P' (U1), while those which begin with the letter 'R' refer to the Micro 'R' (U2).

The two Micros are HITACHI H8/3048, configured in operative mode 7, single chip.

Following are some basic performance features:

- elevated speed in carrying out instruction (2 clock cycles at 16 MHz equal to 125 ns to carry out an addition or subtraction operation between 2 registers);
- 128 KBYTES flash memory;
- 4 KBYTES RAM;
- 70 I/O signals, almost all having more than one function, plus 8 input only.
- All signals are buffered in input or in output:
- the buffer for the T.T.L. signals is 74HC244;
- the driver for the open-collector signals is ULN2803;
- the operational amplifier LM324 is used for the analogic signals.

The following two tables list the 78 Micro 'P' and Micro 'R' signals. For each the following are specified:

- the name of the signal on the pin of the Micro (e.g.: P30)
- the type: I for output. 0 for input;
- the buffer used for the input or the output (e.g.: U7);
- the connector onto which the signal is sent (e.g.: J5);
- the name of the buffer input signal (if in input) or at buffer output (if in output);
- a brief description of functioning (e.g.: Motor A enable).

**TAB.1 Micro ‘P’ Signals**

Micro Signal	I/O	Buffer	Conn.	Signal Connector	Function
P10	I/O			R14	Synchronization 1 micro ‘R’
P11	I/O			R15	Synchronization 2 micro ‘R’
P12	I/O			R16	Synchronization 3 micro ‘R’
P13	I	U21	J9/11	PINP	Generic
P14	I/O			R10	Synchronization 1 micro ‘P’
P15	I/O			R11	Synchronization 2 micro ‘P’
P16	I/O			R12	Synchronization 3 micro ‘P’
P17	O	U17	J9/14	PDAR1	Pump Diluter
P20	I	U15	J7/21	PINP1	Sample rack 1 in place
P21	I	U15	J7/23	PINP2	Sample rack 2 in place
P22	I	U15	J7/25	PINP3	Sample rack 3 in place
P23	I	U15	J17/1	PINP4	Sample rack 4 in place
P24	I	U15	J17/3	PINP5	Sample rack 5 in place
P25	I	U15	J17/5	PINP6	Sample door switch
P26	I	U15	J17/7	PINP7	Stat door switch
P27	I	U15	J17/9	PINP8	Sample button
P30	O	U7	J7/2	PENA	Motor A enable
P31	O	U7	J7/4	PENB	Motor B enable
P32	O	U7	J7/6	PENC	Motor C enable
P33	O	U7	J7/8	PEND	Motor D enable
P34	O	U7	J7/10	PDIRA	Motor A c.w./c.c.w.
P35	O	U7	J7/12	PDIRB	Motor B c.w./c.c.w.
P36	O	U7	J7/14	PDIRC	Motor C c.w./c.c.w.
P37	O	U7	J7/16	PDIRD	Motor D c.w./c.c.w.
P40	O	U9	J7/1	PUT1	ST-BY lamp button (Control Panel)
P41	O	U9	J7/3	PUT2	SAMPLE lamp button (Control Panel)
P42	O	U9	J7/5	PUT3	STAT lamp button (Control Panel)
P43	O	U9	J7/7	PUT4	READY lamp button (Control Panel)
P44	O	U9	J7/9	PUT5	Not used
P45	O	U9	J7/11	PUT6	Not used
P46	O	U9	J7/13	PUT7	Not used
P47	O	U9	J7/15	PUT8	Not used
P50	O	U8	J7/18	PPDA	Motor A power down function
P51	O	U8	J7/20	PPDB	Motor B power down function
P52	O	U8	J7/22	PPDC	Motor C power down function
P53	O	U8	J7/24	PPDD	Motor D power down function
P60	O	U14	J7/17	PUT9	Reserved 1
P61	O	U14	J7/19	PUT10	Reserved 2
P62	O	U21	J9/1	PG1	ADC1 gain A
P63	O	U21	J9/5	PG2	ADC2 gain B
P64	O	U21	J9/23	PSCLK	ADC1-ADC2 serial clock
P65	O	U21	J9/10	PTTL	Not used

P70	I	U23A	J17/10	PAN1	External temperature N.T.C. sense
P71	I	U23B	J9/16	PAN2	Preheater current sense
P72	I	U23C	J9/25	PAN3	Preheater N.T.C. sense
P73	I	U23D	J13/9	PAN4	Reagent N.T.C. sense
P74	I	U16	J9/3	PFOTO1	Optosensor 1 signal
P75	I	U16	J9/7	PFOTO2	Optosensor 2 signal
P76	I	U16	J9/9	PFOTO3	Optosensor 3 signal
P77	I	U16	J9/17	PFOTO4	Optosensor 4 signal
P80	I	U14	J17/11	PINP9	Stat button
P81	I	U14	J17/13	PINP10	Ready button
P82	I	U14	J17/15	PINP11	St-by button
P83	I	U14	J17/16	PINP12	Liquid 1 alarm
P84	I	U14	J17/12	PINP13	Liquid 3 alarm
P90	O	U25	J1/1	P90OUT	Auxiliary TX
P91	O	U25	J11/3	X91	RS232 TX
P92	I	U25	J1/3	P92IN	Auxiliary RX
P93	I	U25		X93	RS232 RX
P94	I	U16	J9/8	PRDY2	ADC2 data ready
P95	I	U16	J9/15	PRDY1	ADC1 data ready
PA0	I	U6	J13/8	POSR1	Optosensor reagent 1 signal
PA1	I	U6	J13/7	POSR2	Optosensor reagent 2 signal
PA2	O	U8	J7/26	PCLKA	Motor A clock
PA3	I	U6	J13/6	POSR3	Optosensor reagent 3 signal
PA4	O	U8	J17/2	PCLKB	Motor B clock
PA5	I	U6	J13/5	POSR4	Optosensor reagent 4 signal
PA6	O	U21	J9/21	PPWM	Preheater P.W.M.
PA7	O	U17	J13/4	PDAR2	Air pump
PB0	O	U8	J17/6	PCLKD	Motor D clock
PB1	I	U21	J9/6	PLEAK	Sampling leak assy
PB2	O	U8	J17/4	PCLKC	Motor C clock
PB3	I	U16	J9/24	PADC1	ADC1 digital output
PB4	I	U21	J9/4	PSLI	Level sensor
PB5	I	U16	J9/2	PADC2	ADC2 digital output
PB6	I				MICRO P address
PB7	I	U14	J17/14	PINP14	Liquid 2 alarm

**TAB.2 Micro ‘R’ Signals**

Micro Signal	I/O	Buffer	Conn.	Signal Connector	Function
R10	I/O			P14	Synchronization 1 micro ‘P’
R11	I/O			P15	Synchronization 2 micro ‘P’
R12	I/O			P16	Synchronization 3 micro ‘P’
R13	I	U22	J10/11	RINP	Power led opto (3)
R14	I/O			P10	Synchronization 1 micro ‘R’
R15	I/O			P11	Synchronization 2 micro ‘R’
R16	I/O			P12	Synchronization 3 micro ‘R’
R17	O	U17	J10/14	RDAR1	Diluter pump
R20	I	U19	J8/21	RINP1	Leak
R21	I	U19	J8/23	RINP2	Interlock
R22	I	U19	J8/25	RINP3	Not used
R23	I	U19	J18/1	RINP4	Not used
R24	I	U19	J18/3	RINP5	Not used
R25	I	U19	J18/5	RINP6	Not used
R26	I	U19	J18/7	RINP7	Not used
R27	I	U19	J18/9	RINP8	Not used
R30	O	U13	J8/2	RENA	Motor A enable
R31	O	U13	J8/4	RENB	Motor B enable
R32	O	U13	J8/6	RENC	Motor C enable
R33	O	U13	J8/8	REND	Motor D enable
R34	O	U13	J8/10	RDIRA	Motor A c.w./c.c.w.
R35	O	U13	J8/12	RDIRB	Motor B c.w./c.c.w.
R36	O	U13	J8/14	RDIRC	Motor C c.w./c.c.w.
R37	O	U13	J8/16	RDIRD	Motor D c.w./c.c.w.
R40	O	U11	J8/1	RUT1	EV1 valve
R41	O	U11	J8/3	RUT2	EV2 valve
R42	O	U11	J8/5	RUT3	EV3 valve
R43	O	U11	J8/7	RUT4	EV5 valve
R44	O	U11	J8/9	RUT5	Cover lock
R45	O	U11	J8/11	RUT6	MP2 pump
R46	O	U11	J8/13	RUT7	MP3 pump
R47	O	U11	J8/15	RUT8	MP4 pump
R50	O	U12	J8/18	RPDA	Motor A power down function
R51	O	U12	J8/20	RPDB	Motor B power down function
R52	O	U12	J8/22	RPDC	Motor C power down function
R53	O	U12	J8/24	RPDD	Motor D power down function
R60	O	U20	J8/17	RUT9	Reserved 1
R61	O	U20	J8/19	RUT10	Reserved 2
R62	O	U22	J10/1	RG1	ADC1 gain A
R63	O	U22	J10/5	RG2	ADC2 gain B
R64	O	U22	J10/23	RSCLK	ADC1-ADC2 serial Clock
R65	O	U22	J10/10	RTTL	Lamp off



R70	I	U24A	J18/10	RAN1	External temperature N.T.C. sense
R71	I	U24B	J10/16	RAN2	Plateheater current sense
R72	I	U24C	J10/25	RAN3	Plate N.T.C. sense
R73	I	U24D	J14/9	RAN4	Not used
R74	I	U18	J10/3	RFOTO1	Optosensor 1 signal
R75	I	U18	J10/7	RFOTO2	Optosensor 2 signal
R76	I	U18	J10/9	RFOTO3	Optosensor 3 signal
R77	I	U18	J10/17	RFOTO4	Optosensor 4 signal
R80	I	U20	J18/11	RINP9	Not used
R81	I	U20	J18/13	RINP10	Peristaltic pump A sense
R82	I	U20	J18/15	RINP11	Peristaltic pump B sense
R83	I	U20	J18/16	RINP12	Not used
R84	I	U20	J18/14	RINP13	Not used
R90	O	U25	J3/1	R90OUT	Auxiliary TX
R91	O	U25	J11/3	X91	RS232 TX
R92	I	U25	J3/3	R92IN	Auxiliary RX
R93	I	U25		X93	RS232 RX
R94	I	U18	J10/8	RRDY2	ADC2 data ready
R95	I	U18	J10/15	RRDY1	ADC1 data ready
RA0	I	U6	J14/8	ROSR1	Not used
RA1	I	U6	J14/7	ROSR2	Not used
RA2	O	U12	J8/26	RCLKA	Motor A clock
RA3	I	U6	J14/6	ROSR3	Not used
RA4	O	U12	J18/2	RCLKB	Motor B clock
RA5	I	U6	J14/5	ROSR4	Not used
RA6	O	U22	J10/21	RPWM	Plateheater P.W.M.
RA7	O	U17	J14/4	RDAR2	Not used
RB0	O	U12	J18/6	RCLKD	Motor D clock
RB1	I	U22	J10/6	RLEAK	Leak plate assy
RB2	O	U12	J18/4	RCLKC	Motor C clock
RB3	I	U18	J10/24	RADC1	ADC1 digital output
RB4	I	U22	J10/4	RSLI	Not used
RB5	I	U18	J10/2	RADC2	ADC2 digital output
RB6	I				Micro 'R' address
RB7	I	U20	J18/12	RINP14	Not used

## PROGRAMMABLE LOGIC

It is made up of U10 (PAL 16V8) and has the following functions:

- reset generation for the two Micros:

Each Micro is reset: if RESA is present, generated by U5 (DS1233-10), when downloading begins (RESB signal) or for the entire length of program downloading of the other Micro (BOOT and SELP signals, the latter is provided by the first serial connector J2).

The relative equations in PALASM format are:

$$/RESP = /RESA + RESB + BOOT * /SELP$$

$$/RESR = /RESA + RESB + BOOT * SELP$$

- Generation of X91 associated with one of two TX lines (Micro 'P' P91 and Micro 'R' R91) to be directed toward the PC both in the operative phase and in the BOOT phase.

R91 is chosen if SELP is low (downloading phase Micro 'R') or if R10 is low (operative phase in which Micro 'R', after ascertaining that the channel is free, wants to transmit): in all other cases X91 brings back P91.

The corresponding equation in PALASM format is:

$$/X91 = /SELP * /R91 + SELP * (/R91 * /R10 + /P91 * R10)$$

- generation of the signals for programming the Micro 'R' flash

Activate PMD2E in the BOOT phase if SELP is high or in the operative phase if Micro 'P' wants to write in flash (low P66);

Activate PVPPE in the BOOT phase if SELP is high;

Activate RMD2E in the BOOT phase if SELP is low or in the operative phase if Micro 'R' wants to write in flash (low R66);

Activate RVPPE in the BOOT phase if SELP is low.

The relative equations in PALASM format are:

$$\text{PMD2E} = \text{/P66} + \text{BOOT} * \text{SELP}$$

$$\text{PVPPE} = \text{BOOT} * \text{SELP}$$

$$\text{RMD2E} = \text{/R66} + \text{BOOT} * \text{/SELP}$$

$$\text{RVPPE} = \text{BOOT} * \text{/SELP}$$

## FLASH PROGRAMMING VOLTAGE

The PMD2E, PVPPE, RMD2E, RVPPE signals pilot 4 U17 Darlingtion sections, whose outputs (PEMD2, PEVPP, REMD2, REVPP signals) activate, respectively, the PNP Q2, Q3, Q1 and Q4 transistors which bring back the 12V respectively on pin MD2 and pin VPP of Micro 'P' and Micro 'R'.

## SERIAL INTERFACES

The first serial of each Micro (signals P90 and P92 for Micro 'P'; R90 and R92 for Micro 'R') is dedicated to communicating with its own auxiliary device (respectively J1 and J3).

For the reception channel of the second serial:

- Micro 'P' is linked to serial 1 of the PC (P93 RX).
- Micro 'R' can be:
  - linked to serial 2 of the PC (R93 RX2 coming from J5) shifting the tin soldered jumper from pins 2-3 of J12 to pins 2-1
  - linked in parallel with those of Micro 'P' and therefore linked to series 1 of the PC (R93 RX :tin soldered jumper on pins 2-1 f J12).

For the transmission channel of the second serial:

- Micro 'R' is linked to serial 2 of the PC (R91 TX2).
- Micro 'P' can:
  - be linked to serial 1 of the PC (P91 TX) shifting the tin soldered jumper from pins 2-3 of J11 to pins 2-1
  - share the transmission to the PC with the other Micro (tin soldered jumper on pins 2-1 of J11). In this case the arbitration between the two Tx signals is carried out in the BOOT phase via the SELP signal and in the operative phase via the R10 and P1 signals which the two Micros exchange and which respectively receive on signals P14 and R14 everything governed by programmable logic).

The two Micros carry out program downloading via J2: in addition to the TX and RX signals also the following signals are utilized for this function:

- SELP which goes directly on the programmable logic and which if low, signals that it wants to load the program on Micro 'R'.
- RTS and DTR which after the RS232-TTL conversion become BOOT and RESB and are brought back on the programmable logic for reset generation and for enabling the programming voltage for Micros 'P' and 'R'.

The same program is loaded on the two Micros.

In the operative phase, the format of the communication messages with the PC contains an address 'number':

- '1' for those messages which can be processed only by Micro 'P', recognized because the PB6 input is high;
- '2' for those messages which can be elaborated only by Micro 'R', recognized because the RB6 input is low.

The four TTL-RS232 conversions:

- P90-P90OUT: transmission toward Micro 'P's auxiliary device;
- R90-R90OUT: transmission toward Micro 'R's auxiliary device;
- R91-TX2: transmission toward Micro 'R's PC,
- TX5V TX: transmission toward either Micro 'P's or Micro 'R's PC,

and the four RS232-TTL conversions:

- P92IN-P92: reception from the auxiliary device for Micro 'P';
  - R92IN-R92: reception from the auxiliary device for Micro 'R';
  - RX2-R93: reception from the PC for Micro 'R';
  - RX-P93-R93: reception from the PC for Micros 'P' and 'R'
- are carried out through the U25 and U26 (MAX202) devices.

**CLOCK**

The system clock is supplied by quartz at 16MHz X1 connected to the Micro 'P' internal oscillator pins EXTAL and XTAL: output TTL of this oscillator (U1/67) is brought back on input XTAL of the Micro 'R' oscillator.

**LED**

There are two LEDs

- LD1: green LED for VM presence
- LD2: green LED for +12.5V presence.

**BUZZER**

The acoustic signals are channeled to the LS1 buzzer and can be activated via the PDAR2 and/or RDAR2 signals of the two Micros.

## 2.5.2 STEPPER MOTORS DRIVER BOARD

### INTRODUCTION

The "STEPPER MOTORS DRIVER BOARD" is the interface between the control board and the stepper motors.

### GENERAL CHARACTERISTICS

This board has been designed to be used in various systems or sub-groups which require mechanical movement .

It can manage up to four stepper motors and has an open collector (Darlington) interface for 16 signals, whose output can be sent to various devices including electro-valves, continuous current motors.

### CIRCUIT DESCRIPTION

**The following description refers to the [SE]00188-01 Rev.B. electrical wiring diagram**

### POWER SUPPLY

The main power supply, VM (24 V), applied via J3 and indicated by the lit green D22 LED, is used for the stepper driver power portion (TA8435H).

The secondary power supply, 12.5 Volts, also comes through connector J3, is filtered by condenser C13, monitored via the green D21 LED and used to generate the 5 Volts necessary for the stepper driver logic section and for the pull-up resistors.

The 5 volts are obtained from the 12.5 Volts through U5 (LM7805) and the C14 filter condenser.

## STEPPER MOTORS DRIVER

The four drivers are represented by the U1, U2, U3 and U4 components.

Following is a description of the functioning of the first (U1) of the four identical circuits.

Later on in this manual a brief description of the functioning of the NMB SDI-C403 device or TOSHIBA's equivalent TA8435H will be provided.

The driver under examination is made up of:

**A logic section** for decoding signals arriving from the control board, via the J4 and J7 connectors, which uses the following inputs:

DIRA (pin 5): determines the direction of the rotation of the motor;

RES (pin 2): reset;

CLKA (pin 6): controls phase excitation frequency: the pilot mode depends on the M1A (pin 8) and M2A (pin 9) inputs, set via the SW1DIP SWITCH according to the following table:

M2A	M1A	Function
0(on)	0(on)	1 full step
0(on)	1(off)	1/2 step
1(off)	0(on)	1/4 step (default setting)
1(off)	1(off)	1/8 step

Please note that the other clock input (pin 7) is connected to the +5 volts.

Use of the 1/4 and 1/2 steps, with respect to the full step, make it possible to lessen vibration and improve positioning accuracy.

**A power section**, (made up of 2 H-shaped bridges), that generates the voltage to be applied to the two stepper motor windings (AM62, AM62N, BM62, and BM62N respectively on pins 23, 20, 19, and 16) on connector J1.

The latter are supplied when the RES signal is high and the ENA signal (coming from the control board, via connector J4) is at a low logic level.

The state of the ENA signal level, and therefore that of the relative Stepper Motor driver, is indicated by the yellow D17 LED.

The D1, D3, D5 and D7 Schottky diodes are to protect the internal power circuits.

**A power control and monitoring section** which, via the C5 (3.3nF) condenser generates a circa 60 kHz internal oscillation frequency used to control the Pulse Width Modulation (PWM).

Said control is affected by:

the voltage levels, produced by the potential drop on the R1 and R2 power resistors, present on pins 21 and 18;

the PDA (pin 10) signal, arriving from the control board (via connector J4), which conditions the internal switching logic intervention threshold; to 0.8 volts if it is kept at a high logic level and 0.5 volts if at a low logic level.

## **OPEN COLLECTOR DRIVERS**

Two integrated circuits are present on the board, U6 and U7 (ULN2803), and constitute an open collector interface for 16 signals:

10 of which come from the control board (UT1...UT10) via J4 with their outputs connected to J2;

6 gathered via J5 and placed in output on J6.

These outputs are used to pilot various types of power devices: continuous current motors, solenoids, electrovalves, etc.



### 2.5.2.1 TA 8435 H

#### INTRODUCTION

TA 8435 H controls the PWM for setting the current of the Stepper Motor windings to a constant value. The device is an integrated circuit for micro-step piloting, used in order to have an efficient and low vibration control over the Stepper Motors.

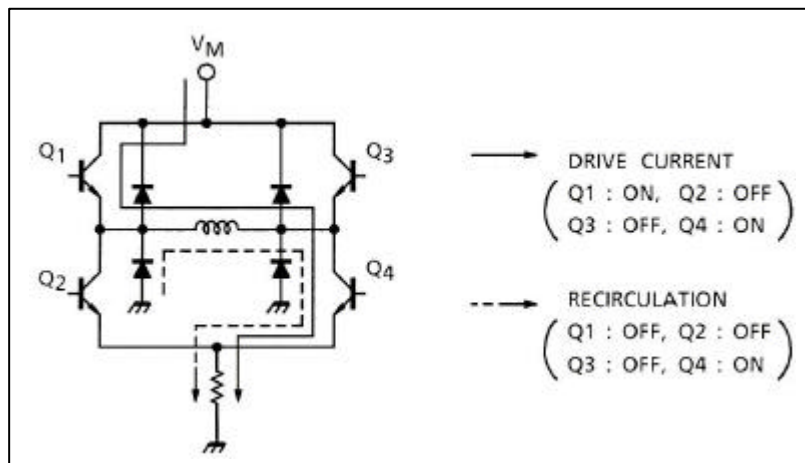
#### MICRO-STEP PILOTING

TA 8435 H controls a Stepper Motor in micro-steps; with a maximum resolution of 1/8 of a step. The current levels of the A and B phases are regulated in micro-steps within the integrated circuit in such a manner that the size and rotation angle of the vector between the micro-steps remains constant. When the clock signals are placed in input (CK1, CK2), the Stepper Motor will rotate in micro-steps.

#### CONTROL OF PULSE WIDTH MODULATION AND OUTPUT CURRENT

Output current pathway flow diagram (PWM control)

TA 8435 H controls PWM, placing the upper transistor in an on/off switching mode. In this case current flow is as illustrated below:



Control of the output current via the detection resistor of the output and input current REF-IN: the motor current (maximum current for micro-step piloting  $I_{\emptyset}$ ), is regulated as shown in the following equation, using the REF-IN in-put and the detection resistor of the external current RNF:

$$I_{\emptyset} = V_{REF} / R_{NF}$$

Where:


for REF-IN = HIGH  $\Rightarrow V_{REF} = 0.8 \text{ V}$

for REF-IN = LOW  $\Rightarrow V_{REF} = 0.5 \text{ V}$

## CONTROL LOGIC

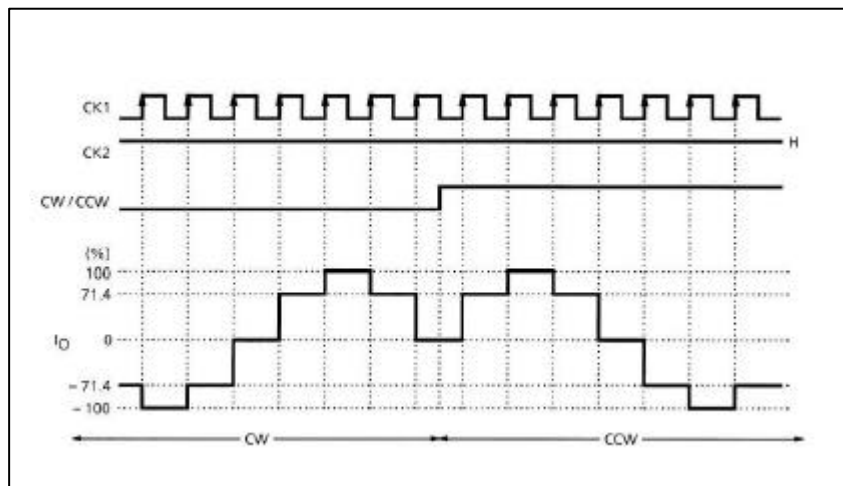Clock input for rotation direction control

To switch the rotation direction from forward to backward there are two clock inputs: CK1 and CK2, and input CW/CCW and the first of the two control methods are shown in the diagram below:



## 1/2 STEP PILOTING

The selection between: full step, 1/2 step, 1/4 step and 1/8 step, is carried out with binary coding on the two M1 and M2 inputs.

### 2.5.3 SAMPLING INTERFACE ASSY

The "Sampling Interface Assy" board controls the arm, the diluter, the pre-heater and wash well liquid leakage.

The Home Sensors (vertical, external arm, internal arm and diluter) are connected to the board via the J9, J10, J11 and J13 connectors.

The Vertical Home, External Arm Home and Internal Arm Home LEDs are placed in series with their power absorption control via the ALFD3 signal, connected to J2 Pin 11 that reports information on the Analytical Control Board (CPU).

Information from the leakage sensor enters on the J12 connector, is processed by U1A, exits on J2 Pin 6 and arrives through J10 to the Analytical Control Board for management.

The N.T.C. resistor detects the temperature of the pre-heater and is connected to J7 on Pins 5 and on Pins 3 and 4 (ground). These Pins are connected to the U2B amplifier, whose output via J2 Pin 25 carries information to the CPU board (Analytical Control Board) on input AN2 (P72) of the Micro R U2.

The pre-heater heating element is connected via J7 Pins 6 and 7 and is powered by Q1 which is controlled by P.W.M. generated by the Micro R of the CPU. Furthermore, it can be excluded by switching the LS1 relay through Q6, conditioned by the OFFPR signal generated by the CPU.

### 2.5.3.1 SAMPLING INTERFACE ASSY

#### ❖ INPUT:

- Vertical Home (J9)
- External Arm Home (J10)
- Internal Arm Home (J11)
- Diluter Home (J13)
- NTC + (J7 Pin 5)
- LEAK (J12)
- PDIL (coming from the CPU) (J2 Pin 14)

#### ❖ OUTPUT:

- Vertical Arm (J3)
- Internal Arm (J4)
- External Arm (J5)
- Diluter (J6)
- Diluter Pump (J8)
- Pre-heater (J7)
- LEAK (J2)

#### **2.5.4 LEVEL SENSOR BOARD 1-2**

The integrated U4 circuit generates a trigger signal for U3 where the sampling probe is connected on Pin 6 via C2.

The capacity variation, produced by the contact between the sampling probe and a liquid, determines a voltage variation on Pin 3 of U3, which ranges from a value close to zero to one which tends towards + 5 V.

This signal is then further processed by U2 and Q1 and is placed in output from the board on J2 Pin 1.

### 2.5.5 REACTION TRAY INTERFACE

This board controls various functions:

J2, J3, J4 and J5 are the input connectors for the opto sensors of the filters wheel, of the washing station, of the reaction plate and of the cuvette holders. This information is transferred to the Analytical Control Board (CPU) using the OS2, OS3, OS4 and OSB signals via J1.

The ALFD3 signal controls power supply to the series of the three emitting diodes contained in the opto sensors connected to J2, J3 and J4. ALFD1 functions in a like manner relative to the opto sensor connected to J5.

Liquid leakage signal enters on J7 and is processed by U1A, whose output (LEAKP) returns via J1 on the CPU.

Plate temperature is detected by the N.T.C. sensor connected to J8 and sent to the Micro P.

The cuvette holder solenoid is turned on by activating the BCU signal produced by a 'Darlington' controlled by the Micro P.

The PWMP signal, that controls the reaction plate heater, is also produced by the Micro P. Its power absorption is controlled via the P71 signal.

Furthermore, the board distributes, toward other modules, input and/or output signals managed by the CPU.

24V and 12.5V arrive to the board. The 5V (U5) used for the amplifiers and the negative voltage (U4) that is sent to both the Sample and Reference PRE-AMPL/ADC are produced here.

### 2.5.5.1 REACTION TRAY INTERFACE

#### ❖ **INPUT:**

- Home Washing Station (J3)
- Home Filters Wheel (J2)
- Home Reaction Tray (J4)
- Home Cuvette Holder (J5)
- Reaction Tray Leak (J7)
- Reaction Tray Temperature Sensor (J8)
- (PWMP (PWM for the Heater) (J1)
- BCU (Activation of the Cuvette Holder Solenoid) (J1)
- OFFPL (Activation of the Lamp, Feedthrough (J1)

#### ❖ **OUTPUT:**

- Cuvette Solenoid Holder (J6)
- Reaction Tray Heater (J10)
- ALFD3 (three light emitting diodes in series) (J1)
- OS3 (Washing Station Home Sensor) (J1)
- OS2 (Filters Wheel Home Sensor) (J1)
- OS4 (RTC Tray Home Sensor) (J1)
- OSB (Cuvette Holder Home Sensor) (J1)
- ALFD1 (Cuvette Holder Sensor Diode Control) (J1)
- P72 (Reaction Plate Temperature) (J1)
- LEAKP (Reaction Plate Leak) (J1)
- P71 (Reaction Tray Heater Power Control) (J1)
- V- (Voltage –12V for PRE-AMPL/ADC) (J1)
- PWRES (Heater Power) (J1)
- OFFPL (Lamp Activation) (J9)

### **2.5.6 REACTION CHAMBER MOTOR INTERFACE**

This board distributes the motors signals - those of the: filters wheel (J2), washing station (J3), and reaction plate (J4), part of the Reaction Group Assy.

The connector J1 interfaces the Stepper Motors Driver Board (30-00188-01) which is managed by the Micro P of the CPU via connector J1. Signal distribution for the three motors is carried out via the three above cited connectors.



### 2.5.7 PRE-AMPL/ADC

The FD1 photodiode detects light energy arriving from the optics. The amplifier (U2) amplifies this signal and sends it to the A/D serial converter ADS 1250.

Both U1 and U4 produce the +5V; the former for the analogic segment and the latter for the digital segment. U5 generates the –5V and U6 the voltage reference VREF. Y1 is the oscillator connected to the U3 clock.

The G0 and G1 signals determine the U3 gain (1, 2, 4 or 8).

The Board interfaces with the Reaction Tray Interface Board (00193-01) via the JP1 connector that sends the signals to the Micro P.

### 2.5.8 PHOTOMETER LAMP PWS

This board powers the 6V 10W halogen lamp.

The main device of this circuit is a type LT 1083CP U1 integrated stabilizer which, under normal operative conditions, is regulated, via the 1 Kohm P1 potentiometer, to 6V for the lamp.

When the system is in Stand-by, the lamp is not turned off. Instead, the voltage is lowered from 6V to 1.2V by the Micro P which acts on the Q1 transistor via the OFFPL signal. This is done in order to not over-stress the lamp by turning it on and off frequently.

### **2.5.9 SAMPLE RACKS IDENTIFICATION BOARD**

This board checks that the samples racks are inserted correctly.

The sole power supply used in the board is Vcc (5V), derived from Vm (24V) that arrives to the board via the J1 connector.

Racks sensing is carried out via OPB848 type opto sensors U1, U2, U3, U4, and U5. These sensors' outputs arrive to connector J7 of the Control Panel (30-00342-00) via J1. They exit through the same J1 and go towards J2 of the Stepper Motors Driver Board (30-00188-01) and from here, via J4 and J7 enter on J7 and J17 of the Analytical Control Board (CPU) (30-00195-01) for the management functions performed by the Micro P (U1).

### **2.5.10 CONTROL PANEL**

This board controls the levels of the Rinse, Water and Cleaning Solution bottles, the closure of the sample rack and reagent rack doors, the buttons – indicators placed on the front panel of the instrument – and the acoustic alarm buzzer.

The St-by and Ready buttons are only lit indicators. The STAT button, instead, actually turns off the acoustic alarms.

Signals regarding the presence of the five sample racks pass through this board.

### **2.5.11 HYDRAULICS INTERFACE BOARD**

This board is used to distribute the following DARLINGTON commands:

EV1, EV3, EV4, EV5, EV6, MP2, MP3, MP4, INTERLOOK, or INTERLOCK.

In addition, this board is used to distribute the peristaltic pump's motor signals.

U2 generates +5 V , U1 is the operator that keeps the leak (D9 on = leak alarm) under control.

The parts connected to J9 and J10 are designed to monitor the fluctuating liquids of the peristaltic pump's two canals.

The microswitch signal transits to J8 which identifies if the cover is open (IL=1) or closed (IL=0).

## CHAPTER 3

---

### - INSTALLATION –

#### INDEX

3	INSTALLATION.....	2
3.1	UNPACKING .....	4
3.2	INSTALLATION	
3.2.1	Installation site characteristics .....	4
3.3	ELECTRICAL POWER SUPPLY .....	5
3.4	COMPUTER - INSTRUMENT CONNECTION.....	5
3.5	LABLES AND SYMBOLS .....	6
3.6	INSTALLATION ACTIVITY.....	7
3.6.1	Software Installation Procedure.....	9
3.7	FIRMWARE AND MASTER HARDWARE UPGRADING.....	11
3.8	ANALYTICAL PERFORMANCE.....	13

## CHAPTER 03

---

### - INSTALLATION -

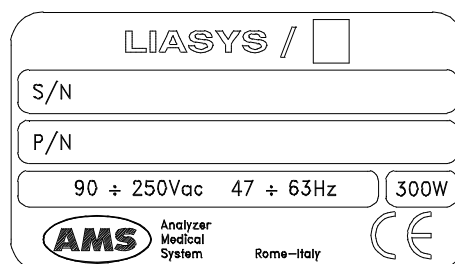
### INDEX

3.1 UNPACKING .....	2
3.2 INSTALLATION.....	4
3.2.1 INSTALLATION SITE SPECIFICATIONS .....	4
3.2.2 ELECTRIC CURRENT REQUIREMENTS.....	4
3.2.3 CONNECTION OF THE ACCESSORIES .....	5
3.2.3.1 POWER SUPPLY .....	5
3.2.3.2 COMPUTER - INSTRUMENT CONNECTION .....	5
3.2.4 ATTENTION .....	6
3.2.5 SYMBOLS.....	6
3.2.6 REGULATORY COMPLIANCE.....	6
3.3 INSTALLATION ACTIVITY .....	7
3.4 SOFTWARE INSTALLATION PROCEDURE .....	9
3.5 FIRMWARE UPGRADING .....	11

### 3.1 UNPACKING


The **ILAB 300 PLUS** is packed and delivered in two separate wooden crates: one contains the analyzer itself and the other the computer, along with its accessories. In the event that the order not include the PC component, packing and delivery will involve one wooden crate plus a corrugated cardboard box. The packing has been expressly studied and designed to insure maximum protection of the contents during shipping and handling. It is therefore extremely important that the crate(s)/box be carefully examined upon delivery in order to ascertain their integrity. Special attention should be dedicated to examining the color of the “Shock Watch” glued to the crates, which must show the color ‘white’. A ‘red’ “Shock Watch” indicates that the crate(s) have experienced some sort of ‘shock’ during handling, transport and/or delivery. This fact must be noted by the courier on the delivery note, as must any and all visible external damage (for example: holes, dents, rips or tears, water marks, etc.) evident at the moment of delivery. This will simplify matters in the event of any future claims for damages.

Upon arrival of the crate(s)/box, take out the delivery note and make sure that all the items on the packing list are included in the crates and are undamaged. Make sure the series number on the delivery note/packing list corresponds to that impressed on the plate on the right side of the instrument.



Open the crate(s)/box from the top and very carefully take out:

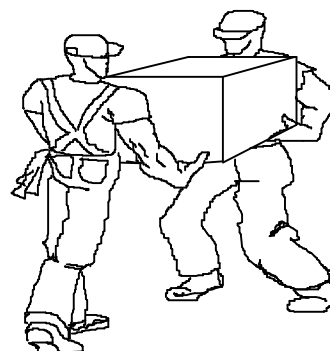- the instrument;
- the computer and accessories.



#### **MAKE SURE THAT THE UNPACKING IS CARRIED OUT BY TWO PEOPLE.**

Do not discard the delivery crate(s)/box or the packing material until the correct functioning of the instrument has been ascertained.

Remove all the items from the crate(s)/box very carefully.

Remove the adhesive tape from the cover of the samples and reagents housing, from the front panels and from the samples and reagents racks.



Before connecting the "ILab 300 Plus", remove the protective packing material placed under the sampling arm and under the wash station group.

**Warning:** in the event that it is necessary to repack any or all of the delivered item(s), the following procedures must be carefully followed:

- Reposition the protective packing material under the sampling arm and under the wash station group.
- Tape down (using masking tape if possible) the cover of the samples and reagents housing, the front panels, and the samples and reagents racks.
- Remove the probe from the sampling arm and place it inside a cuvette. Then cap the cuvette and tape the cap down.
- Be very careful to not bend the wash station cannulas when repositioning the protective packing material.
- Fill the empty spaces around the accessories packed in the crate using "*pluriballs*" or other suitable packing material.

## 3.2 INSTALLATION

The ILAB 300 PLUS must only be installed by a qualified technician who has been authorised and trained to do so. During its installation the system will be checked once again to ensure correct functioning. The persons who are required to operate the ILAB 300 PLUS system must have received the adequate training. This should also include the "know-how" of the normal maintenance for the instrument. A description of the maintenance will be found in Chapter 7 of this manual.

ILAB 300 PLUS is a complex system, and it is therefore extremely important that it is correctly installed in order to fully guarantee fine performance. If the installation and use directions, given in this manual, are not correctly followed and/or security indications are not respected, AMS cannot guarantee correct functioning of the instrument. Apart from this, the security of the operator could be placed at risk.

### 3.2.1 INSTALLATION SITE SPECIFICATIONS

Ascertain that the ILAB 300 PLUS system is not exposed to direct sunlight, draughts, dust or strong magnetic fields. In addition, please take note of the following conditions required for the location of the installation:

<b>USE</b>	In covered and dry place
<b>DEGREE OF POLLUTION</b>	2
<b>INSULATION CLASS</b>	I
<b>INSTALLATION CATEGORY</b>	II
<b>TEMPERATURE</b>	between 18°-32°C
<b>HUMIDITY</b>	20% ÷ 85%
<b>ALTITUDE</b>	Max 3000 m
<b>LOCATION</b>	Shelf or table with a minimum surface of 110x70 cm stable and free of vibration
<b>VENTILATION</b>	Leave a minimum distance of 10 cm around the instrument to permit air circulation . Make sure that the fan (situated under the reagent compartment) is not blocked by any object

### 3.2.2 ELECTRIC CURRENT REQUIREMENTS

The power voltages to which the instrument is adapted are indicated on the left-hand side (see fig. 1). It must be plugged into a plug of the correct voltage.

<b>VOLTAGE</b>	90 ÷ 250 Vac 47/63 Hz ± 10%
<b>FUSES</b>	6.3 Amp/T - 6.3 x 32

**NOTE:** IT IS ADVISABLE TO MAINTAIN THE MAXIMUM STABILITY OF THE ELECTRICAL CURRENT IN THE LABORATORY. WHERE THIS IS NOT POSSIBLE OR ASCERTAINABLE, USE OF THE FOLLOWING SUPPLEMENTARY DEVICES IS RECOMMENDED:

### **ELECTRONIC STABILIZER**

Used to stabilise the electric voltage in the laboratory. Any stabiliser with a power potential greater than 0.5 KW, currently available on the market, can be used.

### **NO-BREAK MODULE UPS - (Uninterrupted Power Supply)**

This module provides two important functions:

- stabilises the main-line power
- supplies current to the instrument in case of a main-line power failure.

## **3.2.3 CONNECTION OF THE ACCESSORIES**

### **3.2.3.1 POWER SUPPLY**



Fig. 1 Plug (use the feeder cable supplied with the instrument).

The sticker indicates the power supply voltage and the values of the fuses.

### **3.2.3.2 COMPUTER - INSTRUMENT CONNECTION**

The instrument and the Personal Computer are connected by one serial RS232 standard cable (P/N. 23935005501), which provides the hardware support for the communication

### 3.2.4 ATTENTION

The following label is found at the rear of the instrument.



**NOTE:** THE REAR PANELS OF THE INSTRUMENT MUST NEVER BE OPENED WITHOUT HAVING FIRST SWITCHED THE INSTRUMENT OFF AND DISCONNECTED THE ELECTRICITY CABLE.

THE MAINTENANCE AND CLEANING PROCEDURES FOUND IN CHAPTER 07 OF THIS MANUAL MUST BE RESPECTED AT ALL TIMES. REMEMBER TO FOLLOW THE DECONTAMINATION PROCEDURE IN CASE OF INSTRUMENT REMOVAL (SEE CHAPTER 07) .

### 3.2.5 SYMBOLS



**ATTENTION: READ THE INSTRUCTIONS IN THE USER MANUAL**



**TERMINAL OF TOTAL MASS PROTECTION (CONDUCTOR)**

### 3.2.6 REGULATORY COMPLIANCE

The ILAB 300 PLUS instrument complies with:

**Safety :**

- EN 61010 - 1 : 1993 + A2 : 1995 therefore meeting the essential requirements of L.V.D. 73/23/EEC and 93/68/EEC

**EMC :**

- EN 61326-1 therefore meeting the essential requirements of Electromagnetic Compatibility Directive 89/336/EEC - 92/31/EEC

### 3.3 INSTALLATION ACTIVITY

Customer \_\_\_\_\_

Instrument S/N \_\_\_\_\_

Service Engineer \_\_\_\_\_

Date \_\_\_\_\_

During the **preinstallation** visit the laboratory area should be inspected in the light of the following requirement:

- Temperature of the laboratory between 18 and 32 °C? ☐ Yes ☐ No
- Main Power supply between 90 e 250 Vac? ☐ Yes ☐ No
- Is it necessary an electronic stabilizer? ☐ Yes ☐ No
- Worktable:
  - ☐ Minimum size 110 x 70 cm plus PC and Printer ☐ Yes ☐ No
  - ☐ Presence of vibration? ☐ Yes ☐ No
  - ☐ Presence of good ventilation (minimum 10 cm) ? ☐ Yes ☐ No
  - ☐ No direct exposure to sunshine ? ☐ Yes ☐ No

Installation Activity	
⇒ Unscrew the screws fixing the crate located in the lower sides.	<input type="checkbox"/>
⇒ Raise the wood crate and take out all the packaging material.	<input type="checkbox"/>
⇒ Check the integrity of all the parts.	<input type="checkbox"/>
⇒ Check the Packing List	<input type="checkbox"/>
⇒ Place the instrument on the workbench and level it by adjusting the front feet.	<input type="checkbox"/>
⇒ Check the ground connections	<input type="checkbox"/>
⇒ Raise manually the washing station	<input type="checkbox"/>
⇒ Remove the top cover and check that the cuvettes are well fitted in their slot.	<input type="checkbox"/>
⇒ Connect the ILab 300 Plus to the P.C (see 3.4)	<input type="checkbox"/>

<b>Activity</b>	
⇒ Install the ILab 300 Plus program provided with the CD Rom. Perform all the P.C settings. (see Installation Procedure point 1 to 11 sect. 3.6.1)	<input type="checkbox"/>
⇒ Keep the arm completely up and turn the instrument , the P.C and the Printer on.	<input type="checkbox"/>
⇒ After removing the protection tube, install the probe. Verify that the probe is perfectly straight. (see Installation Procedure point 12 to 14 sect. 3.6.1)	<input type="checkbox"/>
⇒ Reconstitute the Rinse and cleaning solutions by following the instructions sheet in the packing.	<input type="checkbox"/>
⇒ Place the three bottles in position and fill them with the relevant liquid (H <sub>2</sub> O, Rinse, Cleaning). Close with the covers and insert the tubes.	<input type="checkbox"/>
⇒ Place the two waste tank (one is the biological waste) * Not supplied. Insert the waste tubings. The “RED” connector is the biological outlet.	<input type="checkbox"/>
⇒ Check the tubings for any squashing	<input type="checkbox"/>
⇒ Boot ILab 300 Plus program and run the “Diagnostic”; select “Diluter” and click on the <b>Probe Wash Test</b> command to execute the priming. Perform the “Washing Well Level”. (see Installation Procedure point 15 to 16 sect. 3.6.1)	<input type="checkbox"/>
⇒ Open the Global Setting menu in the “Diagnostic” program. Check and in case adjust the probe in all the positions. (Samples / Reagents / Standards / Controls / Diluent / Washing Pot) (see Installation Procedure point 17 to 18 sect. 3.6.1)	<input type="checkbox"/>
⇒ Check the optic system. (see Sect.6.2.3)	<input type="checkbox"/>
⇒ Check the temperature ( see Sect. 6.6.1)	<input type="checkbox"/>
⇒ Check the Backlash ( see Sect. 6.1.6.5)	<input type="checkbox"/>
⇒ Run two cuvettes washing cycles (see Installation Procedure point 19 sect. 3.6.1)	<input type="checkbox"/>
⇒ Perform the WBL for all the cuvettes and checks the results. (see Installation Procedure point 19 sect. 3.6.1)	<input type="checkbox"/>
⇒ Set the methods, the control serums, standards and reagents.	<input type="checkbox"/>
⇒ Calibrate the methods for testing the ILab 300 Plus. Usually one cinetic and one substrate. Remember to run the Reagent Blank first if necessary	<input type="checkbox"/>

⇒ Run a precision test on the test previously calibrated, with the control serum. (See the analytical performance 3.8)	<input type="checkbox"/>
⇒ Check the results.	<input type="checkbox"/>

### 3.4 SOFTWARE INSTALLATION PROCEDURE

#### Installation Procedure Instrument New

1. Turn on the P.C.
2. Confirm all the steps in order to complete the installation of Windows 98. Type the number of Product Key .(It is glued beside the P.C)
3. Confirm date and time.
4. In order to avoid that, Windows 98 requires a password at the start up, select "Start" "Settings" "Control Panel" "Network". At the "Primary Network Logon" function select "Windows Logon".
5. Select "Start" "Settings" "Control Panel" "International Settings" and select "English(United States)". Click on "Time" and select "HH:mm:ss". Click on Date and select "dd-MMM-yy". Reboot the computer.
6. When required confirm the path to install the driver of the monitor.
7. Set the "Display Resolution" at 800 x 600 Pixels
8. Disable all the options in the "Power Management" properties. Select "Standby" > Never "Turn off the Monitor > Never, "Turn Off Hard Disks > Never". In "Power Schemes" select " Home/Office Desk".
9. Disable the "Screensaver" by selecting "None"
10. Disable the "Background" by selecting "None"
11. In the "BIOS" select APM and disable "Power Management" (Suggested)
12. In the "BIOS" select Peripheral and disable both the Sound and the LAN. (Suggested)
13. Set the maximum number of colors at 65.536 (16 bit). To change the settings, select "Display Properties" > "Settings" > "Colors" → "High Colors".
14. **Select "Start" "Settings" "Control Panel" "Device Manager" "Ports". Select "Communication Port COM1" "Properties" "Port Setting" "Advanced". Adjust both the "buffers" at the minimum level.**
15. Install the printer's driver Epson LX 300 + (the CD is in the packaging of the printer). Select "Settings" "Printer". With the right button of the mouse click on the printer's icon. Select "Properties". "Paper". Set on "A4" "Portrait" and "Tractor" in Paper Source. Select "Graphics" and set the "120 x 144" resolution. Select "Device Options" "Printer Quality" " High Speed ON".

#### New Instruments.

16. Install the **ILab 300 Plus software**. Insert the auto executable CD Rom. In the mask "ILab 300 Plus rel. 1.0.8 Multilanguage Setup" click on "Next". There are three choices: Complete, Upgrade and Service. Select "Complete". Follow the steps

- as required. The system will install the MasterSoftware and the MasterHardware. Select YES when is asked if maintain the existing files. Chose the language.
17. Run again the installation program. Select "Service". Install the program "TEST ILAB" and the program "WINZIP".
  18. Start the ILab 300 Plus program. **P.S : Turn on the instrument first.** The instrument will take a minimum of 10 minutes to reach the correct temperature. *Meantime it is not possible to run analysis but it is possible to enter the different parts of the software.*
  19. Select the diagnostic program.
  20. Reset the arm by selecting the Global Setting folder. Install the probe.
  21. Select the Diluter folder and run several time the "Probe Wash Test" as far as the rinse tube is filled with liquid.
  22. Select the "Configuration" folder. Type the password "1234" and run the "Washing Well Level" test to optimize the level of liquid in the washing pot.
  23. Home all the parts of the instrument by selecting the folders in the diagnostic program. Check all the sensors looking at the color at the bottom of the screen.
  24. Select the "Global Setting" folder and check the probe in all the positions.
  25. Exit the "Diagnostic" program. Select "Start" (green arrow) and run a washing cycle and a "Water Blank Level". (WBL). Check the result.

**NOTE :** when the setting is completed save the new Analyzer.MDB file. To access the Windows program, logout the ILab 300 Plus software by typing "Ctrl Alt 5" (Ctrl Alt Shift – with Software Revision up to 1.0.7.) Respond YES when the system prompts " System Shutdown is Turned Off". Select Shutdown then OK. Make a copy of the Analyzer. MDB file in the Analyzer folder either in a different folder or on a floppy disk.

## Installation Procedure Upgrade.

**Important: With the software version 1.x.x , the system requires at least the "Firmware 5.7".**

1. Exit the ILab 300 Plus software by pressing first "Ctrl Alt Shift –" (Ctrl Alt 5 if already Version 1.0.8). Select "Shutdown".
2. Insert the CD Rom ILab 300 Plus Ver. 1.X.X. Wait the program starts. In the mask "ILab 300 Plus rel. 1.0.8 Multilanguage Setup" click on "Next". There are three choices: Complete, Upgrade and Service. Select "Upgrade". Follow the steps as required. The system will install the new MasterSoftware and the new MasterHardware. Select YES when is asked if maintain the existing files. Choose the language.
3. If necessary run again the installation program. Select "Service". Install the new program "TEST ILAB".
4. Enter the "ILab 300 Plus" software. Check the "Parameters - Option". (Password: SH)

**Important: The file "Analyzer.mdb for software version 1.0.8 is NOT COMPATIBLE with the software version 1.0.7. and previous.**



**NOTE :** when the setting is completed save the new Analyzer.MDB file.

## 3.5 FIRMWARE UPGRADING

### SOFTWARE FIRMWARE UPGRADING PROCEDURE

#### A) Firmware upgrading procedure

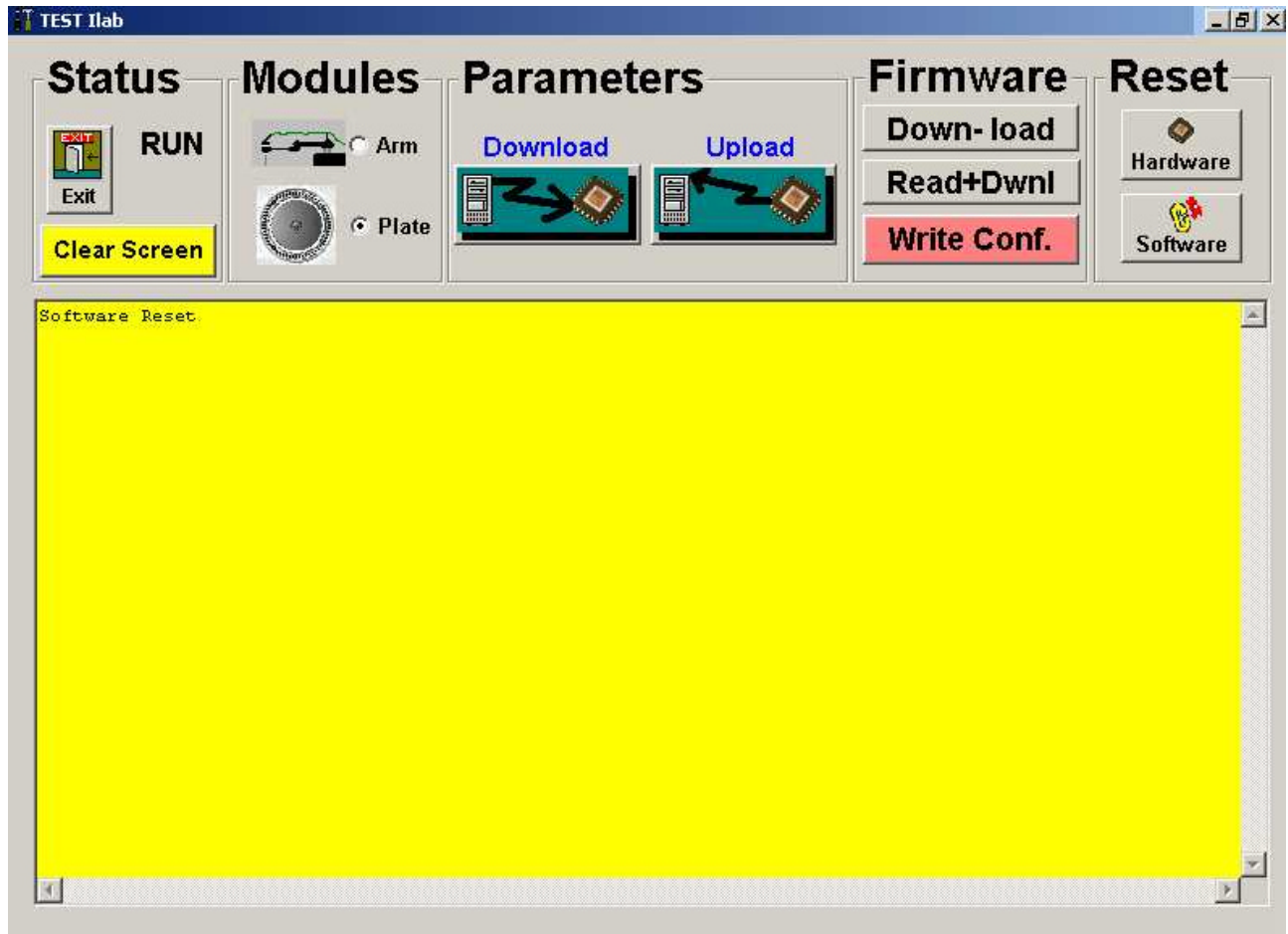
To verify the revision of the firmware, point out the mouse on the figure of the chip at the bottom of the screen in the ILab 300 Plus program. The version of the firmware will appear.

**Open the black rear cover and remove the plastic protection to have access at the Analytical Control Board.**

1. Copy the file ALYSSWXX.A20 in the TestILab folder. (The file ALYSSWXX.A20 is the firmware).
2. Start the "Test ILab" program by double clicking on Test.exe.
3. With the selection on "Plate" (Micro 1), click on the "**Read Down Load**" key.
4. A window named "Write Log" with When the downloading is completed press "Write Confirm" key will appears. Click "OK"
5. It will appears the window "HFD Hitachi Flash Download". Click on "File" key and select "\*.a20" on the "list files of type" by clicking on it.
1. Click on "ALYSSWXX.A20" file (or a:\ ALYSSWXX.A20 file if it is located on a floppy disk) in the folder "TestILab" in order to select it on the "File name" section and click "OK"
2. Click on "Activate" key and wait until the Micro 1 programming operation finishes
3. When finished , a window named "Report" will appear, click "OK"
4. Exit by clicking on "Quit" key and click "OK" to confirm.
5. Click on the "**Write Confirm**" red button to transfer the setting data on the Micro 1. When the "Write Confirm" red button disappears, it means that the data transfer has been completed.
6. Short circuit pins 1 and 2 of the J4 connector on the Analytical Control Board. Select "Arm" (Micro 2) and click on the "**Read Down Load**" key.

7. A window named "Write Log" with When the downloading is completed press "Write Confirm" key will appears. Click "OK"
8. Click on "File" key and select "\*.a20" on the "list files of type" by clicking on it.
9. Click on "ALYSSWXX.A20" file (or a:\ALYSSWXX.A20 file if it is located on a floppy disk) in the folder "TestILab" in order to select it on the "File name" section and click "OK"
10. Click on "Activate" key and wait until the Micro 2 programming operation finishes.
11. When finished, a window named "Report" will appear, click "OK"
12. Exit by clicking on "Quit" key and click "OK" to confirm.
13. Click on the "**Write Confirm**" red button to transfer the setting data on the Micro 2. When the "Write Confirm" red button disappears, this means that the data transfer has been completed.
14. Remove the short circuit from the pins 1 and 2 on the J4 connector on the Analytical Control Board P/N 00195-01
15. Click on the "HW Reset" key and verify a sound signal from the Analytical Control Board to confirm that the firmware upgrading procedure has been successfully completed.
16. Click on the "Exit" key in order to close the application.

## TestILAB new Graphic Interface Vers. 1.0



The functionality of the keys for the updating of the firmware is not changed. Under the box called **FIRMWARE** there are the “Download” “Read + Download” and “Write Conf.” Keys. These buttons maintain the same operative sequence as in the former version.

Following are the improvements:

**Parameters**, introduces two functions:

- **Upload** which enable to read the system setting parameters from the Micro Arm and Micro Plate. These parameters are saved in two different files in the “Analyzer” folder. The two files are: Micro\_arm.ams and Micro\_plate.ams
- **Download** which make possible both to open the files Micro\_Arm.ams and Micro\_Plate.ams and write the system parameters previously saved into the micros This is useful when a replacement of the “Analytical Control Board” is necessary. The files with the system setting parameters can be reloaded. The files Micro\_Arm.ams and Micro\_Plate.ams are automatically update every time the Diagnostic is selected or during the Diagnostic shutdown.

**Modules:** Must be selected before the Upload and/or Download,

ex. If you want to save the Arm parameters, select "ARM" and press UPLOAD.

ex. If you want to save the Plate parameters, select "PLATE" and press UPLOAD.


**From file to Micro** select Arm and press Download (BOX PARAMETERS) then select Plate and pressDownload (BOX PARAMETERS).

At the Micro Arm are linked the arm position parameter, the pre-heater and the diluter.

At the Micro Arm are linked the temperature of the reaction disk and the setting of the temperature adjustment's probe.




COLLEGAMENTI											
J2-1	J1-14		J3-1	J1-14		J4-1	J1-14		J5-1	J1-13	+24V
J2-2	J1-14		J3-2	J1-14		J4-2	J1-13		J5-2	J1-13	+24V
J2-3	J1-11		J3-3	J1-12		J4-3	J1-12		J5-3	J1-10	+12.5V
J2-6	J1-1		J3-6	J1-9		J4-6	J1-7		J5-6	J1-8	GND
J2-5	J1-2		J3-5	J1-9		J4-5	J1-7		J5-5	J1-8	GND
									J6-1	J1-13	
									J6-2	J1-13	
									J6-3	J1-15	
									J6-4	J1-15	

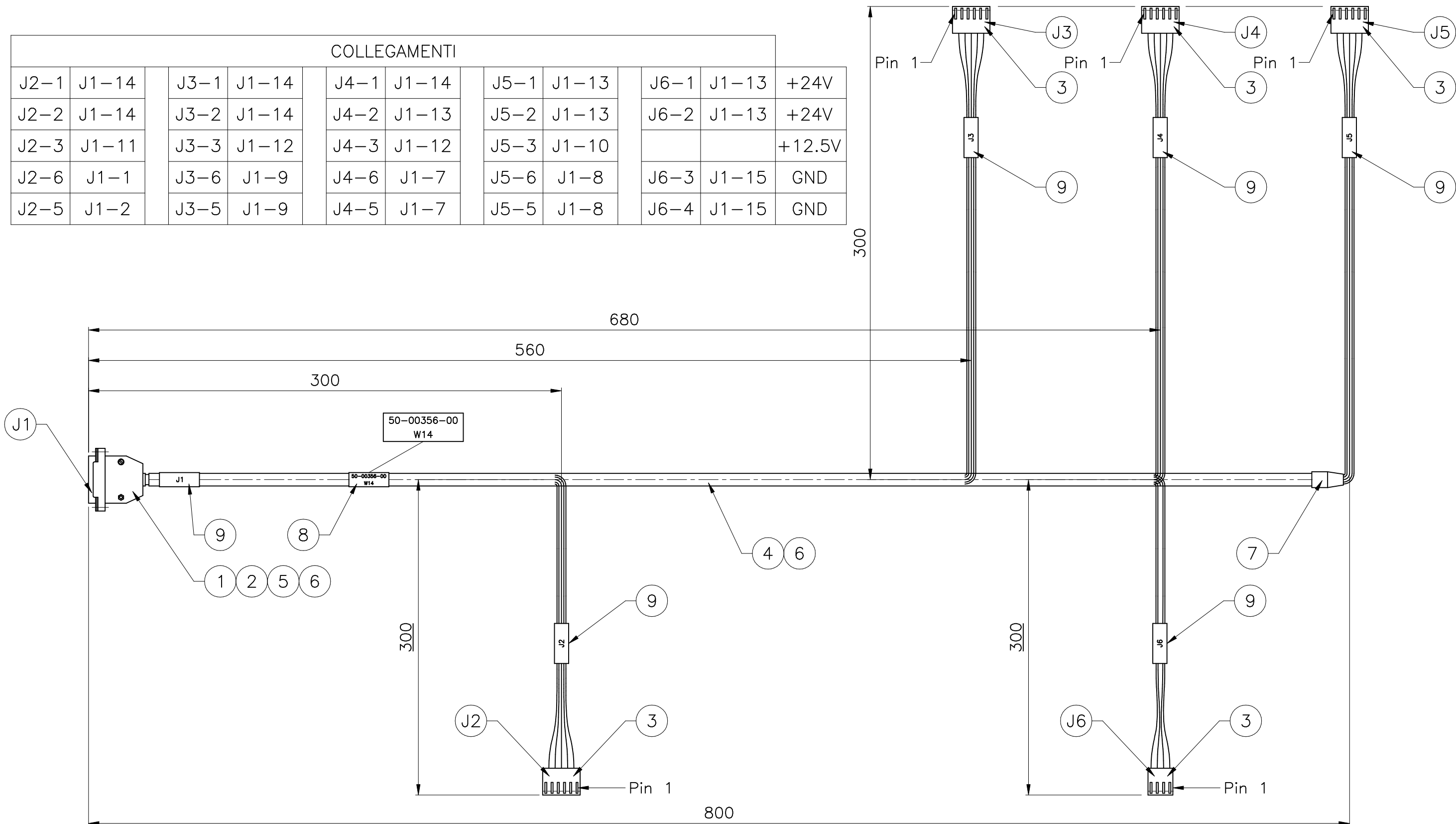


NOTA : -Le targhette Rif.7 e 8 devono essere bianche,  
 e stampate con inchiostro nero e caratteri alti 2mm.

0	Emission				
Rev.	Description	Checked	Date	Approved	Date
<b>AMS</b> Analyzer Medical System Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
		Drawn	G.Pucci	Date	18.07.2001
Description		Scale	1:2.5	Sheet	1 of 1
CAVO ALIMENTAZIONE		Car.	[50]	Rev.	0
		Approved	A.Gagliarducci	Date	18.07.2001

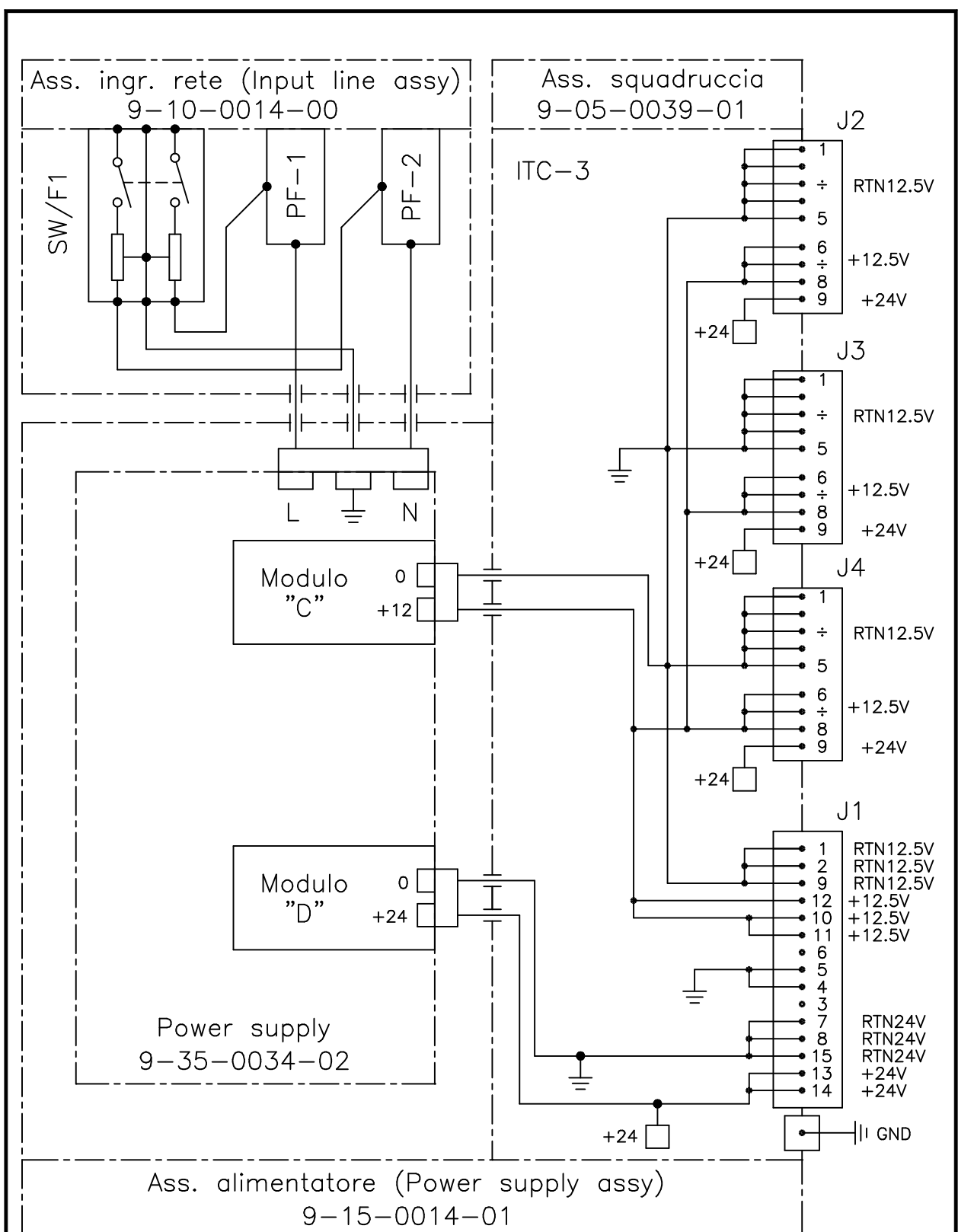

0	Emission				
Rev.	Description	Checked	Date	Approved	Date
<b>AMS</b> Analyzer Medical System Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
Description		Scale	Sheet	Checked	Date
SCHEMA ELETTRICO ALIMEN.			1 of 1	R.Cornacchia	08.01.2002
SCHEMATIC POWER SUPPLY		Drawing	Rev.	Approved	Date
		9-SE-15-0014-01	0	A.Gagliarducci	08.01.2002




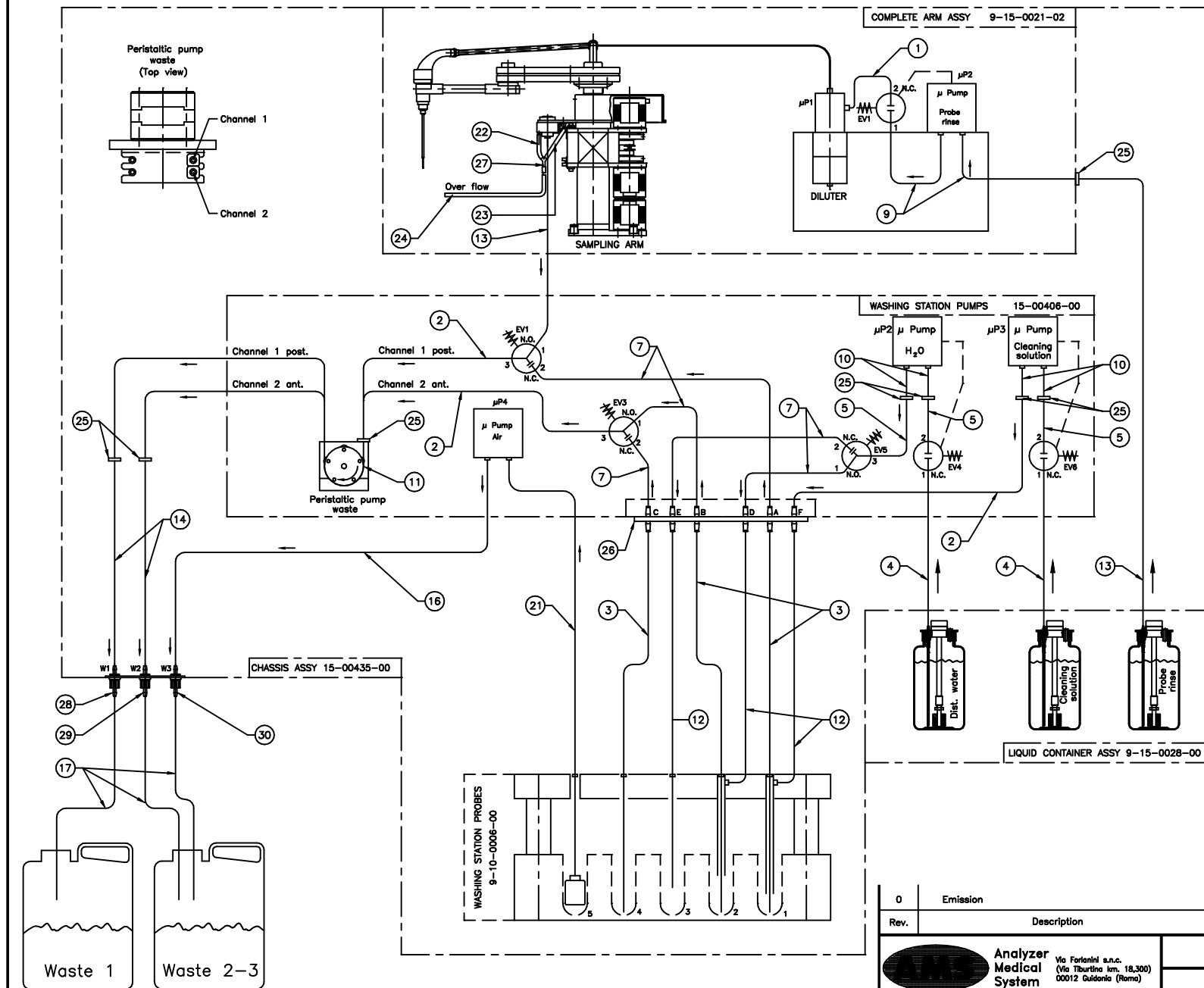



* 30	Kit fitting 9-65-0033-03 (Blue)	1	/
* 29	Kit fitting 9-65-0033-02 (Yellow)	1	/
* 28	Kit fitting 9-65-0033-00 (Red)	1	/
* 27	"Y" fitting 9-35-0053-00	2	/
* 26	Manifold assy 05-00405-00	1	/
* 25	Fitting 9-35-0051-00	9	/
* 24	PVC tubing øi3-øe5	1	350
* 23	PVC tubing øi3-øe5	1	160
* 22	PVC tubing øi3-øe5	1	30
* 21	PVC tubing øi4-øe6	1	850
* 20	PVC tubing øi3-øe5	1	270
* 19	PVC tubing øi3-øe5	1	230
* 18	PVC tubing øi3-øe5	1	70
* 17	Silicon tubing øi4-øe6	3	1200
* 16	Silicon tubing øi4-øe6	1	400
15			
14	Silicon tubing øi2-øe4	2	300
13	Silicon tubing øi2-øe4	2	1100
12	Silicon tubing øi2-øe4	3	600
11	Peristaltic pump tubing (Pharmed)	2	350
10	Silicon tubing øi2-øe4	4	20
9	Silicon tubing øi2-øe4	2	100
8			
7	Silicon tubing øi1.6-øe3.2	5	70
6			
5	Silicon tubing øi1.6-øe3.2	3	95
4	Silicon tubing øi1.6-øe3.2	2	1000
3	Silicon tubing øi1.6-øe3.2	3	600
2	Silicon tubing øi1.6-øe3.2	3	130
1	Teflon tubing AWG16 - VW1	1	160
REF.	DESCRIPTION	Q.ty	L=mm

\* -Not included in the tubes kit



0	Emission				
Rev.	Description	Checked	Date	Approved	Date
<b>Analyzer Medical System</b> Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
Description HYDRAULIC DIAGRAM (LIASYS NEW)		Scale 1:1.5	Sheet 1 of 1	Drawn G.Pucci	Date 12.03.2002
Car. [SI]	Drawing 00457-00	Rev. 0	Checked A.Gagliarducci	Approved A.Gagliarducci	Date 12.03.2002



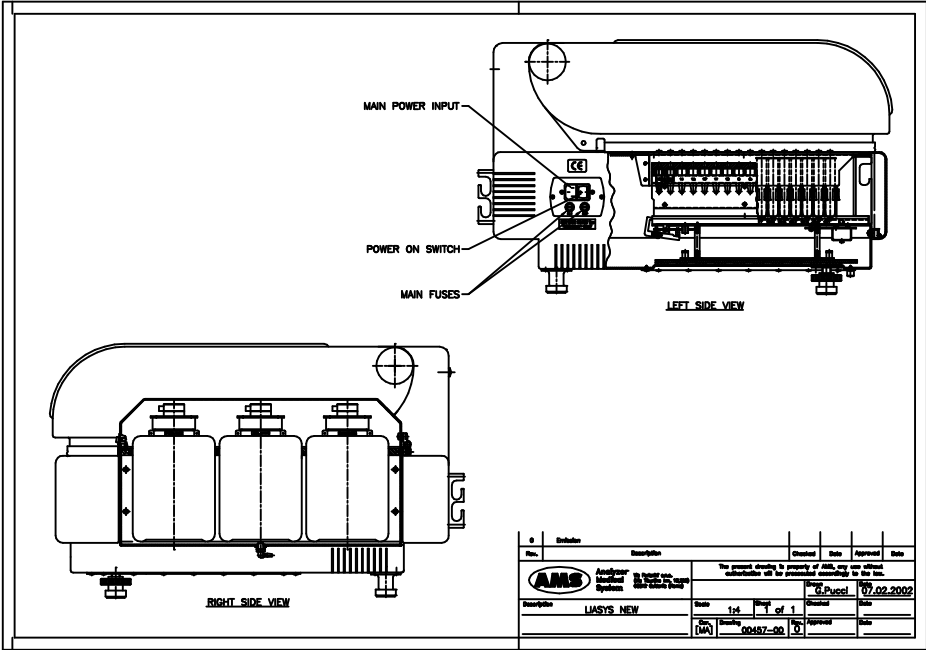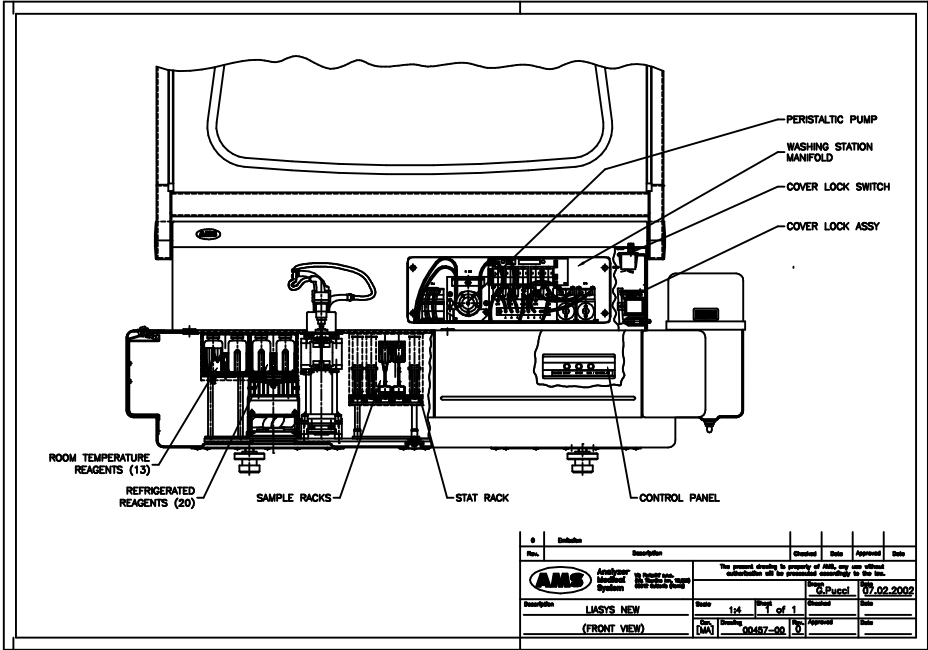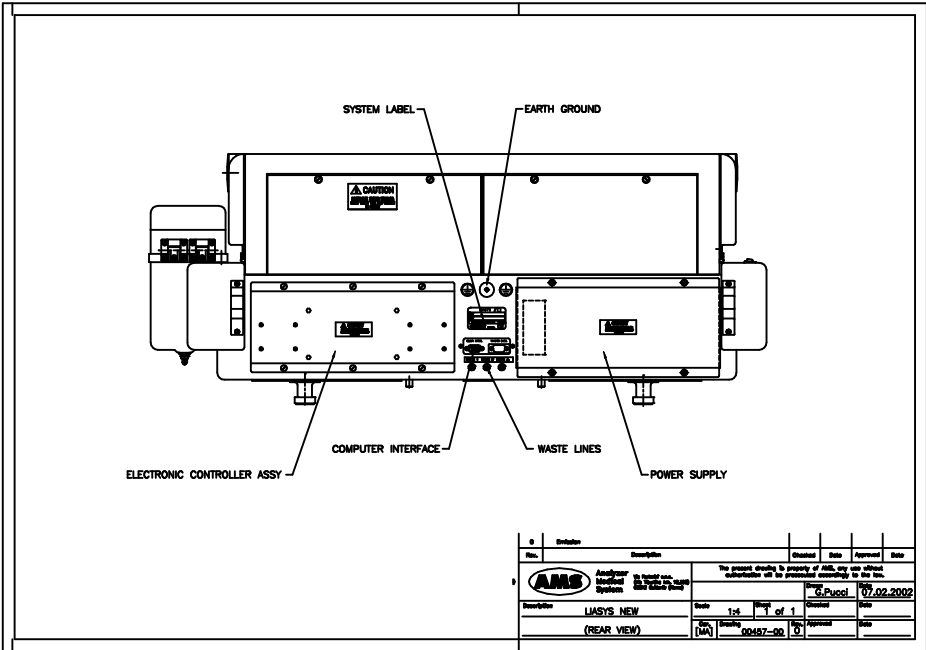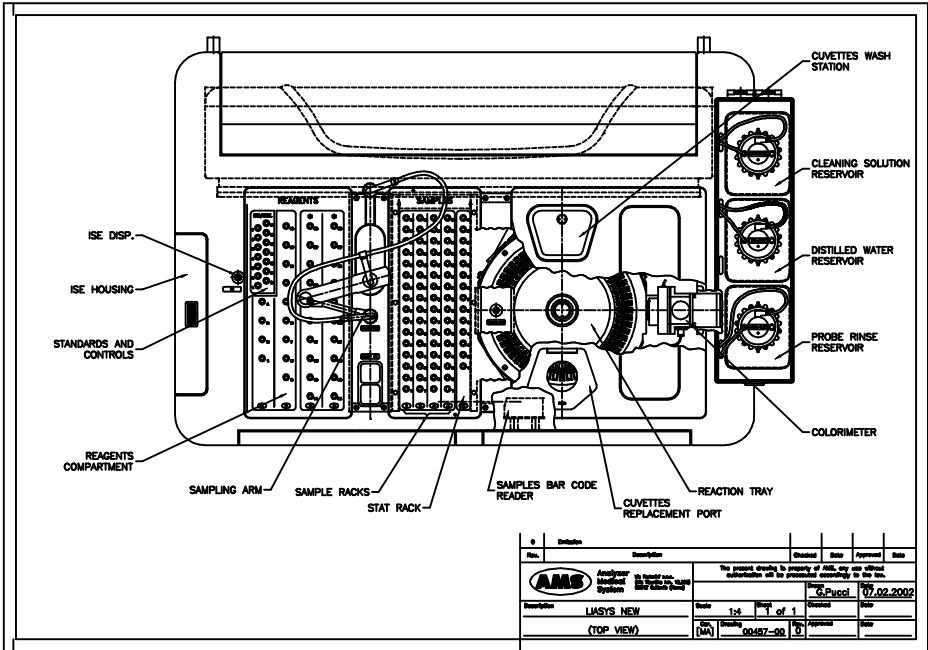


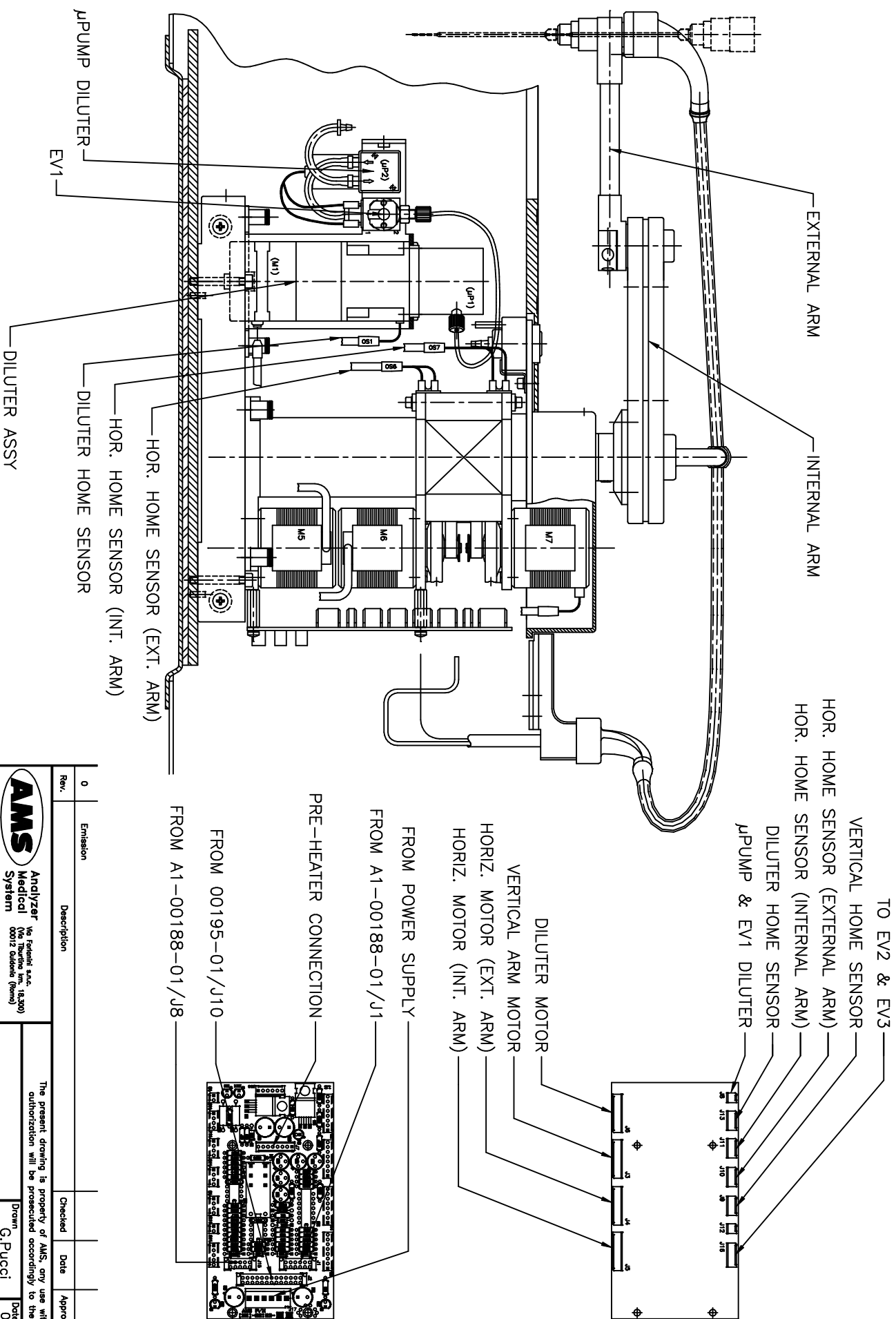


Rev.	0	Emission	Description	Checked	Date	Approved	Date
<p>The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.</p>							
<p><b>AMS</b> Analyzer Medical System  <small>AMS Medical s.p.a. Via Fagnola s.n.c. 16.300 00012 Gattorna (Rovato)</small></p>				Drawn	G. Pucci	Date	03.07.2001
<p>COMPLETE ARM ASSY</p>				Checked	L. Mancini	Date	22.05.2002
<p>Scale 1:2 Sheet 1 of 1</p>				Rev.	0	Approved	A. Gagliarducci
<p>Drawing 9-MA-15-0021-02</p>				Date 22.05.2002			

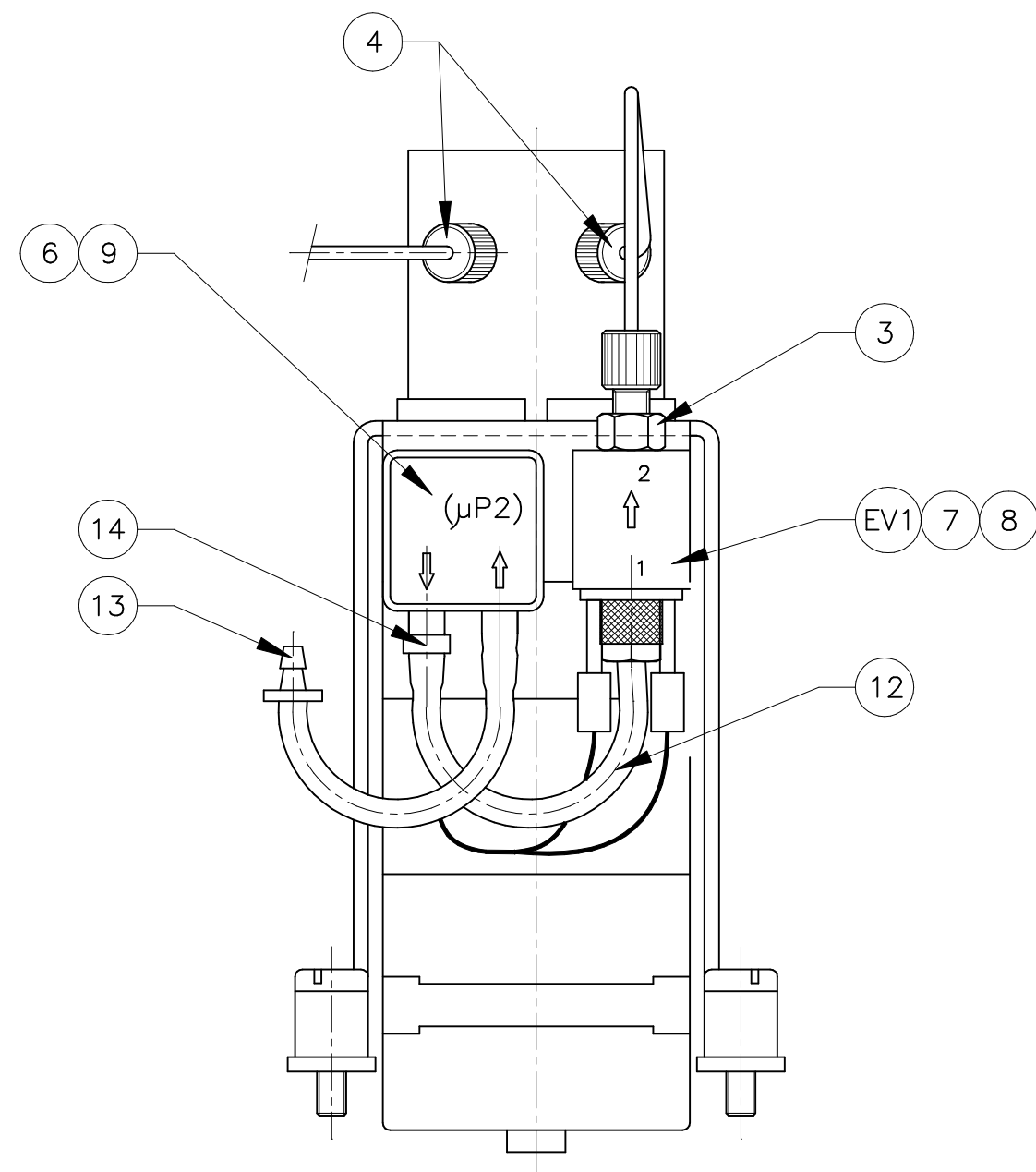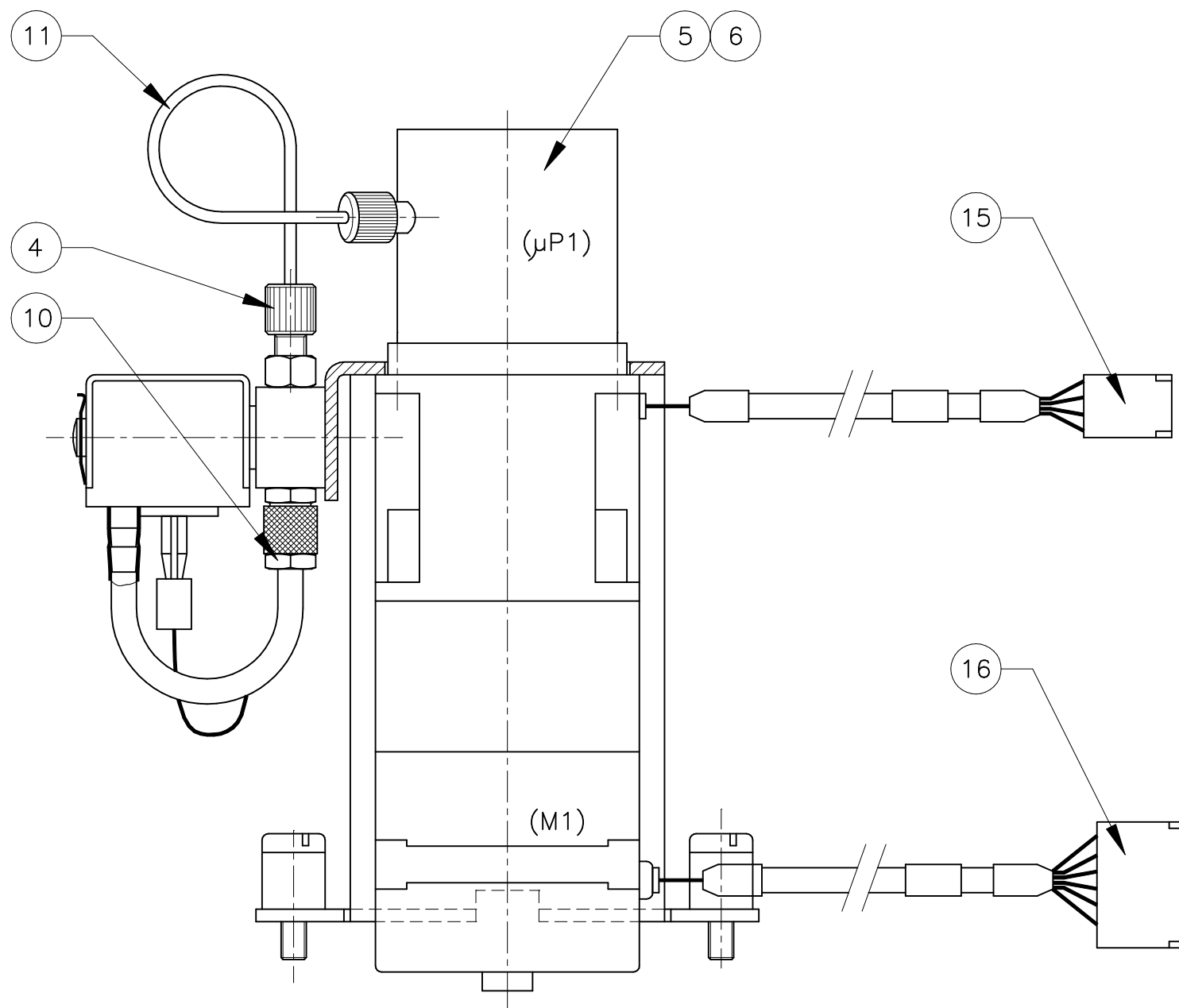


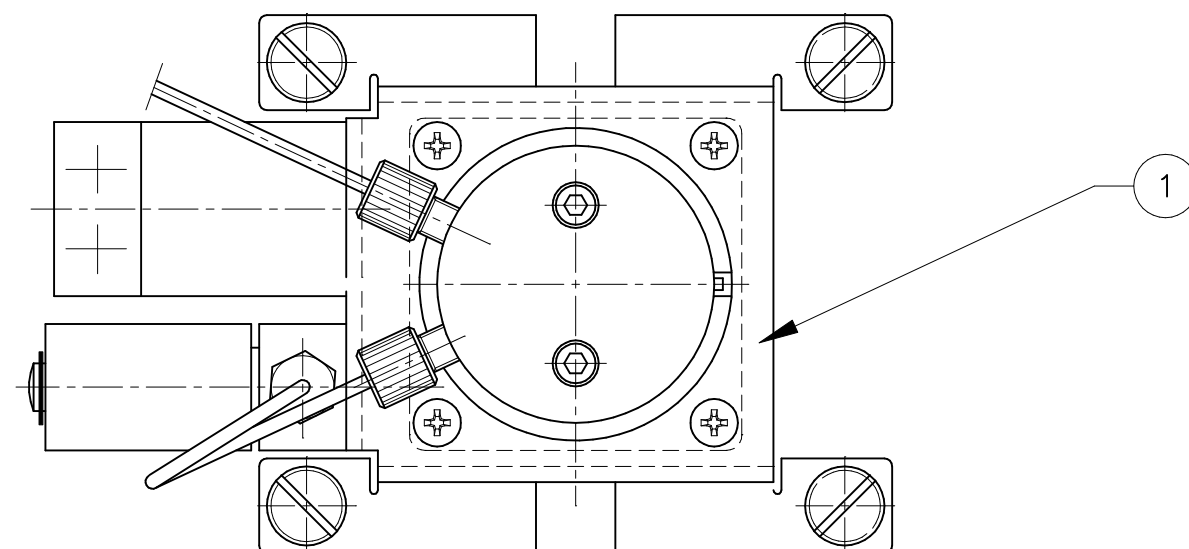
NOTA :—Montare i raccordi pos 3 e 10 con nastro di teflon.




D	Riesame e approvazione (MD14-01 del 06.09.2000)	S.D.S.	06.09.2000	A.G.	06.09.2000
C	Modificato Rif.1 e 5	A.G.	18.02.2000	S.DeSantis	18.02.2000
B	Modificato Rif.4 ed inseriti Rif.9,15,16.	A.G.	16.02.99	S.DeSantis	16.02.99
A	Modificato Rif. 3 e 4; modificato Rif.10 era 9-01-0043-00	A.G.	28.10.98	S.DeSantis	28.10.98
0	Emission				
Rev.	Description	Checked	Date	Approved	Date
Material		<div> <div>AMS</div> <div> <b>Analyzer Medical System</b> </div> <div> <small>Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)</small> </div> </div>		Drawn <b>Carboni</b>	Date <b>12.05.98</b>
				Checked <b>S.D.S.</b>	Date <b>12.05.98</b>
Treatment		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.		Approved <b>A.G.</b>	Date <b>12.05.98</b>
		Description <b>ASSIEME DILUITORE</b>		Scale <b>1:1</b>	Sheet <b>1 of 1</b>
Finish		<b>DILUTER ASSY</b>		Drawing <b>9-15-0007-00</b>	Rev. <b>D</b>

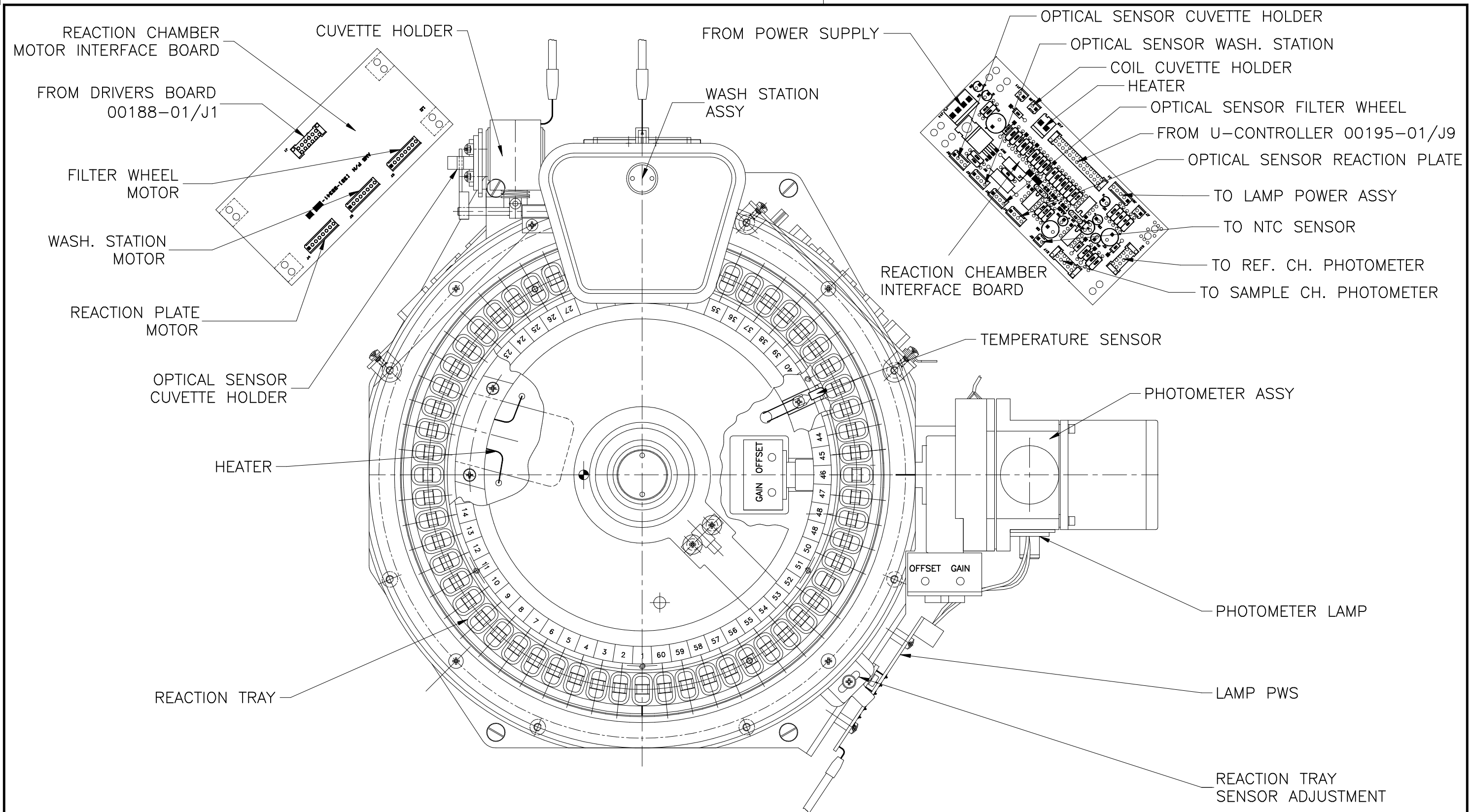



0	Emission				
Rev.	Description	Checked	Date	Approved	Date
<b>AMS</b> Analyzer Medical System <small>Via Forlanini s.n.c.  (Via Tiburtina km. 18,300)  00012 Guidonia (Roma)</small>		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
Description REACTION CHAMBER ASSY		Drawn A.Nicoletti		Date 10.12.2001	
		Scale 1:2	Sheet 1 of 1	Checked A.G.	Date 10.12.2001
		Drawing 9-MA-15-0008-01	Rev. 0	Approved A.G.	Date 10.12.2001

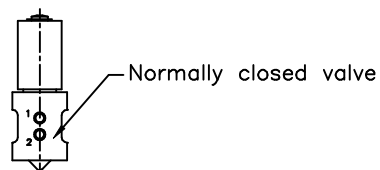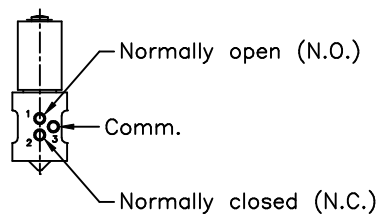

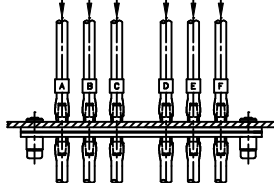
NOZZLE 1 CL. SOL. DISPENSE  
 NOZZLE 3 DIST. WATER DISPENSE  
 NOZZLE 2 DIST. WATER DISPENSE  
 NOZZLE 4 ASPIRATION  
 NOZZLE 2 ASPIRATION  
 NOZZLE 1 ASPIRATION



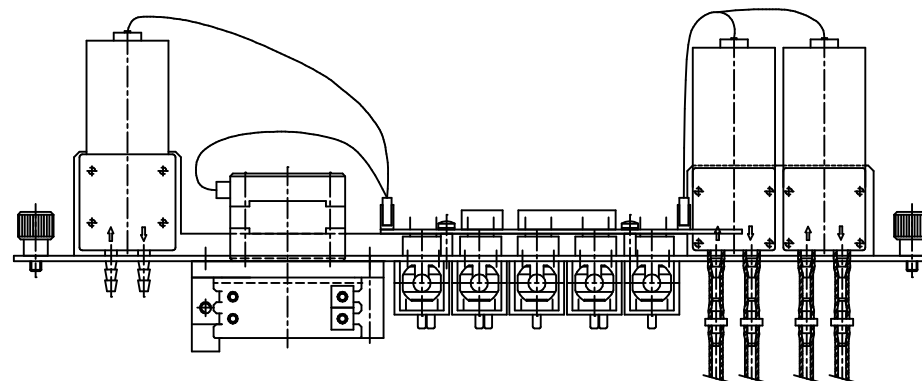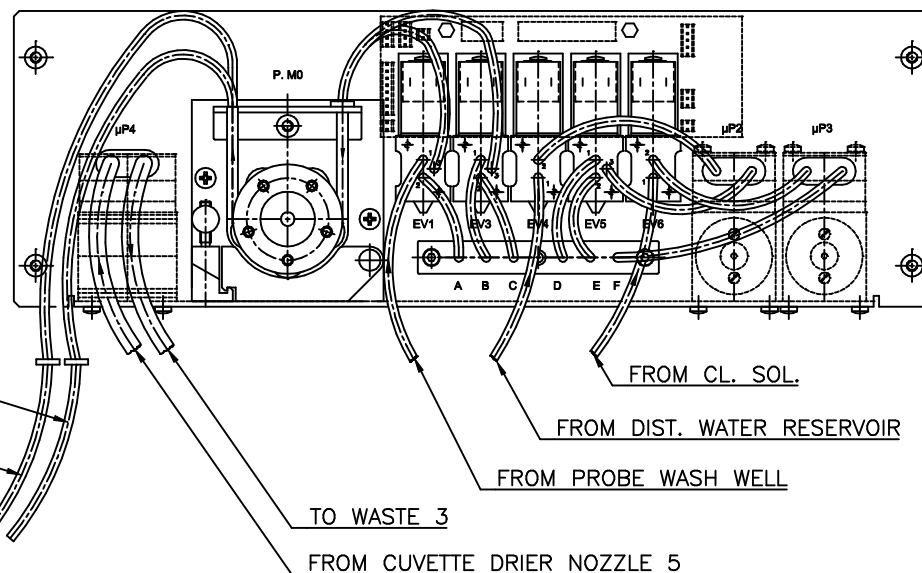TO WASTE 1  
 (1st CH.)  
 TO WASTE 2  
 (2st CH.)



0	Emission				
Rev.	Description	Checked	Date	Approved	Date
<b>AMS</b> Analyzer Medical System <small>Via Forlanini s.n.c.          (Via Tiburtina km. 18,300)          00012 Guidonia (Roma)</small>		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
Description WASHING STATION PUMPS		Scale 1:2		Sheet 1 of 1	
		Car. [MA] Drawing 00406-00		Rev. 0	
		Checked A.Gagliarducci		Date 20.11.2001	
		Approved A.Gagliarducci		Date 20.11.2001	





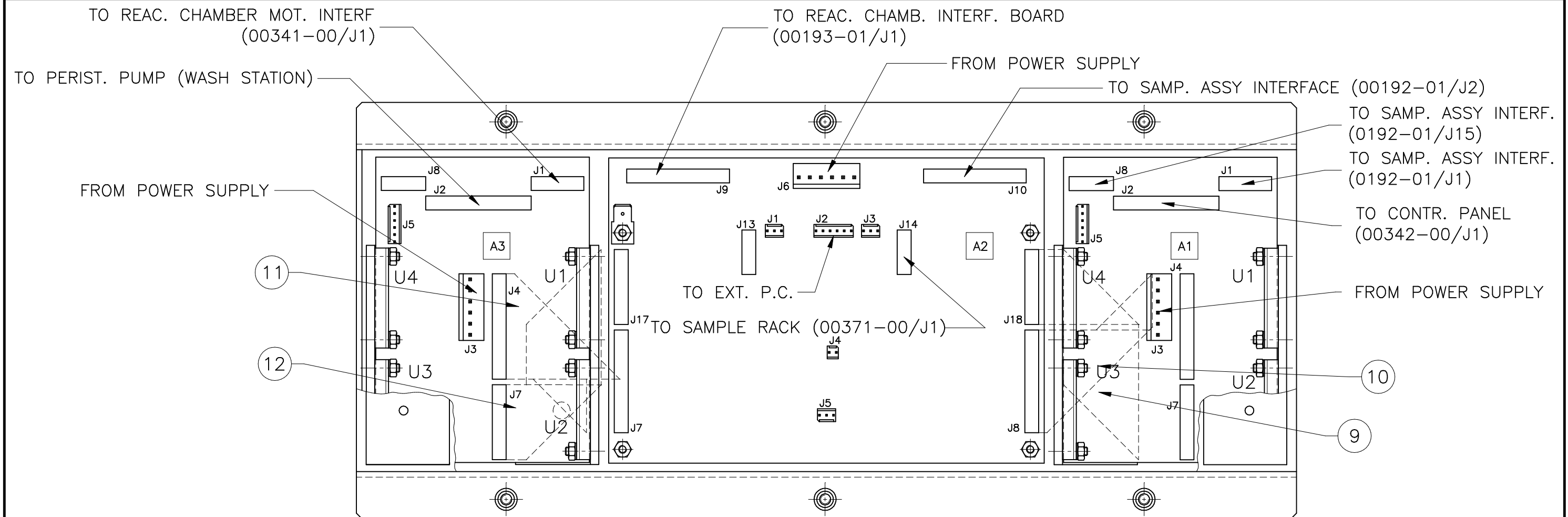


U1	FILTER WHEEL MOT.
U2	WASH STAT. MOT.
U3	REACTION PLATE MOT.
U4	PERISTALTIC PUMP.


U1	HOR. MOT. (EXT. ARM)
U2	HOR. MOT. (INT. ARM)
U3	VERTICAL MOT. (ARM)
U4	DILUTER MOTOR


CONNECTION			
from	to	P/N wire	Ref.
J7-A1	J18-A2	50-00351-00	9
J4-A1	J8-A2	50-00365-00	10
J4-A3	J7-A2	50-00352-00	11
J7-A3	J17-A2	50-00350-00	12

0	Emission (Come da richiesta di modifica N°95)				
Rev.	Description	Checked	Date	Approved	Date
Material	<b>AMS</b> Analyzer Medical System Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)	Drawn	A.Nicoletti	Date	12.12.2001
		Checked	A.G.	Date	12.12.2001
Treatment	The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.	Approved	A.G.	Date	12.12.2001
		Description	ELECTRONIC CONTROLLER ASSY		
Finish		Scale	1:1.5	Sheet	1 of 1
		Car.	[MA]	Drawing	00202-01
		Rev.	0		

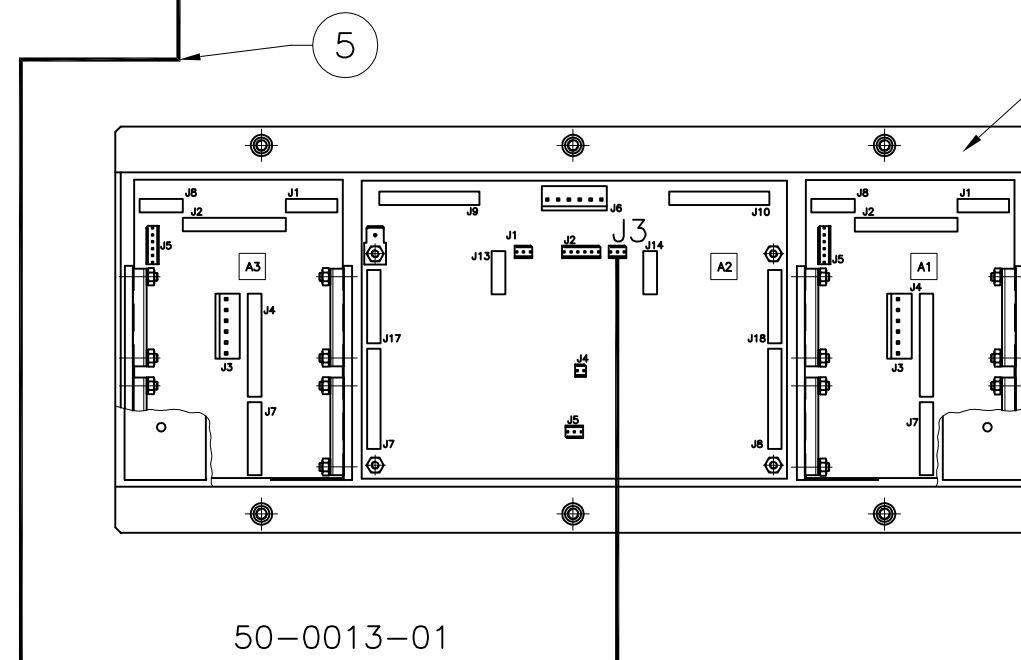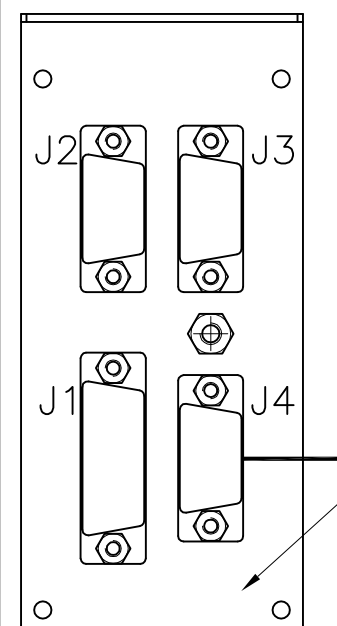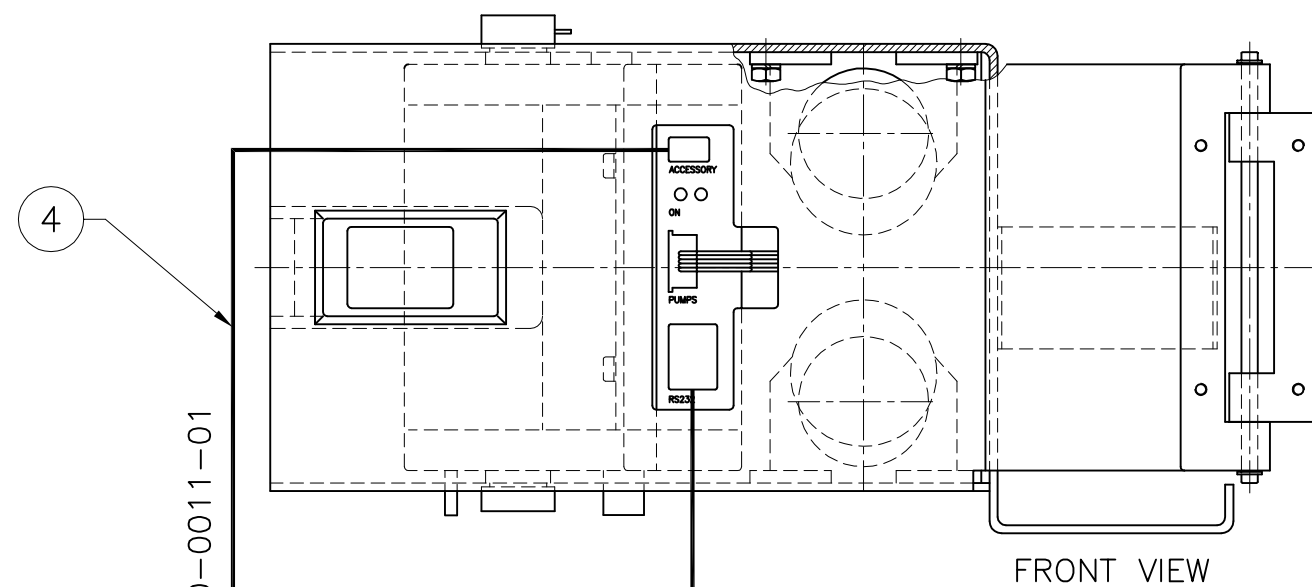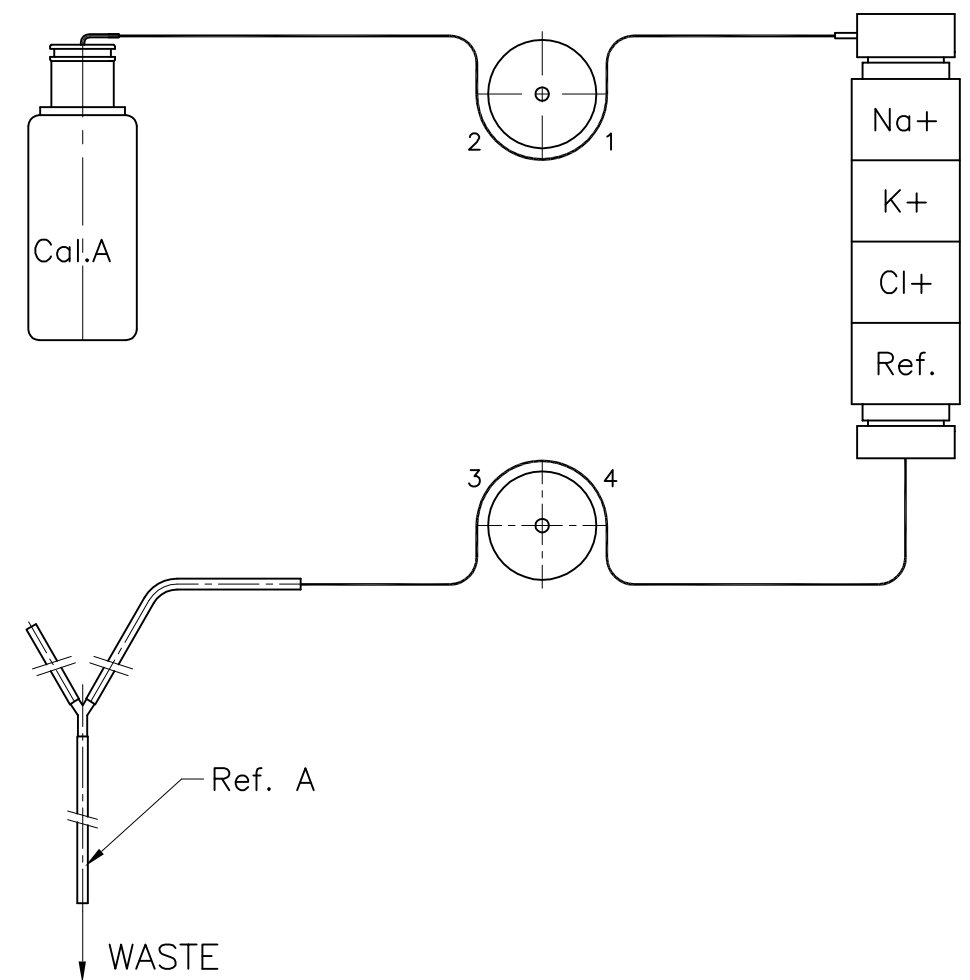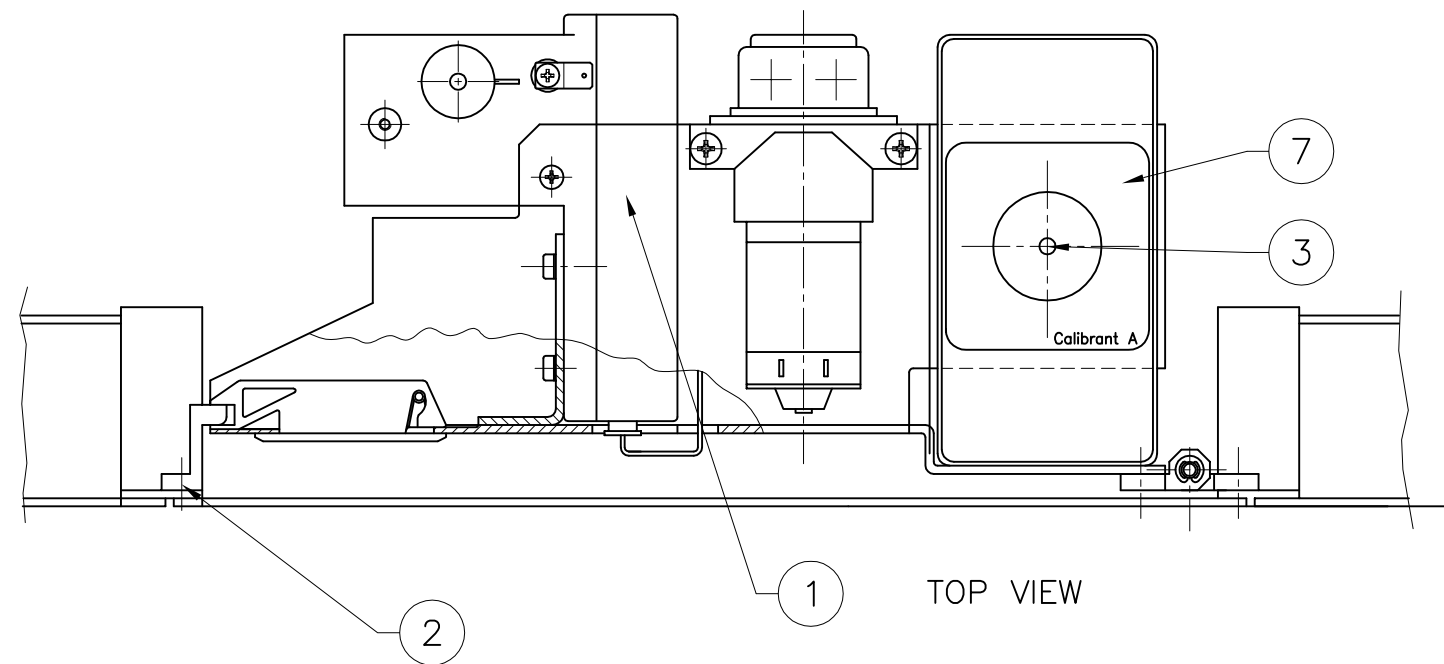


ELECTRONIC CONTROLLER ASSY


0	Emission				
Rev.	Description	Checked	Date	Approved	Date
Material	<div><div></div><div><div>Analyzer Medical System</div><div>Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)</div></div></div>	Drawn G.Pucci	Date 17.03.2003		
		Checked L.Mancini	Date 17.03.2003		
Treatment	The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.	Approved A.Gagliarducci	Date 18.03.2003		
		Description KIT MODULO ISE		Scale 1:1.5	Sheet 1 of 1
Finish	ISE MODULE KIT	Drawing 65-0019-00		Rev. 0	













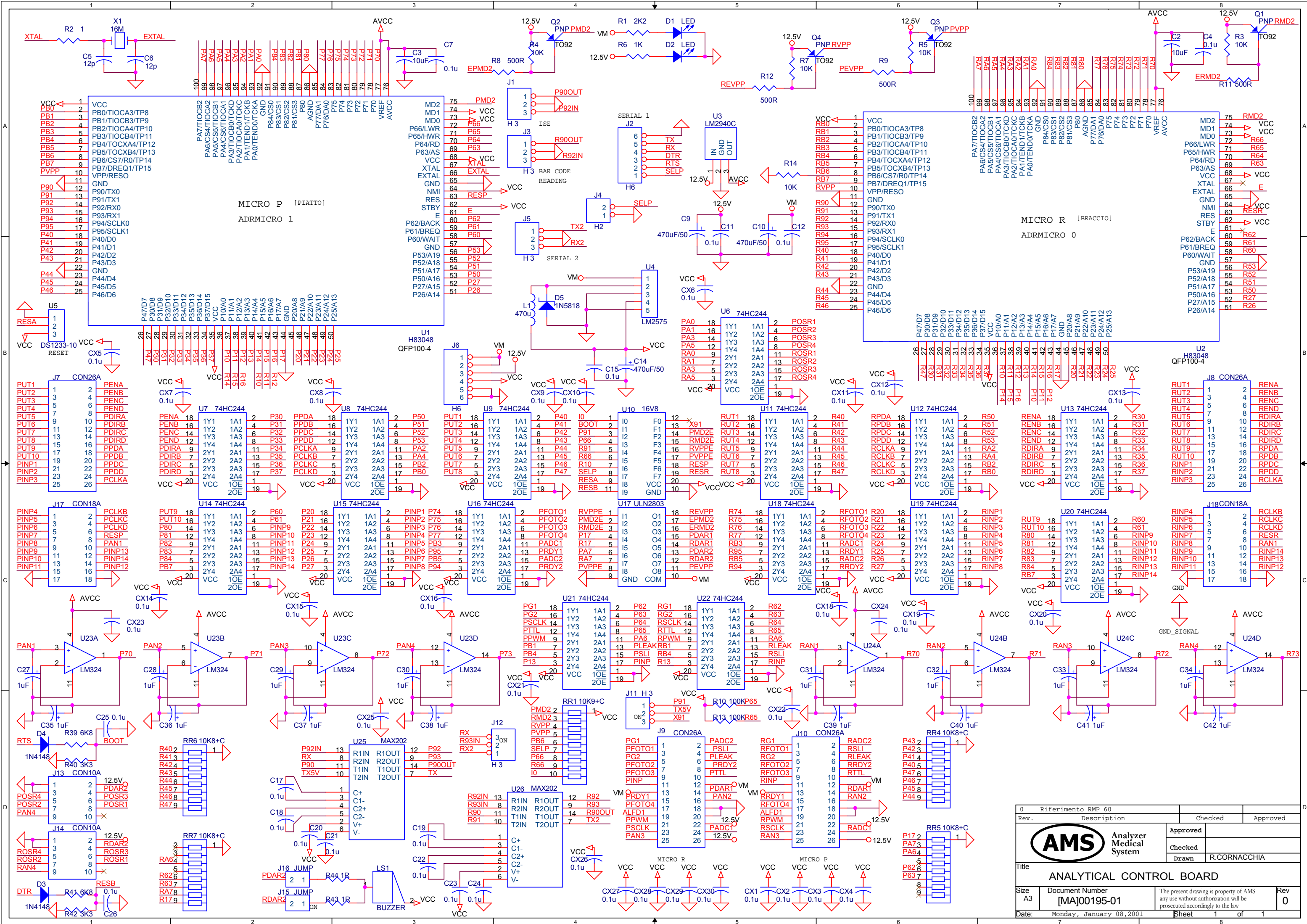



0	Riferimento RMP 60		
Rev.	Description	Checked	Approved
<div><div>AMS</div><div>Analyzer Medical System</div></div>		Approved	
		Checked	
		Drawn	R.CORNACCHIA
Title			
ANALYTICAL CONTROL BOARD			
Size	Document Number	The present drawing is property of AMS any use without authorization will be prosecuted accordingly to the law	
A3	[MA]00195-01		
Date:	Monday, January 08, 2001	Sheet	1 of 1





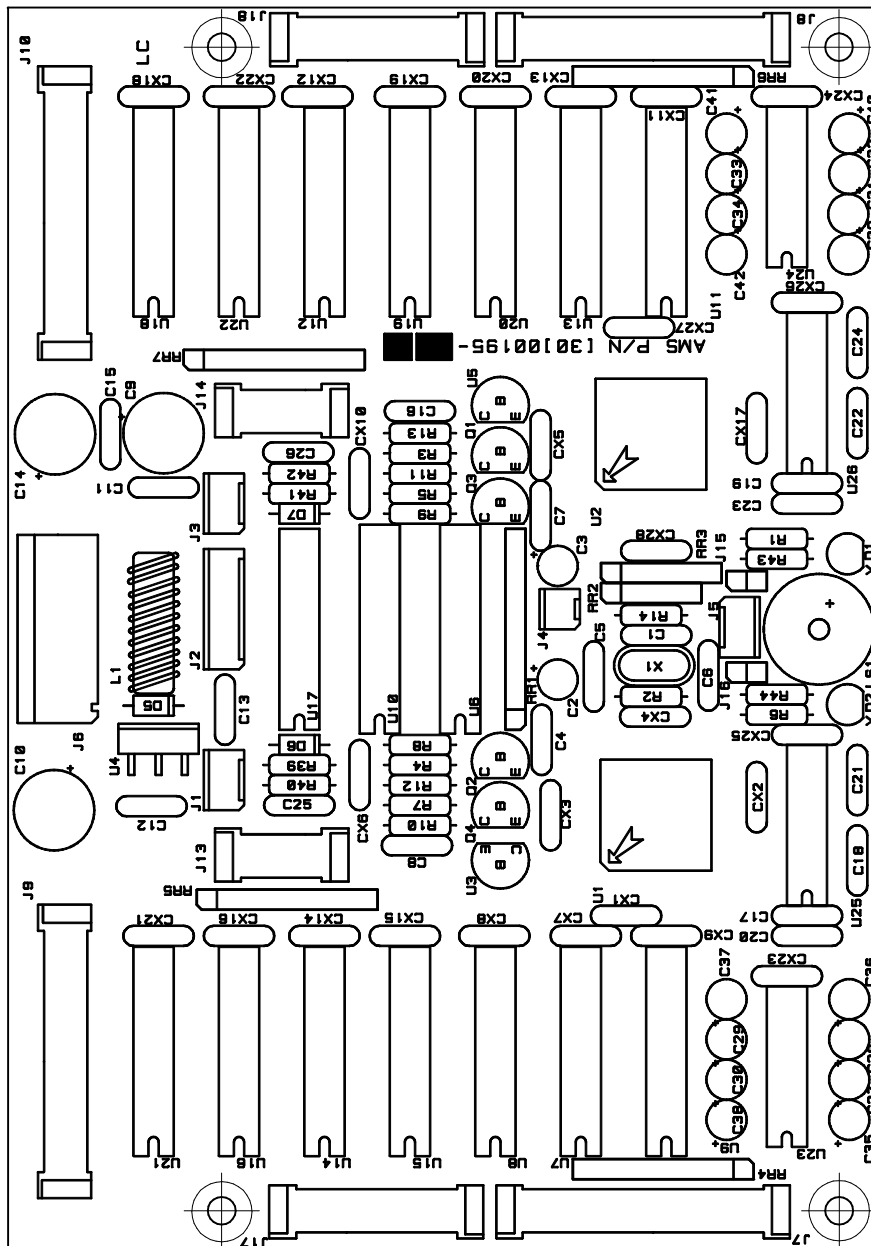

0	Emission (Come da richiesta di modifica N°60)				
Rev.	Description	Checked	Date	Approved	Date
 <b>AMS</b> Analyzer Medical System Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
Description ANALYTICAL CONTROL BOARD		Scale	Sheet 1 of 1	Drawn A.Nicoletti	Date 12.12.2000
		Car. [MA]	Drawing 00195-01	Rev. 0	Approved A.G.

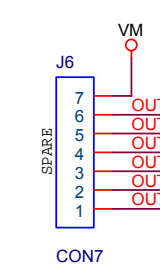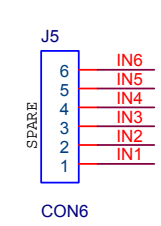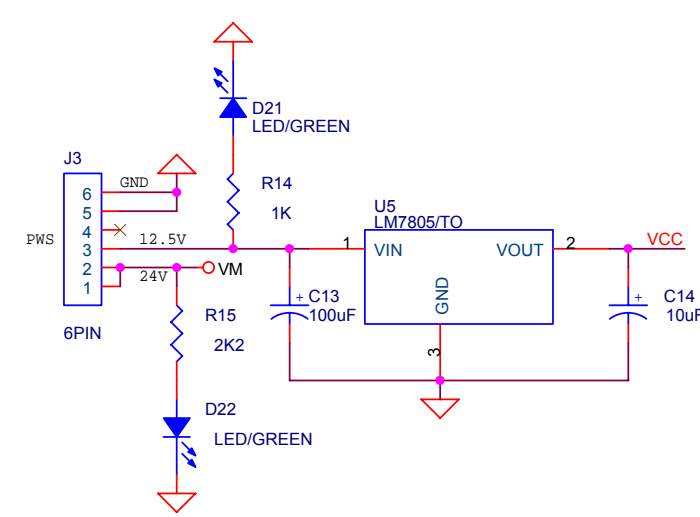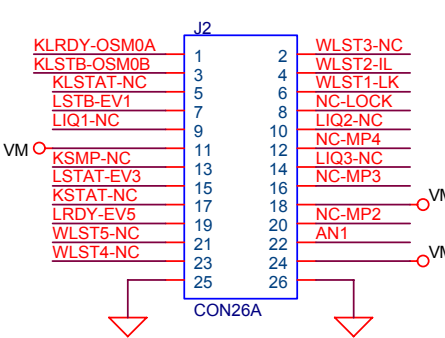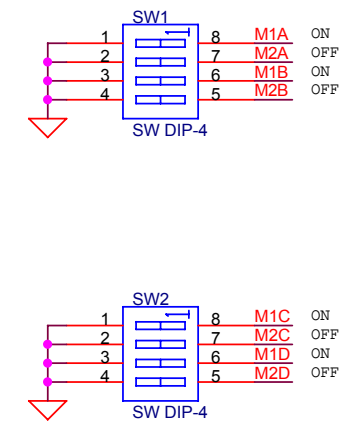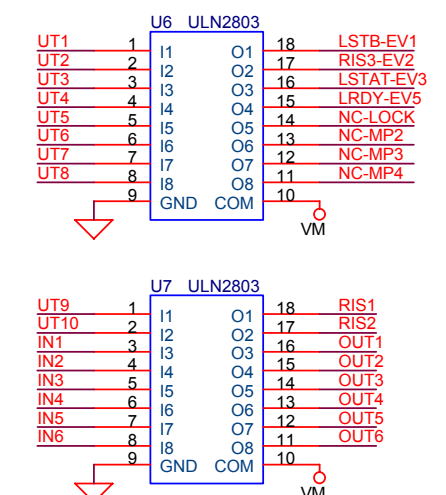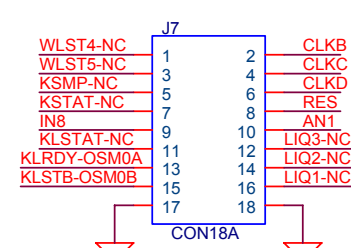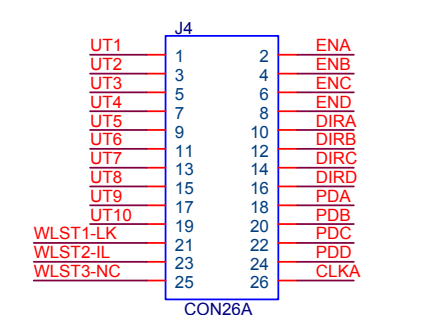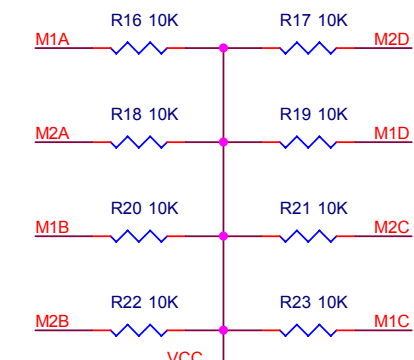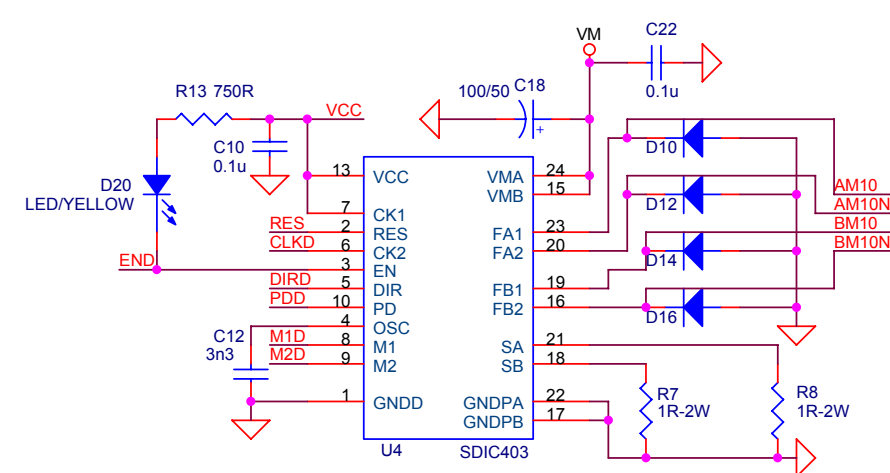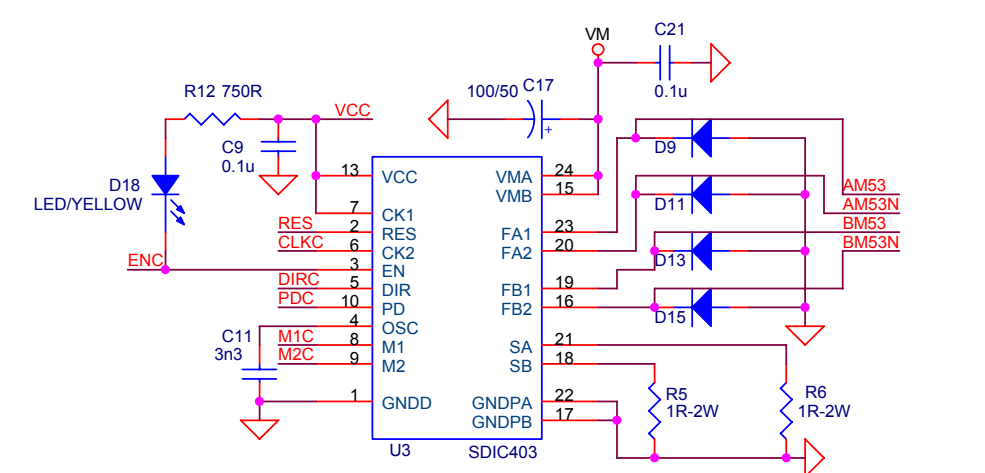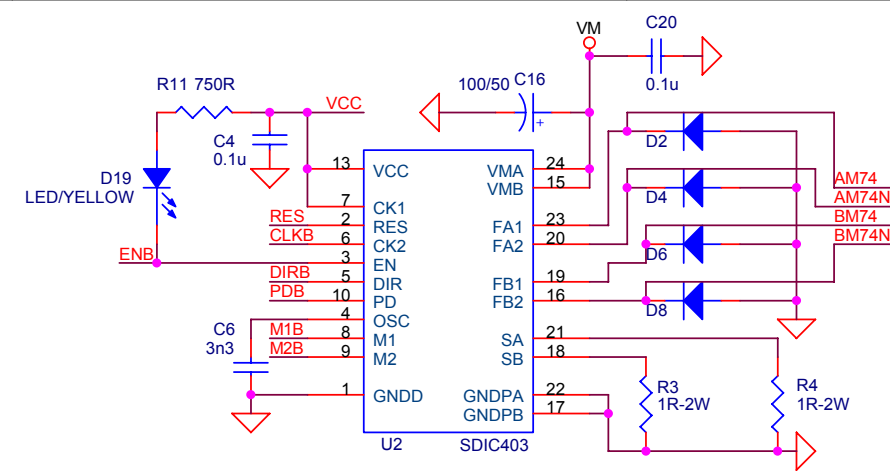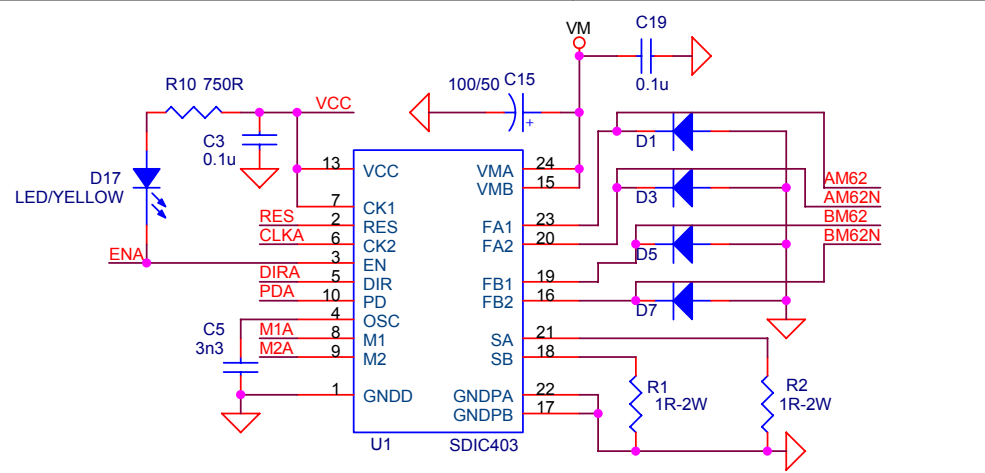


0	Riferimento RMP60		
Rev.	Description	Checked	Approved
 <b>AMS</b> Analyzer Medical System		Approved	
		Checked	
		Drawn	R.CORNACCHIA
Title STEPPER MOTORS DRIVER BOARD			
Size A3	Document Number [MA]00188-01	The present drawing is property of AMS any use without authorization will be prosecuted accordingly to the law	
Date:	Tuesday, January 09, 2001	Sheet 1	of 1

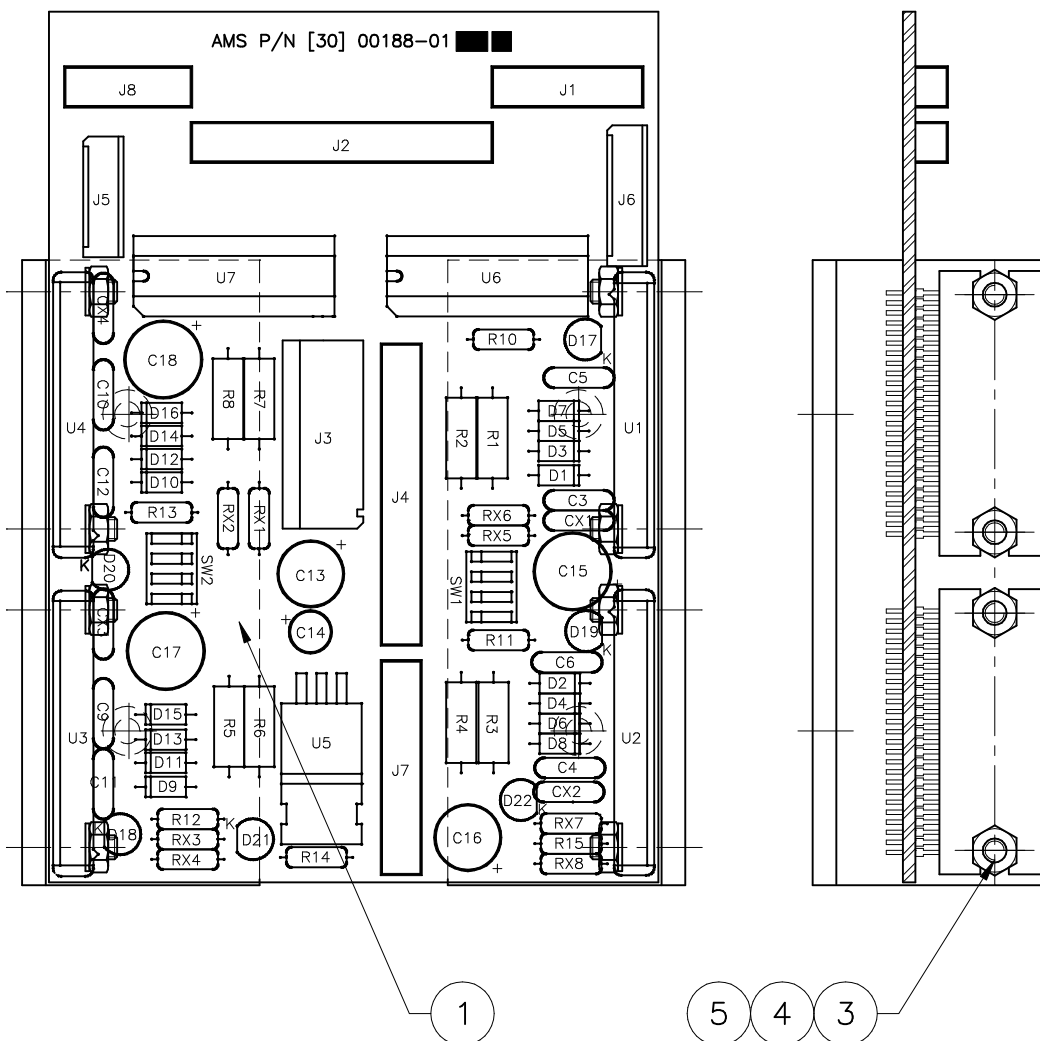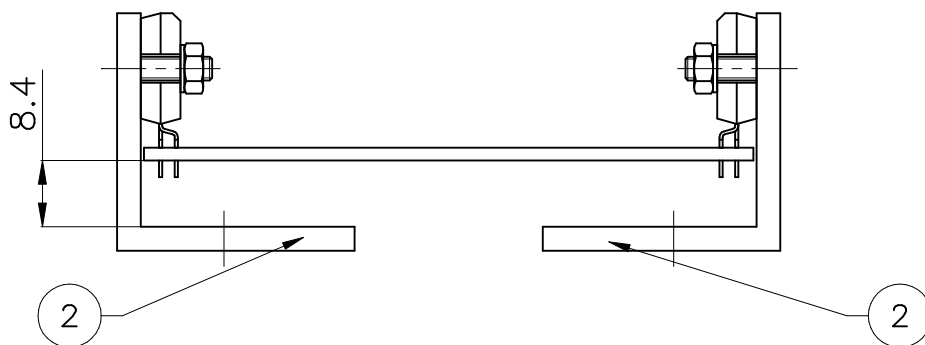



NOTA : -Montare i Rif.2 dopo aver steso, sui componenti da dissipare, la giusta quantità di thermocompound Rif.6  
 -Serrare i Rif. 3,4 e 5 con una coppia pari a 13Nmm.


A	Come da richiesta di modifica N°152 (aggiunte basette dissipatrici)				
0	Emission (Come da richiesta di modifica N°95)	A.N.	06.08.2001	A.G.	06.08.2001
Rev.	Description	Checked	Date	Approved	Date

<b>AMS</b> Analyzer Medical System Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
		Drawn G.Pucci		Date 04.09.2002	
Description STEPPER MOTORS DRIVER BOARD		Scale 1:1	Sheet 1 of 1	Checked L.MANCINI	Date 10.10.2002
		Car. [30]	Drawing 00188-01	Rev. A	Approved A.G. Date 10.10.2002

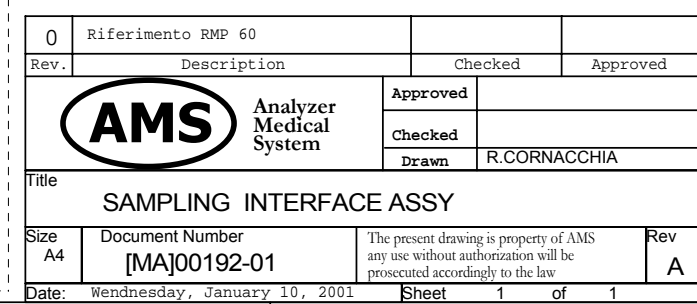


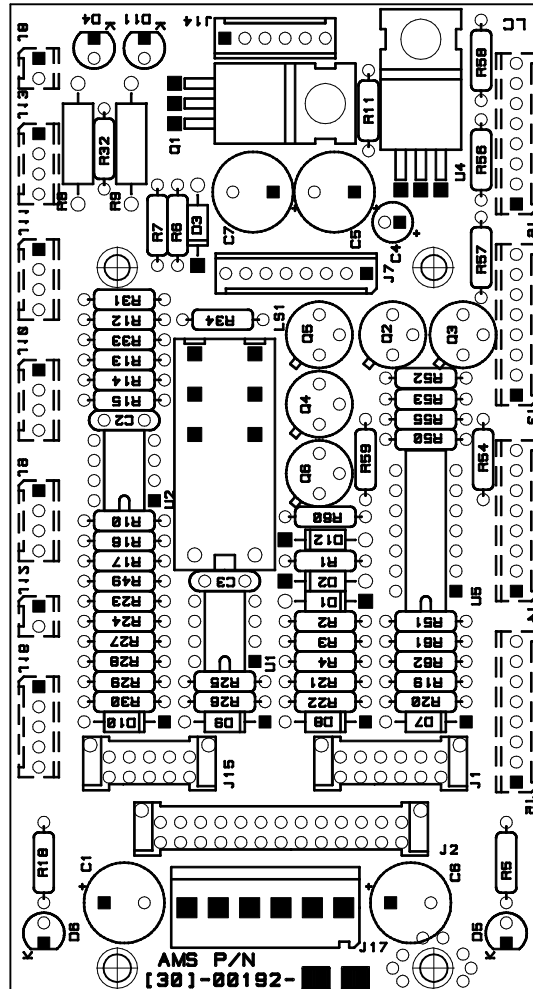


0	Emission				
Rev.	Description	Checked	Date	Approved	Date
 <b>Analyzer Medical System</b> Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
Description SAMPLING INTERFACE ASSY		Scale		Sheet 1 of 1	Drawn A.Nicoletti Date 01.06.2001
		Car. [MA]		Drawing 00192-01	Checked R.C. Date 01.06.2001
		Rev. 0	Approved A.G.	Date 01.06.2001	

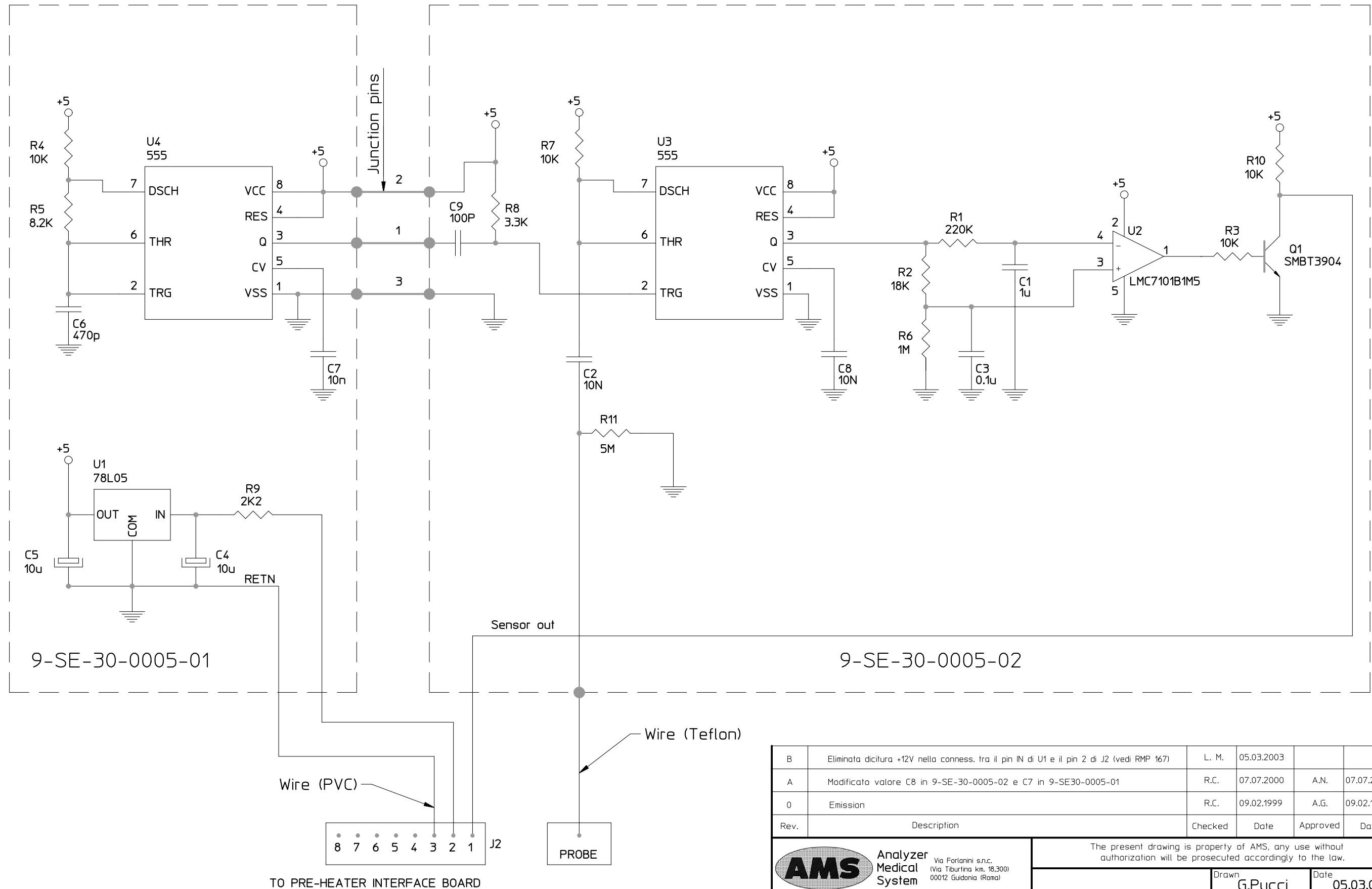


B	Eliminata dicitura +12V nella conness. tra il pin IN di U1 e il pin 2 di J2 (vedi RMP 167)	L. M.	05.03.2003		
A	Modificato valore C8 in 9-SE-30-0005-02 e C7 in 9-SE30-0005-01	R.C.	07.07.2000	A.N.	07.07.2000
0	Emission	R.C.	09.02.1999	A.G.	09.02.1999
Rev.	Description	Checked	Date	Approved	Date
<b>AMS</b> Analyzer Medical System Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
		Drawn <b>G.Pucci</b>		Date <b>05.03.03</b>	
Description <b>SCHEMATIC</b>		Scale <b>1:1</b>	Sheet <b>1 of 1</b>	Checked <b>L. Massenzi</b>	Date <b>05.03.03</b>
Drawing <b>9-SE-30-0005-XX</b>		Rev. <b>B</b>	Approved <b>A. Gagliarducci</b>	Date <b>05.03.03</b>	

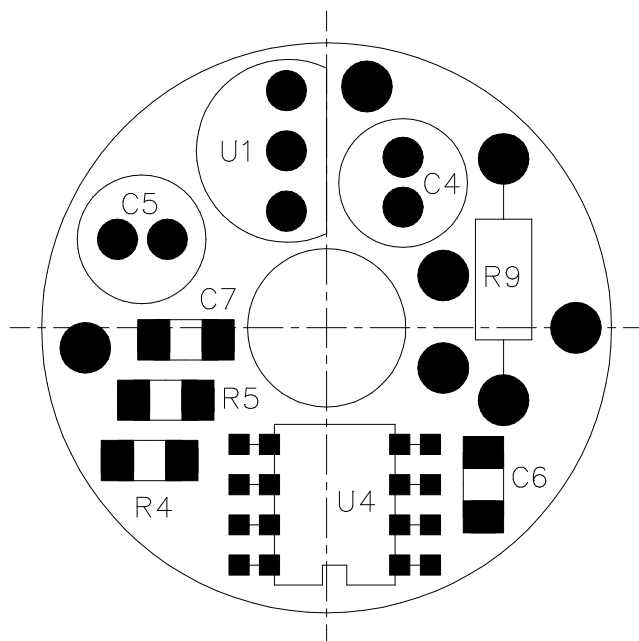


A	Modificato valore C7		07.07.2000		07.07.2000
0	Emission				
Rev.	Description	Checked	Date	Approved	Date


 <b>Analyzer Medical System</b> Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)	The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
			Drawn G.Pucci	Date 09.02.99
Description LEVEL SENSOR BOARD	Scale 4:1	Sheet 1 of 1	Checked R.C.	Date 09.02.99
	Drawing 9-MA-30-0005-01		Rev. A	Approved A.G.
				Date 09.02.99

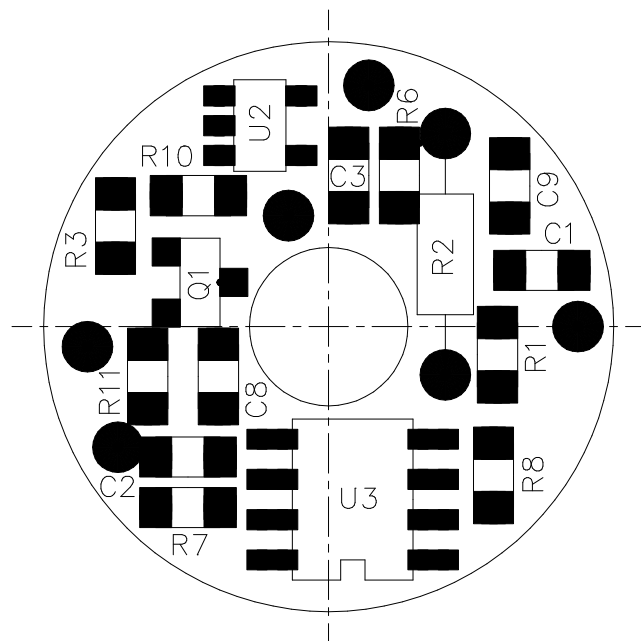



A	Modificato valore C8		07.07.2000		07.07.2000
0	Emission				
Rev.	Description	Checked	Date	Approved	Date


 <b>Analyzer Medical System</b> Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)	The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
			Drawn <b>G.Pucci</b>	Date <b>09.02.99</b>
Description <b>LEVEL SENSOR BOARD</b>	Scale <b>4:1</b>	Sheet <b>1 of 1</b>	Checked <b>R.C.</b>	Date <b>09.02.99</b>
	Drawing <b>9-MA-30-0005-02</b>		Rev. <b>A</b>	Approved <b>A.G.</b>
		Date <b>09.02.99</b>		

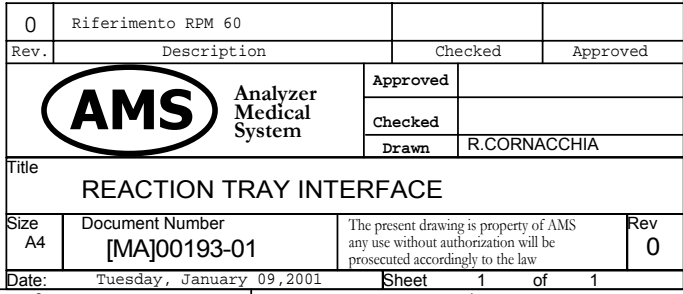

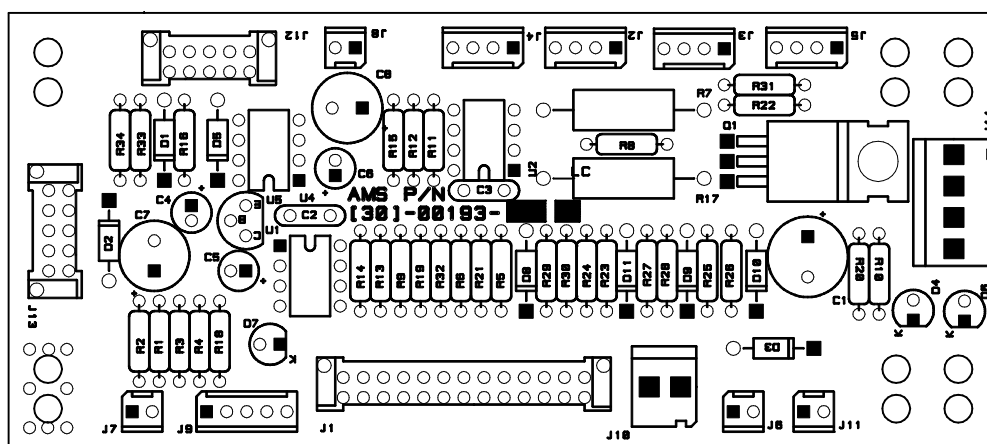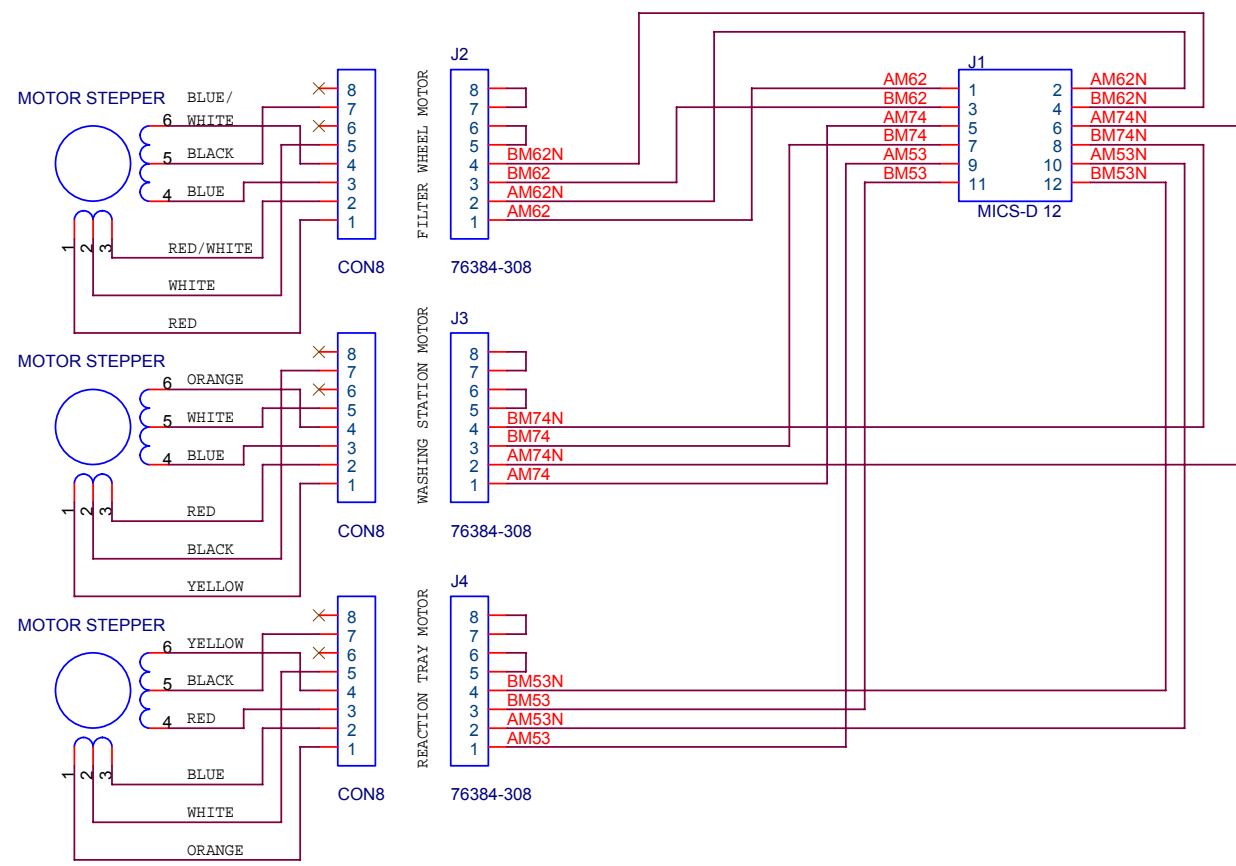
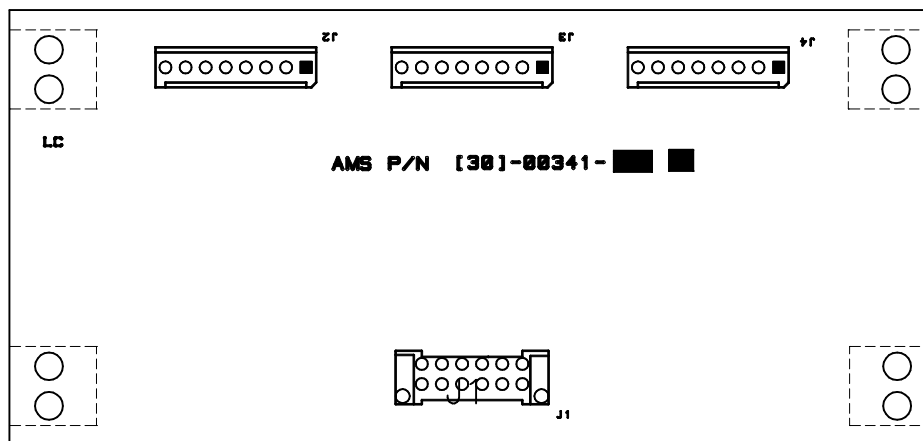





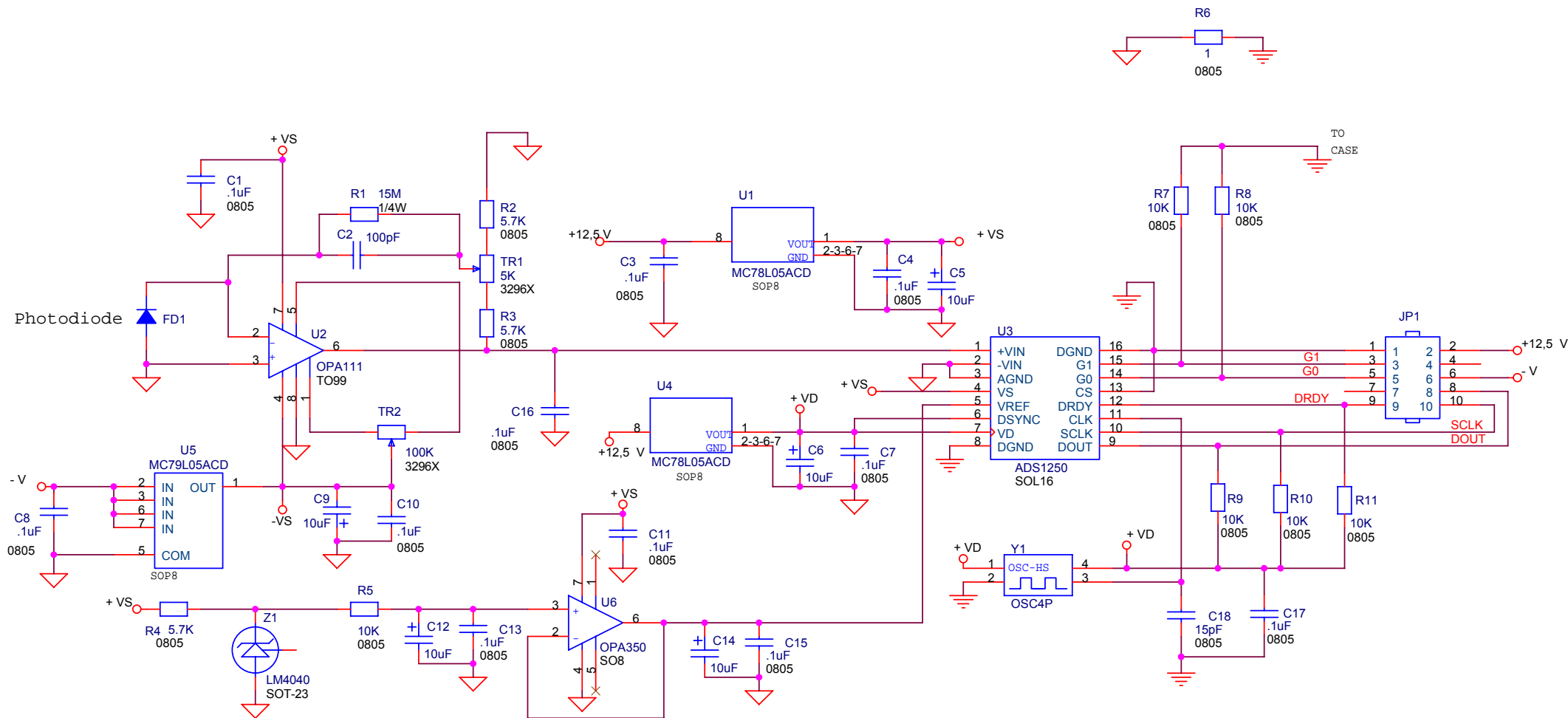
0	Emission (Come da richiesta di modifica N°60)					
Rev.	Description		Checked	Date	Approved	Date
 <b>Analyzer Medical System</b> Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)			The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
			Drawn <u>A.Nicoletti</u>		Date <u>12.12.2000</u>	
Description <u>REACTION TRAY INTERFACE</u>			Scale _____	Sheet <u>1 of 1</u>	Checked <u>R.C.</u>	Date <u>12.12.2000</u>
Car. <u>[MA]</u>			Drawing <u>00193-01</u>	Rev. <u>0</u>	Approved <u>A.G.</u>	Date <u>12.12.2000</u>




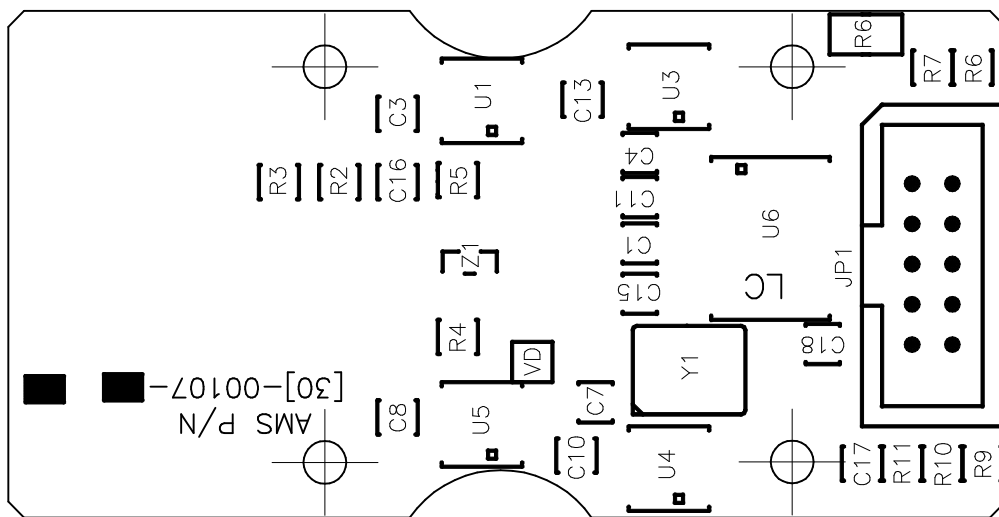
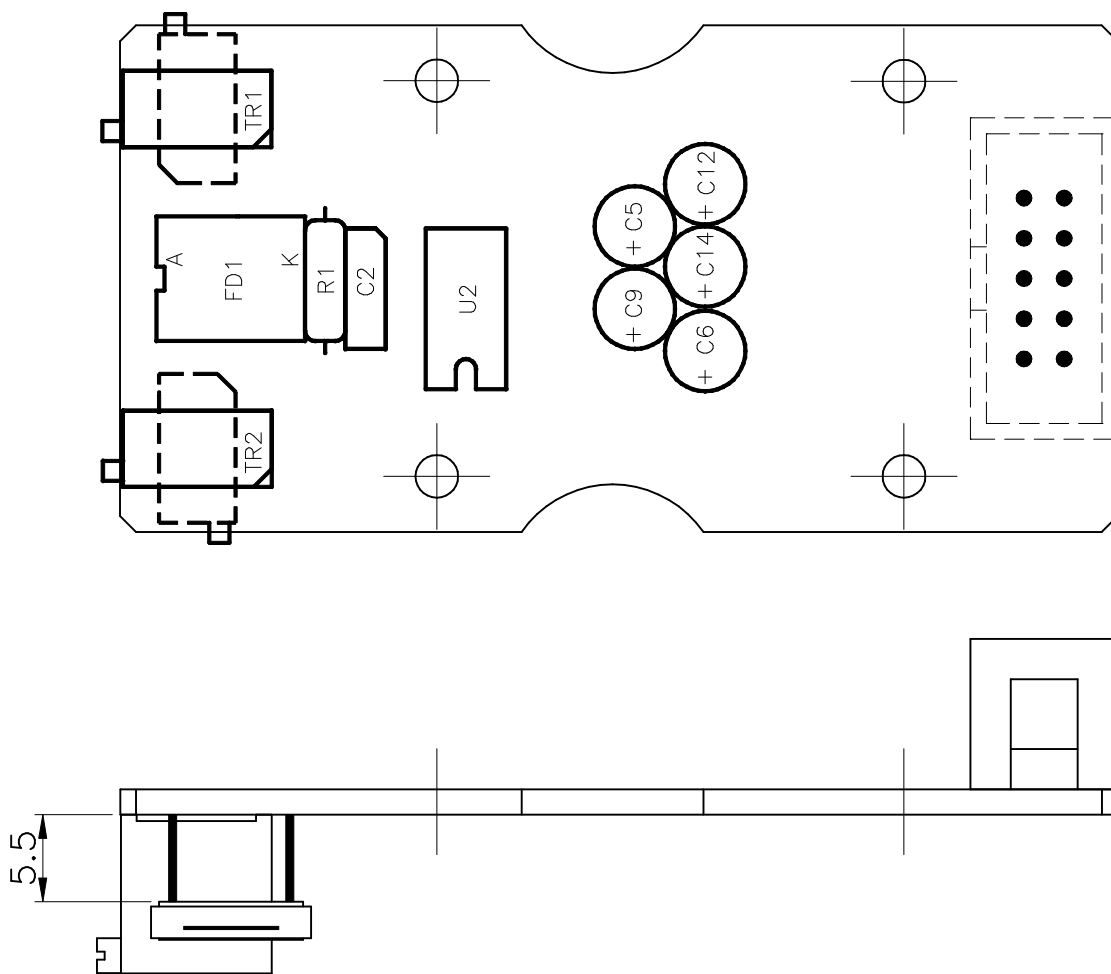
0	Emission			
Rev.	Description		Checked	Approved
 Analyzer Medical System			Approved	
			Checked	
			Drawn	R.CORNACCHIA
Title				
REACTION CHAMBER MOTOR INTERFACE				
Size	Document Number		The present drawing is property of AMS any use without authorization will be prosecuted accordingly to the law	Rev 0
A4	[MA]00341-00			
Date:	Wednesday, May 30, 2001		Sheet	1 of 1




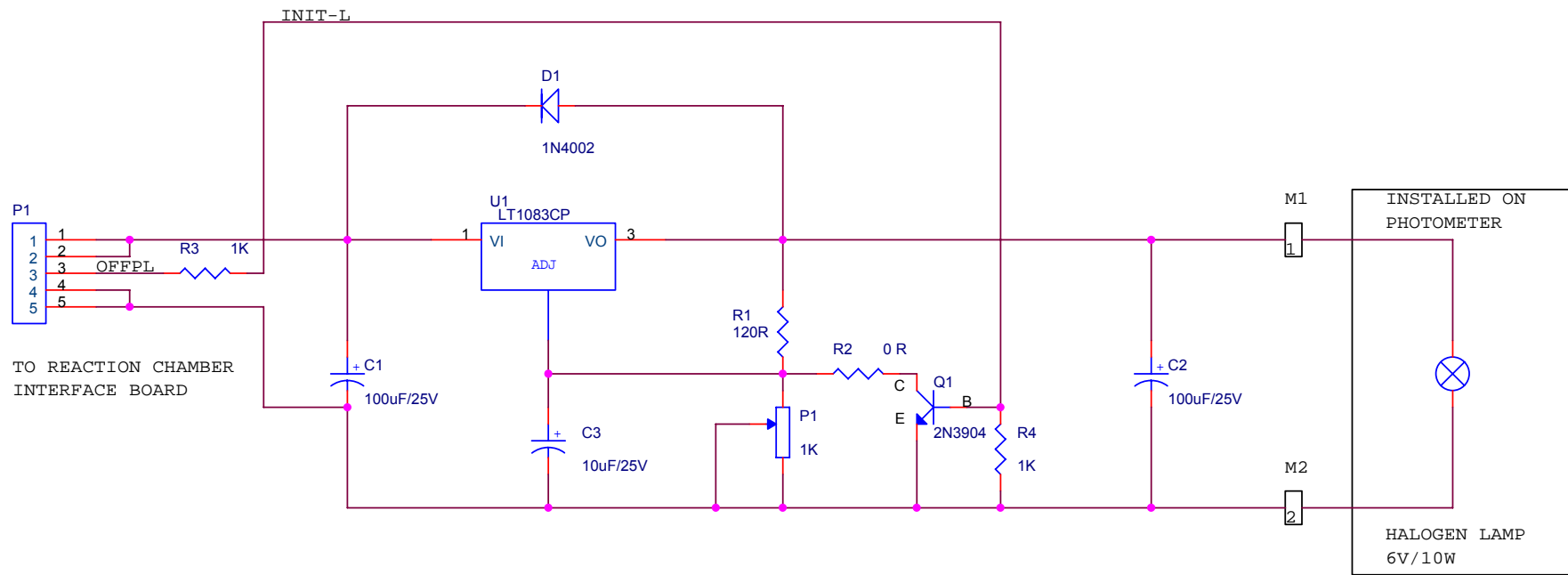
0	Emission				
Rev.	Description	Checked	Date	Approved	Date
 <b>Analyzer Medical System</b> Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
Description <b>REACTION CHAMBER</b> <b>MOTOR INTERFACE</b>		Scale	Sheet 1 of 1	Drawn A.Nicoletti	Date 30.05.2001
		Car. [MA]	Drawing 00341-00	Rev. 0	Approved A.G.




0	Emission		
Rev.	Description	Checked	Approved
 <b>Analyzer Medical System</b>		Approved	
		Checked	
		Drawn	R.CORNACCHIA
Title			
PRE-AMPL/ADC			
Size A4	Document Number [MA]00107-00	The present drawing is property of AMS any use without authorization will be prosecuted accordingly to the law	
Date: Friday, February 01, 2002	Sheet 1 of 1	Rev 0	




0	Emission				
Rev.	Description	Checked	Date	Approved	Date
 <b>Analyzer Medical System</b> Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
		Drawn G.Pucci		Date 07.09.2000	
Description PRE-AMPL/ADC		Scale 2:1	Sheet 1 of 1	Checked R.C.	Date 07.09.2000
		Car. [MA]	Drawing 00107-00	Rev. 0	Approved A.G.
				Date 07.09.2000	




0	Emission			
Rev.	Description		Checked	Approved
			Approved	
			Checked	
			Drawn	R.CORNACCHIA
Title				
PHOTOMETER LAMP PWS				
Size	Document Number		The present drawing is property of AMS any use without authorization will be prosecuted accordingly to the law	Rev
A4	9-[MA]-30-0004-00			0
Date:	Saturday, January 12, 2002		Sheet	1 of 1


CONNECTION	
Board	P1
1	1
2	2
3	3
4	4
5	5


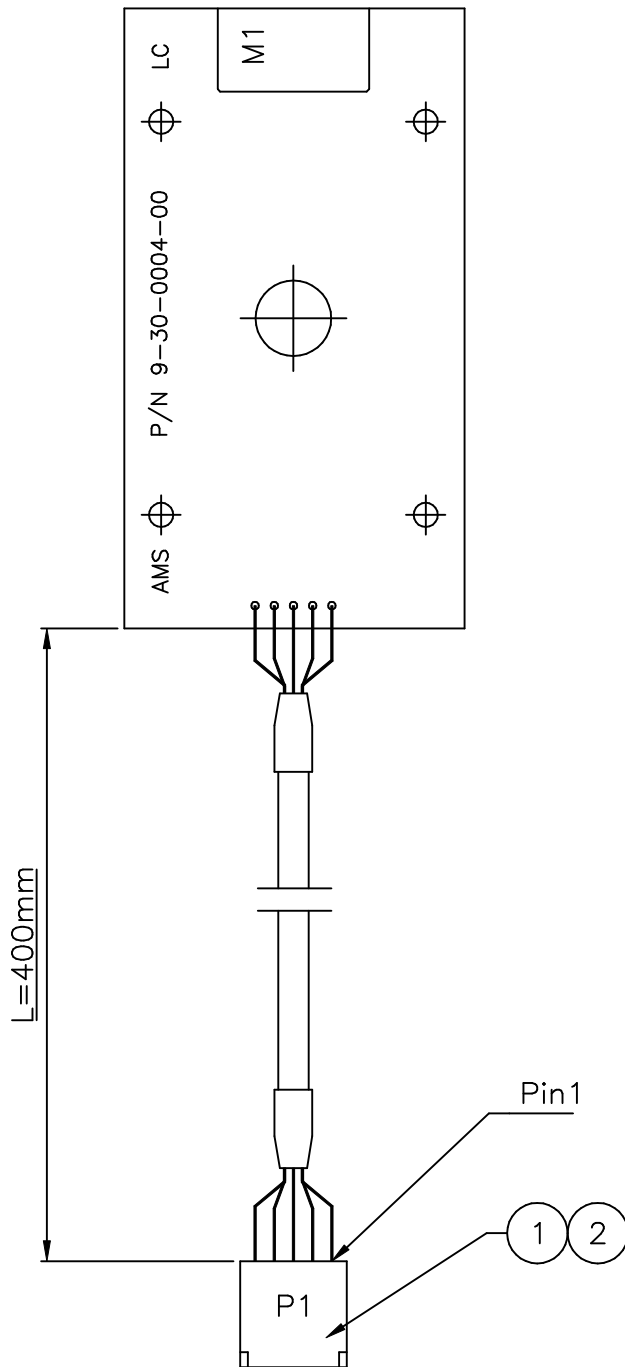



0	Emission				
Rev.	Description	Checked	Date	Approved	Date
<b>AMS</b> Analyzer Medical System Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
Description		Scale	Sheet	Checked	Date
LAMP PWS BOARD		1:1	1 of 1	G.Pucci	04.03.99
		Drawing	Rev.	Approved	Date
		9-ME-30-0004-00	0	A.Gagliarducci	04.03.99



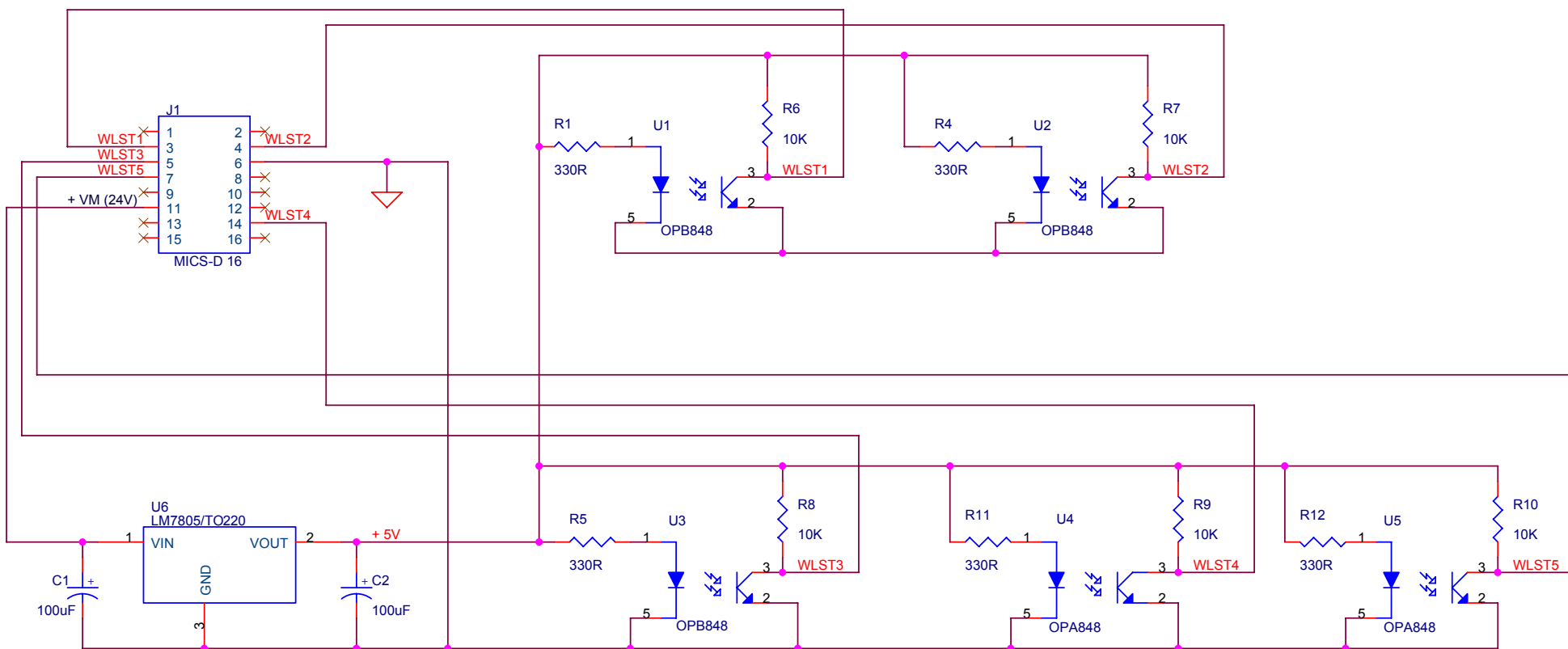


0	Emission		
Rev.	Description	Checked	Approved
 <b>Analyzer Medical System</b>		Approved	
		Checked	
		Drawn	R.CORNACCHIA
Title			
SAMPLE RACKS IDENTIF. BOARD			
Size	Document Number	The present drawing is property of AMS any use without authorization will be prosecuted accordingly to the law	
A4	[MA]00446-00		
Date:	Saturday, January 12, 2002	Sheet	1 of 1
			Rev 0

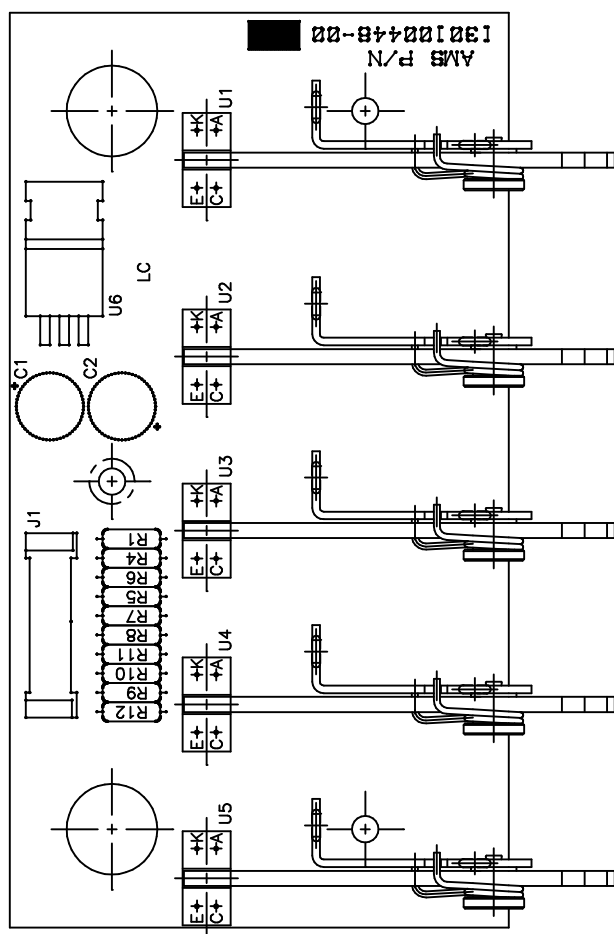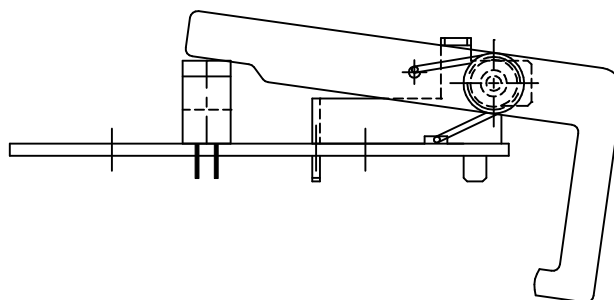


0	Emission				
Rev.	Description	Checked	Date	Approved	Date
 <b>Analyzer Medical System</b> Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
Description		Scale	Sheet	Checked	Date
SAMPLE RACKS			1 of 1	R.Cornacchia	11.01.2002
IDENTIFICATION BOARD		Car. [30]	Drawing 00446-00	Rev. 0	Date 11.01.2002
				Approved A.Gagliarducci	Date 11.01.2002





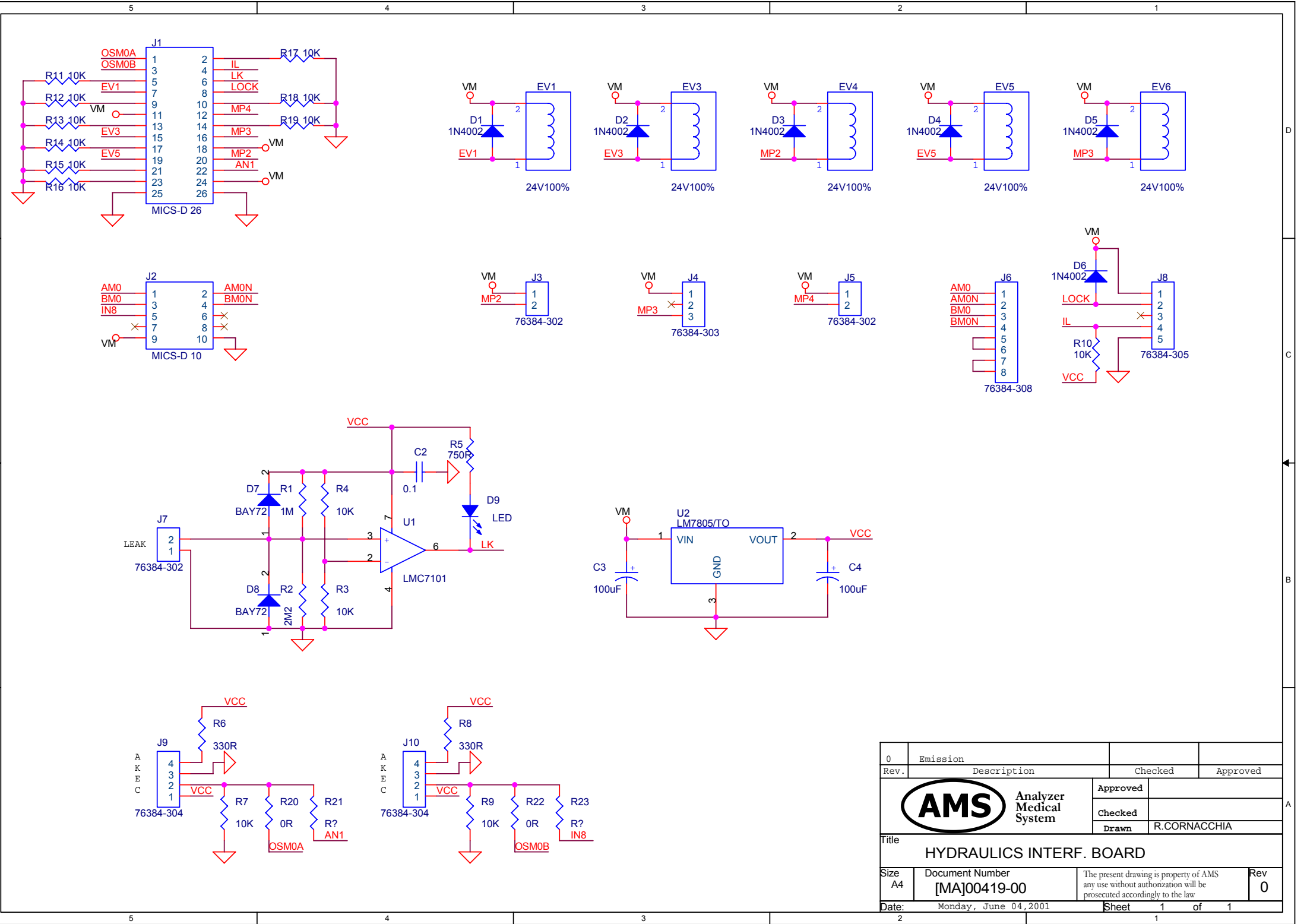

0	Emission		
Rev.	Description	Checked	Approved
<div>AMS</div> <div>Analyzer Medical System</div>		Approved	
		Checked	
		Drawn	R.CORNACCHIA
Title			
HYDRAULICS INTERF. BOARD			
Size	Document Number	The present drawing is property of AMS any use without authorization will be prosecuted accordingly to the law	
A4	[MA]00419-00		
Date:	Monday, June 04, 2001	Sheet	1 of 1

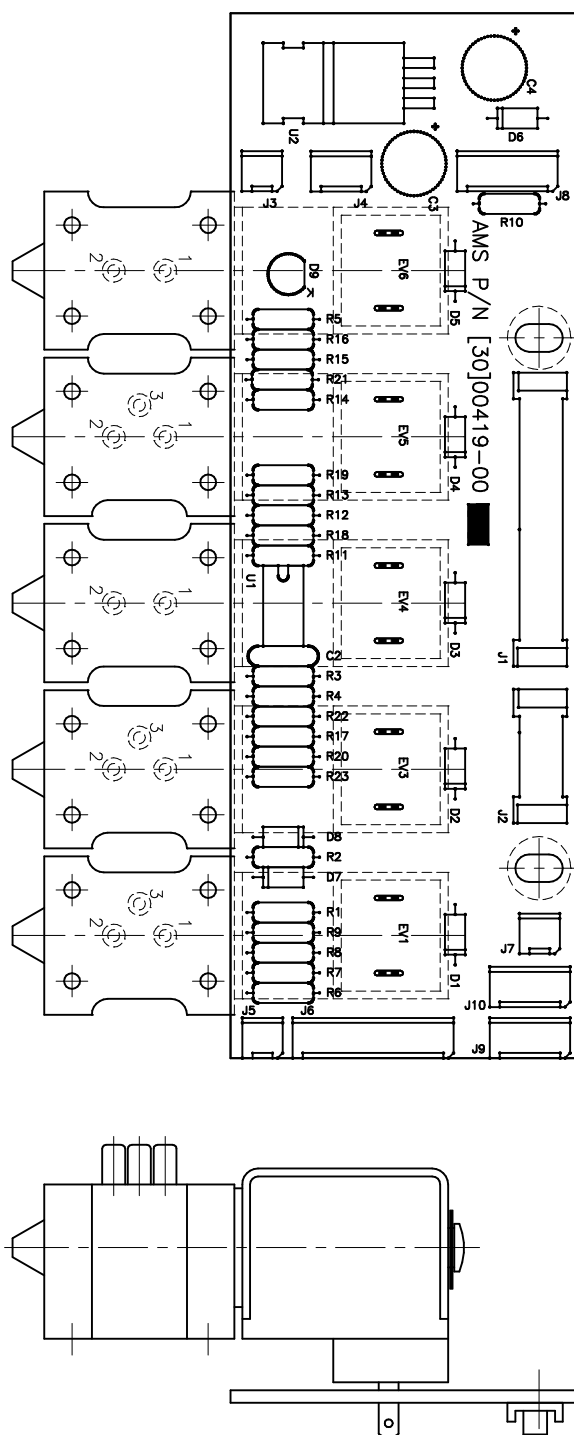

0	Emission				
Rev.	Description	Checked	Date	Approved	Date
<b>AMS</b> Analyzer Medical System Via Forlanini s.n.c. (Via Tiburtina km. 18,300) 00012 Guidonia (Roma)		The present drawing is property of AMS, any use without authorization will be prosecuted accordingly to the law.			
Description HYDRAULICS INTERFACE BOARD		Drawn G.Pucci		Date 23.11.2001	
		Scale 1:1		Sheet 1 of 1	
		Car. [MA]		Rev. 0	
		Drawing 00419-00		Approved A.G.	
				Date 23.11.2001	

## CHAPTER 05

---

### - DIAGNOSTIC PROGRAM -

### INDEX

5 DIAGNOSTIC PROGRAM	3
5.1 STATUS BAR	7
5.2 TEST FOLDER	8
5.3 “ARM” FOLDER	9
5.3.1 VERTICAL COMMAND AREA	10
5.3.2 HORIZONTAL COMMAND AREA	10
5.3.3 STATUS COMMAND AREA	11
5.3.4 ARM SETTING COMMAND AREA	11
5.4 “DILUTER” FOLDER	12
5.4.1 OPERATIONS AREA	13
5.4.2 SENSOR STATE AREA	13
5.4.3 USER AREA	13
5.4.4 STATUS COMMAND AREA	13
5.5 “OPTIC” FOLDER	14
5.5.1 OPERATIONS AREA	15
5.5.2 SENSOR STATE AREA	15
5.5.3 FILTERS AREA	15
5.5.4 DIGITAL CONVERSION AREA	16
5.5.5 STATUS COMMAND AREA	16
5.6 “PLATE” FOLDER	17
5.6.1 OPERATIONS AREA	18
5.6.2 SENSOR STATE AREA	18
5.6.3 USER AREA	18
5.6.4 MOTOR ON/OFF AREA	18
5.6.5 REACTION PLATE AREA	18
5.6.6 STATUS COMMAND AREA	18

5.7 “TEMPERATURES” FOLDER	19
5.7.1 TEMPERATURES AREA	19
5.7.2 OFFSET AREA	20
5.7.3 TEMPERATURE TEST AREA	20
5.7.4 STATUS COMMAND AREA	21
5.8 “WASH” FOLDER	21
5.8.1 OPERATIONS AREA	22
5.8.2 SENSOR STATUS AREA	22
5.8.3 TEST AREA	22
5.8.4 STATUS COMMAND AREA	22
5.8.5 SENSOR STATE	22
5.9 “USER” FOLDER	23
5.9.1 USER TEST MICRO 1 AREA	24
5.9.2 USER TEST MICRO 2 AREA	24
5.9.3 STATUS COMMAND AREA	24
5.10 “MISCELLANEOUS” FOLDER	25
5.10.1 GENERAL RESET AREA	26
5.10.2 SENSOR CHECKING AREA	26
5.10.3 LOG RAW DATA AREA	26
5.10.4 CONTINUOUS TEST RUN	26
5.10.5 SENSOR STATUS AREA	27
5.10.6 STATUS COMMAND AREA	27
5.11 “GLOBAL SETTING” FOLDER	28
5.11.1 TABLE POSITION	29
5.11.2 KEYBOARD	29
5.11.3 DOWN ARM	30
5.11.4 SETTING OF	30
5.12 “MACRO SETTING” FOLDER	32
5.13 “CONFIGURATION” FOLDER	33
5.14 “ISE” FOLDER	37




## 5 DIAGNOSTIC PROGRAM

The diagnostic program enables the operator to perform a complete check of each of the ILab 300 Plus' module functions.

This program has a folder structure, with each folder containing functions pertaining to the specific module. To launch the program the operator has to click with the left side of the mouse, on the **Diagnostic** area located on the lower right side of the screen of the System Monitor (Fig.1).

The ILab 300 Plus must be in the stand-by state to access this area.

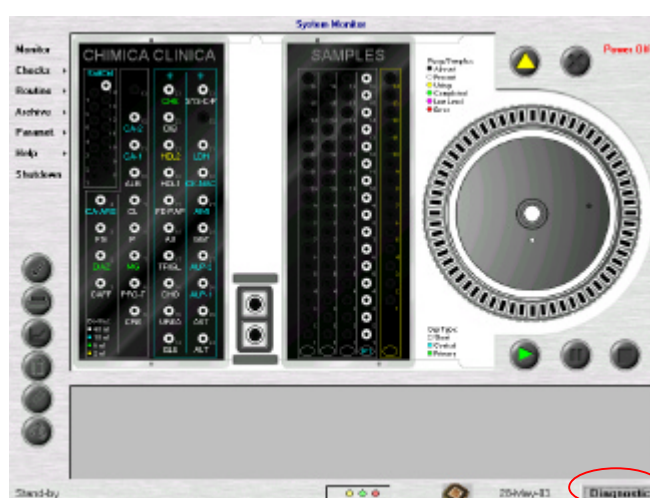
Fig. 1



Click Here

Once the program starts, the following screen appears (Fig. 2)

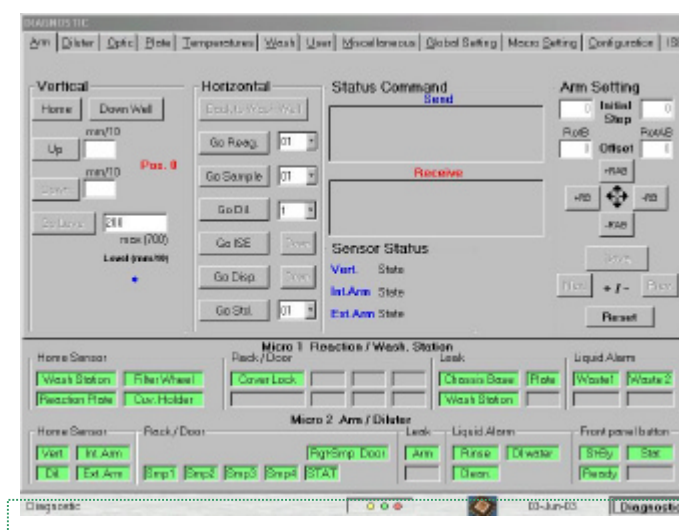Fig. 2



Test Folders

Optical Switch Sensor Status (O.S.S.S.)

Status Bar


The diagnostic program is subdivided into three distinct areas: test folders, optical switch sensor status (O.S.S.S.), status bar. The last two areas are present on every diagnostic window.

The three areas enable the operator to verify multiple functions as specified below:

**Test Folders:** checks the functionality of the ILab 300 Plus' various sub-systems; The individual functions are illustrated later on.

**Optical Switch Sensor Status (O.S.S.S.):** visualizes the sensor status pertaining to micro 1 that controls the reaction plate and the cuvette wash station (Fig. 3), as well as to micro 2 that manages the arm and the dilutor (Fig. 4)

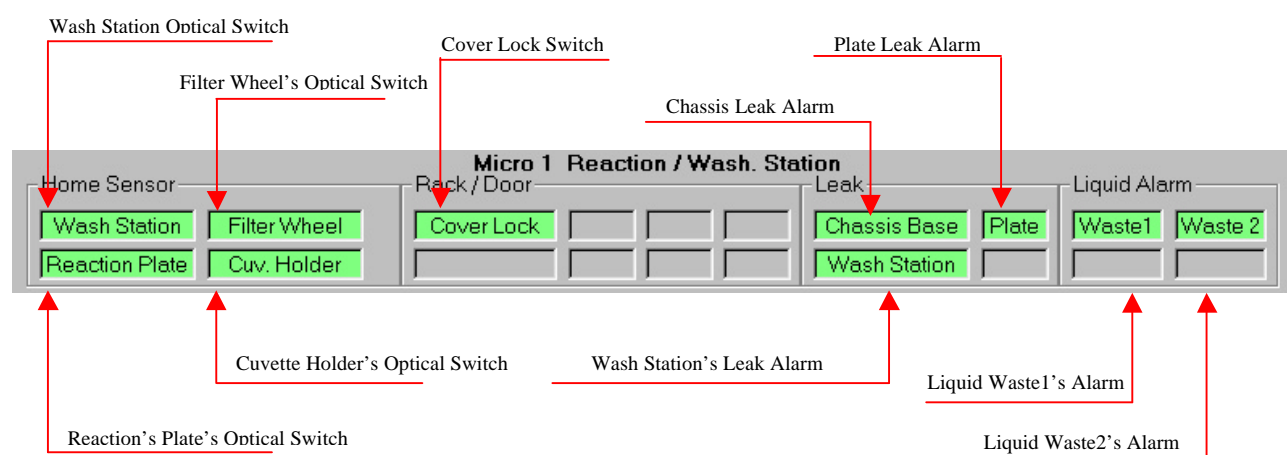Fig. 3



Furthermore, there are three types of unused fields: empty, inactive or those reserved for future use. While the program is running, the fields given in the figure can be either of two distinct colors, red or green.

The fields:

- **Wash Station's Optical Switch (O.S.)**
- **Filter Wheel's O.S.**
- **Reaction Plate's O.S.**
- **Cuvette Holder's O.S.**

turn green when the specific device is in the home position, and turn red when not in the home position.

The field

- **Cover Lock**

turns green when the cover is closed, and turns red when open.

The fields:

- **Chassis Leak Alarm** (Instrument Leaking)
- **Wash Station Leak Alarm** (Wash Station Leaking)
- **Plate Leak Alarm** (Reaction Plate Leaking)


turn green when there is no leaking, and red when there is.

The following two fields signal the alarms for the bottles containing the liquid waste material:

- **Liquid Waste1 Alarm**
- **Liquid Waste2 Alarm**

turning green when the liquid level is not excessive, and red when the waste bottle indicated by the alarm is almost full.

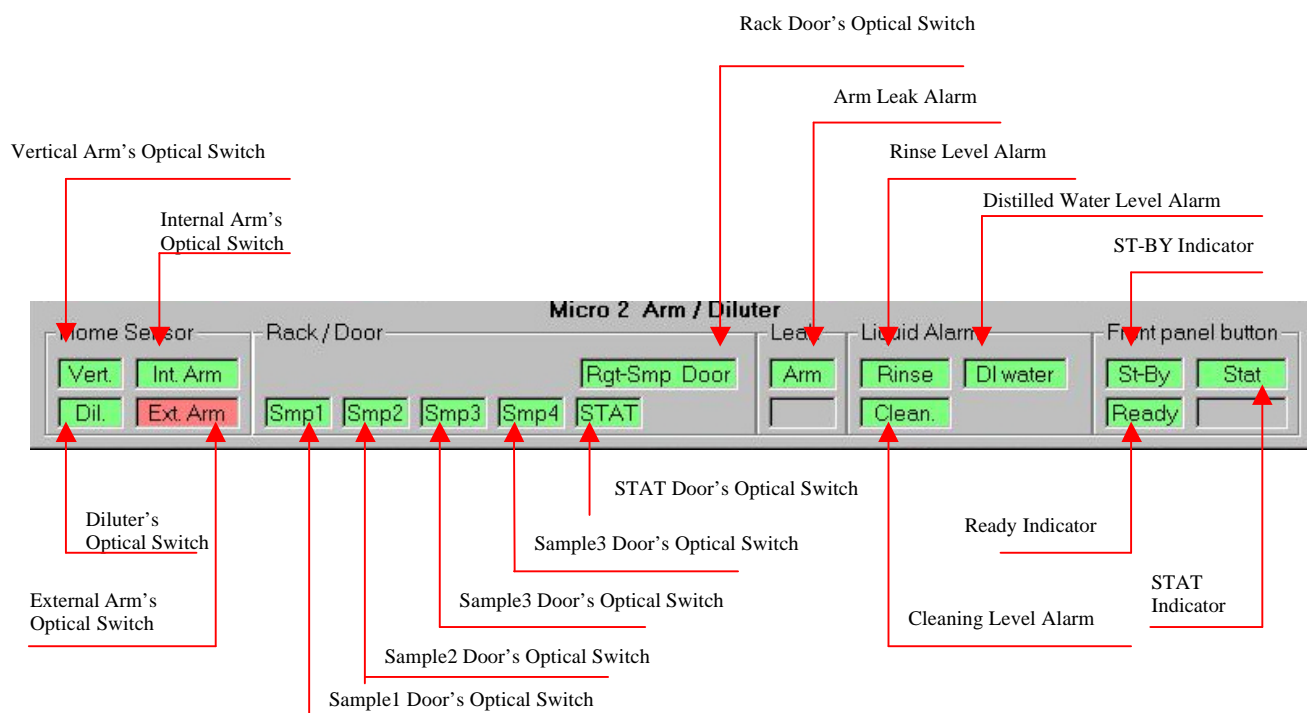Fig. 4



The fields in the figure can be two different colors: either green or red (for the indicators in Home Sensor, Leak, Liquid Alarm and front panel button), and either green or dark gray (for the indicators in Rack/Door except for the field Rgt-Smp Door that can be either green or red).

The fields:

- **Vertical Arm's O.S.**
- **Internal Arm's O.S.**
- **External Arm's O.S.**
- **Diluter O.S.**

turn green when the relative device is in the Home position, and turn red when the relative device is not in the Home position.

The following fields signal the presence of the rack in the appropriate home position (reagents, samples, standard and controls):

- **Sample1 Rack's O.S.**
- **Sample2 Rack's O.S.**
- **Sample3 Rack's O.S.**
- **Sample4 Rack's O.S.**
- **Sample 5 Rack's O.S. (STAT)**

turn green when the rack is in the appropriate home position, and turn dark gray when the rack is not present (not inserted).

The field below signals the general door opening of any of the racks:

- **Rack Door's O.S.**

turns red when the door is open, and green when the door is closed.

The field below signals the ILab 300 Plus's arm leak alarm:

- **Arm Leak Alarm**

turns green if there is no leaking, and red if there is.

The following fields check the liquid level of the distribution/loading bottles:

- **Rinse Level Alarm**
- **Distilled Water Level Alarm**
- **Cleaning Level Alarm**

turning green if the liquid level is above the predetermined minimum limit, and red if the liquid quantity is below the predetermined minimum level (the bottle is almost empty).

The following fields signal the state of the button on the ILab 300 Plus's operator panel:


- **ST-BY Indicator**
- **STAT Indicator**
- **READY Indicator**

turning green when the button is not pushed in, and red if it is.

## 5.1 STATUS BAR

Monitors the instrument's status while the program is running (Fig.5).

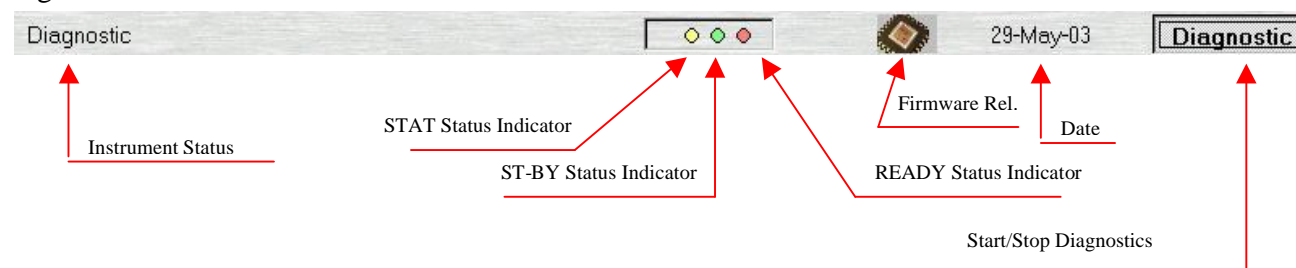Fig. 5



## 5.2 TEST FOLDER

The following is a summary of each individual test folder function.

The ILab 300 Plus's diagnostic program is subdivided into 12 folders:

<b>Arm:</b>	checks if the arm is functioning correctly;
<b>Diluter:</b>	checks if the diluter is functioning correctly;
<b>Optic:</b>	checks if the optical block is functioning correctly;
<b>Plate:</b>	checks if the reaction plate is functioning correctly;
<b>Temperature:</b>	checks the temperature of the pre-heater and reaction plate;
<b>Wash:</b>	checks if the wash station is functioning correctly;
<b>User:</b>	checks if the user functions are working properly;
<b>Miscellaneous:</b>	performs general checks;
<b>Global Setting:</b>	the adjustment of the arm position can be checked;
<b>Macro Setting:</b>	the operator can identify which macros have been loaded;
<b>Configuration:</b>	Through the use of a password, only qualified technical personal can adjust certain parameters (This file is not accessible to the Laboratory operator).
<b>ISE:</b>	Checks if the ISE Module is functioning properly and allows the operator to do maintenance or repair operations.

### 5.3 “ARM” FOLDER

**DIAGNOSTIC**

Arm Diluter Optic Plate Temperatures Wash User Miscellaneous Global Setting Macro Setting Configuration ISE

**Vertical**

Home Down Well

mm/10

Up

mm/10

Down

Go Level 200 max (700)

Level (mm/10)

**Horizontal**

Back to Wash Well

Go Reag. 01

Go Sample 01

Go Dil. 1

Go ISE Down

Go Disp. Down

Go Std. 01

**Status Command**

Send

Receive

**Arm Setting**

Initial Step 0

RotB 0

Offset 0

RotAB 0

+RAB

+RB

-RB

-RAB

Save

Next + / - Prev.

Reset

**Sensor Status**

Vert. State

Int.Arm State

Ext.Arm State

**Micro 1 Reaction / Wash. Station**

Home Sensor

Wash Station Filter Wheel

Reaction Plate Cuv. Holder

Rack / Door

Cover Lock

Leak

Chassis Base Plate

Wash Station

Liquid Alarm

Waste1 Waste2

**Micro 2 Arm / Diluter**

Home Sensor

Vert. Int. Arm

Dil. Ext. Arm

Rack / Door

Smp1 Smp2 Smp3 Smp4 STAT

Rgt-Smp Door

Leak

Arm

Liquid Alarm

Rinse DI water

Clean.

Front panel button

St-By Stat

Ready

Diagnostic

03-Jun-03

Diagnostic

**Warning:** the improper use of the functions described in this folder can damage the sampling probe.

The Arm File is subdivided into four areas:

- ◆ **VERTICAL**
- ◆ **HORIZONTAL**
- ◆ **STATUS COMMAND**
- ◆ **ARM SETTING**

### 5.3.1 VERTICAL COMMAND AREA

(5 commands)

- Home:** Brings the arm vertically to a home position;
- Down Well:** Brings the arm's axis "z" to the height of the well;
- Up:** Raises the arm tenths of a millimeter as preset in the square found immediately on the right;
- Down:** Lowers the arm tenths of a millimeter as pre-determined by the square found immediately on the right;
- Go Level:** Brings the probe to either the reagent level, the sample or the standard and reads the level. Before giving this command, it is necessary to set the parameter for the probe's maximum lowering level in the max area (mm);

### 5.3.2 HORIZONTAL COMMAND AREA

(7 commands)

- Back to wash well:** Brings the arm back to the home position;
- Go Reag.:** Brings the arm to the reagent position specified in the sheet menu next to the command;
- Go Sample:** Brings the arm to the sample position specified in the sheet menu next to the command;
- Go Dil.:** Brings the arm to the dilution position specified in the sheet menu next to the command;
- Go ISE:** Brings the arm to the ISE position;
- Go Disp.:** Brings the arm to the dispense position;
- Go Std.:** Brings the arm to the standard/control position specified in the sheet menu next to the command;



### 5.3.3 STATUS COMMAND AREA

<b>Send:</b>	Visualizes the commands sent from the program to the instrument hardware;
<b>Receive:</b>	Visualizes the answers of the instrument's hardware to the commands sent by the program;
<b>Sensor State:</b>	Indicates the sensor status relative to the arm's vertical axis, the internal arm and the external arm.

### 5.3.4 ARM SETTING COMMAND AREA

(8 commands)

<b>+RAB:</b>	Performs the clockwise rotation of the two arms (to adjust the position);
<b>-RAB:</b>	Performs the counter-clockwise rotation of the two arms (to adjust the position);
<b>+RB:</b>	Performs the clockwise rotation of the external arm (to adjust the position);
<b>-RB:</b>	Performs the counter-clockwise rotation of the external arm (to adjust the position);
<b>Save:</b>	Saves the adjustment data;
<b>Next:</b>	Brings the arm to the next sample, standard, dilution or reagent position
<b>Prev.:</b>	Brings the arm to the previous sample, standard, dilution or reagent position
<b>Reset:</b>	Completely resets the arm

In this area the boxes can be activated when the arm is positioned on a sample, a reagent, a standard or a control:

<b>Original step:</b>	Visualizes the number of base steps.
<b>Offset:</b>	Visualizes the number of offset steps.

## 5.4 “DILUTER” FOLDER

**DIAGNOSTIC**

Arm **Diluter** Optic Plate Temperatures Wash User Miscellaneous Global Setting Macro Setting Configuration ISE

**Operations**

Home

1 ul

Asp Disp

**Sensor Status**

Diluter State

**User**

Probe Wash Test

**Status Command**

Send

Receive

**Micro 1 Reaction / Wash. Station**

Home Sensor: Wash Station, Filter Wheel, Reaction Plate, Cuv. Holder

Rack / Door: Cover Lock

Leak: Chassis Base, Plate, Wash Station

Liquid Alarm: Waste1, Waste2

**Micro 2 Arm / Diluter**

Home Sensor: Vert, Int. Arm, Dil, Ext. Arm

Rack / Door: Smp1, Smp2, Smp3, Smp4, STAT

Leak: Arm

Liquid Alarm: Rinse, DI water, Clean

Front panel button: St-By, Stat, Ready

The Diluter Folder is subdivided into four areas:

- ◆ **OPERATIONS**
- ◆ **SENSOR STATE**
- ◆ **USER**
- ◆ **STATUS COMMAND**

### 5.4.1 OPERATIONS AREA

(3 commands)

- Home:** Brings the diluter to the Home position;
- Asp:** Aspirates the micro-litres as predetermined in the ul. box (may be set manually through + and – or by using the cursor located immediately above);
- Disp.:** Distributes the micro-litres as predetermined in the ul. box (may be set manually through + and – or by using the cursor located immediately above);

### 5.4.2 SENSOR STATE AREA

In the Sensor State area the diluter sensor status is indicated.

### 5.4.3 USER AREA

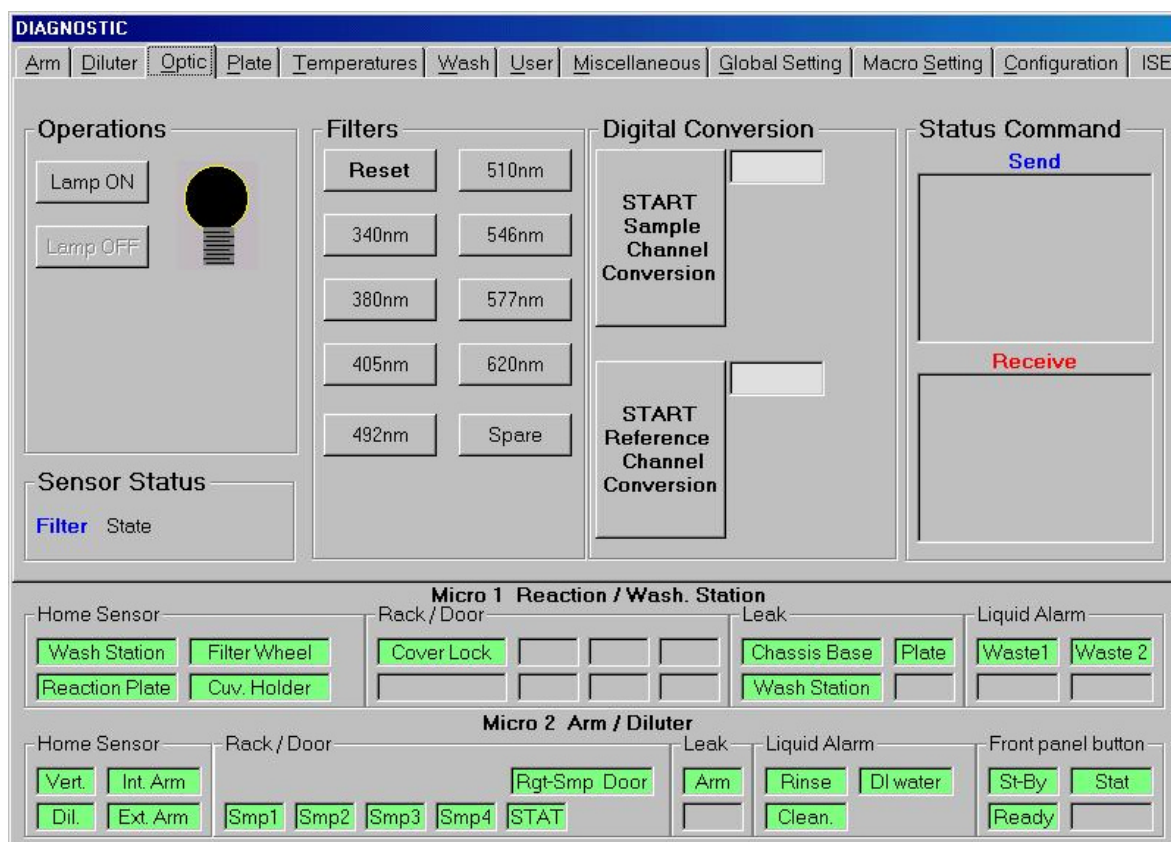
(1 command)

Probe Wash Test: Activates and deactivates the probe wash system

### 5.4.4 STATUS COMMAND AREA

- Send:** Visualizes the commands sent by the program to the instrument hardware;
- Receive:** Visualizes the instrument hardware's response to the commands sent by the program.

## 5.5 “OPTIC” FOLDER



The Optic folder is subdivided into five areas:

- ◆ **OPERATIONS**
- ◆ **SENSOR STATUS**
- ◆ **FILTERS**
- ◆ **DIGITAL CONVERSION**
- ◆ **STATUS COMMAND**

### 5.5.1 OPERATIONS AREA

(2 commands)

**Lamp ON:** Turns on the Optic Lamp;

**Lamp OFF:** Turns off the Optic Lamp;

### 5.5.2 SENSOR STATE AREA

This area indicates the status of the wheel filters.

### 5.5.3 FILTERS AREA

(10 commands)

**Reset :** Positions the wheel filters in home position (dark);

**340 nm:** Positions filter n° 1 in front of the reading sensor;

**380 nm:** Positions filter n° 2 in front of the reading sensor;

**405 nm:** Positions filter n° 3 in front of the reading sensor;;

**492 nm:** Positions filter n° 4 in front of the reading sensor;

**510 nm:** Positions filter n° 5 in front of the reading sensor;

**546 nm:** Positions filter n° 6 in front of the reading sensor;

**577 nm:** Positions filter n° 7 in front of the reading sensor;

**620 nm:** Positions filter n° 8 in front of the reading sensor;

**Spare:** Positions filter n° 9 in front of the reading sensor (Optional);

#### 5.5.4 DIGITAL CONVERSION AREA

(2 commands)

**START/STOP Sample Channel Conversion:** Performs an analogical/digital conversion of the principal channel;

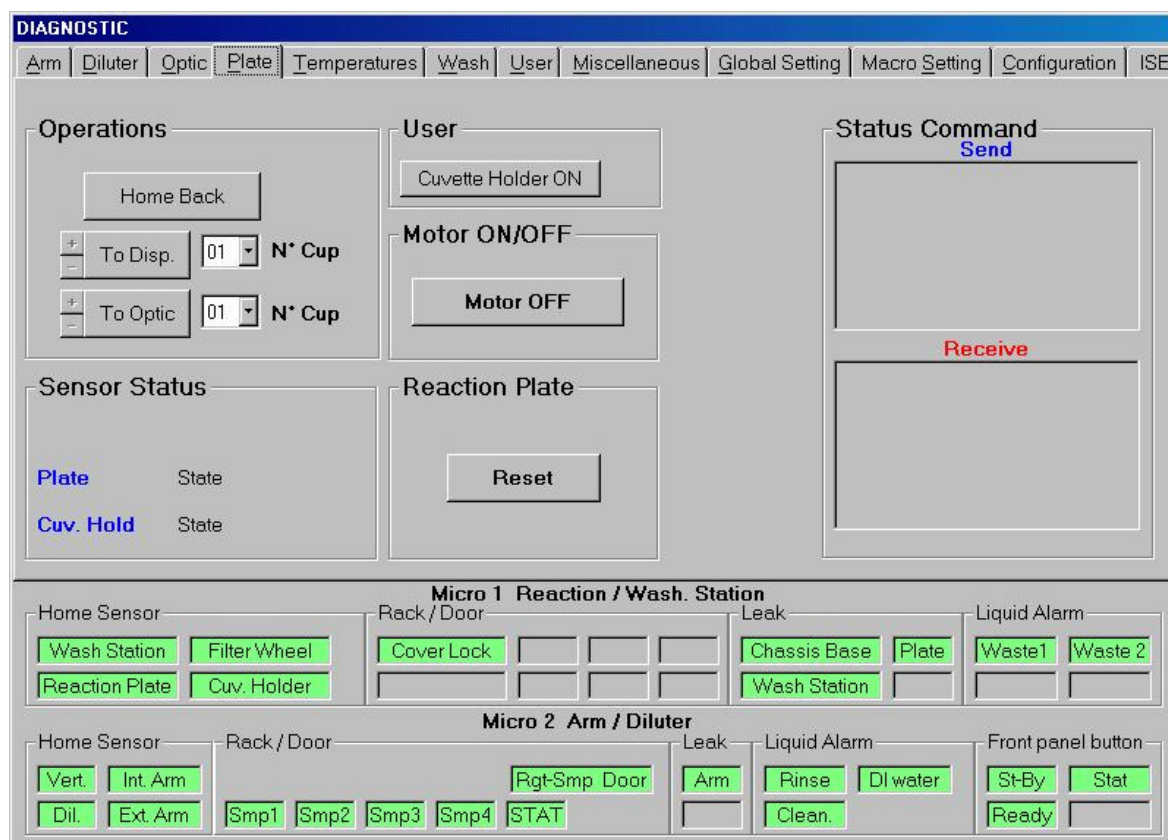
**START/STOP Reference Channel Conversion:** Performs an analogical/digital conversion of the reference channel;

#### 5.5.5 STATUS COMMAND AREA

**Send:** Visualizes the commands sent by the program to the instrument hardware;

**Receive:** Visualizes the instrument hardware's response to the commands sent by the program.

## 5.6 “PLATE” FOLDER



The Plate Folder is subdivided into six areas:

- ◆ **OPERATIONS**
- ◆ **SENSOR STATUS**
- ◆ **USER**
- ◆ **MOTOR ON/OFF**
- ◆ **REACTION PLATE**
- ◆ **STATUS COMMAND**

### 5.6.1 OPERATIONS AREA

(3 commands)

- Home Back:** Performs the same movement made previously, but in the opposite direction
- To Disp.:** Brings the cuvette, indicated on the sheet menu, below the dispense position
- To Optic:** Brings the cuvette, indicated on the sheet menu, to the colorimeter.

### 5.6.2 SENSOR STATE AREA

The status of the plate sensor is indicated in this area.

### 5.6.3 USER AREA

(1 command)

**Cuvette Holder ON/OFF:** Alternately turns ON and OFF the cuvette holder.

### 5.6.4 MOTOR ON/OFF AREA

(1 command)

**Motor ON/OFF:** Engages/Disengages the motors to allow a manual movement of the plate.

### 5.6.5 REACTION PLATE AREA

(1 command)

**Reset:** Performs the resetting of the reaction plate bringing the n°1 cuvette in the distribution position (reaction disp.).

### 5.6.6 STATUS COMMAND AREA

**Send:** Visualizes the commands sent by the program to the instrument hardware;

**Receive:** Visualizes the instrument hardware's response to the commands sent by the program.



## 5.7 “TEMPERATURES” FOLDER

**DIAGNOSTIC**

Arm | Diluter | Optic | Plate | **Temperatures** | Wash | User | Miscellaneous | Global Setting | Macro Setting | Configuration | ISE

### Temperatures

Temp. Test Run

Calib. Sensor

Plate 0 °C

Arm 0 °C

Offset 0 °C

Service Pb. 0 °C

Status Command

Send

Receive

Micro 1 Reaction / Wash. Station

Home Sensor

Wash Station

Filter Wheel

Reaction Plate

Cuv. Holder

Rack / Door

Cover Lock

Leak

Chassis Base

Plate

Wash Station

Liquid Alarm

Waste1

Waste2

Micro 2 Arm / Diluter

Home Sensor

Vert.

Int. Arm

Dil.

Ext. Arm

Rack / Door

Smp1

Smp2

Smp3

Smp4

STAT

Leak

Arm

Rinse

DI water

Clean.

Liquid Alarm

St-By

Stat

Ready

Front panel button

The Temperature Folder is subdivided into three areas:

**Temperatures**

**Offset**

**Status Command**

### 5.7.1 TEMPERATURES AREA

(4 commands)

- Temp. Test Run:** To enter in the Temperature Test folder
- Plate:** Performs a temperature reading in the reaction plate compartment;
- Arm:** Performs a temperature reading in the arm's pre-heater.
- Calib. Sensor:** To carry out the calibration procedure of the thermometric probe included in the kit P/N 23550120400, as illustrated in Chapter Six “Settings and Adjustments”.

### 5.7.2 OFFSET AREA

(2 commands)

To correlate the thermometric probe included in the kit P/N 23550120400 with another reference thermometric probe.

**Service probe:** The offset value that has been previously inserted and the temperature value read by the thermometric probe appear in the closed fields.

**Cursore:** To enter the offset value  $\pm 3\text{ }^{\circ}\text{C}$ . The offset value will be zero after exiting from the folder.

### 5.7.3 TEMPERATURE TEST AREA

(3 commands)

**START:** To perform the test check and to regulate the temperature of the preheater and of the reaction plate by using the kit P/N 23550120400, as illustrated in Chapter 6 “Setting and Adjustments”.

The number of sample cycles are predefined, whereas the reagent and sample positions for test execution can be selected by the operator.

**Service Probe:** the field in which the temperature of the service probe appears.

**WASH:** To run a wash cycle on the cuvettes used for testing.

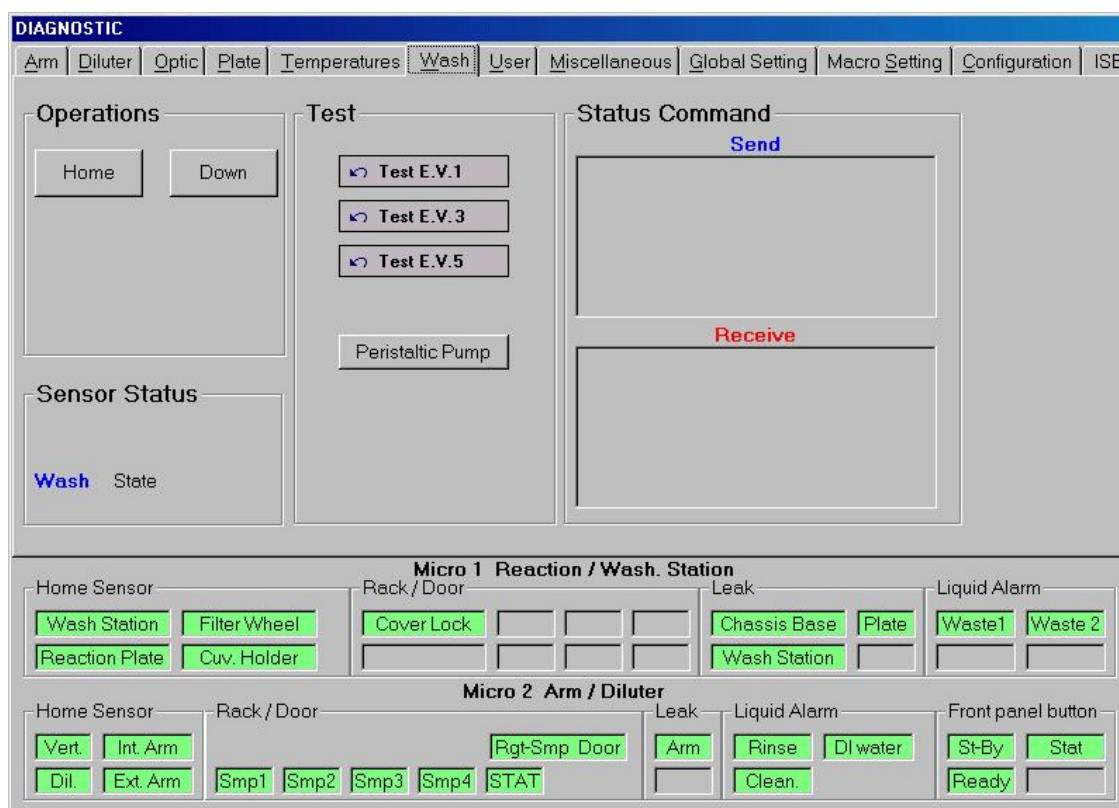
**EXIT:** To exit from the Temperature Test box and to run a wash cycle on the cuvettes if they have been used.

#### 5.7.4 STATUS COMMAND AREA

**Send:** Visualizes the commands sent by the program to the instrument hardware;

**Receive:** Visualizes the instrument hardware's response to the commands sent by the program.

### 5.8 “WASH” FOLDER



The Wash folder is subdivided into four areas:

- ◆ **OPERATIONS**
- ◆ **SENSOR STATUS**
- ◆ **TEST**
- ◆ **STATUS COMMAND**

### 5.8.1 OPERATIONS AREA

(2 commands)

**Home:** Brings the washing device to the Home position;

**Down:** Makes the washing device go down to the washing position.

### 5.8.2 SENSOR STATUS AREA

In this area the sensor state of the wash station is indicated.

### 5.8.3 TEST AREA

(4 commands)

**Test E.V. 1:** Performs the solenoid valve 1 test;

**Test E.V. 3:** Performs the solenoid valve 3 test;

**Test E.V. 5:** Performs the solenoid valve 5 test;

**Peristaltic Pump:** Starts the peristaltic pump, and automatically turns off after a period of tiis predetermined by the user.

### 5.8.4 STATUS COMMAND AREA

**Send:** Visualizes the commands sent by the program to the instrument hardware;

**Receive:** Visualizes the instrument hardware's response to the commands sent by the program.

### 5.8.5 SENSOR STATE

♦ **TEST**

♦ **STATUS COMMAND**

## 5.9 “USER” FOLDER

DIAGNOSTIC											
Arm	Diluter	Optic	Plate	Temperatures	Wash	User	Miscellaneous	Global Setting	Macro Setting	Configuration	ISE
<b>User Test Micro1</b> Main Cover Lock ON			<b>User Test Micro2</b> Alarms ON ● READY ON ● Test BUZZER			<b>Status Command</b> Send Receive					
<b>Micro 1 Reaction / Wash. Station</b>											
Home Sensor		Rack / Door		Leak		Liquid Alarm					
Wash Station	Filter Wheel	Cover Lock				Chassis Base	Plate	Waste1	Waste 2		
Reaction Plate	Cuv. Holder					Wash Station					
<b>Micro 2 Arm / Diluter</b>											
Home Sensor		Rack / Door		Leak		Liquid Alarm		Front panel button			
Vert.	Int. Arm					Rgt-Smp Door	Arm	Rinse	DI water	St-By	Stat
Dil.	Ext. Arm	Smp1	Smp2	Smp3	Smp4	STAT		Clean.		Ready	

The User File is subdivided into three areas:

- ◆ **USER TEST MICRO 1**
- ◆ **USER TEST MICRO 2**
- ◆ **STATUS COMMAND**

### 5.9.1 USER TEST MICRO 1 AREA

(1 command)

**Main Cover Lock ON/OFF:** Alternatively turns on or off the solenoid valve holder and releases the main cover lock.

### 5.9.2 USER TEST MICRO 2 AREA

(3 commands)

**Alarms ON/OFF:** Turns on and off the STAT indicator lamp;

**Ready ON/OFF:** Turns on and off the READY indicator lamp;

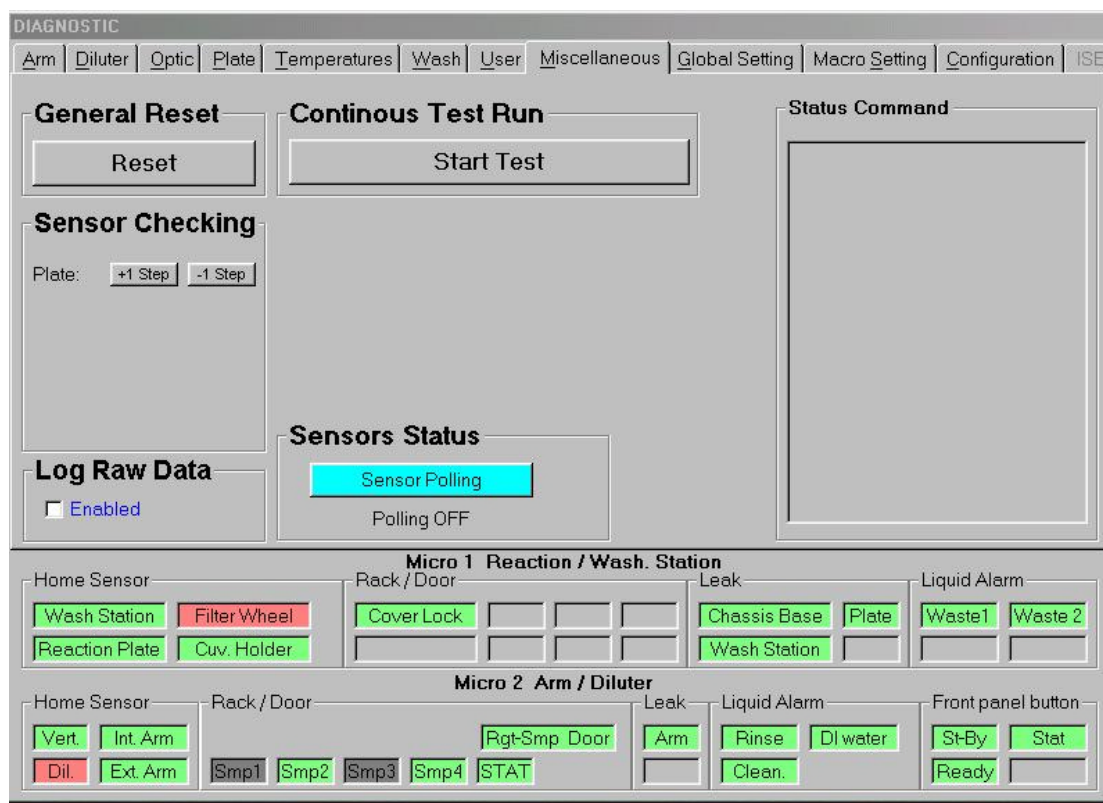
**Test BUZZER:** Performs a complete operational test of the buzzer;

### 5.9.3 STATUS COMMAND AREA

**Send:** Visualizes the commands sent by the program to the instrument hardware;

**Receive:** Visualizes the instrument hardware's response to the commands sent by the program.

## 5.10 “MISCELLANEOUS” FOLDER



The Miscellaneous folder is subdivided into six areas:

- **General Reset**
- **Sensor Checking**
- **Log Raw data**
- **Continuous Test Run**
- **Sensor Status**
- **Status Command**

### 5.10.1 GENERAL RESET AREA

(1 command)

**Reset:** Performs the resetting of the general system.

### 5.10.2 SENSOR CHECKING AREA

(2 commands)

**Plate + 1 Step:** To turn the reaction plate in a counter-clockwise direction to check the home sensor.

**Plate - 1 Step:** To turn the reaction plate in a clockwise direction to check the home sensor.

### 5.10.3 LOG RAW DATA AREA

(1 command)

**Enable:** To save a copy of the files for the technical diagnosis if the instrument malfunctions. The files are memorized in the folder “Log” inside the folder “Analyzer”.

### 5.10.4 CONTINUOUS TEST RUN

(1 command)

**Start Test:** Performs an automatic test during which the modules (Sampling Arm, Reaction plate, Washing Station and Photometer) are moved to all the different positions.

The correct modules positioning is also automatically verified and showed on the following table:





To stop the test push on the “Exit Test” key

### 5.10.5 SENSOR STATUS AREA

(1 command)

**Sensor Polling:** To choose either Polling ON or Polling OFF. When in the Polling ON position, it is possible to check the functionality of the optical filters in various positions.

### 5.10.6 STATUS COMMAND AREA

**Send:** Visualizes the commands sent by the program to the instrument hardware

**Receive:** Visualizes the instrument hardware’s response to the commands.

## 5.11 “GLOBAL SETTING” FOLDER

**DIAGNOSTIC**

Arm | Diluter | Optic | Plate | Temperature | Wash | User | Miscellaneous | **Global Setting** | Macro Setting | Configuration | ISE

Position	Rot A&B	Rot B	Offset A&B	Offset B
Sample 1	-736	1032		
Sample 2	-664	840		
Sample 3	-612	680		
Sample 4	-580	540		
Sample 5	-564	412		
Sample 6	-568	296		
Sample 7	-600	196		
Sample 8	-672	108		
Sample 9	-792	44		
Sample 10	-968	4		
Sample 11	-1184	4		
Sample 12	-1400	40		
Sample 13	-1588	104		
Sample 14	-1740	100		

Prev. ↑ Go ↓ Next

To Use these keys below, this led must be ON

Num Lock Caps Lock Scroll Lock End Exit

Num Lock	/	=	-	Prev.
7	8	9	↑	+
Home	←	5	→	Next
4	↓	6	→	Next
1	End	2000	3	↓
0	Ins	Save	Del	Enter

Setting of

☐ Automatic  
☒ Manual

☒ Sample  
☐ Reagent  
☐ Std/Ctrl  
☐ Dil, Ise, Disp

Refresh

**Micro 1 Reaction / Wash. Station**

Home Sensor: Wash Station, Filter Wheel, Reaction Plate, Cuv. Holder

Rack / Door: Cover Lock

Leak: Chassis Base, Plate, Wash Station

Liquid Alarm: Waste 1, Waste 2

**Micro 2 Arm / Diluter**

Home Sensor: Vert., Int. Arm, Dil., Ext. Arm

Rack / Door: Smp1, Smp2, Smp3, Smp4, STAT, Rgt-Smp Door

Leak: Arm

Liquid Alarm: Rinse, DI water, Clean

Front panel button: St-By, Stat, Ready

**Warning:** the improper use of the functions described in this folder can damage the sampling probe.

The Global Setting Folder is subdivided into four areas:

- ◆ **TABLE POSITION**
- ◆ **KEYBOARD**
- ◆ **DOWN ARM**
- ◆ **SETTING OF**

## 5.11.1 TABLE POSITION

This table has five columns:

<b>Position:</b>	Lists all the possible instrument positions (samples, reagents, standard, controls, diluents, ISE, dispense);
<b>Rot A&amp;B:</b>	Gives the base steps for the complete arm (internal and external) in every position;
<b>Rot B:</b>	Gives the base steps only for the external arm in every position;
<b>Offset A&amp;B:</b>	Gives the corrective steps necessary for the arm's adjustment (internal and external);
<b>Offset B:</b>	Gives the corrective steps necessary for the external arm's adjustment.

Next to the Sample Table there are 3 commands:

<b>Prev:</b>	Returns the arm to the previous position;
<b>Next:</b>	Brings the arm to the successive position;
<b>Go:</b>	Brings the arm to the first position according to the category selected in the area "Setting of" (Sample, reagent, Std/Ctrl, Dil, Ise, Disp).

## 5.11.2 KEYBOARD

The various positions can be improved through the keyboard, saving the operations performed.

The areas visible in Fig.6 may be activated by the mouse or with the help of the numerical buttons on the keyboard. The buttons are listed with their corresponding function:

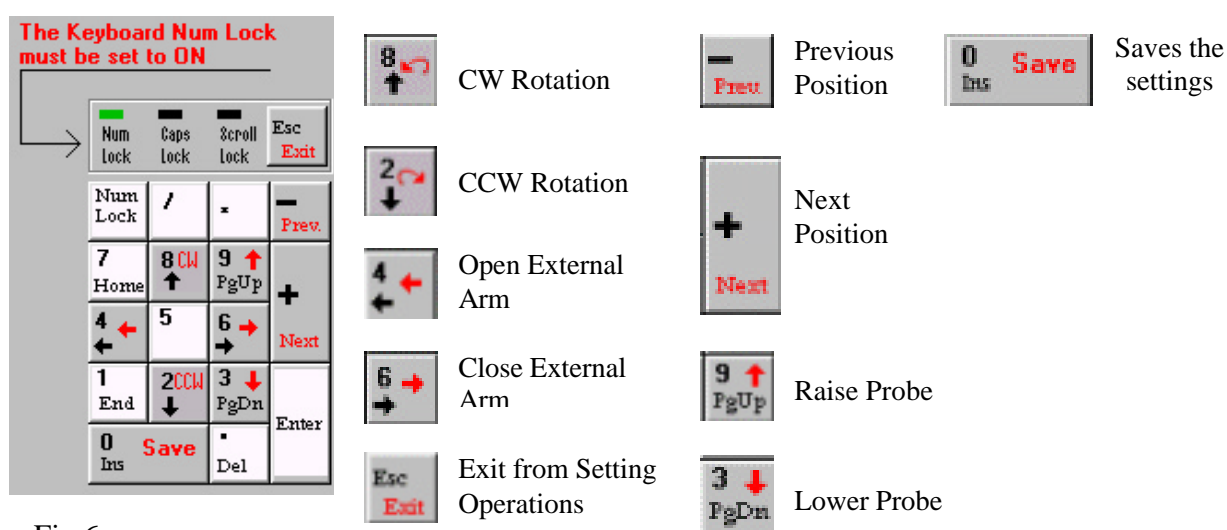


Fig.6

### 5.11.3 DOWN ARM

This area has the option to decide if the arm, once set in a determined position, should automatically command the probe to go down or wait for the manual command.

### 5.11.4 SETTING OF

This area has a selection of location categories that are used to set the functions (Sample, reagent, Std/Ctrl, Dil, Ise, Disp).

### Setting Procedure

- ⇒ Choose the position category in the “Setting of” Area;
- ⇒ Click with the left side of the mouse on the Go area (as soon as the area is clicked on it changes to End);
- ⇒ The ILab 300 Plus resets the arm, and the values of the Offset column appear on the table (Fig.7);

Fig.7

The screenshot displays the LIASYS Diagnostic software interface. At the top, there are tabs for various system components: Arm, Diluter, Optic, Plate, Temperature, Wash, User, Miscellaneous, Global Setting, Macro Setting, and Configuration. The 'Global Setting' tab is currently selected.

Under the 'Global Setting' tab, there are two main sections: 'Sample 1' and 'Setting of'.

The 'Sample 1' section contains a table with the following data:

Position	Rot A&B	Rot B	Offset A&B	Offset B
Sample 1	-736	1032	2	0
Sample 2	-664	840	-1	5
Sample 3	-612	680	-2	6
Sample 4	-580	540	0	7
Sample 5	-564	412	-5	11
Sample 6	-568	296	-3	11
Sample 7	-600	196	1	6
Sample 8	-672	108	-2	11
Sample 9	-792	44	-7	17
Sample 10	-968	4	-10	18
Sample 11	-1184	4	-4	27
Sample 12	-1400	40	11	-6
Sample 13	-1588	104	15	-10
Sample 14	1740	100	14	5

Below the table are 'Prev.' and 'Next' buttons. To the right of the table, there is a red text box that says 'The Keyboard Num Lock must be set to ON'.

The 'Setting of' section contains a numeric keypad with various function keys. The 'Num Lock' key is highlighted in green. Below the keypad, there are two radio buttons: 'Automatic' and 'Manual'. The 'Manual' radio button is selected.

At the bottom of the interface, there are two main sections: 'Micro 1 Reaction / Wash. Station' and 'Micro 2 Arm / Diluter'. Each section contains several status indicators and buttons.

**Micro 1 Reaction / Wash. Station:**

- Home Sensor: Wash Station (green), Filter Wheel (red), Reaction Plate (red), Cuv. Holder (green).
- Rack / Door: Cover Lock (red).
- Leak: Chassis Base (green), Plate (green), Wash Station (green).
- Liquid Alarm: Waste 1 (green), Waste 2 (green).

**Micro 2 Arm / Diluter:**

- Home Sensor: Vert. (green), Int. Arm (green), Dil. (red), Ext. Arm (green).
- Rack / Door: RGT1 (green), RGT2 (green), RGT3 (green), RGT4 (green), Rgt-Smp Door (green), Smp1 (green), Smp2 (green), Smp3 (green), Smp4 (green), STAT (green).
- Leak: Arm (green).
- Liquid Alarm: Rinse (red), DI water (green), Clean. (red).
- Front panel button: St-By (green), Stat (green), Ready (green).

The bottom status bar shows 'Diagnostic' on the left, a row of three colored circles (yellow, green, red) in the center, and '29/11/2001' and 'Diagnostic' on the right.

⇒ The arm centering can be adjusted at the highlighted area by giving commands through the keyboard or mouse:

8 = CW Rotation

2 = CCW Rotation

4 = Open External Arm

6 = Close External Arm

+ = Next Position

- = Previous Position

⇒ Repeat the centering for all the relevant positions;

⇒ Push the Save button to save the operations performed;

⇒ Wait for the arm to automatically reset.

**NOTE:** If the option Manual (default) is chosen in the Down Arm area, to make the arm go up and down it is possible to use either the buttons pgUp and pgDn or click with the mouse on the button with the arm figure in the Down Arm Area.

**WARNING:**

**THE FIRST TIME THAT THE SETTING PROCEDURE IS PERFORMED, IT IS EXTREMELY IMPORTANT TO CHOOSE THE MANUAL OPTION IN THE DOWN ARM AREA.**

**BY CHOOSING THE “AUTOMATIC” OPTION THERE IS A CONSIDERABLE RISK OF BREAKING THE PROBE DUE TO POSITIONS THAT HAVE NOT BEEN PROPERLY CENTERED.**

## 5.12 “MACRO SETTING” FOLDER

**DIAGNOSTIC**

Arm | Diluter | Optic | Plate | Temperatures | Wash | User | Miscellaneous | Global Setting | **Macro Setting** | Configuration | ISE

Read Macro: 01 ☐ Micro1 ☐ Micro2

[Empty text area for macro data]

**Micro 1 Reaction / Wash. Station**

Home Sensor	Rack / Door	Leak	Liquid Alarm
Wash Station	Cover Lock	Chassis Base	Waste1
Filter Wheel		Plate	Waste 2
Reaction Plate		Wash Station	
Cuv. Holder			

**Micro 2 Arm / Diluter**

Home Sensor	Rack / Door	Leak	Liquid Alarm	Front panel button
Vert.		Arm	Rinse	St-By
Int. Arm			DI water	Stat
Dil.	Smp1		Clean.	Ready
Ext. Arm	Smp2			
	Smp3			
	Smp4			
	STAT			

This folder allows the user to read the macro present in the ILab 300 Plus.

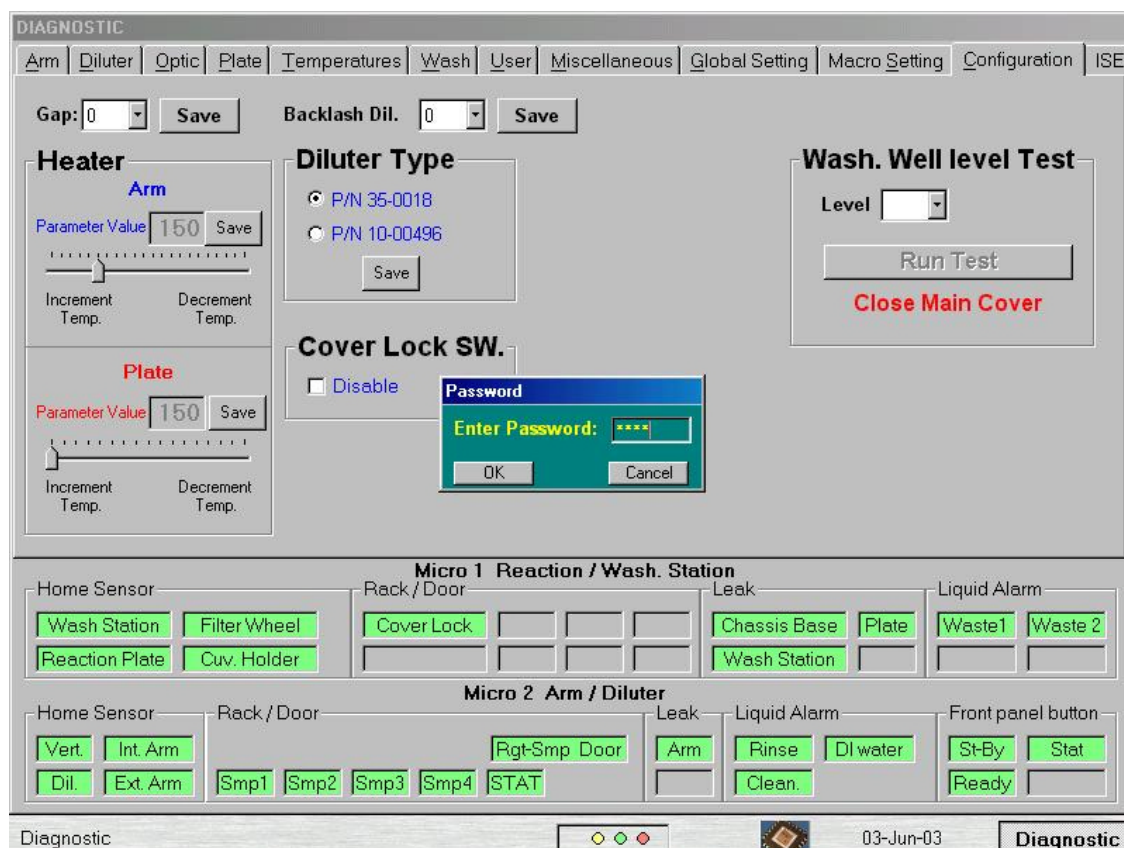
To read the macro follow the below procedure:

- ⇒ Select the micro from which the information will be asked;
- ⇒ Select from the menu sheet “Read Macro” the number of relevant macro;
- ⇒ Click with the left side of the mouse on the Read area;
- ⇒ If the Macro requested exists, the information will appear in the designated area.

### 5.13 “CONFIGURATION” FOLDER

This folder may be accessed only by qualified technical assistant.

It is necessary for the technical assistant to insert a password to access this folder.



The Configuration Folder is subdivided into six areas:

- **Gap**
- **Backlash Dil.**
- **Heater**
- **Diluter Type**
- **Cover Lock SW.**
- **Wash. Well level Test**

**Gap:** defines the volume, in  $\mu\text{L}$ , of the air bubbles that separate the liquids in the sample column; Typical values are between 10 to 15. Big air bubbles produce the best separation among sample, reagent and rinse (Better BIB values). Small air bubbles produce the best sampling precision (Better BIC values). It is necessary to achieve a compromise in order to obtain the best results doing BIB and BIC testing as indicated in the following procedure:

### Precision check of the Analytical plate BIB

Prepare two series of 32 samples each by using **K<sub>2</sub>Cr<sub>2</sub>O<sub>7</sub> 0,04 gr/l** (starting from 4 gr/l concentrated solution and dilute 1:100)

Use the same solution for both samples and reagent

Program a new method with the following parameters:

Method name: **BIB**Type: **End Point**

Sample volume **K<sub>2</sub>Cr<sub>2</sub>O<sub>7</sub>** 0,04 gr/l 3 µL

Reagent volume **K2Cr2O7** 0,04 gr/l 300 µL

Filter: 340 nm

Incubation time: 116 sec

Calculation factor: 1000

## BIB Procedure

1. Select conical cup in **Options**
2. Program the BIB method using the above parameters and configure a reagent rack containing it
3. Enter in the calibration screen and select RBL for BIB method.
4. Put RINSE SOLUTION as reagent and sample (fourteenth position) and Start the Calibration
5. Program 64 BIB ; fill 32 conical samples cup with **K2Cr2O7 0,04 gr/l**
6. Load the racks on sample positions 1 and 2 of the Analyzer.
7. Substitute the RINSE SOLUTION with **K2Cr2O7 0,04 gr/l** on the reagent rack
8. Start the Analyzer. At the end of the last sampling, move the sample racks into positions 3 and 4.



The mean value for each of the 32 sample series must be within **360 ÷ 400**.

The coefficient of variation percentage (CV %) must be lower than **0.7 %** for each series of 32 samples. If there is only one anomalous result that compromises the CV, rerun it and if the result is acceptable, no further operations are required.

### Volume precision check Sampling line **BIC**

Prepare two series of 32 samples each by using **K<sub>2</sub>Cr<sub>2</sub>O<sub>7</sub> 4 gr/l**

Use the same solution for both samples and reagent

Program a new method with the following parameters:

Method name: **BIC**

Type: **End Point**

Sample volume **K<sub>2</sub>Cr<sub>2</sub>O<sub>7</sub> 4 gr/l** 3 µL

Reagent volume RINSE SOLUTION 300 µL

Filter: 340 nm

Incubation time: 116 sec

Calculation factor: 1000

### BIC Procedure

1. Program the BIC method using the above parameters and configure a reagent rack
2. Enter in the calibration screen and select RBL for BIC method.
3. Put RINSE SOLUTION as reagent and sample (fourteenth position) and Start the Calibration
4. Program 64 BIC ; fill 32 conical samples cup with **K<sub>2</sub>Cr<sub>2</sub>O<sub>7</sub> 4 gr/l**
5. Load the racks on sample positions 1 and 2 of the Analyzer.
6. Start the Analyzer. At the end of the last sampling, move the sample racks into positions 3 and 4.

The mean value for each of the 32 sample series must be between **350 ÷ 450**.

The coefficient of variation percentage (CV %) must be lower than **1.6 %** for each series of 32 samples. If there is only one anomalous result that compromises the CV, rerun it and if the result is acceptable, no further operations are required.

**Backlash Dil.:** Due to the mechanical inversion of the motor's rotation it is necessary to mechanically recuperate the Diluter in  $\mu\text{L}$ ; the value to be inserted can be calculated using the following procedure:

1. Enter in the Diagnostic program and make the Sampling Arm "Reset" into the Arm folder.
2. Enter in the **Diluter** folder, press the "HOME" button
3. Push "PROBE WASH TEST" button three times in order to fill the hydraulic sampling line completely
4. Remove the samples protective cover and the sample racks from their housing
5. Enter in the **ARM** folder and push the "GO" button in order to move the probe in the sampling position number 1
6. Put a sample cup containing distilled water under the probe (about 1ml)
7. Enter in the **DILUTER** folder and aspirate 300  $\mu\text{L}$  of distilled water
8. Remove the sample cup containing distilled water and put a plastic sheet under the probe
9. Select 30  $\mu\text{L}$  of distilled water and dispense it on the plastic sheet
10. Select 1  $\mu\text{L}$  of distilled water and aspirate a certain number of times until seeing an aspiration-like effect in the drop on the probe
11. The number of aspirations performed until one sees the aspiration effect represents the BACK LASH value to be inserted into the specific field in the **Configuration** folder
12. Push the **Save** button to memorize the inserted data
13. Push the **Reset** button into the ARM folder
14. Close the Diagnostic program by clicking on the Diagnostic button

**Heater:** to adjust the setting of the sample arm's preheater and reaction plate;

**Diluter Type:** to select the kind of diluter to be used. Each kind of diluter has an individual ratio step per  $\mu\text{L}$ .

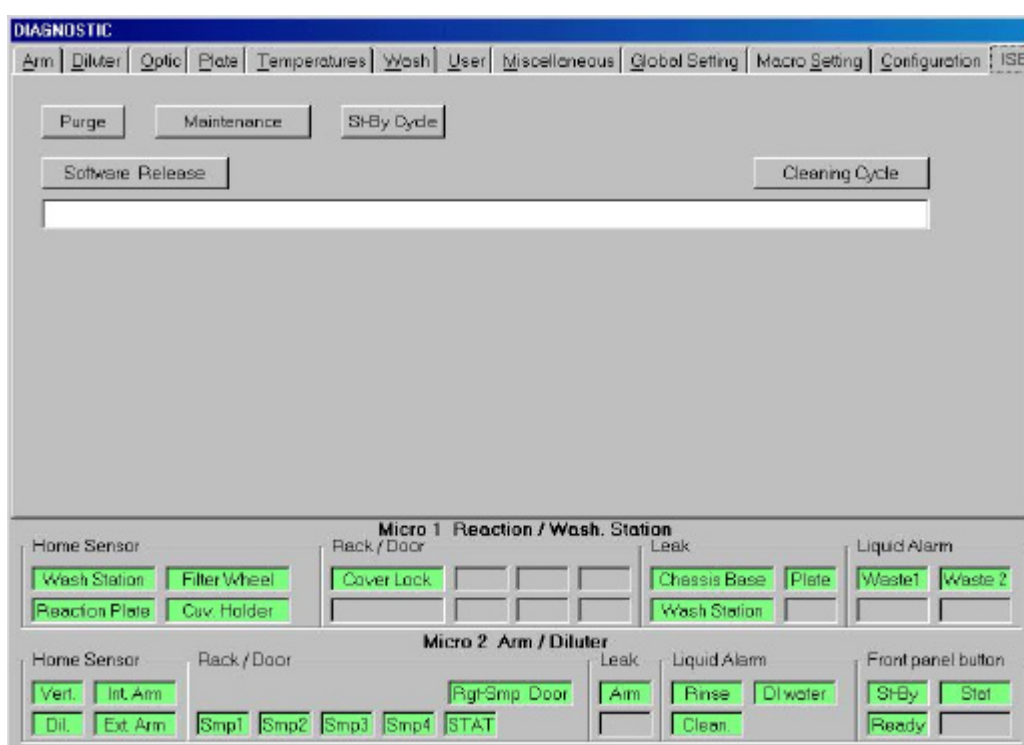
**Cover Lock SW.:** to disable the check for the closure of the ILab 300 Plus cover;

**Wash. Well level Test:** By clicking on the “Run Test” button the system performs an automatic verification and adjustment cycle for the liquid level in the washing well. A number appears in the box “Level” that corresponds to the rotation time of the rinse pump in milliseconds, that is necessary to maintain the correct liquid level in the washing well.

The functioning range is between 400 and 700 milliseconds. If the washing well is not sufficiently filled, a message appears to warn the user to check the hydraulic line.

If the value taken during the test is different from the previous one, it is necessary to click on “Save” in order to save the new data.

## 5.14 “ISE” Folder



This folder permits checks and maintenance activity on the ISE

**Purge:** fills the hydraulic ISE Module line

**Maintenance:** empties the electrodes line for their replacement

**St-By Cycle:** performs a prime cycle to renew the Calibrant A solution in the electrode line

**Software Release:** visualizes the ISE Module Firmware release

**Cleaning Cycle:** performs an electrode washing cycle

## CHAPTER 06

---

### - SETTINGS AND ADJUSTMENTS -

### INDEX

6	SETTINGS AND ADJUSTMENTS .....	4
6.1	SAMPLING ARM.....	4
6.1.1	ALIGNMENT AND ADJUSTMENTS.....	4
6.1.2	SUBSTITUTION OF THE SAMPLING ARM .....	6
6.1.3	SUBSTITUTION OF THE OPTIC SENSORS .....	7
6.1.3.1	SUBSTITUTION OF THE OS7 OPTIC SENSOR.....	7
6.1.3.2	SUBSTITUTION OF THE OS6 AND OS5 OPTO SENSORS.....	9
6.1.3.2a	SUBSTITUTION OF THE OS6 OPTO SENSOR.....	10
6.1.3.2b	SUBSTITUTION OF THE OS5 OPT-COUPLER .....	11
6.1.4	SAMPLING ARM MOTORS .....	12
6.1.4.1	SUBSTITUTION OF THE M7 MOTOR .....	13
6.1.4.2	SUBSTITUTION OF THE M6 MOTOR .....	17
6.1.4.3	SUBSTITUTION OF THE M5 MOTOR .....	19
6.1.5	PRE-HEATER AND LEVEL SENSOR.....	21
6.1.5.1	SUBSTITUTION OF THE PRE-HEATER AND LEVEL SENSOR ASSEMBLY .....	21
6.1.5.2	LEVEL SENSOR FUNCTIONING CHECK .....	25
6.1.6	DILUTER MODULE .....	27
6.1.6.1	SUBSTITUTION OF THE DILUTER MODULE.....	27
6.1.6.2	SUBSTITUTION OF THE DILUTER MICRO-PUMP .....	27

6.1.6.3	SUBSTITUTION OF THE DILUTER ELECTRO-VALVE .....	28
6.1.6.4	DILUTER MODULE - CHECK PROCEDURE FOR MECHANICAL FUNCTIONING .....	28
6.16.5	BACKLASH ADJUSTMENT.....	29
6.2	OPTIC ASSEMBLY.....	31
6.2.1	SUBSTITUTION OPTIC ASSEMBLY .....	31
6.2.2	SUBSTITUTION OF THE PHOTOMETER LAMP.....	33
6.2.3	ELECTRONIC ADJUSTMENT OF THE PHOTOMETER.....	34
6.2.4	SUBSTITUTION OF THE OPTO SENSOR.....	36
6.2.5	SUBSTITUTION OF THE MOTOR.....	37
6.2.6	SUBSTITUTION OF THE OPTIC FILTERS .....	38
6.3	PUMP ASSEMBLY .....	39
6.3.1	SUBSTITUTION OF THE PUMP ASSEMBLY .....	39
6.3.2	SUBSTITUTION OF THE PERISTALTIC PUMP MOTOR.....	40
6.3.2.1	FUNCTION CHECK FOR THE PERISTALTIC PUMP .....	41
6.3.3	SUBSTITUTION OF PART(S) .....	41
6.3.4	MECHANICAL ADJUSTMENT OF THE PERISTALTIC PUMP.....	42
6.3.4.1	HOW TO SUBSTITUTE PERISTALTIC PUMP TUBES .....	43
6.4	REACTION PLATE.....	44
6.4.1	REACTION PLATE - MECHANICAL ADJUSTMENT .....	44
6.4.2	SUBSTITUTION OF THE REACTION PLATE.....	45
6.4.3	SUBSTITUTION OF THE OS3 OPTIC SENSOR.....	46
6.4.4	SUBSTITUTION OF THE M3 MOTOR .....	48
6.5	WASHING STATION ASSEMBLY .....	50
6.5.1	SUBSTITUTION OF THE WASHING STATION ASSEMBLY .....	50

6.5.2	SUBSTITUTION OF THE WASHING STATION OPTIC SENSOR.....	52
6.5.3	SUBSTITUTION OF THE WASHING STATION MOTOR.....	52
6.5.4	WASHING STATION - MECHANICAL ADJUSTMENT AND FUNCTIONAL CHECK.....	53
6.5.4.1	WASHING STATION - HORIZONTAL ALIGNMENT.....	53
6.5.4.2	WASHING STATION - VERTICAL ALIGNMENT .....	54
6.5.4.3	WASHING STATION - FUNCTIONAL CHECK .....	54
6.5.5	WASHING STATION - CUVETTE HOLDER ASSEMBLY .....	55
6.5.5.1	SUBSTITUTION OF THE CUVETTE HOLDER ASSEMBLY MAGNET .....	55
6.5.5.2	SUBSTITUTION OF THE CUVETTE HOLDER ASSEMBLY OPTO SENSOR.....	55
6.5.5.3	ADJUSTMENT AND FUNCTIONAL CHECK OF THE CUVETTE HOLDER MAGNETIC TRIGGER.....	56
6.6	TEMPERATURE CHECK AND ADJUSTMENT .....	58
6.6.1	REACTION CUVETTE TEMPERATURE ADJUSTMENT AND CHECK .....	58
6.6.2	TEMPERATURE CHECK AND ADJUSTMENT USING <u>KIT P/N 55-01204-00</u> .....	59
6.6.3	REFRIGERATION ASSEMBLY TEMPERATURE CHECK.....	68
6.7	RACKS HOLDER ALIGNMENT.....	69
6.7.1	MECHANICAL ALIGNMENT OF THE SAMPLES RACKS HOLDER .....	69
6.7.2	MECHANICAL ALIGNMENT OF THE REAGENTS RACKS HOLDER.....	70


## 6 SETTINGS AND ADJUSTMENTS

### 6.1 SAMPLING ARM

#### 6.1.1 ALIGNMENT AND ADJUSTMENTS

1. Turn on the ILab 300 Plus system (instrument and computer).
2. Launch the Diagnostic program, select "Global Setting". Wait until the instrument has completed the reset procedure.
3. Make sure that all the "home sensors" of the Sampling Arm (Vert., Int. Arm and Ext. Arm) light up green.
4. Make sure that the probe is centered with respect to the wash well. If not, loosen the fastening screws on the external arm and mechanically align the arm on the wash well, paying attention to center the probe with respect to the wash well. If absolutely necessary, also the internal arm can be adjusted - loosen only one of its fastening screws.

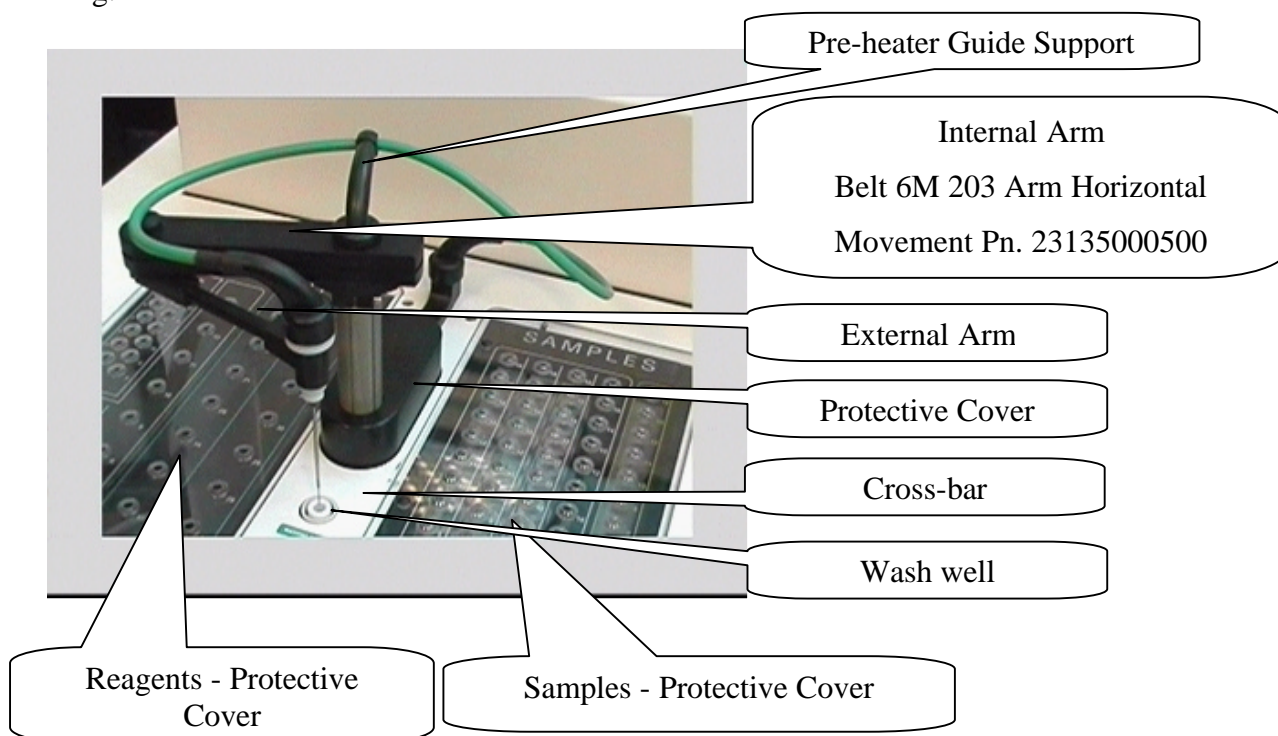Fig.1




5. Tighten the External Arm fastening screws and, if necessary, those of the Internal Arm, centering the probe with respect to the wash well.
6. Perform a complete Arm reset in order to ascertain correct centering of the probe on the wash well. If the centering is not exact, repeat steps 4 and 5.
7. After having mechanically centered the sampling probe on the wash well, adjust all the other positions by using the computer's numeric keyboard as follows:

**N.B.: make sure that the distance between the top edge of the wash well and the tip of the probe is "2.0 +/- 0.5 mm"**

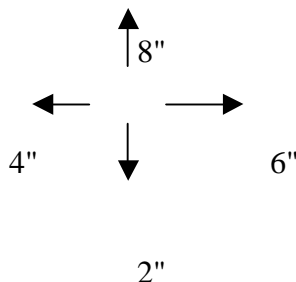**Use the following keys to make the necessary adjustments.**

- To move the arm horizontally



- <Next / + > to move forward one position
- <Prev / - > to move back one position
- <9> to move upward
- <3> to move downward

8. Select "Sample" from "Setting of" and then press the "Go" key. The arm will automatically move itself to position #1 of the Samples rack. Adjust the positioning by using the keys indicated for centering the probe with respect to the hole in the Samples Cover.
9. Move on to the next position by pressing the "Next/+" key and center the probe, as above.
10. Make sure that when adjusting position #16, mechanical end-of-run of the arm does not occur. If this should happen, repeat the mechanical alignment of the sampling probe with respect to the wash well and repeat the procedure from step 4.
11. Press "Save" to memorize the new positions' settings.



12. Select "Reagent" from "Setting of" and then press the "Go" key. The arm will automatically move itself to position #1 of the Reagents rack.  
Adjust the positioning by using the keys indicated for centering the probe with respect to the hole in the Reagents Cover.
13. Move on to the next position by pressing the "Next/+" key and center the probe, as above.
14. Make sure that when adjusting position #33, mechanical end-of-run of the arm does not occur. If this should happen, repeat the mechanical alignment of the sampling probe with respect to the wash well and repeat the procedure from step 4.
15. Press "Save" to memorize the new position settings.
16. Set the positioning of Std/Ctrl and Dil, Ise, Disp. Press "Save" to memorize the new position settings.
17. To exit the diagnostic program, press the "Diagnostic" key
18. To exit the analyzer program press "shutdown" key..

### 6.1.2 SUBSTITUTION OF THE SAMPLING ARM

**N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**

1. Remove the sampling probe from its housing.
2. Remove the samples and reagents racks' protective covers.
3. Remove the samples rack and reagents rack.
4. Unscrew the four cross bar fastening screws. Remove the cross bar by slipping out the Pre-heater guide support (Fig.1).
5. Remove the waste tube from the wash well.
6. Remove the teflon tube from the diluter head by unscrewing the fastening washer, paying special attention to leakage of liquid.
7. Unscrew the four anchored screws indicated in Fig. 5, unplug the J1, J2, J6, J8, J13, J15 and J17 connectors from the Sampling Interface assy. Disconnect the ground wire and following the procedure provided in section 6.1.5.1 remove the Level Sensor and Pre-Heater assembly.

8. Take out the sampling arm and substitute it with the new one.
9. To remount, repeat the above steps in inverse order: from 7 to 1.

**Warning: After having substituted the sampling arm, it is necessary to mechanically align the probe with respect to the wash well and center all the other positions by following the procedures provided in Section 6.1.1 "Alignment and Adjustments".**

### 6.1.3 SUBSTITUTION OF THE OPTIC SENSORS

#### 6.1.3.1 Substitution of the OS7 Optic sensor

1. Remove the sampling probe from its housing.
2. Remove the samples and reagents racks' protective covers.
3. Remove the samples rack and reagents rack.
4. Unscrew the four cross bar fastening screws. Remove the cross bar by slipping out the Pre-heater guide support (Fig.1).
5. Loosen the protective cover's set screw (Fig.2) and lift it.
6. Remove the fastening screws from the optic sensor support, unplug the J11 connector from the Sampling Interface assy and take out the entire assembly (Fig.3).
7. Remove the fastening screws from the optic sensor itself and substitute it with the new one.
8. To remount, repeat above steps in inverse order: from 7 to 1.

**N.B.: when remounting the optic sensor, make sure to center it with respect to the slots, as shown in Fig. 4.**

**N.B.: due to the fact that it may be necessary to further adjust the positioning of the optic sensor, screw down the cross bar and the protective cover only after having made sure that the sampling arm functions correctly by carrying out the procedure described in Section 6.1.1**

**Warning:** Do not tighten the protective cover's set screw as this could result in mechanical malfunction of the sampling arm's vertical movement (Fig. 2).

**Warning:** After having substituted the optic sensor, it is necessary to mechanically align the probe with respect to the wash well and center all the other positions by following the procedures provided in Section 6.1.1 "Alignment and Adjustments".

**N.B.:** the following check is valid only for the OS7 optic sensor

#### **Verification of the correct positioning of the OS7 optic sensor**

Set the Samples position at #16. If it is not possible to obtain an exact centering (due to mechanical end-of-run), move the OS7 optic sensor to the left. After having mechanically aligned the optic sensor, center all the other positions by following the procedures provided in Section 6.1.1

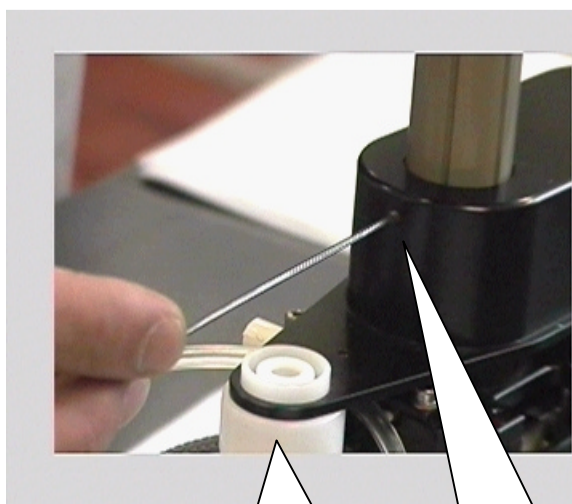


Fig. 2

Washing Well assembly  
PN. 23905006300

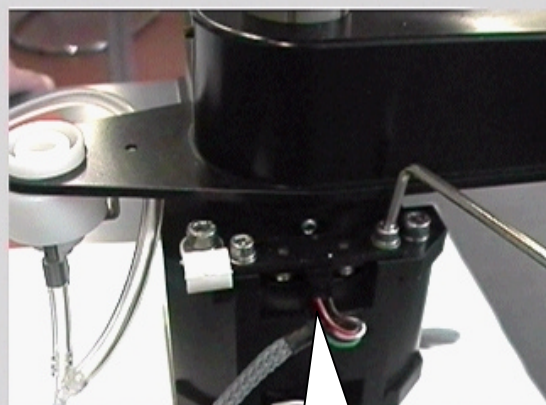


Fig. 3

OS7 Opto coupler assy  
PN. 23910002300

Protective cover set screw

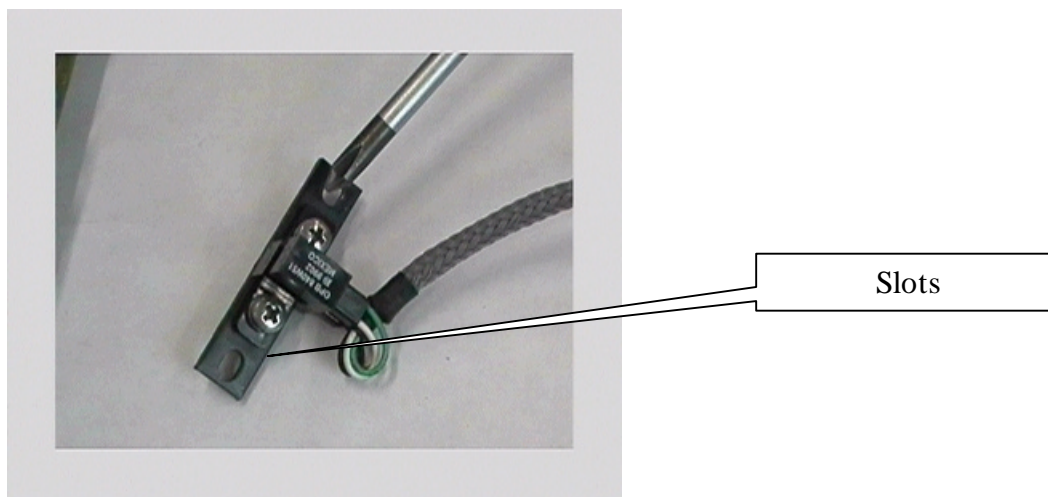


Fig. 4

### 6.1.3.2 Substitution of the OS6 and OS5 Opto sensors

**N.B.:** Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.

In order to arrive at the OS6 and OS5 Opto sensors, it is necessary to remove the arm assembly from its housing. Unplug connectors J9 and J10 from the Sampling Interface assy; disconnect the wash well hydraulic tube and the teflon tube on the diluter. Loosen the four anchored screws. When removing the sampling arm, be extremely careful to not damage the pre-heater.

The following steps, from 1 to 6, are to be carried out regardless of which opto sensor is being substituted (OS6 or OS5).

1. Remove the sampling probe from its housing.
2. Remove the samples and reagents racks' protective covers (Fig.1).
3. Remove the samples rack and reagents rack.
4. Unscrew the four cross bar fastening screws. Remove the cross bar by slipping out the Pre-heater guide support (Fig.1).
5. Remove the waste tube from the wash well and the teflon tube on the diluter.
6. Unscrew the four anchored screws indicated in (Fig.5), unplug the J1, J2, J15 and J17 connectors from the Sampling Interface assy (Fig.5). Disconnect the ground wire and remove the sampling arm.

### 6.1.3.2 a Substitution of the OS6 Opto sensor

- a) Remove the screws which fasten the opto sensor support and take out the entire assembly (Fig.6).
- b) Unplug the J10 connector from the Sampling Interface assy, remove the opto sensor fastening screws and substitute it with the new one.
- c) To remount, repeat the above steps in inverse order: from b) to a) and then from 6 to 1 in Section 6.1.3.2

**Warning:** Do not tighten the protective cover's set screw as this could result in mechanical malfunction of the sampling arm's vertical movement (Fig.2).

**N.B.:** when remounting the opto sensor, make sure to center it with respect to the slots, as shown in Fig. 4.

**Warning:** After having substituted the opto sensor, it is necessary to mechanically align the probe with respect to the wash well and center all the other positions by following the procedures provided in Section 6.1.1 "Alignment and Adjustments".

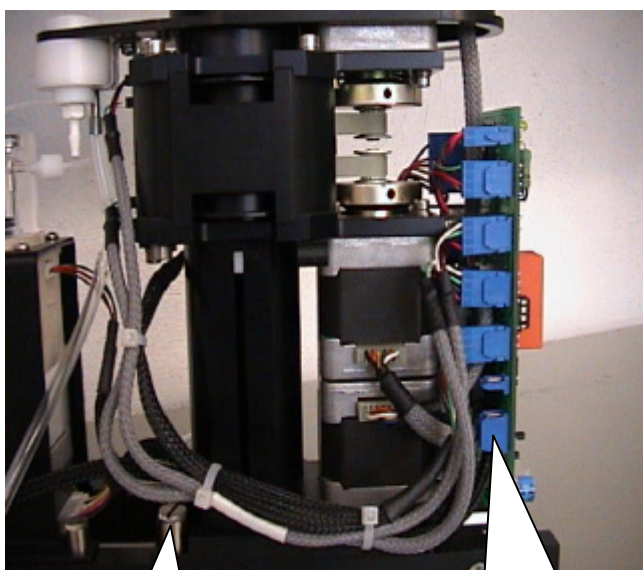


Fig. 5

Anchored screws

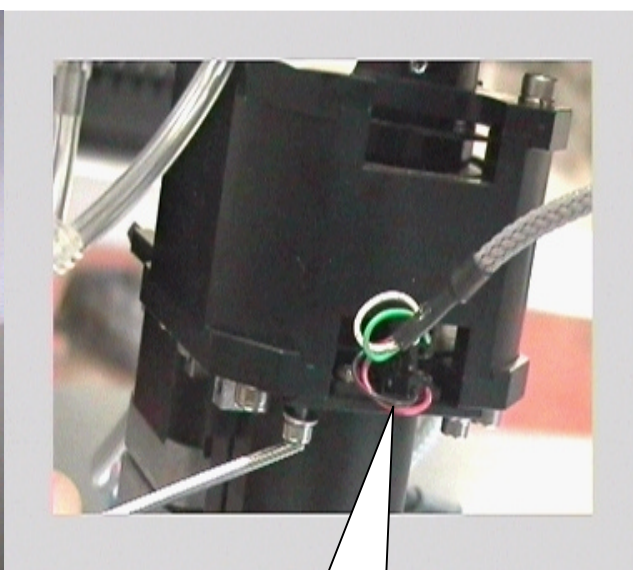


Fig. 6

Sampling Interface Assy  
PN.23300019201OS6 Opto coupler assy  
PN. 23910002300

### 6.1.3.2 b Substitution of the OS5 Opt-coupler

- a) Carry out Steps 1 through 6 as indicated in Section 6.1.3 "Substitution of the OS6 and OS5 Opto sensors.
- b) Remove the screws which fasten the opto sensor support and take out the entire assembly (Fig.7).
- c) Unplug the J9 connector from the Sampling Interface assy, remove the optic sensor fastening screws and substitute it with the new one.
- d) Pay particular attention to the spacer washers placed between the opto sensor and the support (Fig.8)
- e) To remount, repeat the above steps in reverse order: from c) to a) and then from 6 to 1 in Section 6.1.3.2

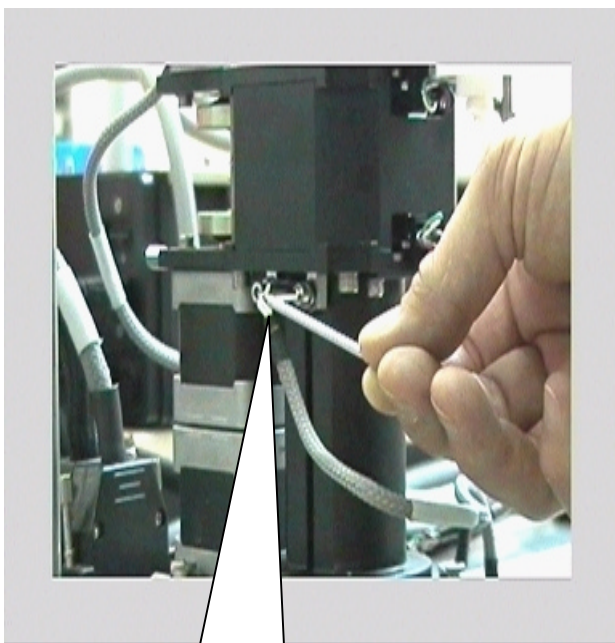


Fig. 7

OS5 Opto coupler assy  
PN. 23910002400

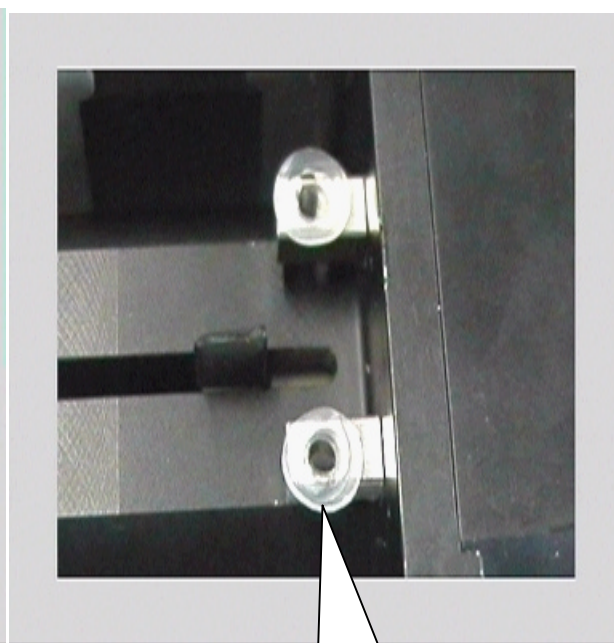


Fig. 8

Spacer washers

**Warning:** Do not tighten the protective cover's set screw as this could result in mechanical malfunction of the sampling arm's vertical movement (Fig. 2).

**N.B.:** make sure that the distance between the top edge of the wash well and the tip of the probe is "2.0 +/- 0.5 mm"

**Warning:** After having substituted the opto sensor, it is necessary to mechanically align the probe with respect to the wash well and center all the other positions by following the procedures provided in Section 6.1.1 "Alignment and Adjustments".

## 6.1.4 SAMPLING ARM MOTORS

**N.B.:** Make sure the instrument (ILab 300 Plus) is turned off before performing this procedure.

**Important:** pay special attention to the Pre-heater when removing or servicing the sampling arm.

In order to arrive at the M7, M6 and M5 motors, it is necessary to remove the sampling arm assembly from its housing; unplug the J1, J2, J3, J4, J5, J15 and J17 connectors from the Sampling Interface assy; disconnect the hydraulic tube from the wash well and the teflon tube on the diluter; remove the four anchored fastening screws from the assembly.



### 6.1.4.1 Substitution of the M7 motor

1. Remove the sampling probe
2. Remove the samples and reagents racks' protective covers.
3. Remove the samples rack and reagents rack.
4. Unscrew the four cross bar fastening screws. Remove the cross bar by slipping out the Pre-heater guide support (Fig.1).
5. Remove the waste tube from the wash well.
6. Remove the teflon tube on the diluter head by unscrewing the fastening washer - pay particular attention to leakage of liquid.
7. Unscrew the four anchored screws indicated in Fig. 3; unplug the J1, J2, J4, J15 and J17 connectors from the Sampling Interface assy; unplug the ground wire; remove the sampling arm.
8. Loosen the fastening screw and lift the sampling arm's protective cover (Fig.4), in order to arrive at the M7 motor.
9. Remove the M7 motor fastening screws (Fig 5).
10. Loosen the motor pulley screws (Fig.6).
11. Take out the M7 motor and substitute it with the new one.(Fig.7).
12. Lift the pulley to approximately 0.5mm from its upper limit and tighten its screws(Fig.8).
13. Push the M7 motor outward, as shown in (Fig.9), and holding it in this position, so as to maintain the belt in tension, tighten the motor fastening screws.
14. To remount, repeat the above steps in inverse order: from 8 to 1.

<p><b>Warning:</b> Do not tighten the protective cover's set screw as this could result in mechanical malfunction of the sampling arm's vertical movement (Fig. 2).</p>
---



**Warning: After having substituted the motor, it is necessary to mechanically align the probe with respect to the wash well and center all the other positions by following the procedures provided in Section 6.1.1 "Alignment and Adjustments".**

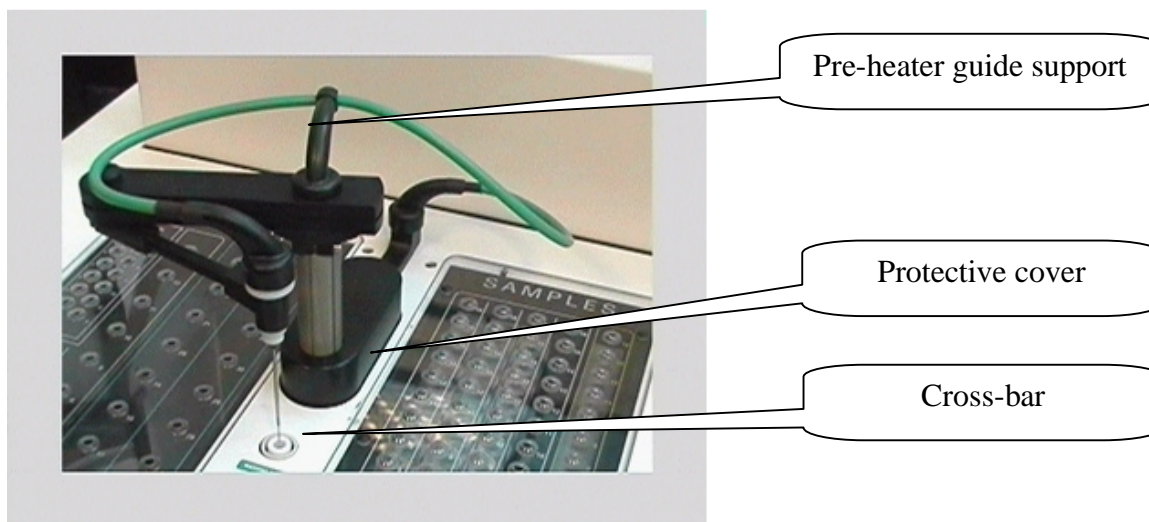


Fig.1

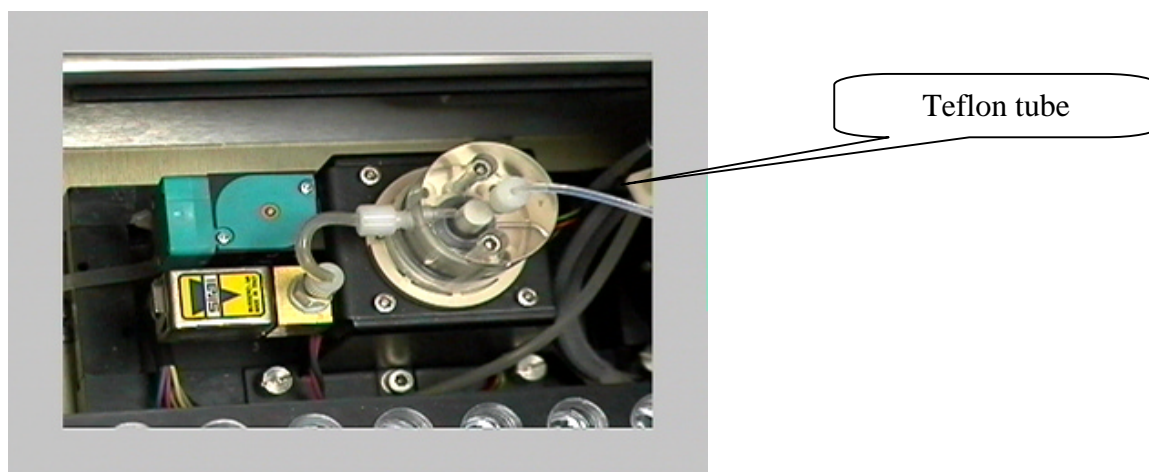


Fig. 2 Dilutor Module



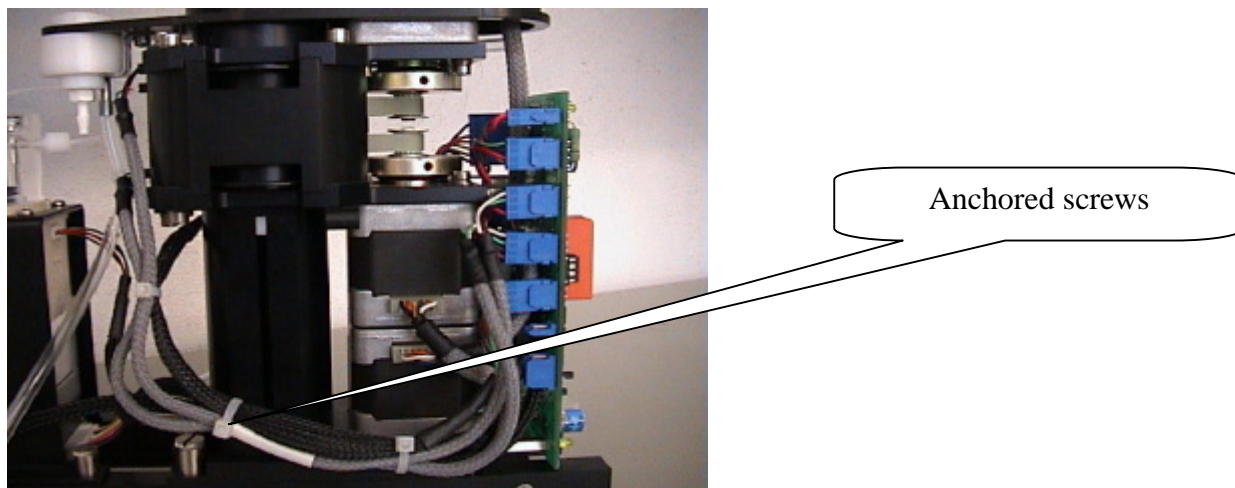


Fig. 3 Arm assembly without base and Sampling Interface Board PN. 23915002000

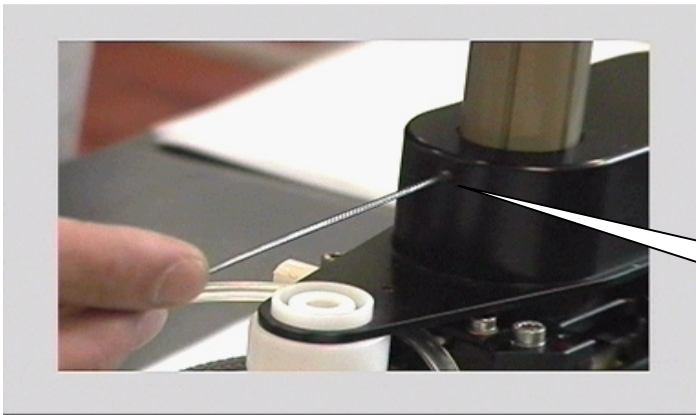


Protective cover set screw

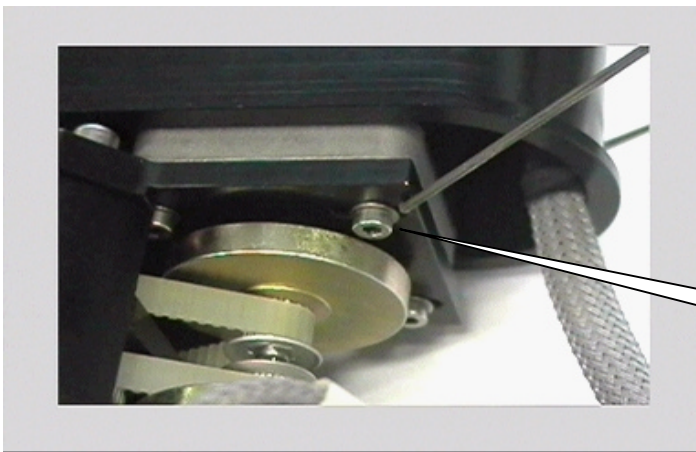Fig. 4




M7 motor fastening screws

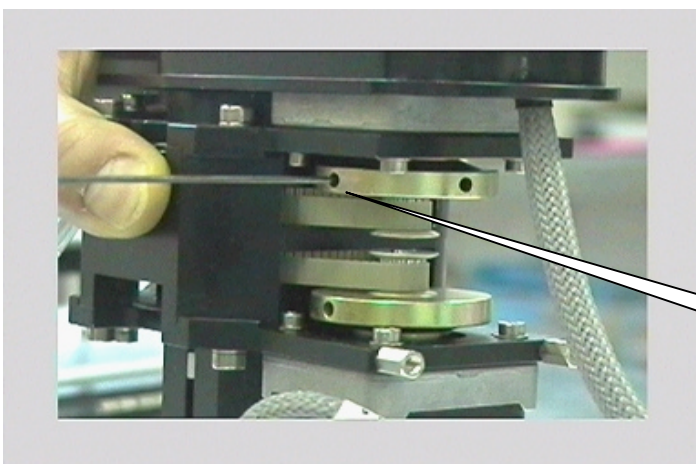Fig. 5




Pulley screws

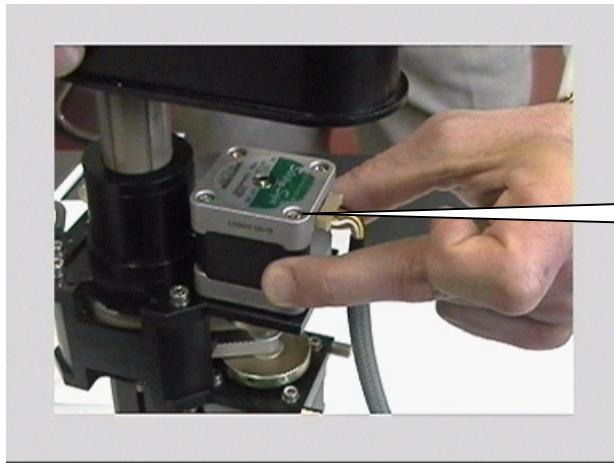Fig. 6



M7 Motor assy PN. 23100039900

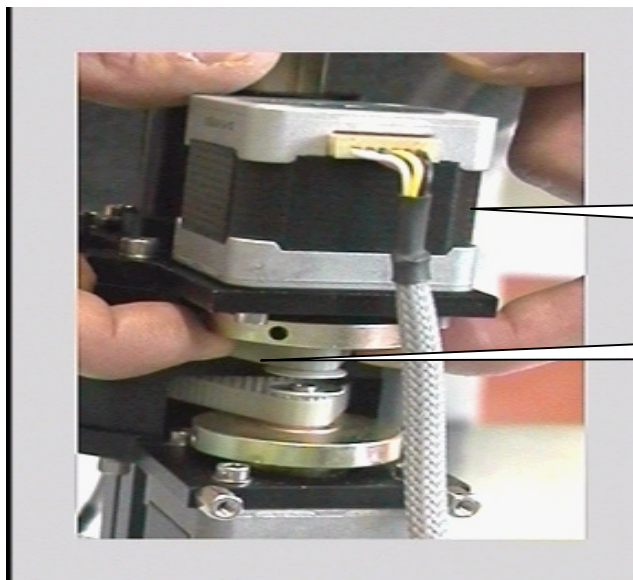
Fig. 7



M7 Motor assy PN. 23100039900

M7 Motor's belt PN. 23935003100

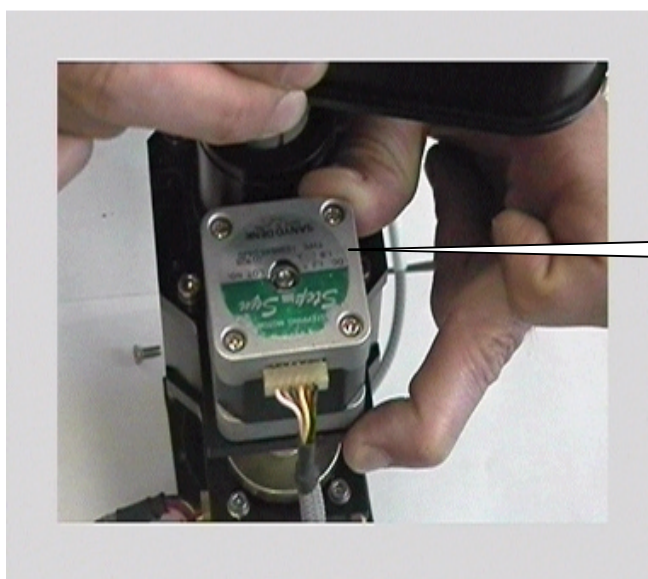Fig. 8



M7 Motor assy PN. 23100039900

Fig. 9




### 6.1.4.2 Substitution of the M6 Motor

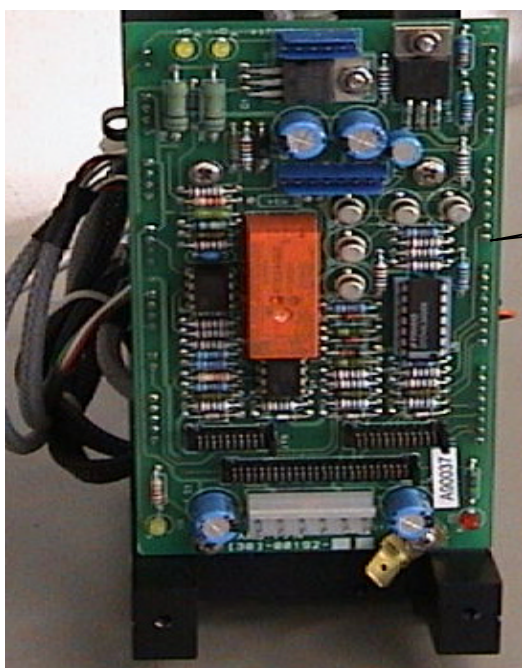**N.B.:** Make sure the instrument (ILab 300 Plus) is turned off before performing this procedure.

**Important:** pay special attention to the Pre-heater when removing or servicing the sampling arm.

1. Remove the sampling probe.
2. Remove the samples and reagents racks' protective covers.
3. Remove the samples rack and reagents rack.
4. Unscrew the four cross bar fastening screws. Remove the cross bar by slipping out the Pre-heater guide support (Fig.1).
5. Remove the waste tube from the wash well.
6. Remove the teflon tube on the diluter head by unscrewing the fastening washer (Fig.2), paying special attention to leakage of liquid.
7. Unscrew the four anchored screws indicated in (Fig.3); unplug the J1, J2, J5, J15 and J17 connectors from the Sampling Interface assy; unplug the ground wire; remove the sampling arm.
8. Take the Sampling Interface assy out by unplugging the connectors and unscrewing the fastening screws (Fig.10).



Sampling Interface assy PN.23300019201

Fig.10

**N.B.: In order to be able to remove the M6 motor from its support, it is necessary to first partially disassemble the module as shown in Fig.12.**

9. Unscrew the three fastening screws (Fig.11) and remove the base as shown in (Fig.12).
10. Remove the motor fastening screws and loosen its pulley screws (Fig.13).
11. Take out the M6 motor and substitute it with the new one
12. Lift the pulley to approximately 0.5mm from its upper limit and tighten its screws(Fig.13)
13. Push the motor outward, as shown in (Fig.9), and holding it in this position so as to maintain the belt in tension, tighten the motor fastening screws.
14. Tighten the base support screws (Fig.11).
15. To remount, repeat the above steps in inverse order: from 8 to 1.

**Warning:** Do not tighten the protective cover's set screw as this could result in mechanical malfunction of the sampling arm's vertical movement (Fig.4).

**Warning:** After having substituted the motor, it is necessary to mechanically align the probe with respect to the wash well and center all the other positions by following the procedures provided in Section 6.1.1 "Alignment and Adjustments".

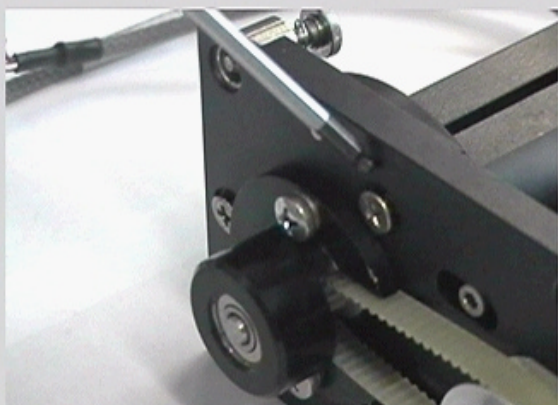


Fig. 11


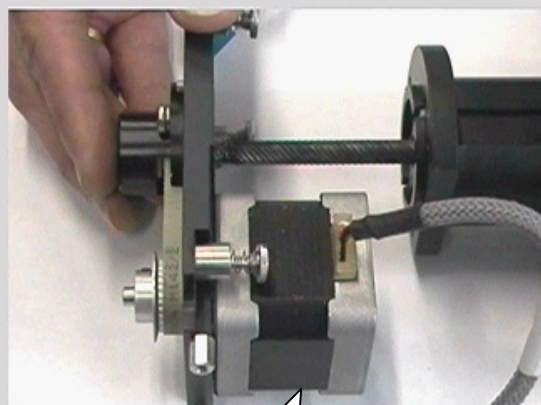


Fig. 12

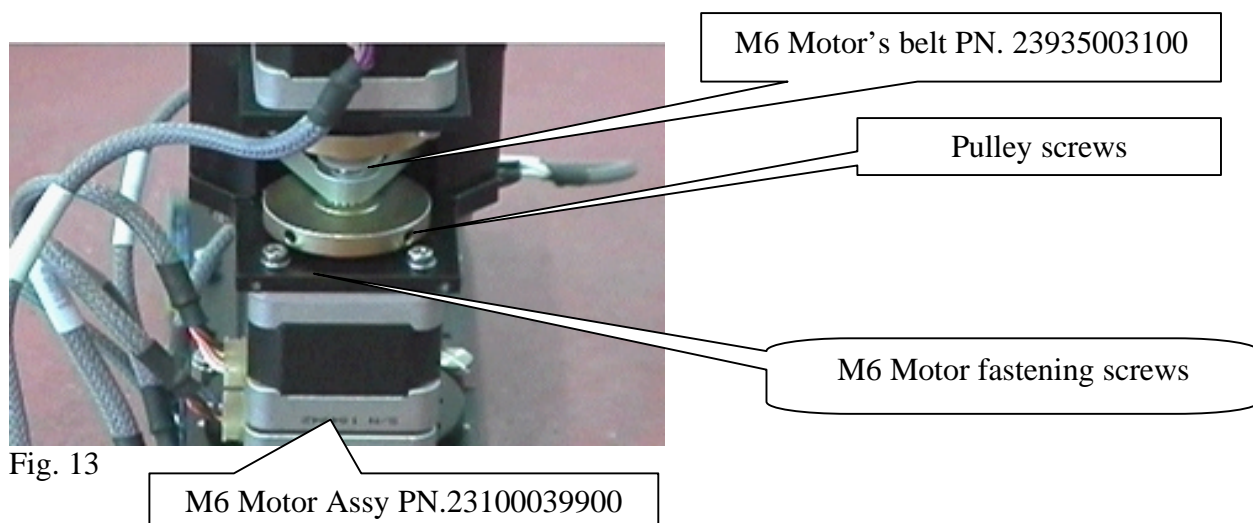M5 Motor Assy (Vertical)  
PN. 23100039900



#### 6.1.4.3 Substitution of the M5 motor

**N.B.:** Make sure the instrument (ILab 300 Plus) is turned off before performing this procedure.

**Important:** pay special attention to the Pre-heater when removing or servicing the sampling arm.

1. Remove the sampling probe.
2. Remove the samples and reagents racks' protective covers.
3. Remove the samples rack and reagents rack.
4. Unscrew the four cross bar fastening screws. Remove the cross bar by slipping out the Pre-heater guide support (Fig.1).
5. Remove the waste tube from the wash well.
6. Remove the teflon tube on the diluter head by unscrewing the fastening washer (Fig.2), paying special attention to leakage of liquid.
7. Unscrew the four anchored screws indicated in (Fig.3); unplug the J1, J2, J3, J15 and J17 connectors from the Sampling Interface assy; unplug the ground wire; remove the sampling arm.
8. Take the Sampling Interface assy out by unplugging the connectors and unscrewing the fastening screws (Fig.10).

**N.B.: In order to be able to remove the M5 motor from its support, it is necessary to first partially disassemble the module as shown in Fig.12.**

9. Unscrew the three fastening screws (Fig.11) and remove the base as shown in (Fig.12).
10. Remove the motor fastening screws and loosen its pulley screws (Fig.14).
11. Take out the M5 motor and substitute it with the new one
12. Lift the pulley to approximately 0.5mm from its upper limit and tighten its screws (Fig.14).
13. Push the motor outward, as shown in (Fig.9), and holding it in this position so as to maintain the belt in tension, tighten the motor fastening screws.
14. Tighten the base support screws (Fig.11).
15. To remount, repeat the above steps in inverse order: from 8 to 1.

**Warning:** Do not tighten the protective cover's set screw as this could result in mechanical malfunction of the sampling arm's vertical movement (Fig. 4).

**Warning:** After having substituted the motor, it is necessary to mechanically align the probe with respect to the wash well and center all the other positions by following the procedures provided in Section 6.1.1 "Alignment and Adjustments".

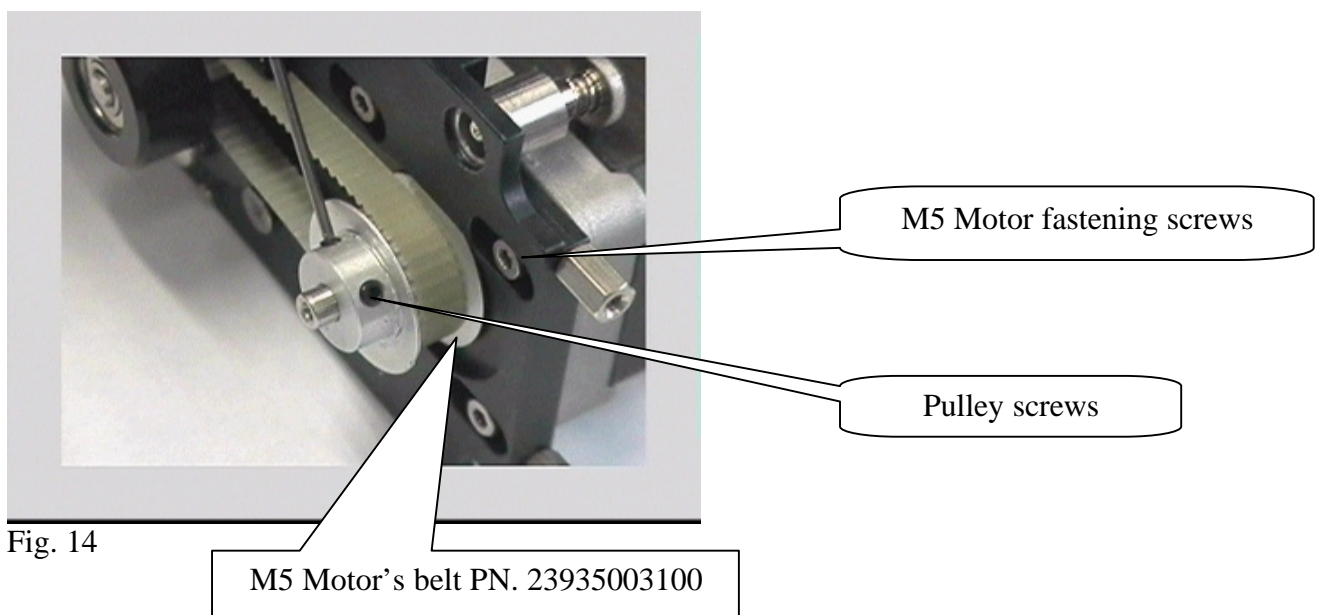


Fig. 14

M5 Motor's belt PN. 23935003100

## 6.1.5 PRE-HEATER AND LEVEL SENSOR

### 6.1.5.1 Substitution of the pre-heater and level sensor assembly

**N.B.:** Make sure the instrument (ILab 300 Plus) is turned off before performing this procedure.

1. Remove the sampling probe.
2. Remove the samples and reagents racks' protective covers.
3. Remove the samples rack and reagents rack.
4. Unscrew the four cross bar fastening screws. Remove the cross bar by slipping out the rear Pre-heater guide support (Fig.2).
5. Remove the waste tube from the wash well.
6. Remove the teflon tube on the diluter head by unscrewing the fastening washer (Fig.3) pay particular attention to leakage of liquid.
7. Unplug connector J7 from the Sampling Interface assy.
8. Lift the teflon washer and expose the half-ring; remove the half-ring (Fig.4) and (Fig.5).
9. Remove the teflon washer and the spring (Fig.6).
10. Remove the Pre-heater assembly (Fig.7) and (Fig.8), making sure that the blocking pin remains in place.
11. Remount the new assembly, repeating the above steps in inverse order: from step 10 to step 1.

**Note 1:** the length of the Pre-heater supply cable from the top of the external arm to the intermediate guide must be approximately 23 cm (Fig.2).

**Note 2:** when remounting the teflon tube in the diluter, make sure that it is inserted correctly.

12. Turn on the ILab 300 Plus system (instrument and computer).
13. Launch the diagnostic program. Select the "Diluter" function and press the "Probe Rinse Cycle" several times in order to fill the hydraulic circuit with liquid.



14. Make sure that no liquid leaks from the probe assembly and/or the dilutor. Make sure that the liquid flows from the probe freely.
15. Select the "Arm" key and then "Go Sample". Set Samples position #16. Make sure that the Pre-heater cable does not hinder the arm while it is in movement. If this situation should occur, turn the Pre-heater cable counter-clockwise in order to optimize the assembly.
16. Press the "Reset" key
17. To exit the diagnostic program, press the "Diagnostic" key
18. To exit the analyzer program, press the "Shutdown" key.

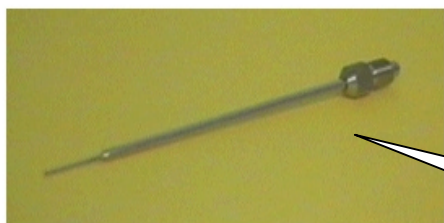


Fig. 1 Probe Assembly PN. 23910006201

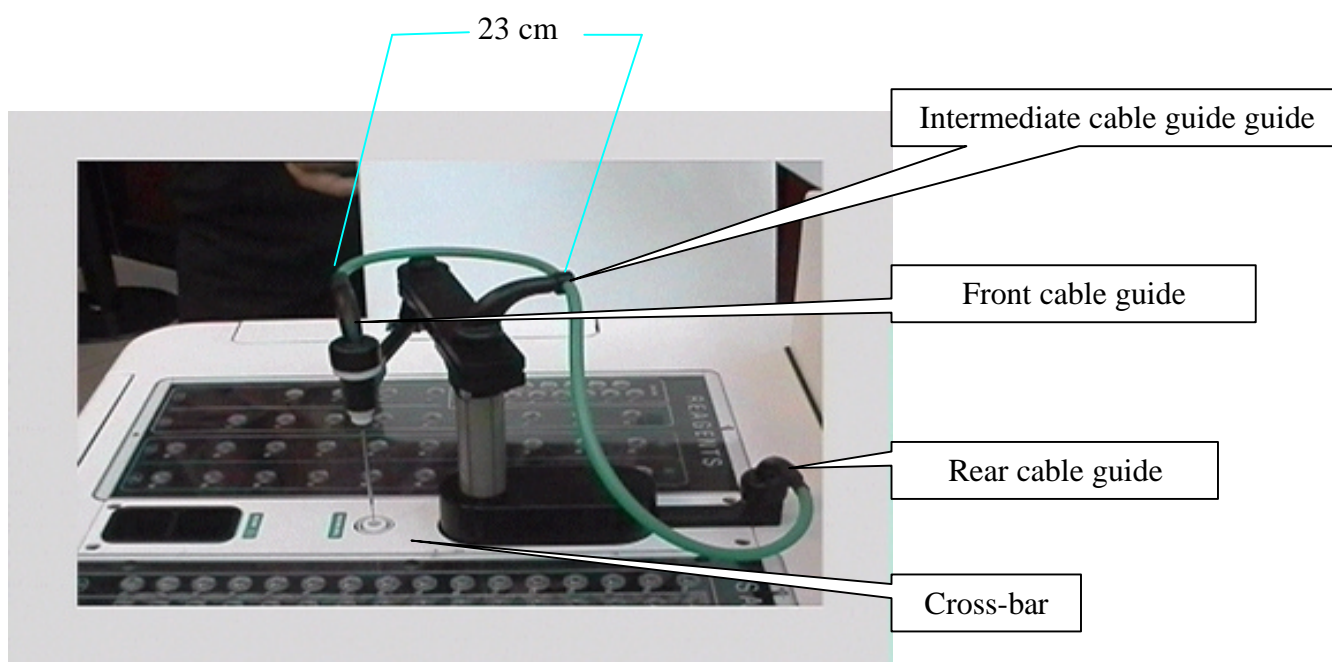


Fig. 2

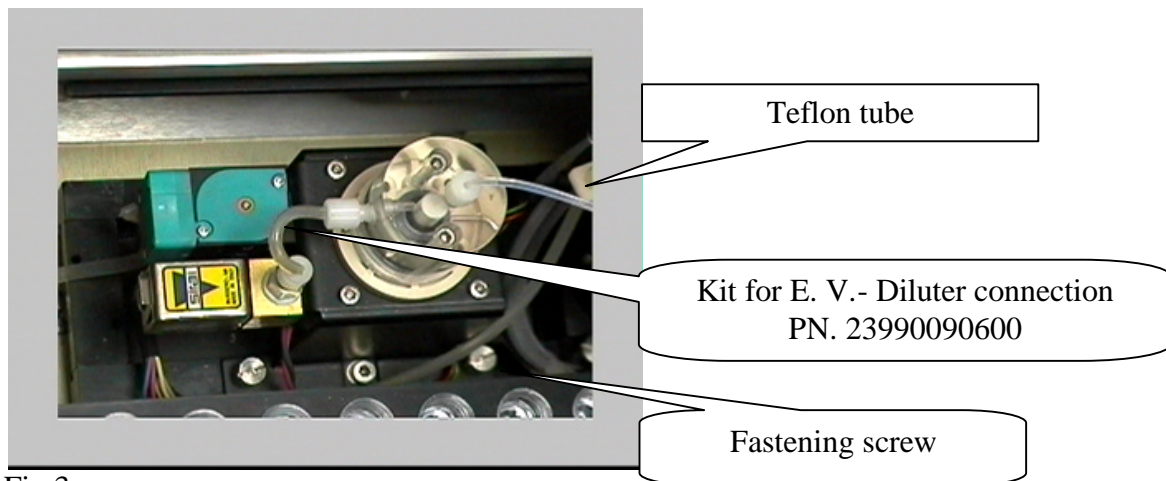


Fig.3

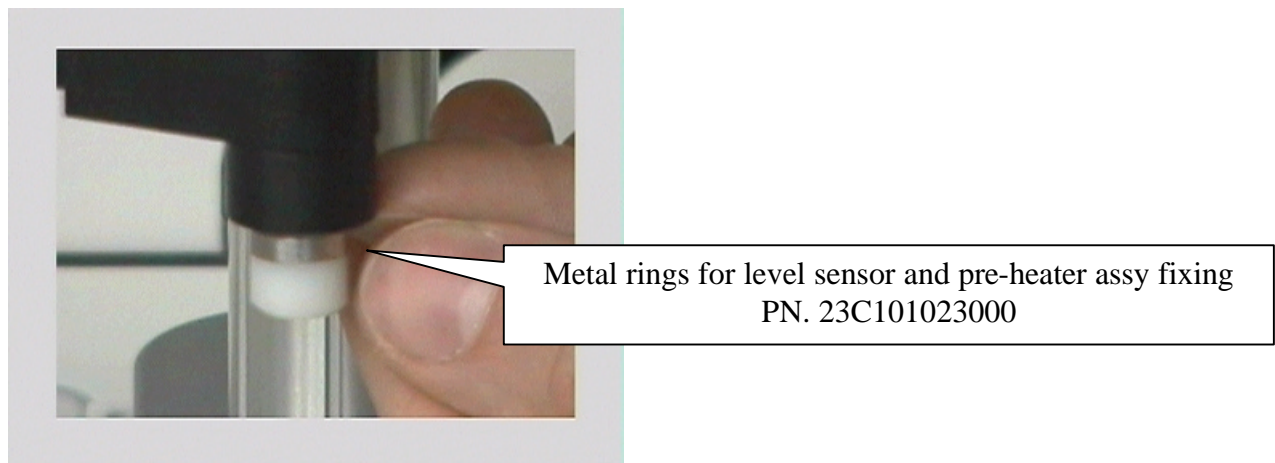


Fig. 4

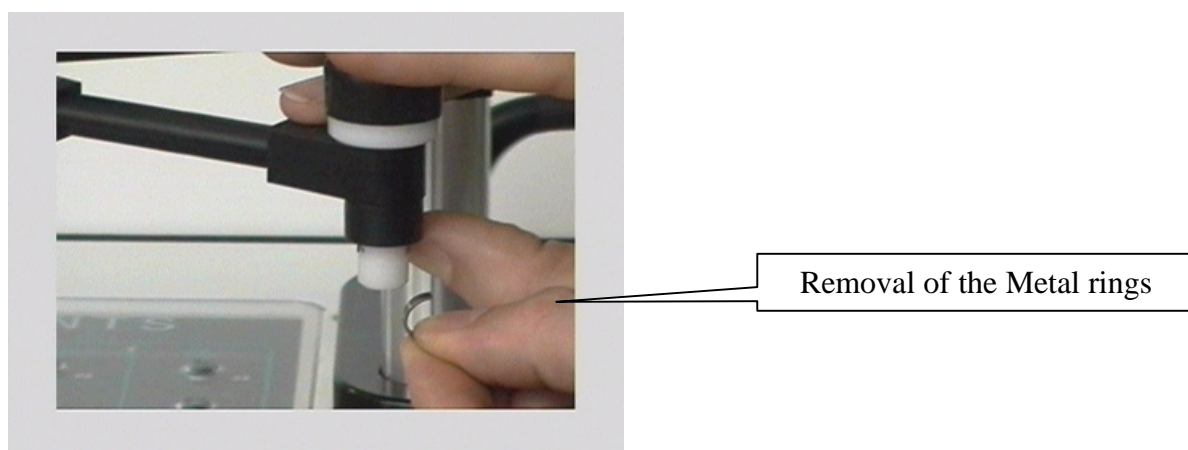


Fig. 5

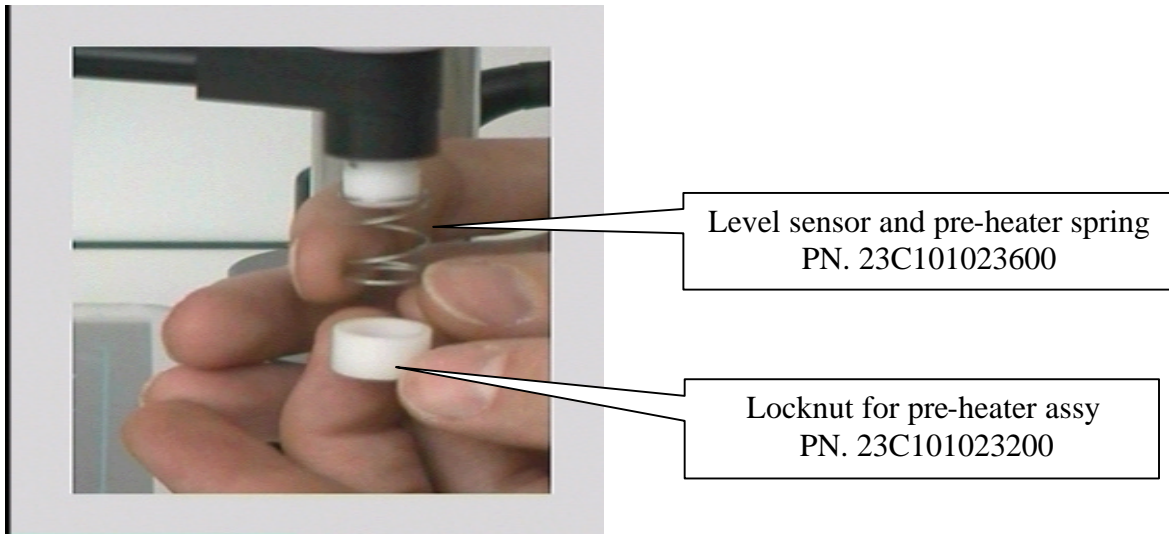


Fig. 6

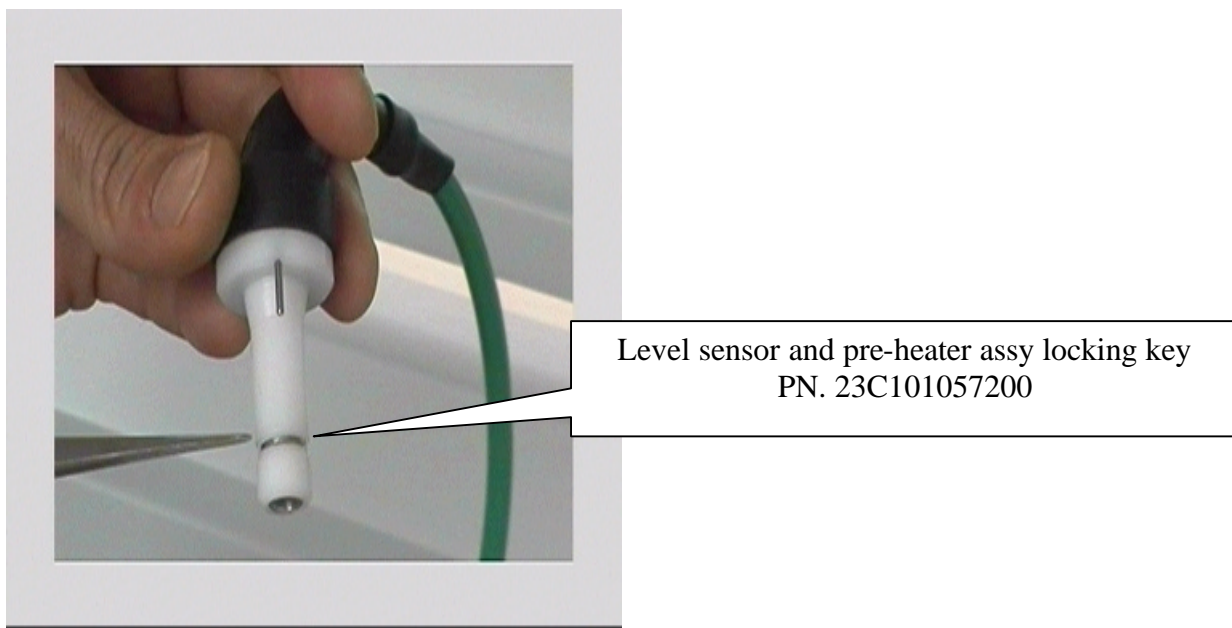


Fig. 7

### 6.1.5.2 Level Sensor functioning Check

1. Turn on the ILab 300 Plus system (instrument and computer).
2. Remove the protective cover of the samples rack.
3. Touch the sampling probe with your finger and make sure that the level sensor LED on the Sampling Interface assy turns on approximately one second, when the probe is touched (Fig. 9).

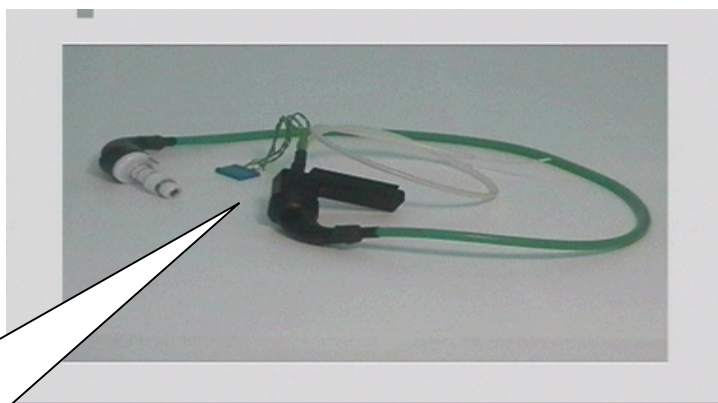


Fig. 8

Level Sensor and Pre-Heater assembly  
PN 23910006300

Sensor Level LED

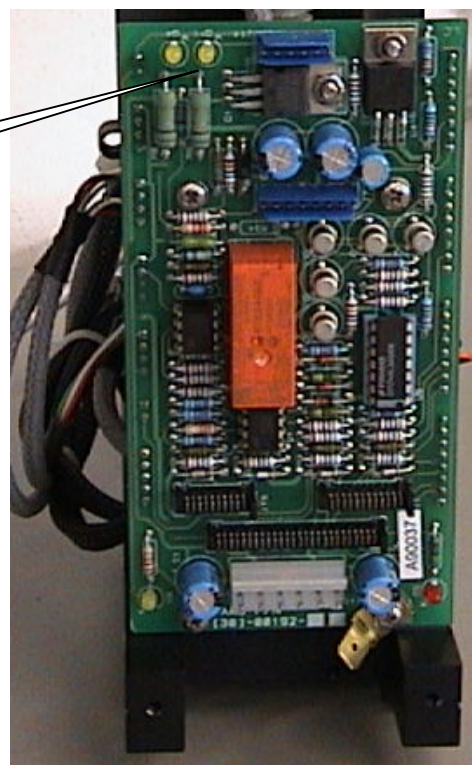


Fig. 9 Sampling Interface assy  
PN. 23300019201

**N.B.:during automatic functioning (routine analysis), every time the sampling probe dips into any liquid whatsoever, the "Level sensor" LED on the Sampling Interface assy will turn on for approximately one second. We advise that you never allow the level sensor to operate when the solution quantity is less than 290 ml.**

4. In order to more specifically verify the correct level-sensing functioning, proceed as follows:
  - a) Place a cup containing liquid in the first samples rack, position nr 1 (Volume  $\geq 290 \mu\text{l}$ ).
  - b) Launch the diagnostic program. Select the "Arm" function and press the "Go Sample" key in order to allow the probe to move into the position corresponding to the sample cup.
  - c) In the window next to the "Go Level", enter the maximum descent of the probe (e.g.500).
  - d) Select "Go Level" and make sure that the probe drops down into the cup until it senses the liquid. In the window to the right of the "Go Level" key the level reached by the probe will appear, expressed in tenths of millimeters
  - e) Repeat step d) several times, each time returning to the Vertical "Home" position. Make sure that the level indicated remains constant = 10 (1 mm).
  - f) Repeat steps b) to e) for positions 8 and 16 of the first rack using the same sample cup. Make sure that the difference between all the positions is = 1mm.
  - g) Repeat steps b) to f) for samples racks three and five (Stat positions 1, 8, 14).
  - h) Press the "Back to Wash Well" key and then select "Go Sample" in order to allow the probe to return to a position in the rack, which does not contain any cup.
  - i) Select "Go Level" and make sure that the level reached is the maximum indicated in the given window is 9999.
  - j) Select the "Reset" key, remove the cups from the racks and exit the program by pressing the "Diagnostic" key.
  - k) To exit the analyzer program press “shutdown” key.

## 6.1.6 DILUTER MODULE

### 6.1.6.1 Substitution of the diluter module

**N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this procedure.**

1. Remove the samples and reagents protective covers (Fig.1).
2. Remove the samples rack and reagents rack.
3. Remove the teflon tube on the diluter head by unscrewing the fastening washer (Fig.2), paying special attention to leakage of liquid.
4. Unscrew the four anchored screws indicated in (Fig.2); unplug the J6, J8, and J13 connectors from the Sampling Interface assy (Fig.3). Remove the diluter and substitute it with the new one.
5. To remount, repeat the above steps in reverse order: from 4 to 1.
6. After mounting the new module, carry out the check procedure for mechanical functioning as described in Section 6.1.6.4

### 6.1.6.2 Substitution of the diluter micro-pump

**N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this procedure.**

1. Remove the samples and reagents protective covers (Fig.1).
2. Remove the samples rack and reagents rack.
3. Remove the teflon tube on the diluter head by unscrewing the fastening washer (Fig.2) paying particular attention to leakage of liquid.
4. Unscrew the four anchored screws indicated in (Fig.2); unplug the J6, J8, and J13 connectors from the Sampling Interface assy (Fig.3).
5. Remove the diluter module from the instrument.
6. Unscrew the two fastening screws on the micro-pump, take out the pump and substitute it with the new one (Fig.2)
7. To remount, repeat the above steps in reverse order: from 6 to 1.
8. After mounting the new micro-pump, carry out the check procedure for mechanical functioning as described in Section 6.1.6.4.



### 6.1.6.3 Substitution of the diluter Electro-valve

**N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this procedure.**

1. Remove the samples and reagents protective covers (Fig.1).
2. Remove the samples rack and reagents rack.
3. Remove the teflon tube on the diluter head by unscrewing the fastening washer (Fig.2) paying particular attention to leakage of liquid.
4. Unscrew the four anchored screws indicated in (Fig.2); unplug the J6, J8, and J13 connectors from the Sampling Interface assy (Fig.3).
5. Remove the diluter module from the instrument.
6. Unscrew the two fastening screws on the electro-valve, take it out and substitute it with the new one (Fig.2)
7. To remount, repeat the above steps in reverse order: from 6 to 1.
8. After mounting the new electro valve, carry out the check procedure for mechanical functioning as described in Section 6.1.6.4.

### 6.1.6.4 Diluter module - Check procedure for mechanical functioning

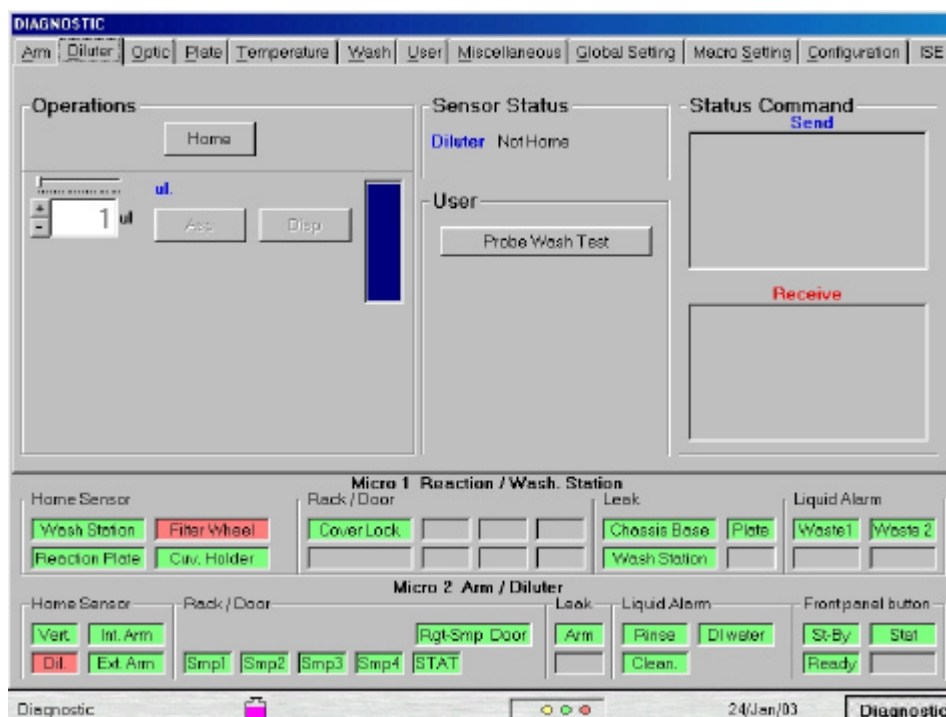
1. Turn on the ILab 300 Plus system (instrument and computer).
2. Remove the protective cover of the samples and of the reagents (Fig.1).
3. Launch the Diagnostic program, select the "Diluter" function" and click on "Home" field. Make sure that the "Home Sensor" Dil. shows green.
4. Select aspiration and then dispense volumes between 3 to 500 microliters clicking on the appropriate windows.
5. Make sure that the diluter both aspirates and dispenses as requested and that there is no leakage of liquid from the hydraulic connections.
6. Exit the Diagnostic program by pressing the "Diagnostic" key.
7. Exit the analyzer program press "shutdown" key.

### 6.1.6.5 BACKLASH

## BACKLASH ADJUSTMENT

By adjusting the “Backlash”, we are able to compensate the backlash of the diluter.

1. Enter the ILab 300 Plus software.
2. Select “Diagnostic”.
3. Select the “Global Setting” folder to reset the Arm
4. Select the “Diluter” folder
5. Run the “Probe Wash Test” twice
6. Select the “Arm” folder
7. Move the probe in the Diluter position. (GO DIL)
8. Return back to the “Diluter” folder



9. Home the Diluter
10. Prepare a sample cup filled with deionized water.
11. Set up the cursor at 300 MicroL
12. By keeping the cup below the probe, (the probe must be submerged into the cup so that it can aspirate the 300 MicrL of liquid), click on “Aspirate” (ASP)
13. Remove the cup. Be careful do not touch the tip of the probe. Any movement of the probe can affect on the adjustment.
14. Set up the cursor at 1 MicroL.
15. Look carefully at the tip of the probe while selecting “DISP”
16. Count the number of click with the DISP button before a drop of liquid appears. (A small bubble).
17. The number of clicks minus 1 is the value that must be set in “Backlash DIL “ in the
18. “Configuration” folder.
19. Save the new value.



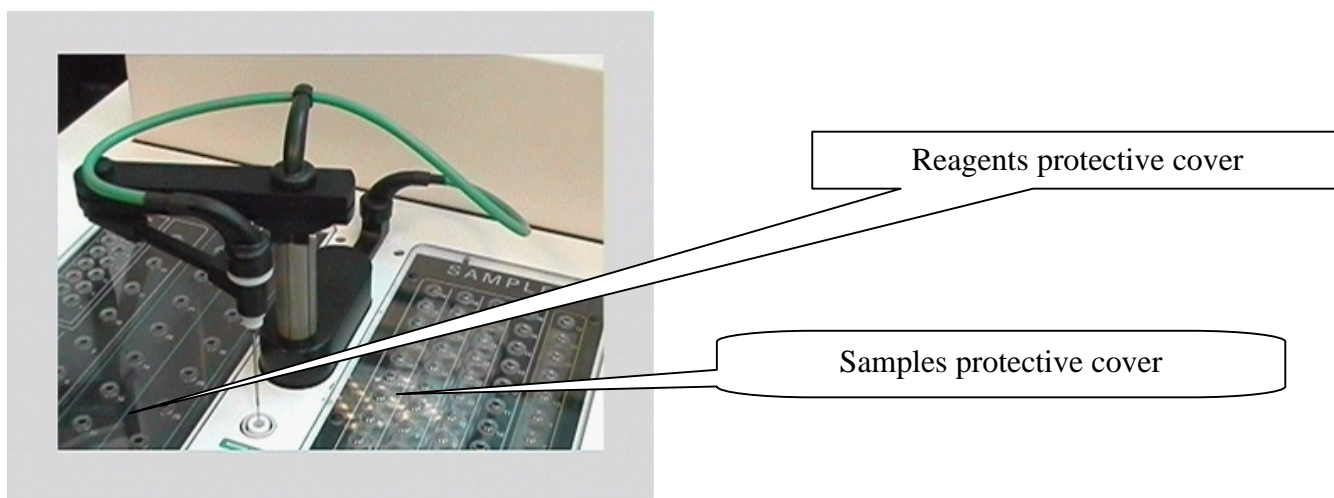


Fig.1

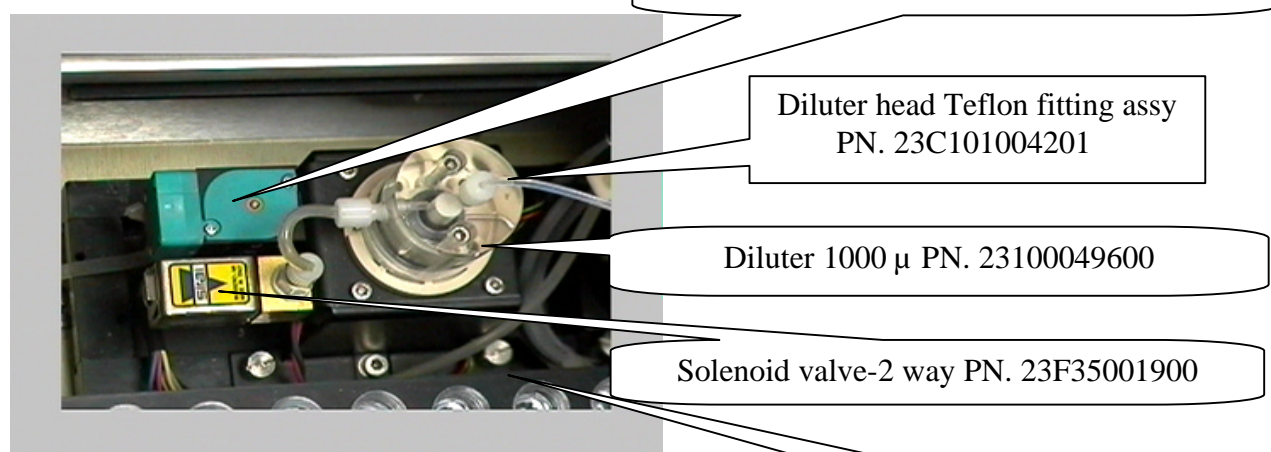


Fig. 2 Diluter Module

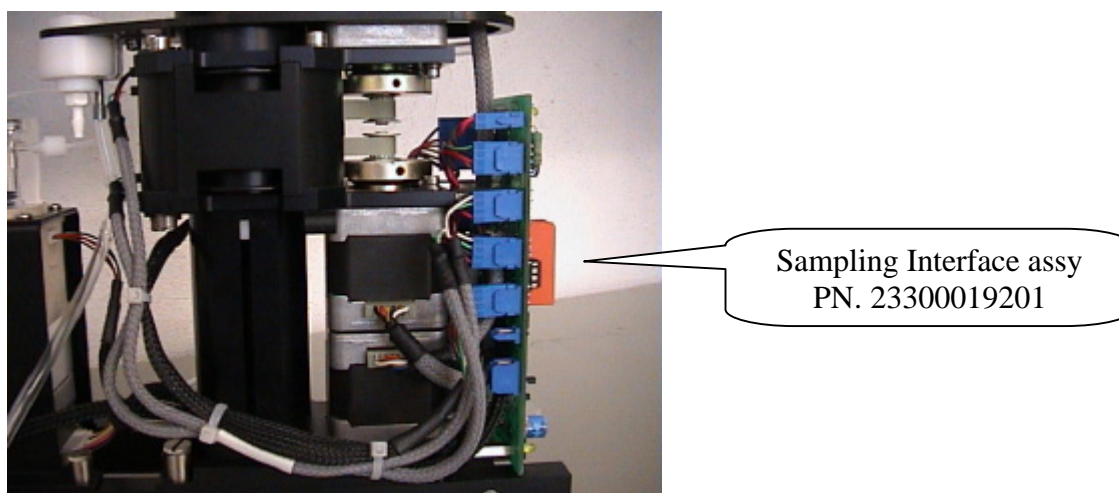


Fig. 3

## 6.2 OPTIC ASSEMBLY

### 6.2.1 SUBSTITUTION OPTIC ASSEMBLY

- **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**

1. Remove the cover panel of the reaction plate (Fig.1).
2. Unscrew the four cross-bar fastening screws and remove the cross-bar (Fig.1).
3. Unscrew the 4 fastening screws on the reaction plate assembly (Fig.2).
4. Remove the J1, J2 connectors and the ground wire from the Reaction Chamber Motor Interface board (Fig.2).
5. Remove the J1, J2, J14 connectors and the ground wire from the Reaction Plate Interface board (Fig.2).
6. Open the peristaltic pump assembly panel cover. Unscrew the two fastening screws to make free the support that join the seven tubes coming from the washing station probes.
7. Disconnect the six tubes coming from the washing station probes and the air tube inserted on the micro pump P4 coming from the porous pad on fifth probe of the washing station.
8. Remove the reaction plate from the instrument.
9. Remove the lamp's power supply wires connected to the Lamp Regulator board by loosening the two screws on the M1 connector (Fig. 4).
10. Remove the J1 connectors from the Preamp/ADC (Fig.5) and the ground wire from the optic group.
11. Unscrew the four fastening screws, remove the optic assembly and substitute it with the new one (Fig. 3).
12. Remount the assembly, repeating the above steps, 11 through 1, in reverse order.
13. After remounting, reseal the lamp and carry out the electronic adjustment check procedures as described in Sections 6.2.2 and 6.2.3

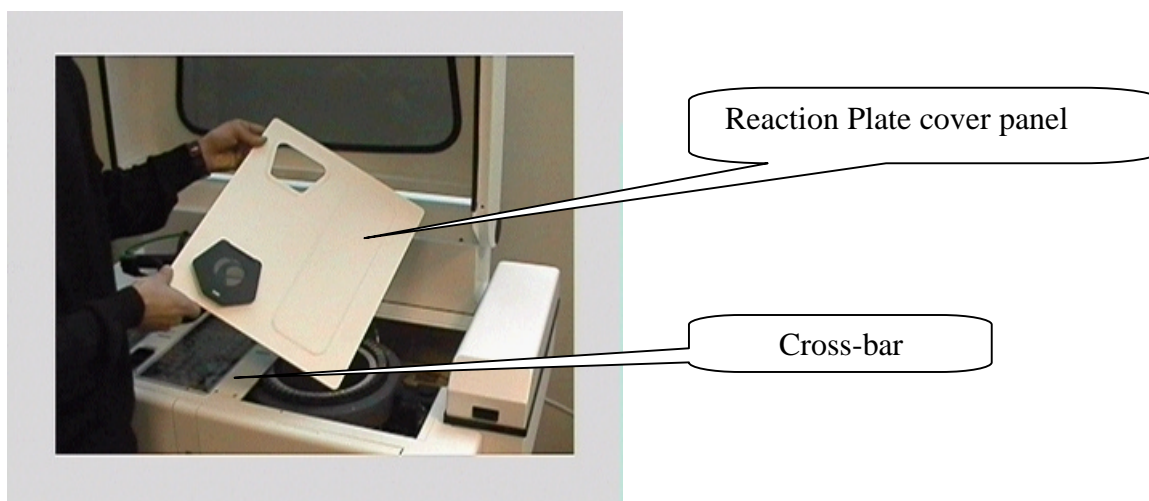


Fig. 1

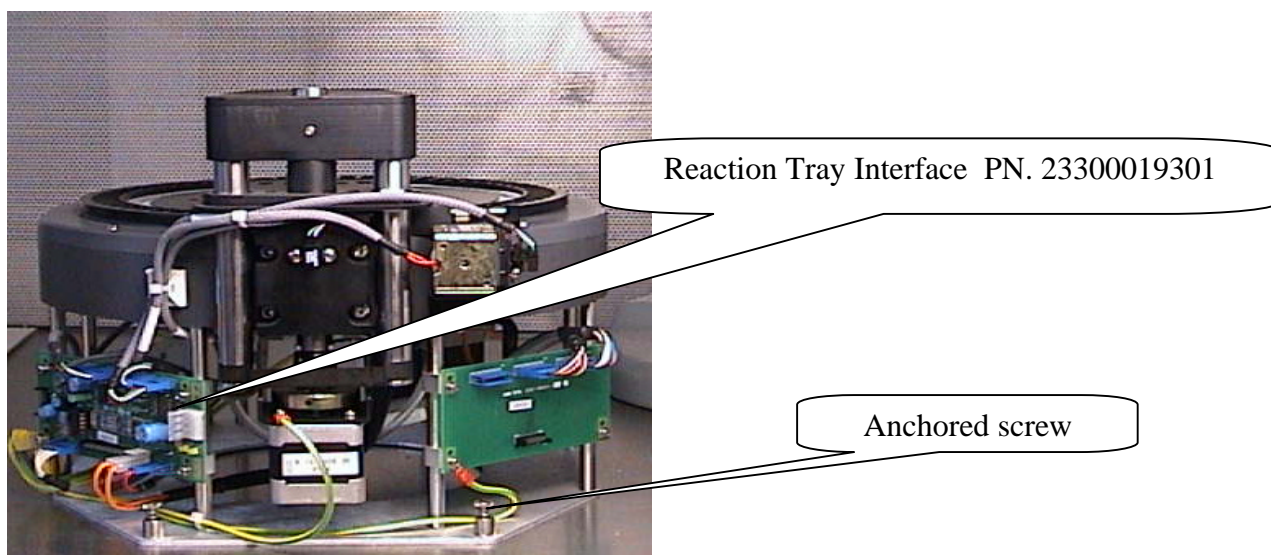


Fig. 2 Reaction Plate assembly PN. 23910006101 \*

\*Reaction Plate assembly is supplied without Washing station and Photometer

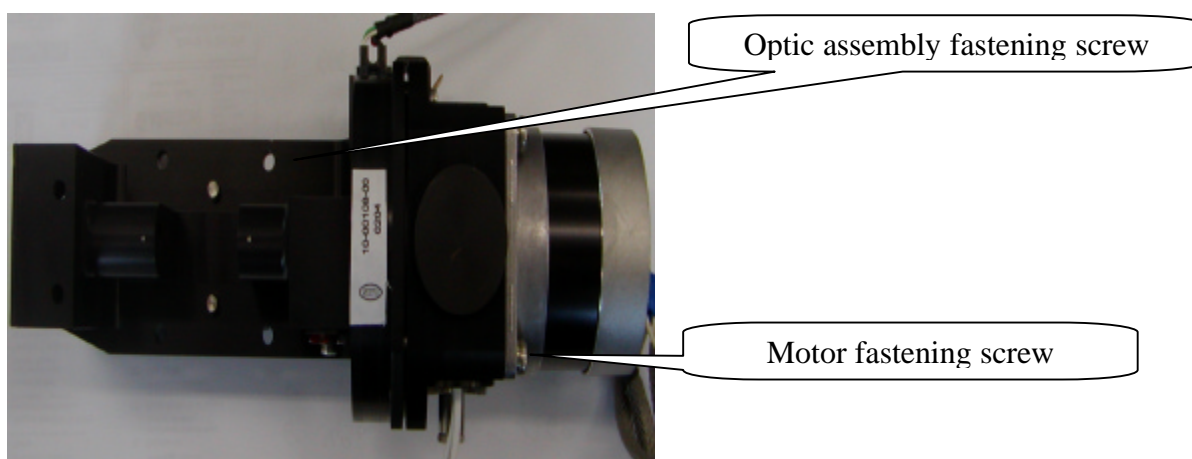



Fig. 3 New photometer assy (2optic Channel)  
PN. 23100010800

## 6.2.2 SUBSTITUTION OF THE PHOTOMETER LAMP

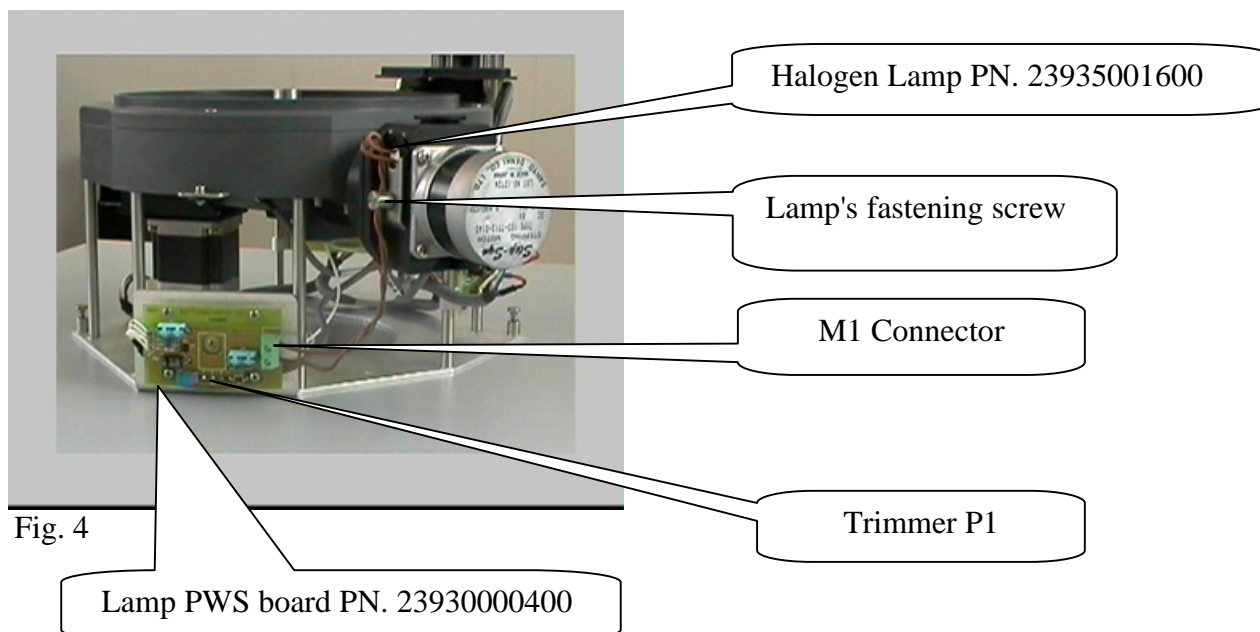➤ **N.B.:** Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.

1. Remove the cover panel of the reaction plate (Fig.1).
2. Remove the lamp's electrical wires connected to the Lamp Regulator board by loosening the two screws on the M1 connector (Fig.4).
3. Unscrew the lamp's fastening screws and remove it from its housing.
4. Remount the new lamp, repeating the above steps, 3 through 1, in reverse order.
5. Carry out the electronic adjustment of the photometer as described in Section 6D.2.1.

**Warning:** Do not touch the glass portions of the lamp.



### 6.2.3 ELECTRONIC ADJUSTMENT OF THE PHOTOMETER

**Warning:** make sure that the photometer lamp and the reaction cuvettes are brand-new before applying the following procedure.

1. Turn on the "ILab 300 Plus" system (first the instrument and then the computer). Remove the reaction plate cover (Fig.1).
2. Launch the diagnostic program, select the "Plate" function and reset by clicking on the appropriate key.
3. Make sure that the "Home Sensor" Reaction Plate window lights up green.
4. Place 500 µL of distilled water in cuvette #31 of the reaction plate.
5. Select the "Optic" function and turn the lamp "On" by clicking on the appropriate field.
6. Reset the filter wheel by clicking on the appropriate key.
7. Make sure that "Home Sensor" Filter Wheel window lights up green.
8. Make sure that the Regulator board for Connector M1 shows 6.0 volts  $\pm$  0.1. If not, turn the P1 trimmer (Fig.4).

**N.B.:** The reaction cuvette must be perfectly clean. If not, carry out a "Cuvette Wash" cycle (Start Work)

**Warning:** Room lighting and daylight can influence the reading. To avoid this interference, replace the reaction plate cover panel before carrying out any check procedures.

9. Click on "Start Sample Channel Conversion" key and make sure that the adjacent window shows a count value from 50 to 100. If not, adjust the Preamp/ADC board "Offset" trimmer inside the analysis plate in order to set the desired value (Fig. 5).
10. After regulating of the main channel offset, it is necessary to click on " Stop Sample Channel Conversion" key.
11. Click on "Start Reference Channel Conversion" and make sure that the adjacent window shows a count value from 50 to 100. If not, adjust the Preamp/ADC board "Offset" trimmer of the reference channel located external to the analysis plate in order to set the desired value (Fig. 5).
12. After regulating of the reference channel offset, it is necessary to click on " Stop Reference Channel Conversion" key.




13. Select the filter positions from 1 to 8 (from 340 nm to 620 nm) and click on "Start/Stop Sample Channel Conversion" key and determine **which filter transmit the highest signal.**
14. Verify in the adjacent window a count value of  $55.000 \pm 1000$ . If not, adjust the Preamp/ADC board "Gain" trimmer inside the analysis plate in order to set the desired value (Fig. 5).
15. After regulating the Gain, select all the other filters and clicking on "Start/Stop Sample Channel Conversion" verify count values between 29000 and 56000

**Warning:** Room lighting and daylight can influence the reading. To avoid this interference, replace the reaction plate cover panel before carrying out any check procedures.

16. Select the highest signal filter, identified in point 13. Click on "Start Reference Channel Conversion" and make sure that the adjacent window shows a count value of  $40.000 \pm 1000$ . If not, adjust the Preamp/ADC board "Gain" trimmer external to the analysis plate in order to set the desired value (Fig. 5).
17. After regulating the Gain, select all the other filters and clicking on "Start/Stop Reference Channel Conversion" verify count values between 29000 and 41000.
18. After concluding these operations, empty cuvette #31 and exit the diagnostic program by clicking on the "Diagnostic" key.
19. Exit the analyzer program press "shutdown" key.

Fig. 5 New photometer assy (2optic Channel)

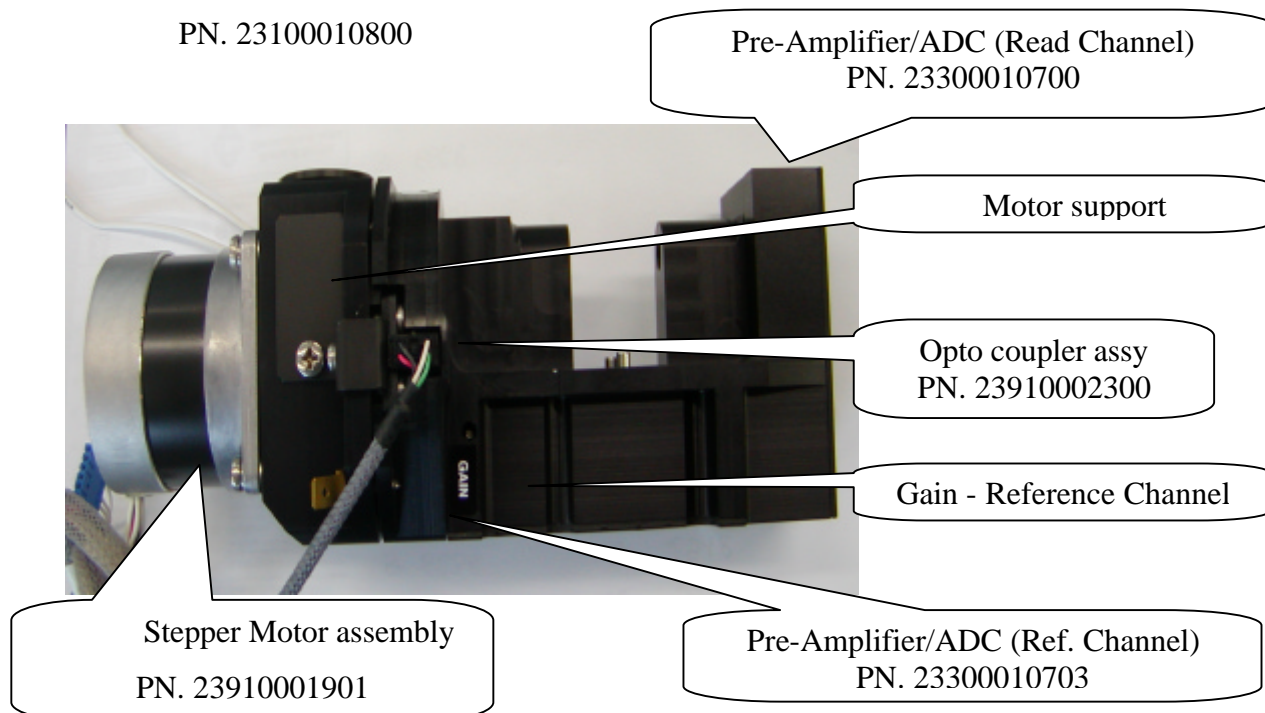PN. 23100010800



#### 6.2.4 SUBSTITUTION OF THE OPTO SENSOR

➤ **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**

1. Remove the cover panel of the reaction plate (Fig.1).
2. Unscrew the four cross-bar fastening screws and remove the cross-bar (Fig.1).
3. Unscrew the 4 fastening screws on the reaction plate assembly (Fig.2).
4. Remove the J1, J2 connectors and the ground wire from the Reaction Chamber Motor Interface board (Fig.2).
5. Remove the J1, J2, J14 connectors and the ground wire from the Reaction Plate Interface board (Fig.2).
6. Open the peristaltic pump assembly panel cover. Unscrew the two fastening screws to make free the support that join the six tubes coming from the washing station probes.
7. Disconnect the seven tubes coming from the washing station probes and the air tube inserted on the micro pump P4 coming from the porous pad on fifth probe of the washing station.
8. Remove the reaction plate from the instrument.
9. Remove the lamp's power supply wires connected to the Lamp Regulator board by loosening the two screws on the M1 connector (Fig. 4).
10. Unscrew the lamp's fastening screws and remove it from its housing.
11. Unscrew the four fastening screws on the motor support in order to access the filter wheel.
12. Position the filter wheel in such a manner as to make it possible to see and unscrew the internal opto sensor support screw. Unscrew the external opto sensor support screw.
13. Unscrew the optic sensor fastening screws from its support and substitute with the new one.
14. Remount the assembly, repeating the above steps, 12 through 1, in reverse order.
15. After remounting, reseal the lamp and carry out the electronic adjustment procedures as described in Sections 6.2.3

### 6.2.5 SUBSTITUTION OF THE MOTOR

➤ **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**

1. Remove the cover panel of the reaction plate (Fig.1).
2. Unscrew the four cross-bar fastening screws and remove the cross-bar (Fig.1).
3. Unscrew the 4 fastening screws on the reaction plate assembly (Fig.2).
4. Remove the J1, J2 connectors and the ground wire from the Reaction Chamber Motor Interface board (Fig.2).
5. Remove the J1, J14 connectors and the ground wire from the Reaction Plate Interface board (Fig.2).
6. Open the peristaltic pump assembly panel cover. Unscrew the two fastening screws to make free the support that join the six tubes coming from the washing station probes. Disconnect the seven tubes coming from the washing station probes and the air tube inserted on the micro pump P4 coming from the porous pad on fifth probe of the washing station.
7. Remove the reaction plate from the instrument.
8. Remove the lamp's power supply wires connected to the Lamp Regulator board by loosening the two screws on the M1 connector (Fig. 4).
9. Unscrew the lamp's fastening screws and remove it from its housing.
10. Unscrew the four fastening screws on the motor support in order to access the pulley.
11. Unscrew the four fastening screws on the motor; loosen the pulley fastening screws (Fig.3).
12. Remove the motor and substitute it with the new one.
13. Remount the assembly, repeating the above steps, 12 through 1, in reverse order.

**Warning:** when tightening the motor pulley fastening screws, make sure that the filter wheel is centered with respect to the opto sensor.

14. After remounting, reseal the lamp and carry out the electronic adjustment procedures as described in Sections 6.2.2 and 6.2.3



## SUBSTITUTION OF THE OPTIC FILTERS

**Interferential filters Kit PN. 23965002900**

➤ **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**

1. Remove the cover panel of the reaction plate (Fig.1).
2. Unscrew the four cross-bar fastening screws and remove the cross-bar (Fig.1).
3. Unscrew the 4 fastening screws on the reaction plate assembly (Fig.2).
4. Remove the J1 connector and the ground wire from the Reaction Chamber Motor Interface board (Fig.2).
5. Remove the J1, J14 connectors and the ground wire from the Reaction Plate Interface board (Fig.2).
6. Open the peristaltic pump assembly panel cover. Unscrew the two fastening screws to make free the support that join the six tubes coming from the washing station probes. Disconnect the seven tubes coming from the washing station probes and the air tube inserted on the micro pump P4 coming from the porous pad on fifth probe of the washing station.
7. Remove the reaction plate from the instrument.
8. Remove the lamp's power supply wires connected to the Lamp Regulator board by loosening the two screws on the M1 connector (Fig. 4).
9. Unscrew the lamp's fastening screws and remove it from its housing.
10. Unscrew the four fastening screws on the motor support in order to access the filter wheel.
11. Unscrew the fastening screw that holds the elastic band that holds the filters.
12. Remove the filter from its housing and substitute it with the new one.
13. Remount the assembly, repeating the above steps, 11 through 1, in reverse order.
14. After remounting, reseal the lamp and carry out the electronic adjustment procedures as described in Sections 6.2.3.

**Warning:** Do not touch the glass parts of the filter.

**Warning:** when substituting an interferential filter due to malfunctioning or functional anomalies, it is always advisable to substitute the entire set of filters mounted on the instrument.


## 6.2 PUMP ASSEMBLY

### 6.3.1 SUBSTITUTION OF THE PUMP ASSEMBLY

**N.B.:** Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.

1. Unscrew the two knobs to remove the door located on the front elevation of the instrument (Fig. 1) and the two noticeable screws located on the rear panel (Fig. 2) to access the washing pump assembly.

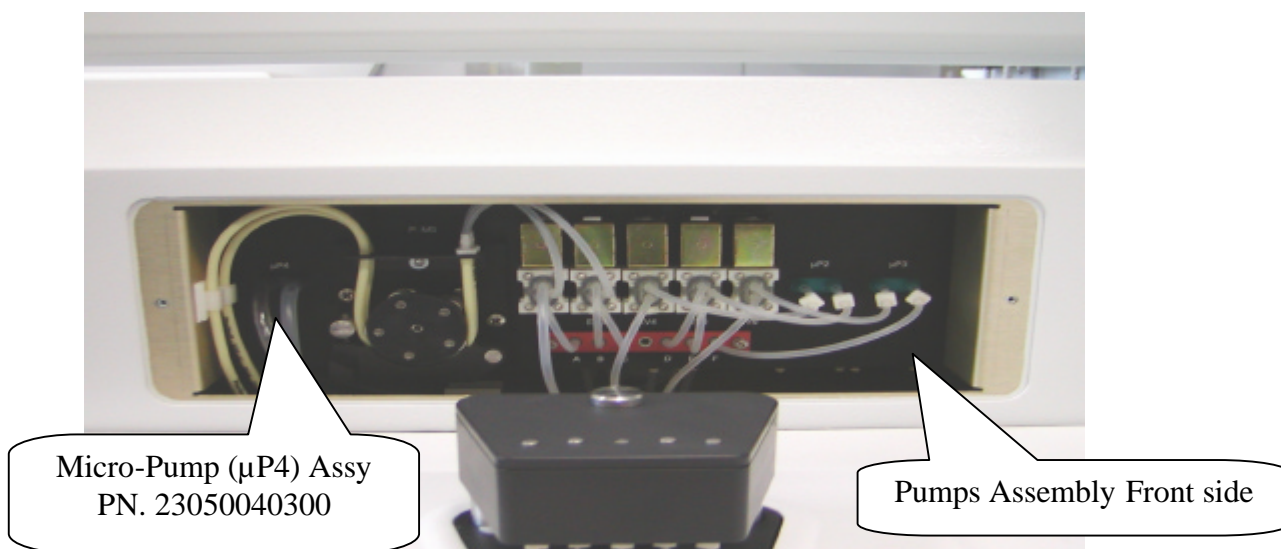Fig. 1



2. Disconnect connectors J1, J2 e J8 from the Hydraulics Interface board (rear side Fig.2)
3. Disconnect the air tubes from  $\mu$ P4 and the waste tubes coming from the peristaltic pump (front side, Fig. 1)
4. Disconnect the three tubes inserted on the EV1, EV4 ed EV6 valves, respectively coming from the washing well, from the distilled water bottle and from the cleaning solution bottle (front side, Fig. 1)
5. Unscrew the four screws which hold the Pumps Assembly. Remove it from its site (rear side Fig.2) and replace the Assembly with the new one.
6. Remount, repeating the above steps, 5 through 1, in reverse order.

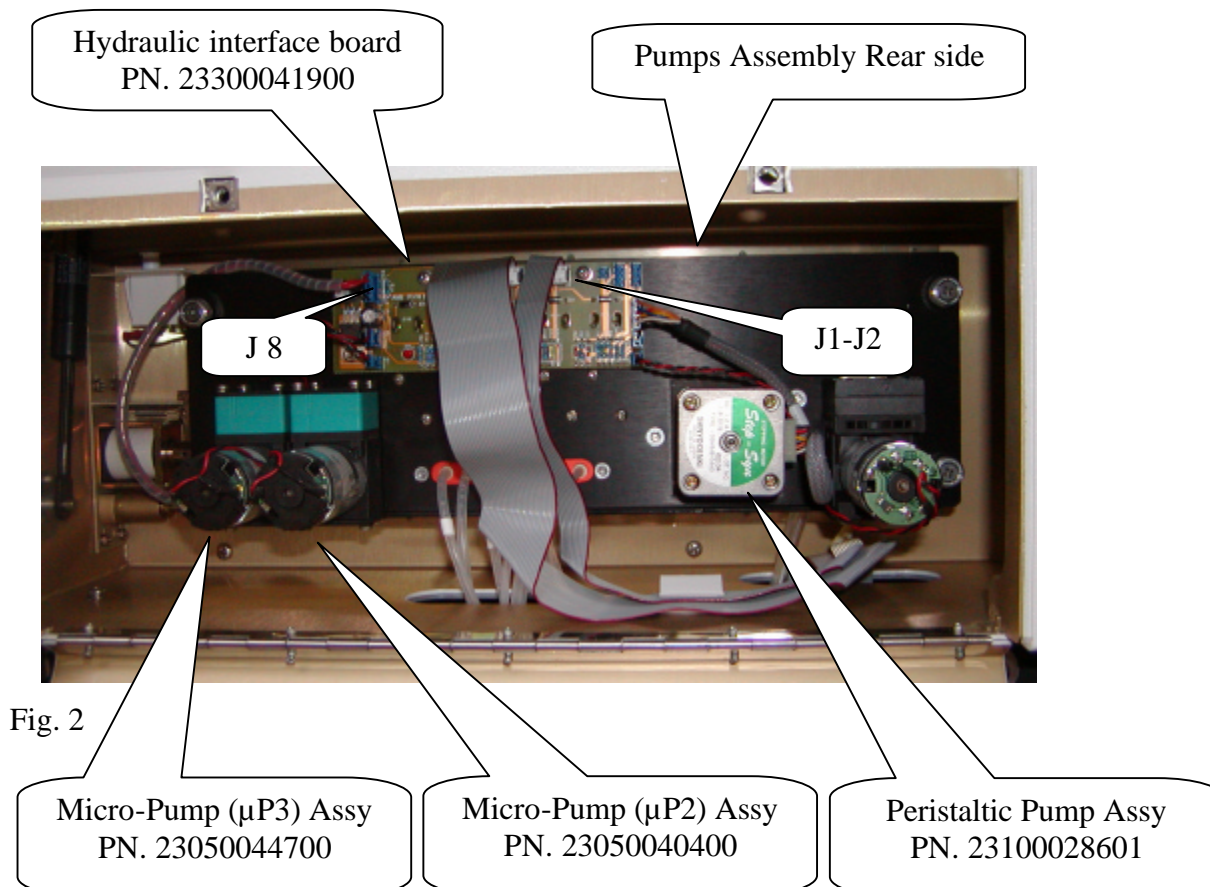


Fig. 2

### 6.3.2 SUBSTITUTION OF THE PERISTALTIC PUMP MOTOR

- **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**
1. Unscrew the two knobs to remove the door located on the front elevation of the instrument (Fig. 1) and the two noticeable screws located on the rear panel (Fig. 2) to access the washing pump assembly.
  2. Unlock the hinged tube guide by lowering and rotating the guide lock to the left in order to free the tubes (Fig. 3).
  3. Disconnect the connector J6 from the Hydraulics Interface Board (Rear side Fig.2)
  4. Unscrew the two motor fastening screws and remove the peristaltic pump assembly from its site (Front side Fig.3)
  5. Unscrew the four screws which hold the motor and loosen the two grains which hold the peristaltic wheel.
  6. Remove the motor and substitute it with the new one.
  7. Remount, repeating the above steps, 6 through 1, in reverse order.

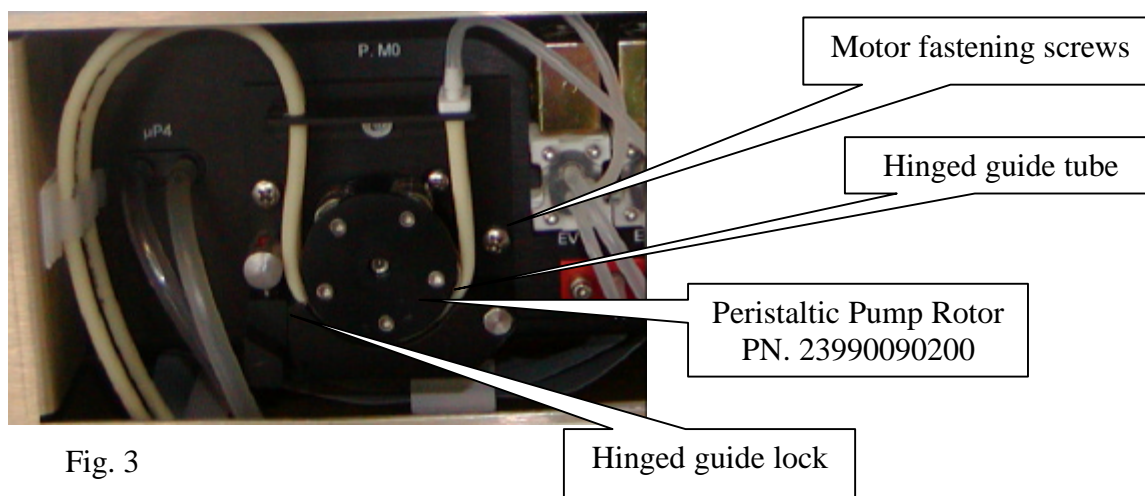


Fig. 3

### 6.3.2.1 Function check for the peristaltic pump

1. Turn on the ILab 300 Plus system and launch the “Analyzer” program. Select "Start Work" and perform one or more cuvette washing cycles. Make sure that the sampling probe washing well empties completely without any leakage of liquid. If leakage should occur, adjust the screw indicated in Fig. 5 using the Peristaltic Pump Setting Gage as illustrated in figure 6 (See Service News N° 9 dated 21.05.2003)
2. Close the "ILab 300 Plus" program by clicking on "Shutdown".

### 6.3.3 SUBSTITUTION OF PART(S)

➤ **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**

1. Unscrew the two knobs to remove the door located on the front elevation of the instrument (Fig. 1) and the two noticeable screws located on the rear panel (Fig. 2) to access the washing pump assembly.
2. For each part to be replaced, if necessary, disconnect the relative electrical power connection, the hydraulic tubes and the fastening screws.
3. If necessary, unscrew the four anchored screws of the assembly and turn it ahead (Fig.2).
4. Remount, repeating the above steps, 3 through 1, in reverse order.

Any necessary adjustments (in order to optimize pump aspiration in relation to the newly substituted tube) can be made by turning the relative screw (Fig. 5) using the Peristaltic Pump Setting Gage as illustrated in figure.

**N.B.: If power loss is experienced in the filling or emptying of liquid (after a period of instrument activity), perform the following operations before substituting or servicing any part(s):**

- Disconnect the liquid input tube from the malfunctioning pump.
- Insert a small amount of liquid using a syringe from the disconnected ends in order to wet the membranes (Fig. 4).

Remove the syringe and reconnect the tube(s); perform functioning check as described in Section 6.3.2.1.

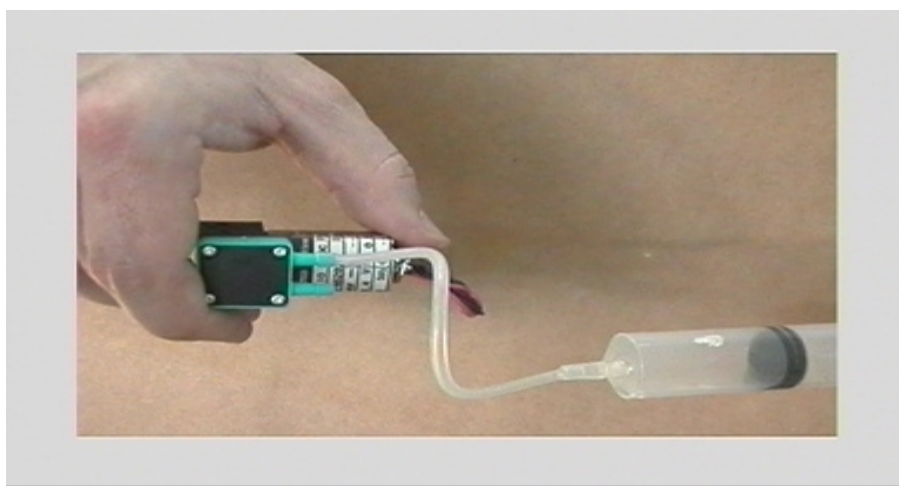


Fig. 4

#### 6.3.4 MECHANICAL ADJUSTMENT OF THE PERISTALTIC PUMP

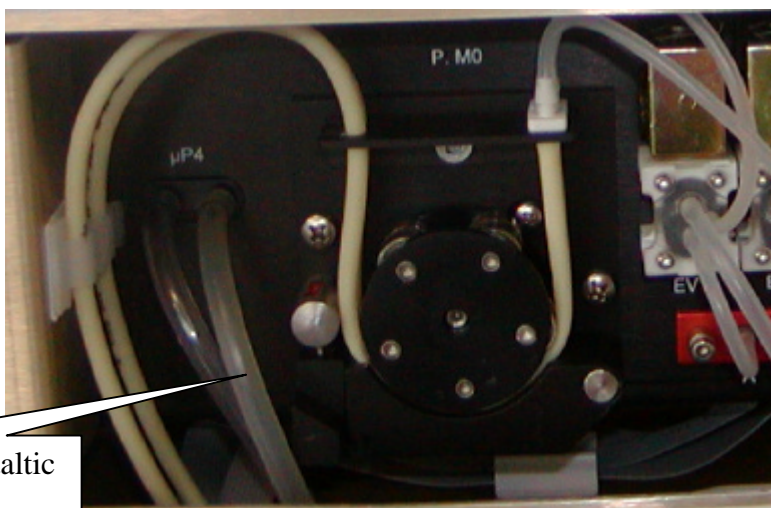
- **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this procedure.**

The peristaltic pump tubes must be substituted every three to six months, depending on instrument use.

### 6.3.4.1 How to substitute peristaltic pump tubes

1. Unscrew the two knobs on the cover, located on the instrument 'headboard' behind the washing station, in order to arrive at the pumps assembly (Fig. 1).
  2. Unhook the hinged guide by lowering and rotating the guide's blocking bracket to the left thereby freeing the tubes (Fig.3).
  3. Pull the tubes off their relative nipples and substitute with new ones.
  4. Remount, repeating the above steps, 2 through 1, in reverse order.
- Always use original replacement parts; never lubricate the peristaltic pump roller bearings.
  - After prolonged instrument inactivity: verify the efficiency of the peristaltic pump tubes.

**Please note:** if the instrument is not able to completely empty the washing well, check the condition of the tube and make sure that the space between the hinged guide and the rotor is approximately 1.5 mm. Any necessary adjustments (in order to optimize pump aspiration in relation to the newly substituted tube) can be made by turning the relative screw (Fig. 5) using the Peristaltic Pump Setting Gage as illustrated in figure 6 (See Service News N° 9 dated 21.05.2003)



Screw for regulating peristaltic pump aspiration

Fig. 5






Peristaltic Pump Setting Gage  
PN. 23010120100

Fig. 6

## 6.4 REACTION PLATE
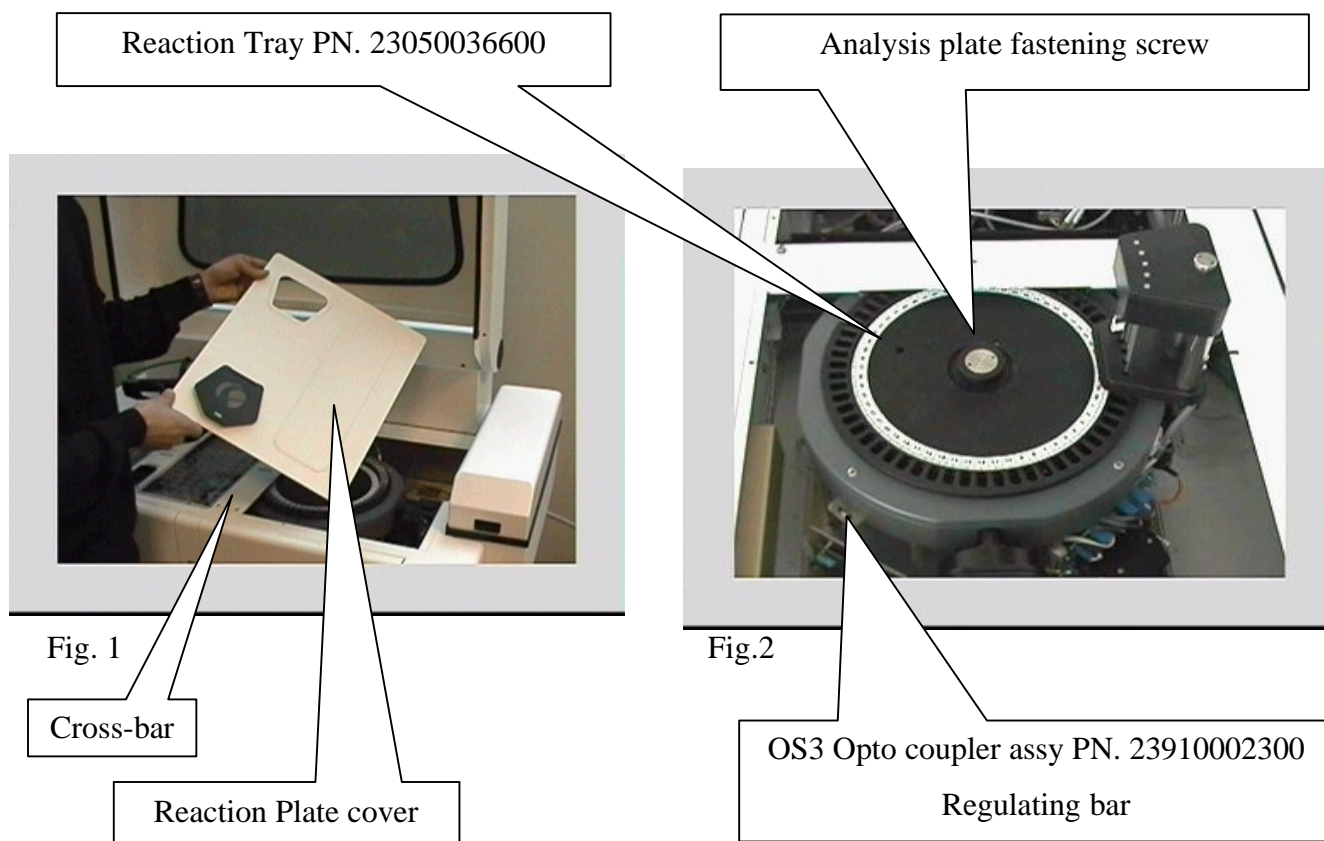
### 6.4.1 REACTION PLATE - MECHANICAL ADJUSTMENT

1. Turn on the ILab 300 Plus system and launch the “Analyzer” program. Remove the reaction plate cover (Fig. 1).
2. Launch the diagnostic program, select the "Optic" function and turn on the light by clicking on the designated key.
3. Reset the filter wheel by clicking on the appropriate key.
4. Select key “Spare” to bring filter position # 9 to the 'reading' position.
5. Select the "Plate" function. Reset by clicking on the appropriate key. Make sure that the "Home Sensor" Reaction Plate window lights up green.
6. Make sure that the light beam inside the cuvette is centered with respect to the sides of the cuvette. If not, adjust the OS3 optic sensor regulating bar (Fig. 2).
7. To exit the diagnostic program press the "Diagnostic" key
8. To exit the analyzer program press “shutdown” key



#### 6.4.2 SUBSTITUTION OF THE REACTION PLATE

➤ **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure**

1. Remove the reaction plate cover (Fig. 1).
2. Unscrew the four cross-bas fastening screws and remove the cross-bar (Fig. 1).
3. Unscrew the four anchored screws of the reaction plate module (Fig. 3).
4. Remove the J1 connector and the ground wire from the Reaction Chamber Motor Interface.
5. Remove the J1, J14 connectors and the ground wire from the Reaction Plate Interface.
6. Remove the top cover of the washing station using the special wrench supplied by the manufacturer.
7. Disconnect the seven tubes attached to the probes on the washing station.
8. Remove the reaction plate from the instrument and substitute it with the new one.
9. Remount, repeating the above steps, 7 through 1, in reverse order.



**N.B.: After substituting the reaction plate, carry out the mechanical adjustment procedure of the reaction plate as described in Section 6.4.1**

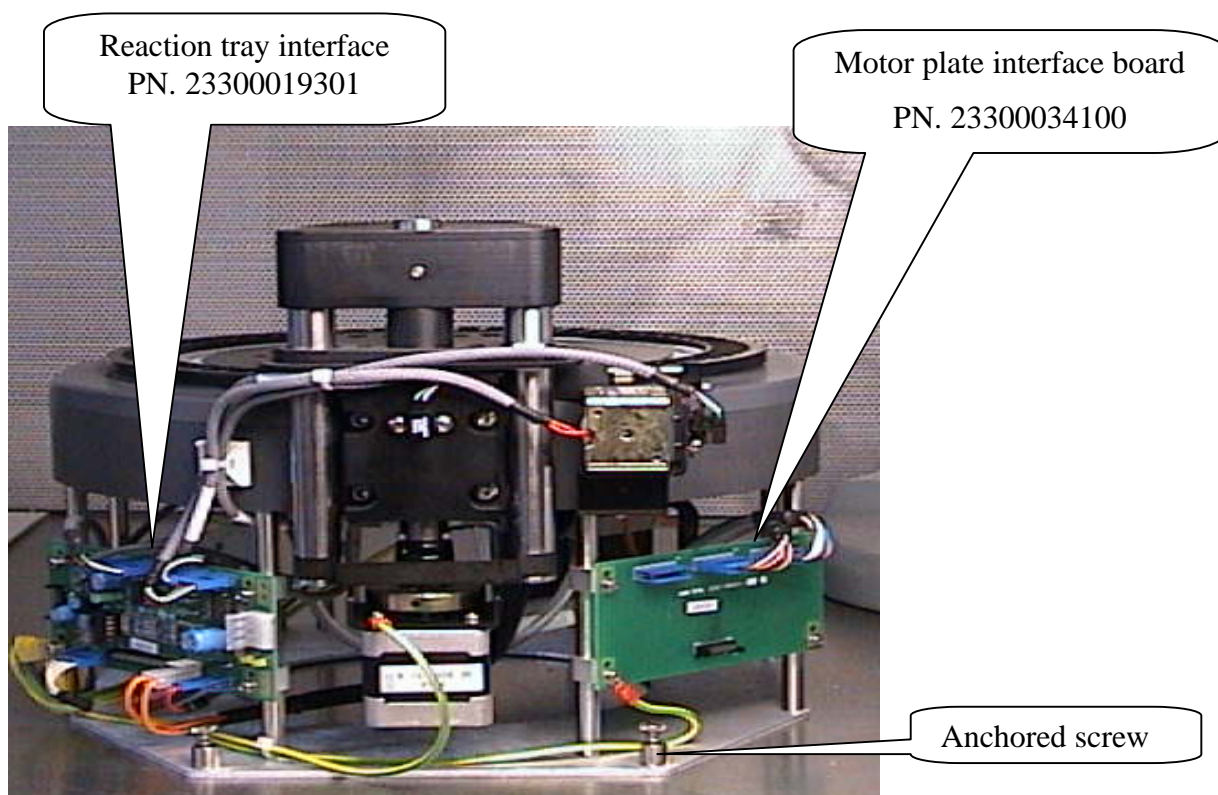


Fig. 3 Reaction Plate assembly PN. 23910006101 \*

\* Reaction Plate assembly is supplied without Washing Station and Photometer

### 6.4.3 SUBSTITUTION OF THE OS3 OPTIC SENSOR

➤ **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**

1. Remove the reaction plate cover (Fig. 1).
2. Unscrew the four cross-bas fastening screws and remove the cross-bar (Fig. 1).
3. Unscrew the four anchored screws of the reaction plate module (Fig. 3).
4. Remove the J1 connector and the ground wire from the Reaction Chamber Motor Interface(Fig. 3).
5. Remove the J1, J4 and J14 connectors and the ground wire from the Reaction Plate Interface (Fig.3).
6. Open the cover panel of the peristaltic pump assembly and disconnect the seven tubes attached to the probes on the washing station.
7. Remove the reaction plate from the instrument.

8. Remove the fastening screw from the analysis plate using the special tool supplied by the manufacturer (Fig. 2).
9. Unscrew the two fastening screws from the bottom cover panel of the washing station and lift it. (Fig. 4).
10. Unscrew the two fastening screws on the cuvette holder assembly and lift it. This will make it possible to remove the analysis plate (Fig. 5).
11. Unscrew the two fastening screws on the optic sensor and substitute it with the new one.
12. Remount, repeating the above steps, 11 through 1, in reverse order.

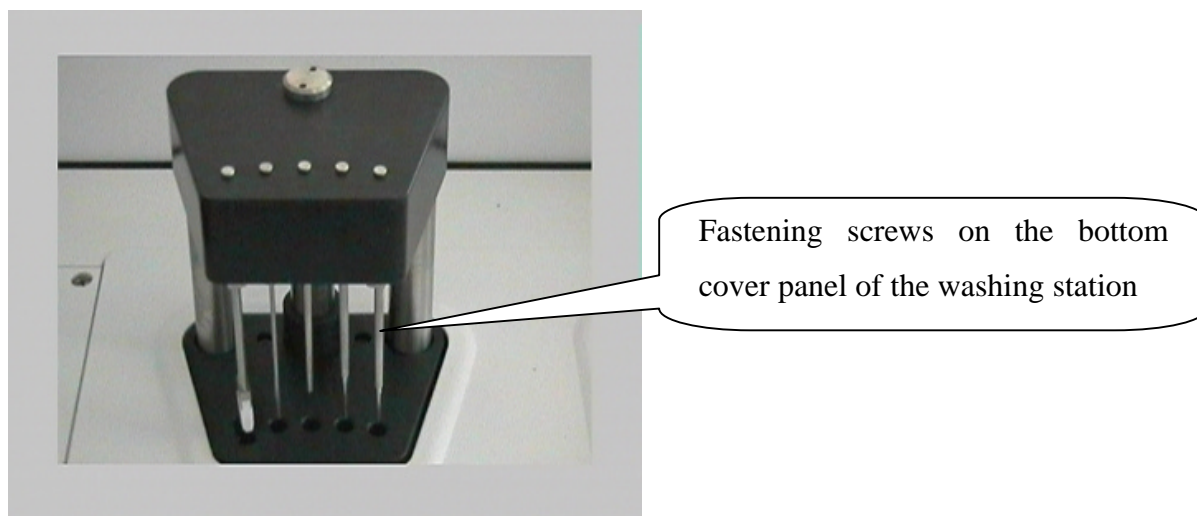


Fig. 4 Washing Station probes assembly PN. 23910000600

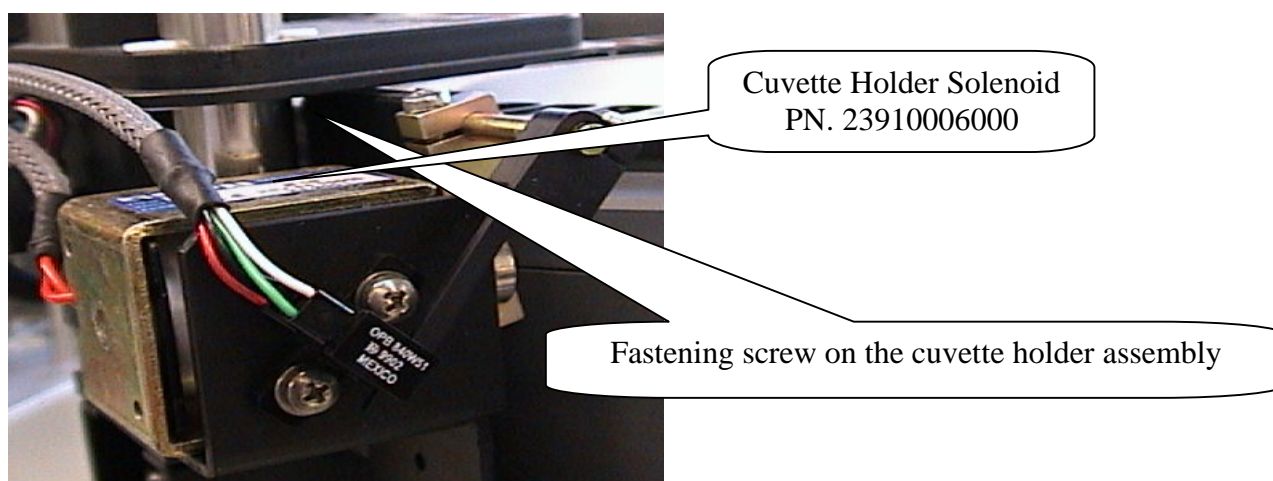


Fig. 5 Cuvette Holder

**N.B.: After substituting the optic sensor, carry out the mechanical adjustment procedure of the reaction plate as described in Section 6.4.1**

#### 6.4.4 SUBSTITUTION OF THE M3 MOTOR

➤ **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**

1. Remove the reaction plate cover (Fig. 1).
2. Unscrew the four cross-bas fastening screws and remove the cross-bar (Fig. 1).
3. Unscrew the four anchored screws of the reaction plate module (Fig. 3).
4. Remove the J1 and J4 connectors and the ground wire from the Reaction Chamber Motor Interface (Fig. 3).
5. Remove the J1, J14 connectors and the ground wire from the Reaction Plate Interface (Fig. 3).
6. Open the cover panel of the peristaltic pump assembly and disconnect the seven tubes attached to the probes on the washing station.
7. Remove the reaction plate from the instrument.
8. Unscrew the four fastening screws that hold the motor support (Fig. 7).
9. Remove the four fastening screws on the motor (Fig. 8).
10. Remove the motor and substitute it with the new one.
11. Remount, repeating the above steps, 9 through 1, in reverse order.

**Warning:** Belt tension must be regulated by turning the M4 screw (shown in Fig. 7) clockwise by 1.3 Kg/cm.


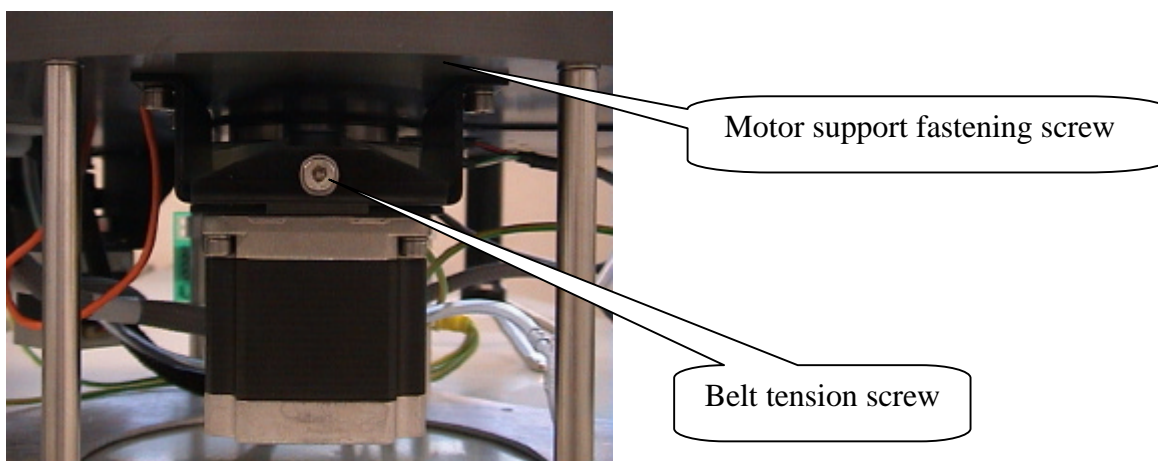


Fig. 7

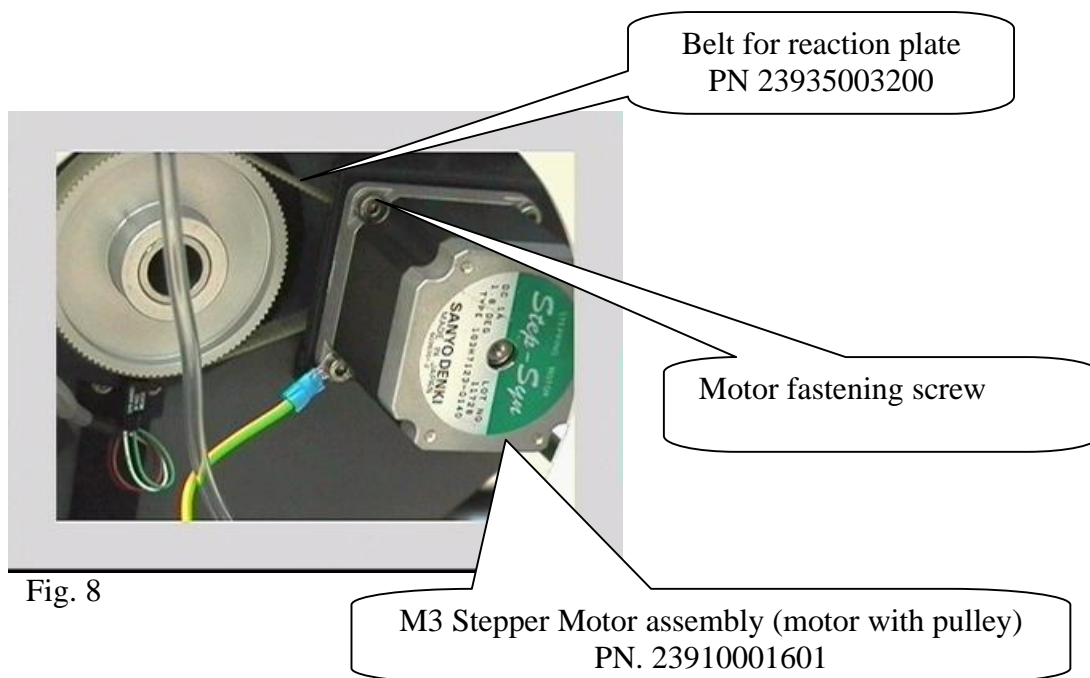

**N.B.: After substituting the motor, carry out the mechanical adjustment procedure of the reaction plate as described in Section 6.4.1**

## 6.5 WASHING STATION ASSEMBLY

### 6.5.1 SUBSTITUTION OF THE WASHING STATION ASSEMBLY

➤ **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**

1. Remove the reaction plate cover (Fig. 1).
2. Unscrew the four cross-bas fastening screws and remove the cross-bar (Fig. 1).
3. Unscrew the four anchored screws of the reaction plate module (Fig. 2).
4. Remove the J1 and J3 connectors and the ground wire from the Reaction Chamber Motor Interface (Fig. 2).
5. Remove the J1, J3, J5, J6, J14 connectors and the ground wire from the Reaction Plate Interface (Fig. 2).
6. Remove the top cover of the washing station using the special tool supplied by the manufacturer (Fig. 5).
7. Disconnect the seven tubes attached to the probes on the washing station (Fig. 6).
8. Remove the reaction plate from the instrument.
9. Unscrew the four fastening screws on the washing station assembly (Fig. 3).
10. Remove the washing station assembly and replace it with the new one.
11. Remount, repeating the above steps, 10 through 1, in reverse order.
12. After substituting and remounting, carry out the mechanical adjustment and functional check procedures as described in section 6.5.4

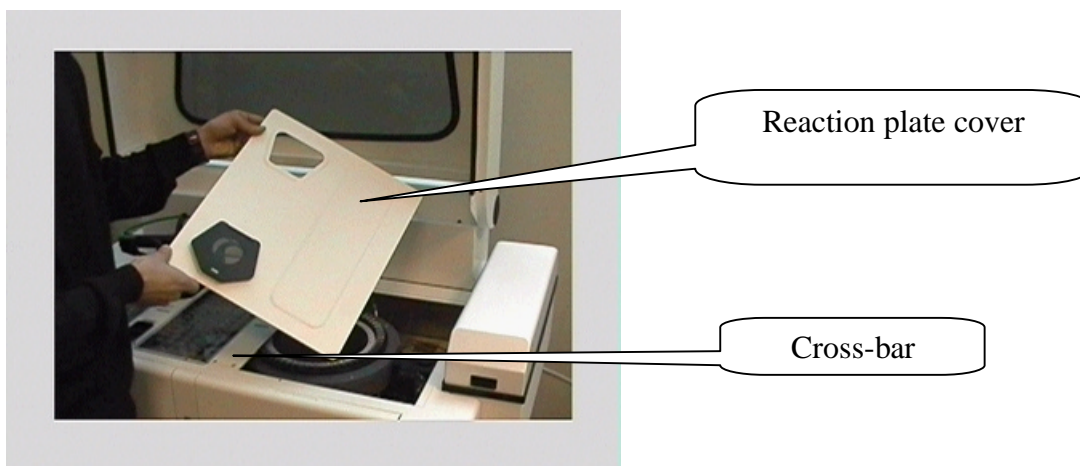


Fig. 1



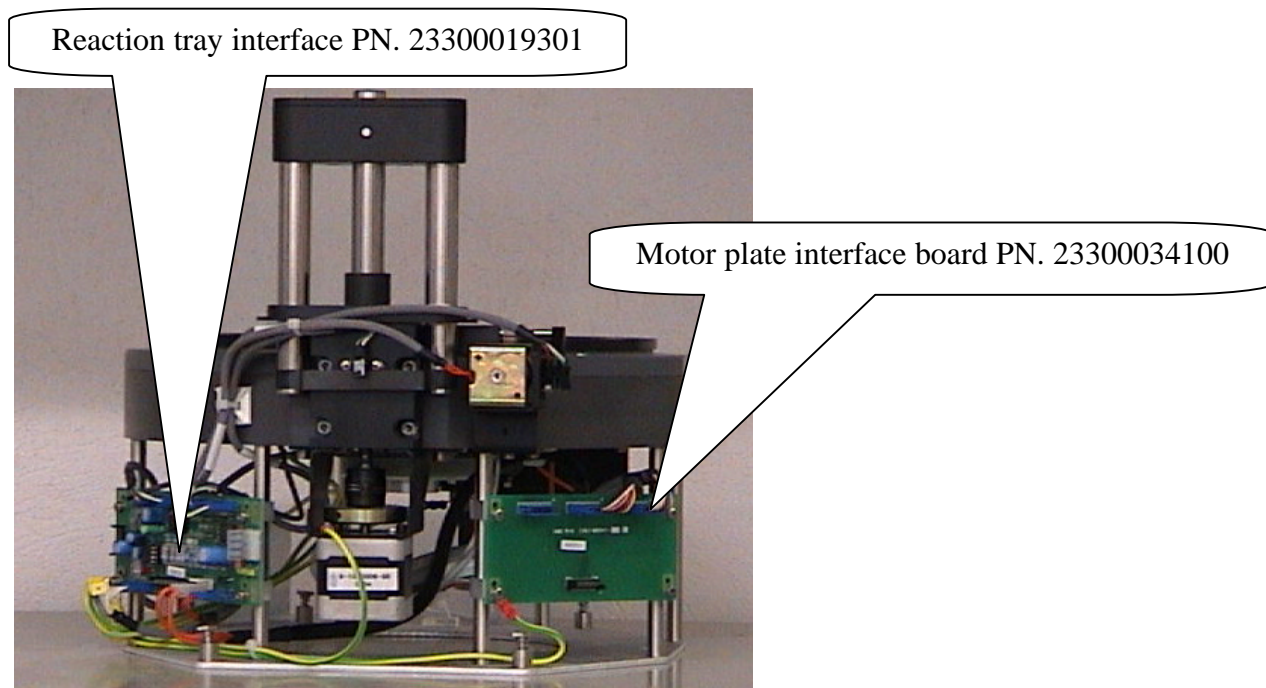


Fig. 2 Reaction Plate assembly PN. 23910006101 \*

\* Reaction Plate assembly is supplied without Washing Station and Photometer

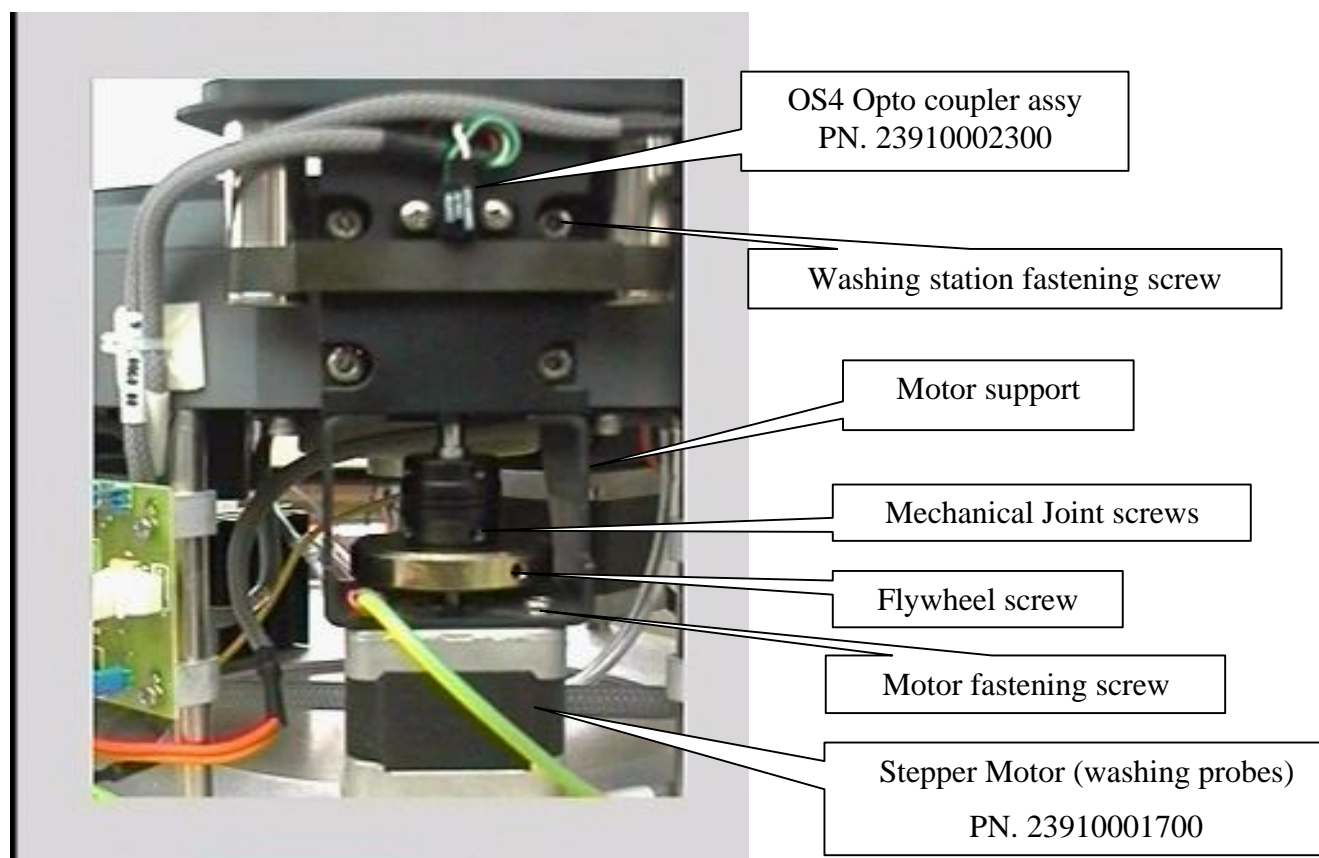


Fig. 3 Washing Station probes assembly PN. 23910000600

### 6.5.2 SUBSTITUTION OF THE WASHING STATION OPTIC SENSOR

- **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**

1. Remove the reaction plate cover (Fig. 1).
2. Unscrew the four cross-bas fastening screws and remove the cross-bar (Fig. 1).
3. Unscrew the four anchored screws of the reaction plate module (Fig. 2).
4. Place the reaction plate vertically inside the instrument.
5. Remove the J3 connector from the Reaction Plate Interface (Fig.2).
6. Unscrew the two fastening screws on the optic sensor and substitute it with the new one (Fig. 3).
7. Remount, repeating the above steps, 6 through 1, in reverse order.
8. After substituting and remounting, carry out the mechanical adjustment and functional check procedures as described in section 6.5.4

### 6.5.3 SUBSTITUTION OF THE WASHING STATION MOTOR

- **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure..**

1. Remove the reaction plate cover (Fig. 1).
2. Unscrew the four cross-bas fastening screws and remove the cross-bar (Fig. 1).
3. Unscrew the four anchored screws of the reaction plate module (Fig. 2).
4. Place the reaction plate vertically inside the instrument.
5. Remove the J3 connector from the Reaction Plate Interface (Fig.2).
6. Loosen the two flywheel and motor joint screws (Fig. 3).
7. Unscrew the four fastening screws on the motor support and remove it (Fig. 3).
8. Unscrew the four fastening screws on the motor, substitute it with the new one (Fig. 3).
9. Remount, repeating the above steps, 8 through 1, in reverse order.
10. After substituting and remounting, carry out the mechanical adjustment and functional check procedures as described in section 6.5.4


## 6.5.4 WASHING STATION - MECHANICAL ADJUSTMENT AND FUNCTIONAL CHECK

### 6.5.4.1 Washing Station - Horizontal Alignment

1. Turn on the "ILab 300 Plus" system and launch the "Analyzer" program
2. Launch the diagnostic program and select the "Plate" function.  
Reset by clicking on the appropriate key. Make sure the "Home Sensor" Reaction Plate window lights up green.
3. Select the "Washing" function. Reset by clicking on the appropriate key. Make sure the "Home Sensor" Washing Station window shows a green light.
4. Select several times "Home" and "Down" by clicking on the appropriate keys. Make sure the washing station probes enter the reaction cuvettes freely and totally unhindered. Pay particular attention to the drying pad placed on probe number five. If the drying pad does experience any friction whatsoever, loosen the two fastening screws on the bottom cover panel of the washing station and mechanically align the pad better (Fig. 5).
5. Upon completion of this procedure, tighten the fastening screws on the bottom cover of the washing station and exit the diagnostic program by clicking on "Diagnostic".
6. Carry out the functional check of the washing station as described in section 6.5.4.3

Fig. 5

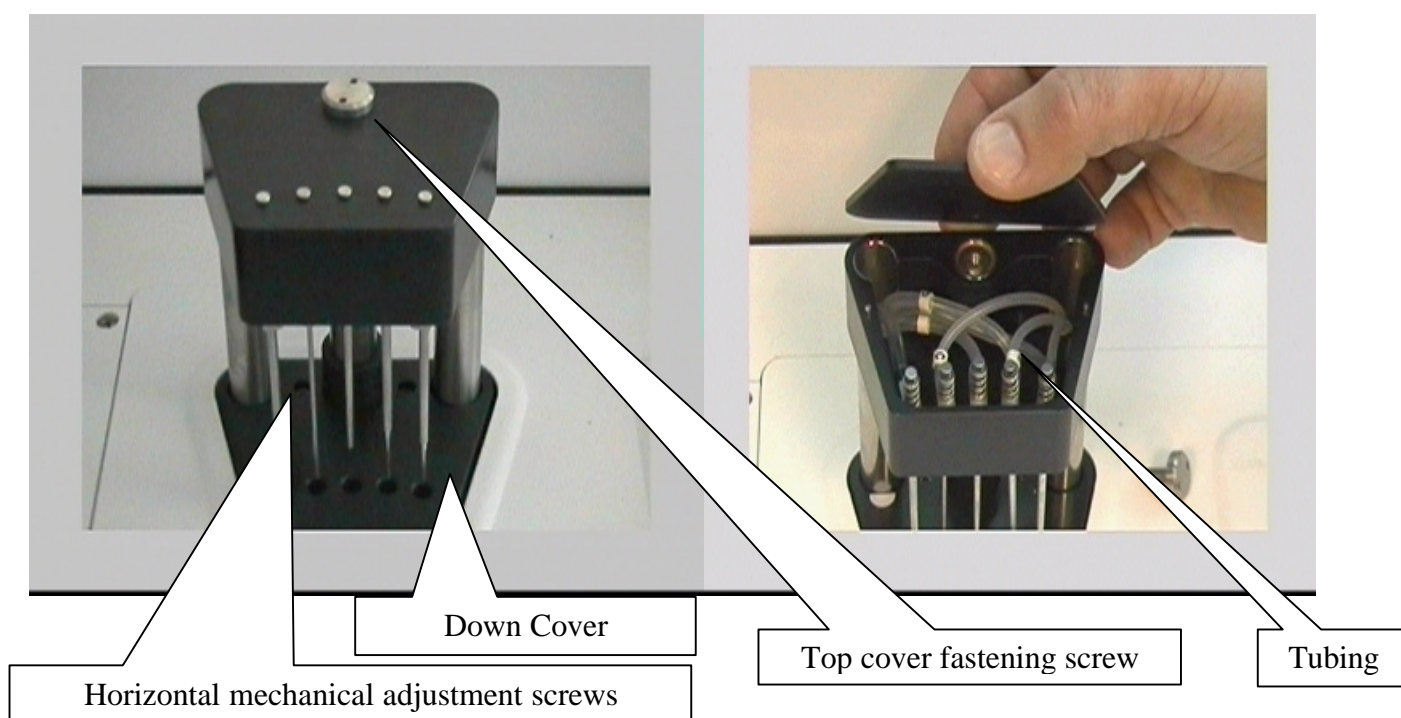Fig. 6





#### **6.5.4.2 WASHING STATION - VERTICAL ALIGNMENT**

1. Turn on the "ILab 300 Plus" system (first the instrument and then the computer).
2. Launch the diagnostic program and select the "Plate" function.
3. Reset by clicking on the appropriate key. Make sure the "Home Sensor" Reaction Plate window lights up green.
4. Select the "Wash" function. Reset by clicking on the appropriate key.
5. Make sure the "Home Sensor" Washing Station window shows a green light.
6. Select several times "Home" and "Down" by clicking on the appropriate keys. Make sure the washing station probes enter the reaction cuvettes freely and totally unhindered. Pay particular attention to the drying pad placed on probe number five.
7. The upper limit must be reached before the mechanical end of run. The lower limit must allow the drying pad to touch the bottom of the cuvette.
8. The washing station does not normally need vertical mechanical adjustment. If, however, it does, loosen the OS4 optic sensor fastening screws (Fig.3) and move it either up or down in order to correctly position the washing station. Upon completion of this procedure, tighten the optic sensor fastening screws and exit the diagnostic program by clicking on "Diagnostic".
9. Carry out the functional check of the washing station as described in section 6.5.4.3.

#### **6.5.4.3 WASHING STATION - FUNCTIONAL CHECK**

1. Turn on the "ILab 300 Plus" system (first the instrument and then the computer).
2. Launch the "ILab 300 Plus" program and, under System Monitor, carry out a wash cycle and reading cycle of the cuvettes using water (WBL).
3. Check system efficiency by comparing the obtained WBL values against the range of acceptance of WBL. Make sure there is no leakage of liquid inside the analyzer.

## 6.5.5 WASHING STATION - CUVETTE HOLDER ASSEMBLY

### 6.5.5.1 SUBSTITUTION OF THE CUVETTE HOLDER ASSEMBLY MAGNET

➤ **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure..**

1. Remove the reaction plate cover (Fig. 1).
2. Unscrew the four cross-bas fastening screws and remove the cross-bar (Fig. 1).
3. Unscrew the four anchored screws of the reaction plate module (Fig. 2).
4. Place the reaction plate vertically inside the instrument.
5. Remove the magnet's two fastening screws (Fig. 7).
6. Remove connector J6 from the Reaction Plate Interface board (Fig. 2).
7. Remove the magnet and substitute it with the new one.
8. Remount, repeating the above steps, 6 through 1, in reverse order.

### 6.5.5.2 SUBSTITUTION OF THE CUVETTE HOLDER ASSEMBLY OPTO SENSOR

➤ **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**

1. Remove the reaction plate cover (Fig. 1).
2. Remove connector J5 from the Reaction Plate Interface board (Fig. 2).
3. Remove the optic sensor's two fastening screws and substitute it with the new one (Fig. 7).
4. Remount, repeating the above steps, 3 through 1, in reverse order.
5. After substituting and remounting, carry out the mechanical adjustment procedures as described in section 6.5.5.3

### 6.5.5.3 ADJUSTMENT AND FUNCTIONAL CHECK OF THE CUVETTE HOLDER MAGNETIC TRIGGER

➤ **N.B.: Make sure the instrument (ILab 300 Plus) is turned off before performing this substitution procedure.**

1. Loosen the screw on the cuvette holder trigger (Fig. 7).
2. Rotate the reaction plate in such a manner as to position the trigger between two cuvettes (Fig. 8).
3. Manually push the pin toward the center of the solenoid - simulating magnetic attraction between the two.
4. Position the trigger in such a manner that the top edge is at the same height as the top edge of the cuvette. Tighten the cuvette holder trigger fastening screw (Fig. 8).
5. Turn on the "ILab 300 Plus" system and launch the "Analyzer" program.
6. Launch the diagnostic program and select the "Plate" function.
7. Reset by clicking on the appropriate key. Make sure the "Home Sensor" Reaction Plate window lights up green.
8. Select the "Wash" function and reset by clicking on the appropriate key.
9. Make sure the "Home Sensor" Washing Station window lights up green.
10. Select "Home" and "Down" several times by clicking on the appropriate keys. Make sure the "Home Sensor" Cuv. Holder window lights up green when the washing station is in the "Home" position and lights up red when the washing station is in the "Down" position. If not, loosen the fastening screw on the opto sensor flag and adjust to obtain the correct signal (Fig. 7).
11. To exit the diagnostic program press the "Diagnostic" key.
12. To exit the analyzer program press "Shutdown" key.

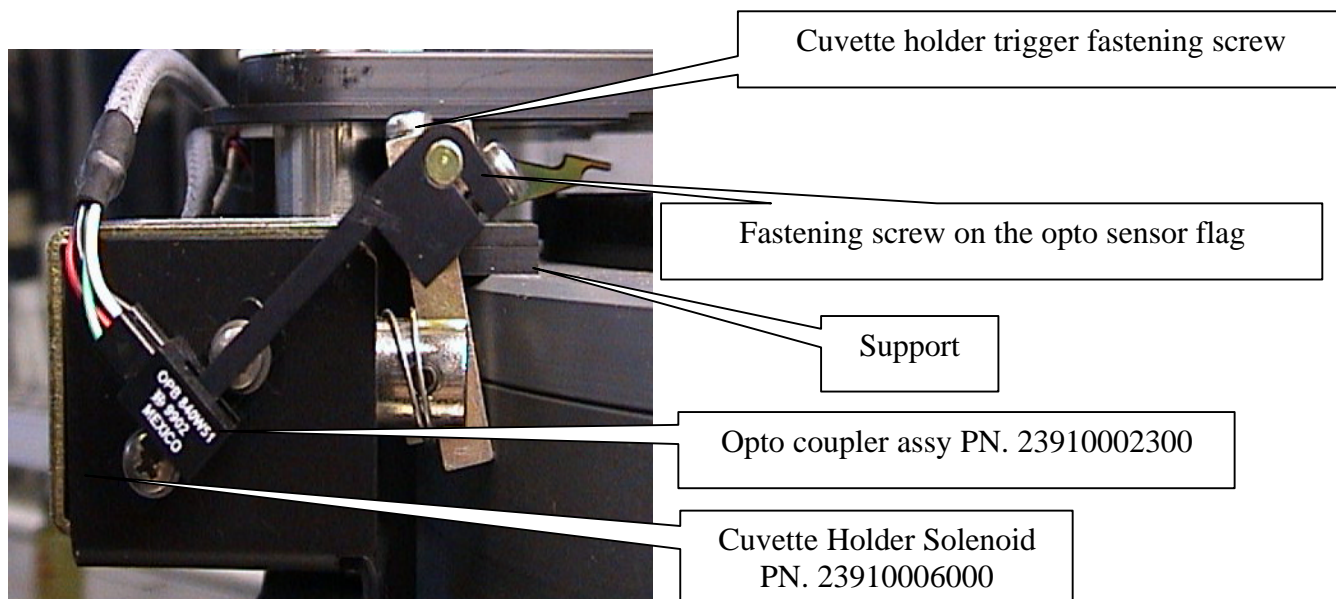


Fig. 7 Cuvette Holder

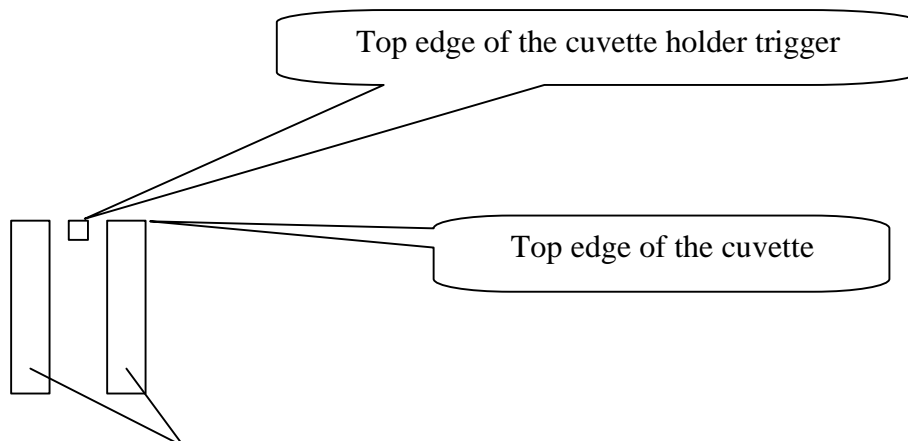


Fig. 8 Reaction cuvettes

## 6.6 TEMPERATURE CHECK AND ADJUSTMENT

### 6.6.1 REACTION CUVETTE TEMPERATURE ADJUSTMENT AND CHECK

1. Turn on the "ILab 300 Plus" system and launch the "Analyzer" program. Wait approximately one hour in order to allow the instrument to reach temperature balance (at room temperature 20°). The warm-up time may vary from 30 minutes (at about 21° C room temperature) up to 100 minutes (at about 18° C room temperature).
2. Launch the diagnostic program and select the "Plate" function.
3. Carry out the cuvette # 16 in dispensation position.
4. Open the door and fill the cuvette # 1 with 50 microliters of distilled water.
5. Select "Configuration" function.
6. Digit the password (1234) in the appropriate window and press "OK".
7. Insert the temperature probe (Model type Chemitec S 501 with probe Ø max 0.9 mm) inside cuvette #1, from the door, until it touches the bottom of the cuvette.
8. Make sure the temperature read, after 5 – 10 minutes, is  $37.0^{\circ}\text{C} - 0.2 + 0.3$ . If not, move the cursor within the "Plate" window to select the desired temperature and then click on "Save" to memorize said temperature.

**N.B.:** counts increase correspond to temperature decrease.  
 counts decrease correspond to temperature increase  
seven counts correspond to about one degree centigrade.

9. Remove the temperature probe from the cuvette and empty it.
10. Repeat the procedure from point 3 to 9 using in sequence cuvette # 16, 31, and 46.
11. To exit the diagnostic program press the "Diagnostic" key.
12. Carry out the "Glucose" test on seven samples using distilled water for sample and reagent.
13. Press, "Start" and once the seventh sample has been dispensed, press "Stop". Wait the sound signal from the analyzer, that means the actual stop.
14. Insert the temperature probe all the way to the bottom of the cuvette in which the sampling probe is dispensing. Make sure the temperature registered is  $37.0^{\circ}\text{C} - 0.2 + 0.3$ .
15. If not, launch the "Diagnostic" program and select the "Configuration" function.
16. Digit the password (1234) in the appropriate window and press "OK".

17. Move the cursor within the "Heater" window to select the desired temperature and then click on "Save" to memorize said temperature.
18. Remove the temperature probe from the cuvette and exit from the "Diagnostic" program.
19. Carry out a washing cycle from "Start Work" mask in order to clean the cuvette.
20. **To exit the "Analyzer" program press "Shutdown" key.**

**N.B.:** counts increase correspond to temperature decrease.  
 counts decrease correspond to temperature increase  
seven counts correspond to about one degree centigrade.

## 6.6.2 TEMPERATURE CHECK AND ADJUSTMENT USING KIT P/N 23550120400

### **IMPORTANT NOTICE**

**BEFORE APPLYING THE FOLLOWING PROCEDURE (PO 55-01204-00).**

CHECK WHETHER THE INTEGRATED **CHIP "TL 071xxx"** IS INSTALLED ON THE **"CONTROL PANEL BD. – POS. U2** (see fig. 1)" AND REPLACE IT WITH THE INTEGRATED **CHIP "LMC 7101 A or B or TVL 2461 cp."**

### **Temperature check and adjustment**

**This document describes how to adjust and check the pre-heater and plate temperature of the Instrument System.**

### **2. Summary:**

**CHECK WHETHER THE INTEGRATED CHIP "TL 071xxx" IS INSTALLED ON THE "CONTROL PANEL BD. – POS. U2 (see fig. 1)" AND REPLACE IT WITH THE INTEGRATED CHIP "LMC 7101 A or B or TLV 2461 cp." BEFORE APPLYING THE FOLLOWING PROCEDURE.**

- Turn on the system and wait in order to allow the instrument working temperature settlement. The warm-up time may vary from 30 minutes (at about 21° C room temperature) up to 100 minutes (at about 18° C room temperature).



- Calibrate the temperature probe by using the A and B resistors



- Put 300 µL (approx.) of distilled water PRE-WARMED at  $37.0 \div 37.3^{\circ}\text{C}$  into the cuvettes # 45 and 46



- Measure the temperature into the cuvettes # 45 and 46 to get an **indication** on the plate temperature



- Launch the temperatures test from the “Temperature” folder inside the “Diagnostic” program in order to **measure** and/or **slightly adjust them**

### 3. Special tools:

**KIT FOR TEMPERATURE ADJUSTMENT P/N 23550120400** composed by the following items:

- Temperature probe
- Calibrator Resistor A (10 KO)
- Calibrator Resistor B (5,76 KO)

### 4. Procedure description:

With a thermostatic bath, heat **250 ÷ 500 mL** of distilled water up to **37-37.5°C** (see page #3 if a thermostatic should be not not available).

**CHECK WHETHER THE INTEGRATED CHIP “TL O71xxx” IS INSTALLED ON THE “CONTROL PANEL BD. – POS. U2 (see fig. 1)” AND REPLACE IT WITH THE INTEGRATED CHIP “LMC 7101 A or B or TLV 2461 cp.” BEFORE APPLYING THE FOLLOWING PROCEDURE.**

- Turn on the system and wait in order to allow the instrument working temperature settlement. The warm-up may vary from 30 minutes (at about 21° C room temperature) up to 100 minutes (at about 18° C room temperature).
- During warm-up time, remove the protection cover of the Control Panel Board (Front Side Light Indicators) (See fig. 1).
- During warm-up time, launch the “Diagnostic” program and then select the “Temperature” Folder in the “Diagnostic” Program
- During warm-up time, select Calibration Sensor and carry out the three steps by following the indications showed on the screen.
- Replace the calibration resistor with the temperature probe.

**Note 1:** Now the temperature probe and the service probe function can be used as a thermometer.

Then route the temperature probe tip through the small door for cuvettes replacement. (See fig. 2)

To avoid the temperature probe damage be sure that its wires do not interfere with the reaction plate rotation.

**By executing the reaction plate Reset, the cuvettes # 45 and 46 will be accessible through the small door so allowing the following operations:**

- a. Pre-warm the micropipette tip by repeatedly aspirating and dispensing the liquid from/into the 250-500 ml container. Put 300 µL (approx.) of distilled water (pre-warmed at 37.0÷37.3°C into the cuvettes # 45 and 46.
- b. Wait 2-3 minutes then enter in the “Temperature” Folder and select “Temp. Test Run”.
- c. Measure the temperature by dipping and gently shaking the temperature probe into the liquid in the cuvettes.
- d. Take note about the measured temperatures.

**Note 2:** The reaction plate temperature must be set at 37.0 ÷ 37.2 °C



- e. If necessary, modify the value in the PLATE parameter field, (Configuration Folder in the Diagnostic Program) by taking into account that:

### **7 COUNTS CORRESPOND TO ABOUT 1°C**

**If the difference should be higher than 2°C, wait a few more minutes and repeat steps a.-b.-c.-d.-e.**

- f. Close the main cover and after having checked the correct setting of the reaction plate temperature, select Temp. Test Run in the Temperature Folder inside the Diagnostic program.
- g. Verify that the temperature probe or others tools are correctly positioned so to avoid their interference with the moving parts of the instrument. During this procedure step in fact, the Sampling arm, the washing Station and the Reaction Plate will become active.
- h. By pressing START, the system will sample seven reaction cuvettes in sequence
- i. At the end of the sampling cycles, the main cover will be unlocked.
- j. Open the small door for cuvettes replacement. Measure and take note of the cuvettes temperature starting from the seventh up to the first cuvette by dipping and gently shaking the temperature probe into the liquid.
- k. At the end of measurements, press WASH. The system is ready for a new cycle.

**Note 3:** The measure obtained from the seventh cuvette is the best indication of the pre-heater temperature. The measure obtained on the first cuvette is the nearest to the reaction plate temperature (previously adjusted to  $37.0 \div 37.2$  °C). The objective is to obtain the same temperatures into all the seven cuvettes (37° - 37.2°C).

- l. If necessary, modify the value in the ARM parameter field, (Configuration Folder in the Diagnostic Program) by taking into account that:

### **7 COUNTS CORRESPOND TO ABOUT 1°C**

- m. After any modification in the Parameter values repeat the steps from “f” through “k”.

- n. At the end of the adjustment procedure, get out from the Diagnostic Program and from the Instrument software, execute a WASH cycle then followed by the WBL.

**Further information:**

The reaction plate has been studied in order **to maintain the dispensed liquids at the set temperature and not to heat them up.**

Consequently it is really important to adjust the pre-heater temperature as much as possible near to the reaction plate temperature value ( $37.0 \div 37.2$  °C).

If a thermostatic bath should be not available, then apply the following:

**CHECK WHETHER THE INTEGRATED CHIP “TL 071xxx” IS INSTALLED ON THE “CONTROL PANEL BD. – POS. U2 (see fig. 1)” AND REPLACE IT WITH THE INTEGRATED CHIP “LMC 7101 A or B or TLV 2461 cp.” BEFORE APPLYING THE FOLLOWING PROCEDURE.**

- Turn on the system and wait in order to allow the instrument working temperature settlement. The warm-up may vary from 30 minutes (at about 21° C room temperature) up to 100 minutes (at about 18° C room temperature).
- During warm-up time remove the protection cover of the Control Panel Board (Front Side Light Indicators). (See fig. 1)
- During warm-up time, launch the diagnostic program and select the Temperatures Folder in the Diagnostic Program
- During warm-up time, select Calibration Sensor and carry out the three steps by following the indications showed on the screen
- Replace the calibrator resistor with the temperature probe.

**Note 1:** Now the temperature probe and the service probe function can be used as a thermometer. Route the temperature probe tip through the small door for cuvettes replacement. (See fig 2)

To avoid the temperature probe damage, to be sure its wires do not interfere with the reaction plate rotation.

1. Verify that the temperature probe or others tools are correctly positioned so to avoid their interference with the moving parts of the instrument. During this procedure step in fact, the Sampling arm, the washing Station and the Reaction Plate will become active.
2. Select Temp. Test Run in the Temperature Folder inside the Diagnostic program.
3. By pressing START, the system will sample seven reaction cuvettes in sequence
4. At the end of the sampling cycles, the main cover will be unlocked.
5. **Wait 10 minutes**, open the main cover and the small door for cuvettes replacement. Measure and take note of the cuvettes temperature starting from the third up to the first cuvette by dipping and gently shaking the temperature probe into the liquid.
6. At the end of the measurements, press WASH. The system is ready for a new cycle.

**Note 2:** In this case, cuvettes 1 through 3 can give a good indication of the reaction plate temperature.

7. If necessary, modify the value in the PLATE parameter field, (Configuration Folder in the Diagnostic Program) by taking into account that:

## 7 COUNTS CORRESPOND TO ABOUT 1°C

8. After any modification in the Parameter values repeat the steps from step “1” through “7”.
9. Close the main cover and after having checked the correct setting of the reaction plate temperature, select Temp. Test Run in the Temperature Folder inside the Diagnostic program.
10. Verify that the temperature probe or others tools are correctly positioned so to avoid their interference with the moving parts of the instrument. During this procedure step in fact, the Sampling arm, the washing Station and the Reaction Plate will become active.
11. By pressing START, the system will sample seven reaction cuvettes in sequence
12. At the end of the sampling cycles, the main cover will be unlocked.
13. Open the small door for cuvettes replacement. Measure and take note of the cuvettes temperature starting from the seventh up to the first cuvette by dipping and gently shaking the temperature probe into the liquid.
14. At the end of measurements, press WASH. The system is ready for a new cycle.

**Note 3:** The measure obtained from the seventh cuvette is the best indication of the pre-heater temperature. The measure obtained on the first cuvette is the nearest to the reaction plate temperature (previously adjusted to  $37.0 \div 37.2$  °C). The objective is to obtain the same temperatures into all the seven cuvettes ( $37^{\circ} - 37.2^{\circ}\text{C}$ ).

15. If necessary, modify the value in the ARM parameter field, (Configuration Folder in the Diagnostic Program) by taking into account that:

### **7 COUNTS CORRESPOND TO ABOUT 1°C**

16. After any modification in the Parameter values repeat the steps from “7” through “14”.

17. At the end of the adjustment procedure, get out from the Diagnostic Program and from the Instrument software, execute a WASH cycle then followed by the WBL.

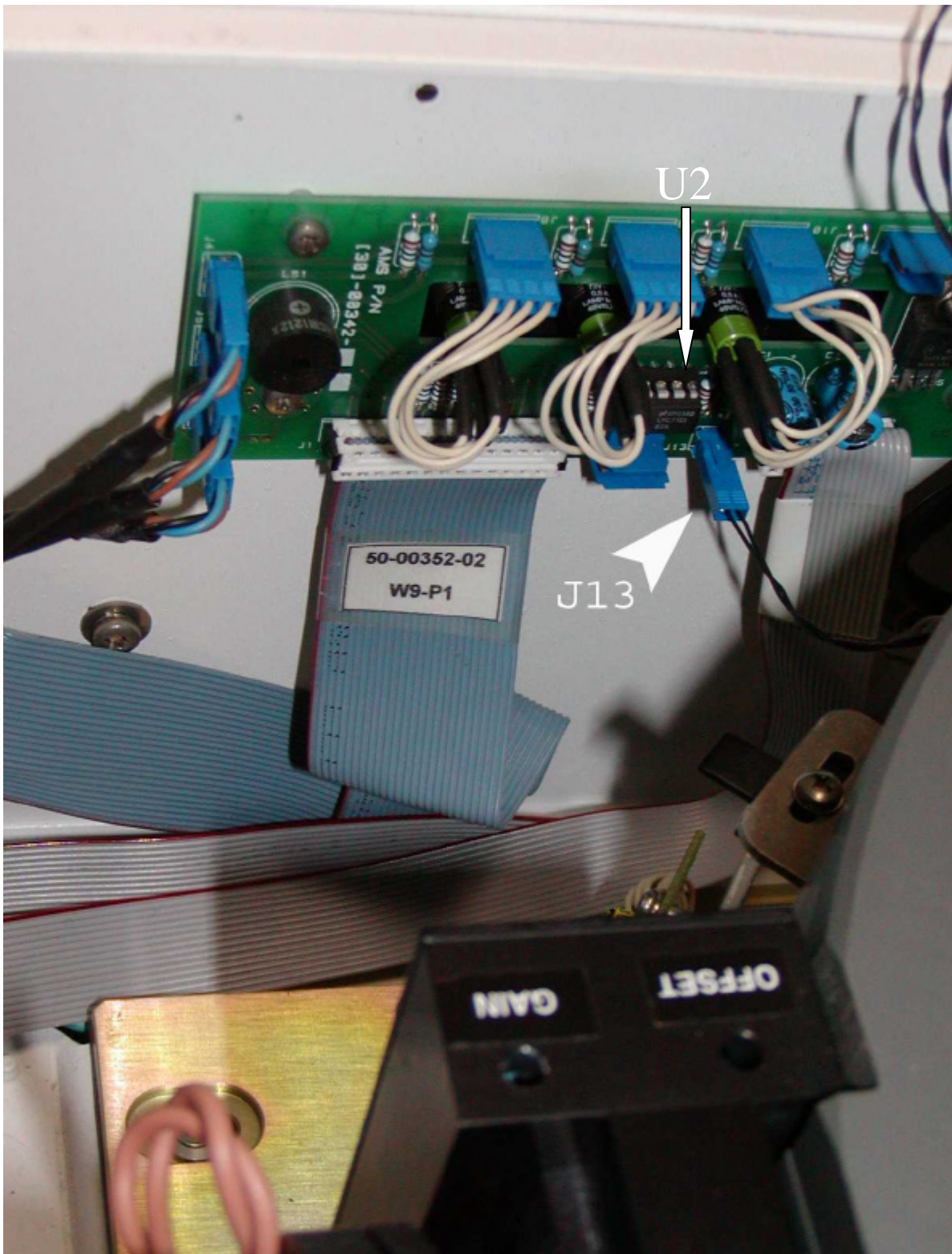


Fig.1

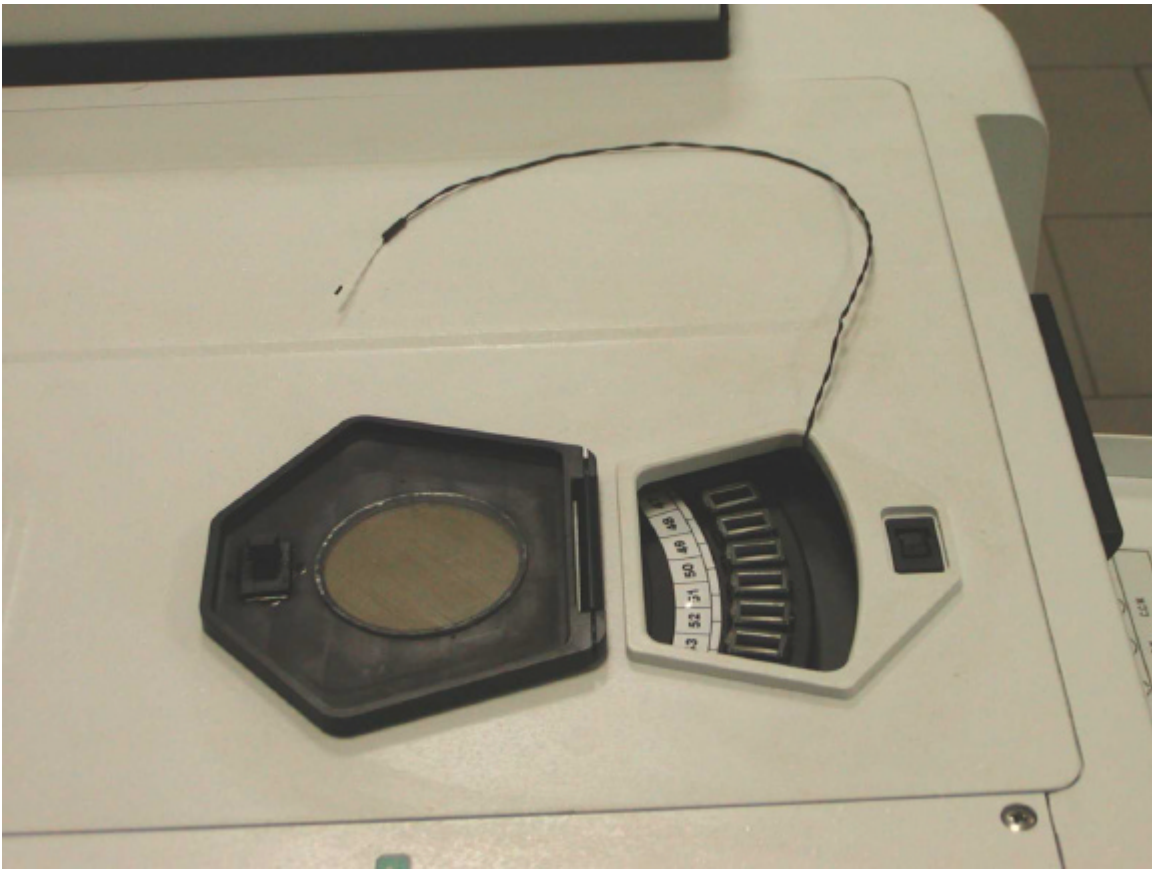



Fig. 2

### 6.6.3 REFRIGERATION ASSEMBLY TEMPERATURE CHECK

**Warning:** Keep the cooling fan protection grid, placed at the bottom of the refrigerated rack, unencumbered and free of dust.

1. Place a reagent container containing distilled water in a position of the refrigerated reagents rack.
2. Turn on the "ILab 300 Plus" system (instrument and computer).  
Wait approximately one hour in order to allow the instrument to reach temperature balance (at room temperature 20° C).
3. Place the temperature probe in the reagent container and check that the temperature is approximately 10°C lower than room temperature. If not, substitute the reagents racks holder (Fig. 1).
4. Remove the reagent container containing distilled water and exit from the analyzer program press "Shutdown" key.

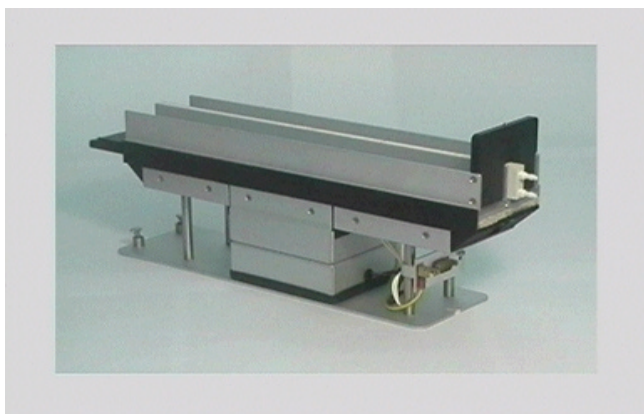Fig. 1



Refrigerated Reagent Rack assembly



## 6.7 RACKS HOLDER ALIGNMENT

### 6.7.1 MECHANICAL ALIGNMENT OF THE SAMPLES RACKS HOLDER

**Warning:** A check of the mechanical alignment of the samples racks must be performed every time a rack is removed (either for servicing or substitution).

**N.B.:** Make sure the instrument (ILab 300 Plus) is turned off before performing this procedure.

1. Open the two front panels and remove the samples and STAT racks.
2. Remove the top cover of the samples (Fig.1).
3. Loosen the four fastening screws on the samples racks holder (Fig. 2).
4. Replace the top cover of the samples. Place empty samples cups in the first and last positions of the racks and reinsert the racks.
5. Align the holes of the caps with those of the top cover of the samples.
6. Once again, remove the top cover of the samples, remove the racks and tighten the fastening screws of the samples racks holder.
7. Reposition the racks and replace the top cover of the samples.
8. Close the two front panels.



Fig. 1

top cover of the samples



Fig. 2

Samples racks holder



### 6.7.2 MECHANICAL ALIGNMENT OF THE REAGENTS RACKS HOLDER

**Warning:** A check of the mechanical alignment of the reagents racks must be performed every time a rack is removed (either for servicing or substitution).

**N.B.:** Make sure the instrument (ILab 300 Plus) is turned off before performing this procedure.

1. Open the front panel and remove the reagents racks.
2. Remove the top cover of the reagents (Fig. 1).
3. Loosen the four fastening screws on the reagents racks holder (Fig. 2).
4. Replace the top cover of the reagents. Place empty reagents cups in the first and last position of the racks and reinsert the racks.
5. Align the holes of the reagents containers with those of the top cover of the reagents.
6. Once again, remove the top cover of the reagents, remove the racks and tighten the fastening screws of the reagents racks holder.
7. Reposition the racks and replace the top cover of the reagents .
8. Close the front panel.

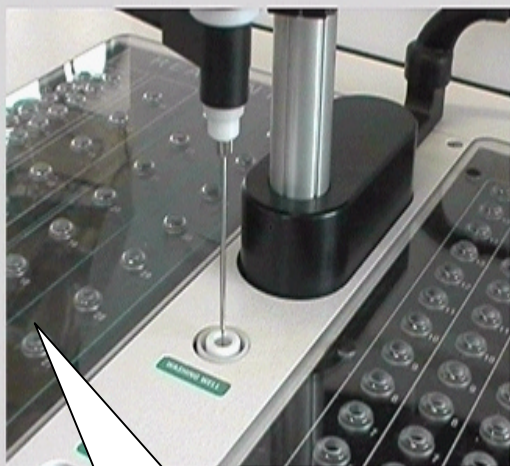


Fig. 1

top cover of the reagents



Fig. 2

Reagents racks holder

## CHAPTER 07

### - MAINTENANCE -

#### INDEX

7.	MAINTENANCE.....	2
7.1	PREVENTIVE MAINTENANCE.....	2
7.2	LIST OF PARTS SUBJECT TO WEAR AND USAGE.....	4
7.3	SAMPLING PROBE - CLEANING PROCEDURE.....	5
7.4	WASH STATION CANNULAS - CLEANING PROCEDURE.....	5
7.5	WASH SOLUTION BOTTLES - CLEANING.....	6
7.5.1	WASH SOLUTION BOTTLES - CLEANING PROCEDURE.....	6
7.6	POROUS PAD WASHING .....	7
7.6.1	POROUS PAD - CLEANING PROCEDURE.....	7
7.6.2	PROCEDURE FOR REPLACING THE POROUS PAD.....	11
7.7	HYDRAULIC CIRCUIT WASHING .....	12
7.7.1	HYDRAULIC CIRCUIT - CLEANING PROCEDURE.....	12
7.8	CHANGING THE PERISTALTIC PUMP TUBES.....	13
7.8.1	PROCEDURE FOR REPLACING THE PERISTALTIC PUMP TUBES.....	13
7.9	CHANGING THE PHOTOMETER LAMP.....	15
7.9.1	PROCEDURE FOR REPLACING THE PHOTOMETER LAMP.....	15
7.10	CHANGING THE REACTION CUVETTES.....	16
7.10.1	PROCEDURE FOR REPLACING THE REACTION CUVETTES.....	16
7.11	CHANGING THE TUBES KIT.....	17
7.11.1	PROCEDURE FOR REPLACING THE TUBES KIT.....	17
7.12	TABLE C – LIST OF THOSE MAINTENANCE PROCEDURES THAT CAN BE PERFORMED BY THE USER AND/OR BY THE MAINTENANCE TECHNICIAN .....	18
7.13	DECONTAMINATION PROCEDURE.....	20
7.14	ANALYTICAL VERIFICATION.....	20
7.15	SWITCH OFF THE INSTRUMENT.....	21

## 7 MAINTENANCE

This chapter contains all those routine operations, which concern instrument maintenance. Said procedures, listed and described below, should be carefully and scrupulously followed in order to guarantee the manufacturer's specifications and the perfect working order of the instrument over time.

### 7.1 PREVENTIVE MAINTENANCE

#### MAINTENANCE SCHEDULE

Table A, illustrated below, lists all those procedures to be carried out by the user/operator and the relative frequency schedule. Strict adherence to said schedule will guarantee the optimal operative efficiency of the instrument.

**TABLE A – MAINTENANCE SCHEDULE**

<b>FREQUENCY</b>	<b>PROCEDURE</b>	<b>NOTES</b>
DAILY – Before launching “Start Work”	Check the levels of all the wash solutions (Rinse, Water and Cleaning)	
DAILY – Before launching “Start Work”	Check the levels of Reagents, Standards and Controls	
DAILY – Before launching “Start Work”	Check the levels of the Waste Bottles and, if necessary, empty them	
DAILY – Before launching “Start Work”	Carry out a WBL cycle	See the User's Manual, Chapter 03– Description of Instrument Software
DAILY – After Shutdown	Clean the Sampling Probe using either paper toweling or gauze	See Procedure 7.3

**TABLE A – MAINTENANCE SCHEDULE**

<b>FREQUENCY</b>	<b>PROCEDURE</b>	<b>NOTES</b>
EVERY TWO WEEKS	Clean the Wash Station's four cannulas	See Procedure 7.4
EVERY TWO WEEKS	Clean the wash solution bottles (Rinse, Water and Cleaning)	See Procedure 7.5
ONCE A MONTH	Clean the Tip	See Procedure 7.6
ONCE A MONTH	Clean the Hydraulic Circuit	See Procedure 7.7
EVERY TWO MONTHS	Change the Reaction Cuvettes	See Procedure 7.10
EVERY SIX MONTHS	Change the Peristaltic Pump Tubes	See Procedure 7.8
EVERY SIX MONTHS	Change the Porous Pad	See Procedure 7.6
ONCE A YEAR	Change all the tubes (Tube Kit)	See Procedure 7.11
YEARLY OR EVERY 2000 HOURS.	Change the Photometer lamp	See Procedure 7.9

- N. B.: the above-described maintenance schedule refers to that situation in which the workload of the *Analyzer* is approximately 500 tests per day. The interval frequency may vary according to the individual instrument's daily workload.

## 7.2 LIST OF PARTS SUBJECT TO WEAR AND USAGE

Description	Type	Quantity	Code
Kit tubing peristaltic pump		2 “	23965004000
Reaction cuvettes		60 “	23965003100
Drying Pad		1 “	23901003800
Halogen Lamp		1 “	23935001600
Tubes Kit – complete		1 “	23965002701
Tygon tubing	1 mt.	1 “	23900125300
Kit for E.V.Diluter connection		1 “	23990090600
Probe		1 “	23905006400
Probe new coating		1 “	23905006401
MicroPump Assembly		1 “	23910002801
Complete Probe Assembly		1 “	23910006200
Comp Probe Assembly New C		1 “	23910006201
Diluter head Teflon Fitting		10 “	23C101004201
Diluter fitting		10 “	23C1010122200
E.V. Rinse fitting		10 “	23C1900125400
Solenoid Valve		1 “	23F35001900
MicroPump (up4) Assy		1 “	23050040300
MicroPump (up2) Assy		1 “	23050040400
MicroPump (up3) Assy		1 “	23050044700
Manifold Assembly		1 “	23050040500
Inlet/Outlet fitting for rinse...		1 “	23010122400
Lamps kit		5 “	23965003500
Na electrode		1 “	23350081400
K electrode		1 “	23350081500
Cl electrode		1 “	23350081600
Reference electrode		1 “	23350081700
Peristaltic Pump Tubing-head		1 “	23350083000
ISE inlet tubing connection	6 mt	1 “	23350083100

### 7.3 SAMPLING PROBE – CLEANING PROCEDURE

- 1) Turn off the *Analyzer*
- 2) Use only lint-free paper toweling or gauze
- 3) Dampen the gauze or paper toweling with distilled water and clean the outside of the sampling probe. Wipe the probe from the top downwards only! This is to avoid that any bits of cloth, paper or lint fibers accidentally enter the probe itself.
- 4) The manufacturer suggests that once weekly the above-described cleaning procedure be performed using, instead of only simple distilled water, a 5% sodium hypochlorite solution to dampen the gauze and then be repeated using distilled water.

### 7.4 WASH STATION CANNULAS – CLEANING PROCEDURE

- 1) Turn off the *Analyzer*.
- 2) Place a sheet of paper under the wash station cannulas in order to keep any extraneous material from accidentally falling into the cuvettes.
- 3) Use only lint-free paper toweling or gauze.
- 4) Dampen the gauze or paper toweling with distilled water and clean the outside of the cannulas. Wipe the cannulas from the top downwards only! This is to avoid that any bits of cloth, paper or lint fibers accidentally enter the cannulas themselves.
- 5) The manufacturer suggests that once weekly the above-described cleaning procedure be performed using, instead of only simple distilled water, a 5% sodium hypochlorite solution to dampen the gauze and then be repeated using distilled water.

## 7.5 WASHING SOLUTION BOTTLES – CLEANING

During normal use and over time, mold and dust can build up inside the wash solution bottles. For this reason, it is extremely important that they be periodically washed. Said cleaning must be thorough and meticulous in order to insure that every trace of mold or residue be removed.

How often the bottles must be cleaned depends on their use and on the quality of the distilled water used in that particular laboratory. However, the manufacturer recommends thorough washing at least once every two weeks.

It is extremely important that the user not underestimate the risks associated with mold and dust particles. They are to be regarded as a serious hazard as they can be the cause of instrument malfunction.

The wash solution bottles are located on the right side of the analyzer.

### 7.5.1 WASHING SOLUTION BOTTLES - CLEANING PROCEDURE

- 1) Turn off the *Analyzer*.

Pull the level sensor connectors out from the bottle caps.

- 2) Take the caps off the bottles and empty them.
- 3) Fill each bottle with a 5% sodium hypochlorite solution.
- 4) Clean the inside of each bottle using a bottlebrush in order to remove all traces of mold and/or residue.
- 5) Leave the sodium hypochlorite solution stand in the bottles for at least ten minutes.
- 6) Empty the bottles, rinse them repeatedly and well with tap water, and then twice more using distilled water.
- 7) Dry the bottles.
- 8) Fill the bottles with their proper solutions.
- 9) Replace the bottles in their respective housings.
- 10) Close the bottles and reconnect the level sensors.
- 11) Carry out two 'Wash Cuvettes' cycles and two 'WBL' cycles. Check system efficiency by comparing the WBL values and the test results obtained against the laboratory's quality control values.

## 7.6 POROUS PAD WASHING

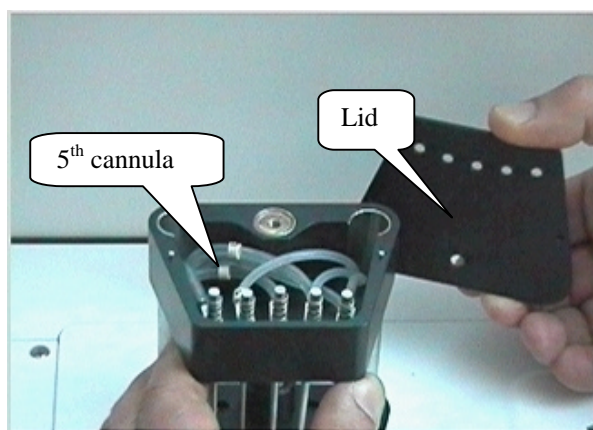
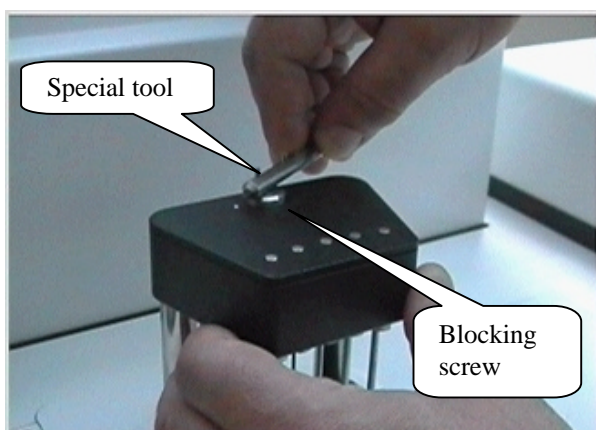
The Tip is used to dry the cuvettes after they have been washed. This drying process is carried out via aspiration and therefore, over time, the Tip will necessarily absorb various contaminating particles.

The manufacturer suggests that the Tip be replaced every six months. Said frequency may vary depending on the workload of the individual laboratory and the operating conditions/environment of the single instrument.



The Tip must be washed regularly in order to guarantee proper functioning and must be replaced with a new Tip as necessary.

### 7.6.1 POROUS PAD - CLEANING PROCEDURE

1. Turn off the *Analyzer*.
2. Remove the top cover (lid) of the Wash Station using the specific tool included among the instrument's accessories.






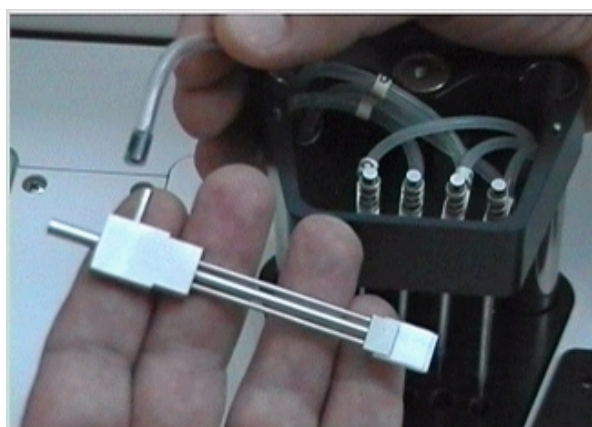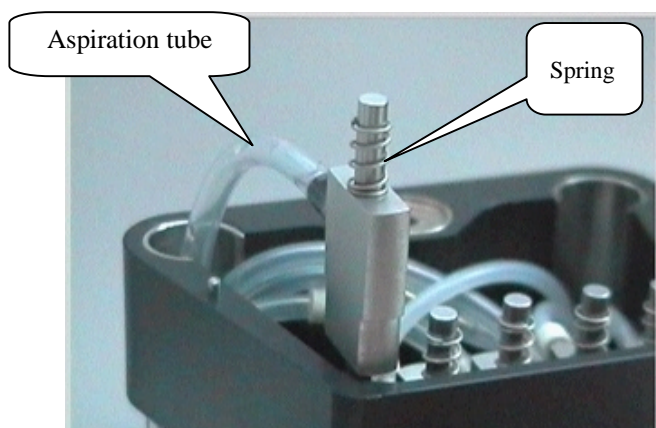3. Remove the cannula that contains the Tip (the 5<sup>th</sup> cannula) and disconnect the aspiration tube
4. Remove the spring






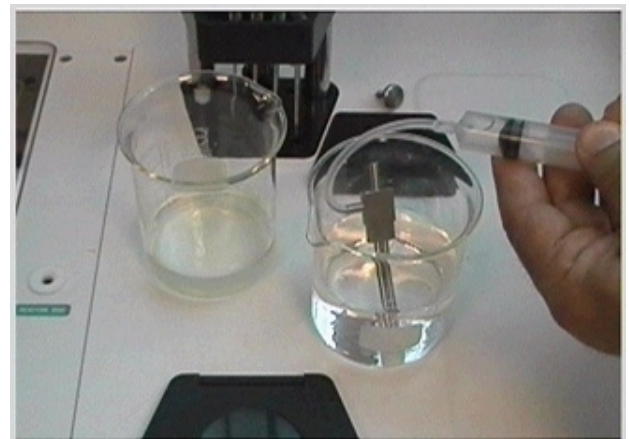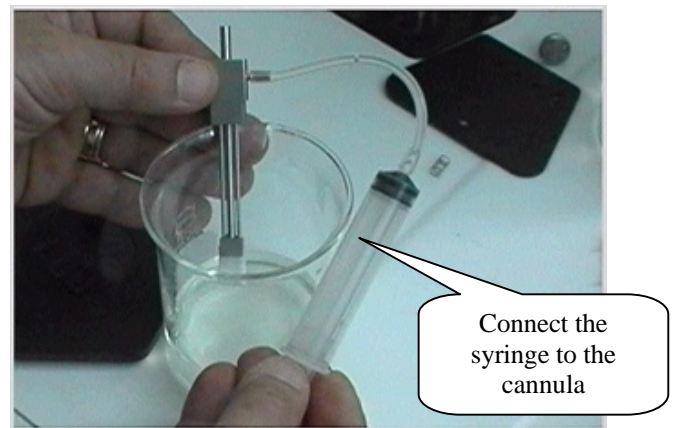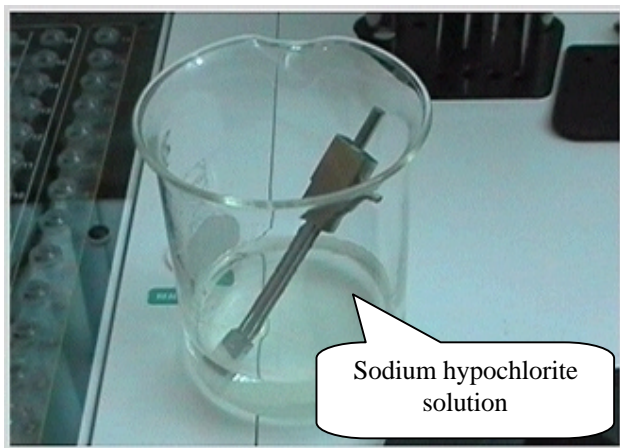5. Once the cannula has been removed, immerse it in a 5% solution of sodium hypochlorite for at least 15 minutes.
6. Attach a 10 ml syringe to the cannula as illustrated in Figure 6.
7. Aspirate and dispense the hypochlorite solution through the cannula (and the Tip) until the latter is completely clean. This aspirating and dispensing forces the liquid through the Tip fibers in both directions. Please see Figure 7.



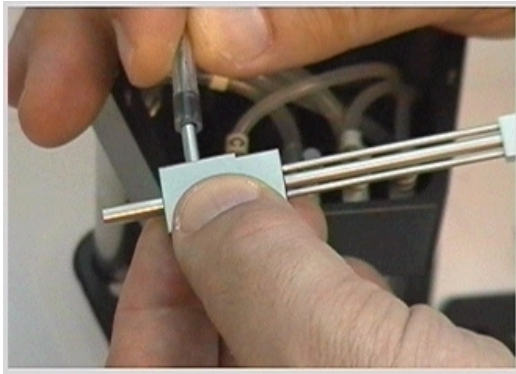8. Once the Tip is clean, repeat the aspirating and dispensing cycle 10 more times using distilled water, then disconnect the syringe.




9. Re-connect the aspiration tube to the cannula




10. Re-position the cannula in its housing being careful to not kink the aspiration tube



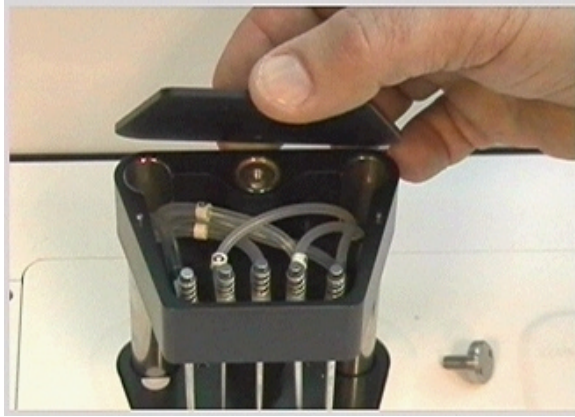11. Re-position the spring on the cannula's cylindrical axel



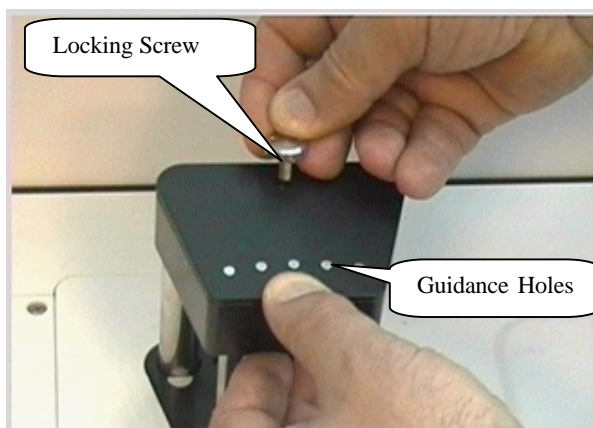12. Close the lid, making sure that the cylindrical tops of the cannulas fit into the corresponding holes on the lid






13. Insert the blocking screw into its setting and screw it down using the special tool included among the instrument's accessories

**IMPORTANT:** Before starting the instrument back up, manually move the wash station upwards as far as it will go.

## 7.6.2 PROCEDURE FOR REPLACING THE POROUS PAD

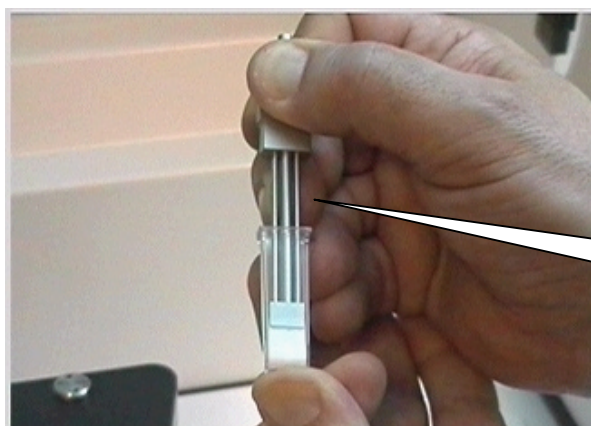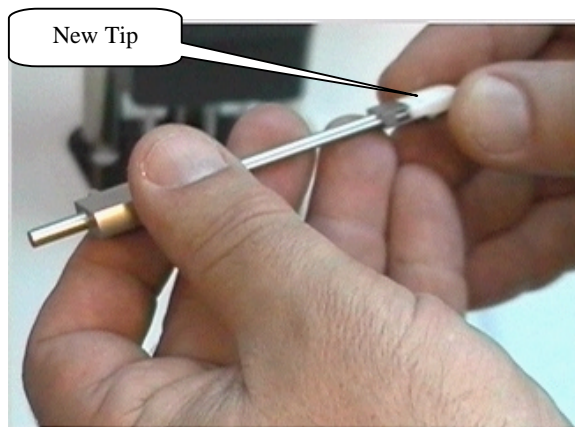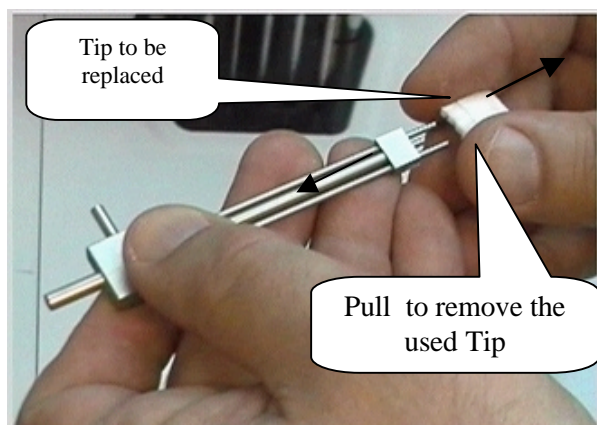1. Turn off the *Analyzer*.
2. Unscrew the screw that fastens the black top cover of the wash station, using the special tool included among the instrument's accessories.
3. Remove the wash station cover lid and take out the cannula containing the Tip. Be very careful to not lose the spring.
4. Disconnect the aspiration tube, remove the used Tip, and replace it with a new one. Press lightly to push the new Tip into place – be careful to not press too hard as this could deform the Tip.
5. Insert the cannula containing the new Tip into a cuvette, pushing it down inside until it takes on the shape of the inside of the cuvette.
6. Remove the cuvette, connect the aspiration tube to the cannula and reposition the cannula back into its proper housing in the wash station. Be careful to not kink the aspiration tube while doing so.
7. Reposition the spring in its proper housing. Remount the wash station coverlid making sure that the cylindrical tops of the cannulas fit into the corresponding holes on the wash station coverlid.

**IMPORTANT:** Before starting the instrument back up, manually move the wash station upwards as far as it will go.





## 7.7 HYDRAULIC CIRCUIT WASHING

During normal use and over time, mold and dust can build up inside the wash bottles and can have a negative effect on the hydraulic circuit, compromising the correct functioning of the micro-pump and valves. This, in turn, can lead to inefficiency in the sampling probe and cuvette washing system.

For this reason, it is extremely important that the hydraulic circuit be periodically washed. Said cleaning must be thorough and meticulous in order to assure that every trace of mold or residue be removed.

How often the hydraulic circuit must be washed depends on the operating conditions/environment of the single instrument and the quality of the distilled water used in that particular laboratory. The manufacturer recommends thorough washing at least once a month.

It is extremely important that the user not underestimate the risks associated with mold and dust particles. They are to be regarded as a serious hazard as they can be the cause of instrument malfunction.

The hydraulic circuit input cannulas are located inside the bottles on the right side of the *Analyzer*.

### 7.7.1 HYDRAULIC CIRCUIT - CLEANING PROCEDURE

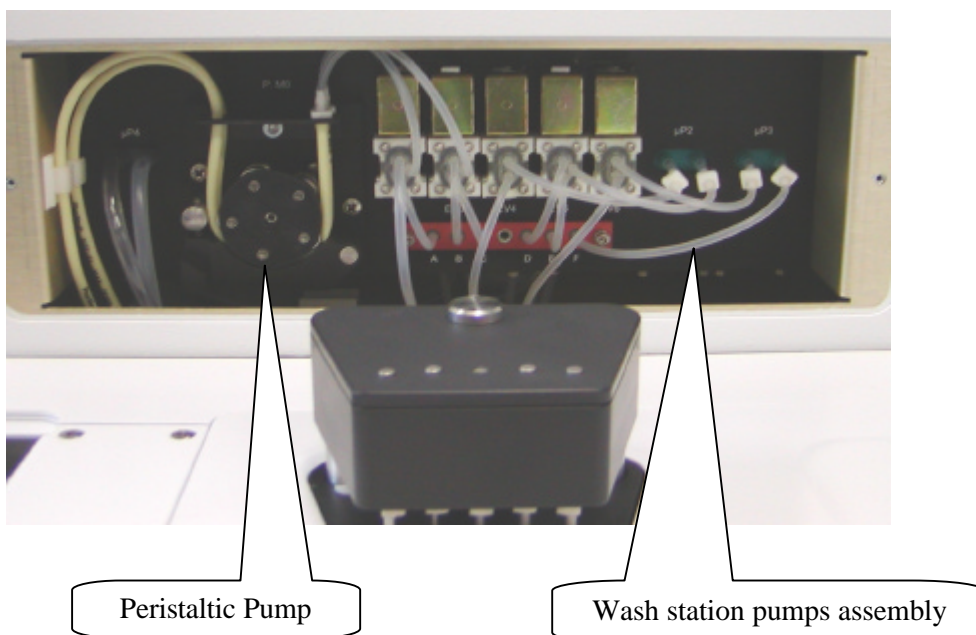
- 1) Turn on the *Analyzer*.
- 2) Prepare a bottle containing 500 ml of a 5% sodium hypochlorite solution
- 3) Insert the three aspiration cannulas, located inside the liquids bottles, into the bottle containing the sodium hypochlorite solution.
- 4) Carry out a 'Wash cuvettes' cycle and then a 'WBL' cycle cycle.
- 5) Wait for fifteen minutes. Clean and dry the three cannulas and then insert them into a bottle containing distilled water.
- 6) Carry out a 'Wash cuvettes' cycle and then a 'WBL' cycle. Repeat.
- 7) Insert the three aspiration cannulas back into their respective bottles. Said bottles should have, in the meantime, been cleaned and filled with a fresh supply of the required solution.
- 8) Carry out two 'Wash cuvettes' cycle and then a 'WBL' cycle. Check system efficiency by comparing the WBL values and the test results obtained against the laboratory's quality control values.

## 7.8 CHANGING THE PERISTALTIC PUMP TUBES

The manufacturer recommends that the Peristaltic Pump tubes be replaced every six months. Said frequency may vary depending on the workload of the individual laboratory. The quality and reliability of these tubes is fundamental to a correct emptying of the cuvettes.

### 7.8.1 PROCEDURE FOR REPLACING THE PERISTALTIC PUMP TUBES


- 1) Turn off the *Analyzer*.
- 2) Remove the panel located behind the wash station (Fig. 1).
- 3) Unhook the hinged guide by lowering and rotating the guide's blocking bracket to the left (Fig. 2).
- 4) Once the hinged guide has been removed, pull the tubes out of their relative nipples.



- 5) Insert the new tubes into their relative nipples
- 6) Position the tubes around the peristaltic pump rotor.
- 7) Close back the hinged guide.

- 8) Manually rotate the peristaltic pump rotor clockwise and check to make sure that the tubes are correctly positioned inside the guide.
- 9) Turn on the *Analyzer* and wait until the instrument has reached its proper operating temperature.
- 10) Have the instrument carry out a 'Wash cuvettes' cycle and then a 'WBL' cycle. Check system efficiency by comparing the WBL values and the test results obtained against the laboratory's quality control values.
- 11) Make sure that there is no leakage and then close back the panel.

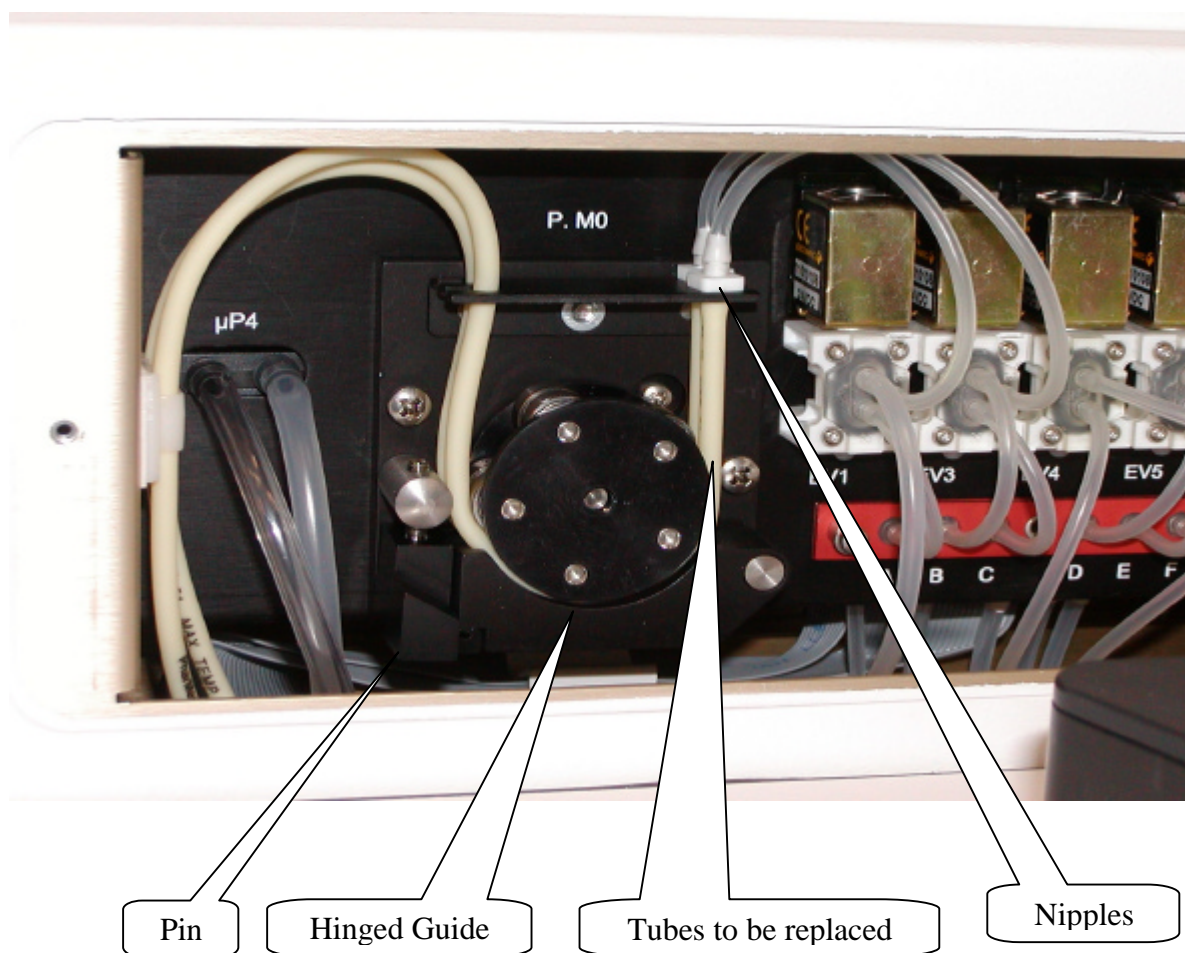Fig.2



## 7.9 CHANGING THE PHOTOMETER LAMP

The manufacturer suggests that the lamp be replaced after approximately 2000 hours of use or one year.

Figure 3 illustrates the photometer lamp, its support base and its power supply wires. There is a small hole on the lamp base useful for its mechanical alignment. The lamp is mounted on the photometer, which is located on the right-hand side of the reaction plate.

### 7.9.1 PROCEDURE FOR REPLACING THE PHOTOMETER LAMP

- 1) Turn off the *Analyzer*.
- 2) Remove the reaction plate cover.
- 1) Disconnect the power supply wires from the Lamp Regulator Board by loosening the clamp screws on the M1 connector (Fig. 4).
- 4) Unscrew the lamp's fastening screw and remove the lamp from its housing (Fig. 4)
- 5) Replace the old lamp with a new one making sure that the pin is correctly inserted in the alignment hole (Fig. 4). Remount the assembly by repeating the above steps 4 through 1 in inverse order.
- 6) Turn the *Analyzer* on and wait until the instrument has reached its proper operating temperature.
- 7) Carry out a 'WBL' cycle. Check system efficiency by comparing the WBL values and the test results obtained against the laboratory's quality control values.

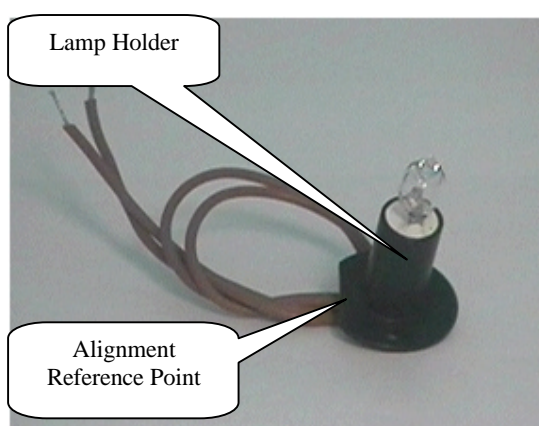


Fig.3

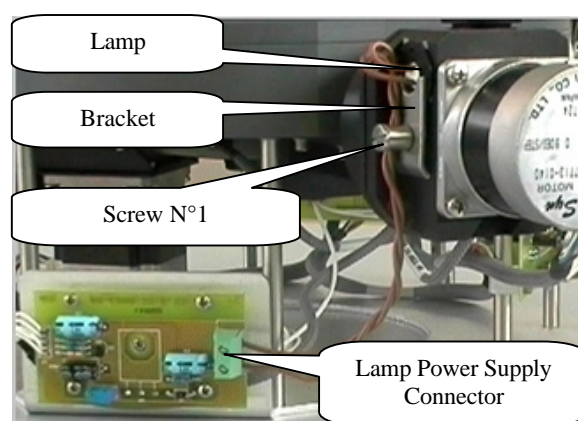


Fig.4

**WARNING!** DO NOT TOUCH THE GLASS PART OF THE LAMP WITH YOUR FINGERS!

IF NECESSARY, USE A CLEAN CLOTH TO REMOVE DUST, OR ALCOHOL TO REMOVE MORE STUBBORN DIRT.



## 7.10 CHANGING THE REACTION CUVETTES

Over time and through normal use the perfect transparency of the cuvettes diminishes. This less-than-perfect transparency has a negative impact on the quality of the optical readings. The manufacturer suggests that the cuvettes be replaced after two months of use. The cuvettes are located inside the reaction plate.

### 7.10.1 PROCEDURE FOR REPLACING THE REACTION CUVETTES

1. Turn on the *Analyzer* and wait until the instrument has reached its proper operating temperature.
2. Select Checks → WBL and click on the cell that must be replaced. Click on **“Change”**, located under the work program (Please see Chapter 03 Operator Manual– Description of Instrument Software).
3. A pull-down menu will appear. Click on **“Plate – First Half”** and then on **OK**.
4. A dialog box will appear asking: **“Do you want to change the cuvettes requested?”** Click on **OK** to confirm. Remove the reaction plate coverlid.
5. Remove the cuvettes numbered 1 through 30 by **VERTICALLY** lifting them out from their housing and then replace them with new cuvettes. Make sure to reinsert the new cuvettes **VERTICALLY**. Moreover, be especially careful to not touch the external surface of the cuvettes dedicated to photometric reading.
6. Select **“Plate – Second Half”** from the pull-down menu and then click on **OK**.
7. A dialog box will appear asking: **“Do you want to change the cuvettes requested?”** Click on **OK** to confirm.
8. Remove the cuvettes numbered 31 through 60 by **VERTICALLY** lifting them out from their housing and then replace them with new cuvettes. Make sure to reinsert the new cuvettes **VERTICALLY**. Moreover, be especially careful to not touch the external surface of the cuvettes dedicated to photometric reading.
9. Carry out a ‘WBL’ cycle. Check system efficiency by comparing the WBL values and the test results obtained against the laboratory’s quality control values.

## 7.11 CHANGING THE TUBES KIT

Over time and through normal use, the tubes become worn.

The manufacturer suggests that the tubes be replaced at least once a year. These tubes are located behind the panel situated behind the wash station.

### 7.11.1 PROCEDURE FOR REPLACING THE TUBES KIT

- 1) Turn off the *Analyzer*.
- 2) Remove the panel located behind the wash station.
- 3) Replace the tubes following the indications provided in the hydraulic diagram (SI 00457-00).
- 4) Turn on the *Analyzer* and wait until the instrument has reached its proper operating temperature.
- 5) Carry out a 'Wash cuvettes' cycle and then a 'WBL' cycle. Make sure that there is no leakage.
- 6) Check system efficiency by comparing the WBL values and the test results obtained against the laboratory's quality control values.

**7.12 TABLE C – LIST OF THOSE MAINTENANCE PROCEDURES THAT CAN BE PERFORMED BY THE USER AND/OR BY THE MAINTENANCE TECHNICIAN**

PROCEDURE	FREQUENCY (*)	EFFECTUATED BY
CLEANING THE SAMPLING PROBE	DAILY	OPERATOR
CLEANING THE FOUR WASH STATION CANNULAS	EVERY TWO WEEKS	OPERATOR
CLEANING THE WASH SOLUTION BOTTLES	EVERY TWO WEEKS	OPERATOR
REPLACING REACTION CUVETTES	EVERY TWO MONTHS	OPERATOR
CLEANING THE POROUS PAD	ONCE A MONTH	OPERATOR
REPLACING THE PROUS PAD	AS NEEDED	OPERATOR
ALIGNING AND ADJUSTMENT OF THE SAMPLING ARM (Std&Ctr/Reagents/Samples/Dispensing)	AS NEEDED (e.g.: after replacing any mechanical part)	SERVICE ENG.
REPLACING THE PERISTALTIC PUMP TUBES	EVERY SIX MONTHS	OPERATOR
REPLACING THE PHOTOMETER BULB	AFTER 2000 HOURS OF USE	SERVICE ENG.
REPLACING THE SAMPLING PROBE	AS NEEDED	OPERATOR
REPLACING THE TUBES (Tubes Kit)	ONCE A YEAR	SERVICE ENG.

<b>PROCEDURE</b>	<b>FREQUENCY (*)</b>	<b>EFFECTUATED BY</b>
REPLACING THE PRE-HEATER AND THE SENSOR LEVEL	AS NEEDED	SERVICE ENG.
MECHANICAL ALIGNMENT OF THE SAMPLING ARM	AS NEEDED	SERVICE ENG.
REPLACING OR ADJUSTING THE OPTIC SENSORS	AS NEEDED	SERVICE ENG.
REPLACING THE BELT	AS NEEDED	SERVICE ENG.
REPLACING A MOTOR	AS NEEDED	SERVICE ENG.
ALIGNING THE WASH STATION/REACTION PLATE	AS NEEDED	SERVICE ENG.
ALIGNING THE CUVETTE BLOCK	AS NEEDED	SERVICE ENG.
REPLACING THE PHOTOMETER	AS NEEDED	SERVICE ENG.
ADJUSTING THE PHOTOMETER	AS NEEDED	SERVICE ENG.
REPLACING THE ELECTRONIC BOARDS/CARDS AND THE MECHANICAL MODULES	AS NEEDED	SERVICE ENG.

N. B.: the above-indicated frequency intervals may vary according to the individual instrument's daily workload.

### 7.13 DECONTAMINATION PROCEDURE

Before replacing any instrument parts, repairing any defective items or performing any instrument maintenance procedure(s), the operator or maintenance technician must carry out the below-described decontamination procedure of the instrument part(s) involved in the operation(s).

This procedure can be performed on:

- the entire *Analyzer*
- those part(s) of the instrument subject to possible contamination

#### Material necessary

- an *ESOFENOL* solution diluted to 6% (60 cc in one liter of distilled water). *ESOFENOL* is an antibacterial and antiviral substance.
- Rubber gloves
- Mask
- Lab coat

#### External surfaces and individual parts

- ➔ Spray the solution all over the instrument, paying particular attention to wetting:
  - the sampling arm
  - the reaction plate (including the cuvettes)
  - the racks
  - the instrument chassis
- ➔ Allow the solution to stand for approximately 30 minutes
- ➔ Wipe the solution off the instrument and the various components using a sponge dampened with distilled water

Carry out a decontamination of the internal hydraulic circuit.

### 7.14 ANALYTICAL VERIFICATION.

**Verify the performance of the iLab 300 Plus every six months, either after the decontamination of the hydraulic circuit or after the replacement of the peristaltic pump tube.**

To verify the instrument performance use two tests. (ISE tests will be included if available). One enzymatic and the other one colorimetric. Use fresh reagents, calibrator and controls. Run the reagent blank and the calibration if required. Perform 20 replicates of two control levels.

Verify that average values are within the limits reported by the insert sheet of the control serum. Verify the C.V %. The coefficient of variation (C.V) should not be greater than what is reported into the "Precision" section of the relevant method insert sheet.

## 7.15 SWITCH OFF THE SYSTEM

### STANDBY MODE.

Select the Shutdown function from the Monitor and then confirm with “Yes”. If requested, the instrument will perform ISE cleaning, and will empty the fluidics. The computer will turn off. Do not power off the instrument. In this mode all temperature regulation (of reagents and reaction disk) are ON.

### SWITCHING OFF

From the standby mode turn the instrument off by pressing the main switch located on the left side of the instrument. Remove the reagents and place them in the refrigerator.

### SHUTDOWN PROCEDURE

This procedure must be carried out whenever the ILab 300 Plus is to be switched off for several days.

- Take out the tubing from the deionized water bottles, rinse and clean the tubes and insert them in a bottle filled with deionized water.
- Click the “Start Button”. Select “Wash all cuvettes”, “Water Blank Level” and then “OK”. Repeat this procedure three times.
- Take the tubing from the bottles and leave them in air.
- Click the “Start Button” and select “Wash all cuvettes”, “Water Blank Level” and “OK”.
- Turn the instrument off.
- Remove the sensor caps and empty the three bottles. Rinse them with deionized water.
- To restart the system place the three tubes in their respective bottles.
- Click the “Start Button”, Select “Wash all cuvettes”, “Water Blank Level” then “OK”. Repeat this procedure twice.

### For ILab 300 Plus with the ISE module:

- Move down the Ise door. Disengage the Ise module to have access at the upper inlet port of the assembly.
- Dispense 500 uL of Calibrator A into the upper port.
- Seal the inlet port with “Parafilm”
- Empty the Calibrator A bottle. Rinse the bottle with distilled water.

## CHAPTER 08

---

### – HOST COMMUNICATION –

#### INDEX

8.1 COMMUNICATION WITH THE HOST COMPUTER.....	2
8.1.1 COMMUNICATION PARAMETERS .....	2
8.2 PROTOCOL SPECIFICATIONS.....	3
8.3 HEADER RECORD (H) .....	4
8.4 PATIENT RECORD (P) .....	5
8.5 TEST ORDER RECORD (O).....	7
8.6 RESULTS RECORD (R) .....	10
8.7 MESSAGE TERMINATOR RECORD (L) .....	12
8.8 REQUEST INFORMATION RECORD (Q): .....	13

## 8.1 COMMUNICATION WITH THE HOST COMPUTER

The *ILab 300 Plus* can be connected to a host computer for the purpose of facilitating results print-out and patient management.

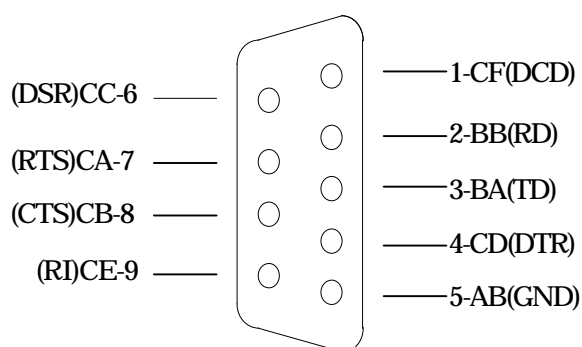
In order to enable communication between the *ILab 300 Plus* and the host computer, select the **Host Link** field under **Options** in the **Parameters** menu.

To activate communication between the *ILab 300 Plus* and the host computer, select **Host-Tx** (please see the software description in Chapter 03 of the User's Manual)

### 8.1.1 COMMUNICATION PARAMETERS

The *ILab 300 Plus* is linked to the managing computer using an RS-232C serial connector having the following specifications:

- Transmission method : Asynchronous, half duplex
- Baud Rate : 9600 Bit/sec.
- Data bits: : 8
- Parity : None
- Stop bit : 1
- Connector : 9 pin type D (male output from the ILab 300 Plus)



Serial connector



## 8.2 PROTOCOL SPECIFICATIONS

This part of Help (**Protocol Specifications**) contains information for the laboratory computer and analyzer. This exchange of data follows specific **ASTM** protocols:

**E 1381-95** Standard Specification for Low-Level Protocol to Transfer Messages between Clinical Laboratory Instruments and Computer Systems;

**E 1394-97** Standard Specification for Transferring Information between Clinical Instruments and Computer Systems.

**ASTM** uses a number of different terms to indicate the way it groups data.

- **Field:** an individual piece of data often referred to as a data field or a data element.
- **Record:** a number of logically related data fields grouped together to form one part of a complete message.
- **Repeat field:** a data field of the same type as the one immediately preceding it. A delimiter separates one instance of a repeat field from the next.
- **Component field:** part of data field that might contain more than one piece of data.

The default communication configuration for the *Analyzer* is the following: "**9600,N, 8,1**".

**ASTM** uses record types that are common and familiar to all laboratory personnel. It uses the following record types:

- **Header Record (H):** contains identifying information about the sending station, conventions that the device uses for field recognition, and the date and time of send station transmission.
- **Patient Record (P):** contains patient information and identification number.
- **Test Order Record (O):** contains information about the assay or requests themselves and includes other data.
- **Result Record (R):** contains information about the outcome of individual tests for an individual patient and follows a sample program record. The results contain the actual measurements derived from the test and a comparison of the individual result to certain ranges specified as norms for the laboratory.

- **Message Terminator Record (L)**: although the ASTM protocol supports three additional record types - a Request for Information Record, a Scientific Record and a Manufacturer's Information Record - the *Analyzer* is not implementing these in the first release and will ignore them.
- **Request Information Record (Q)**: is used by either clinical instruments or computer systems for a remote request for information from its reciprocal system.

The instrument does not send or accept comment records.

### 8.3 HEADER RECORD (H)

Field	Field Title	Down Load	Up Load	Max Len	Description and Valid Values								
1	Record Type ID	<b>R</b>	<b>A</b>	1	This is a required field that contains an “H” identifying it as a header record.								
2	Delimiters	<b>I</b>	<b>A</b>	4	<p>The <i>Analyzer</i> System uses only the four default values shown here. Delimiters may not be duplicated. The field delimiter follows the escape character to separate the delimiter specification from a subsequent field in the header record. Using default values, the first six characters of the header record will appear using the following characters:</p> <p>H I\^&amp;I</p> <table><tr><td>Field Delimiter</td><td>I</td></tr><tr><td>Repeat Delimiter</td><td>\</td></tr><tr><td>Component Delimiter</td><td>^</td></tr><tr><td>Escape Delimiter</td><td>&amp;</td></tr></table>	Field Delimiter	I	Repeat Delimiter	\	Component Delimiter	^	Escape Delimiter	&
Field Delimiter	I												
Repeat Delimiter	\												
Component Delimiter	^												
Escape Delimiter	&												
3	Message Control ID	<b>I</b>	<b>N</b>										
4	Access Password	<b>I</b>	<b>N</b>										
5	Sender Name or ID	<b>I</b>	<b>A</b>	<b>10</b>	‘SHAnalyzer’: This is the name of the device that is sending the data.								
6	Sender Street Address	<b>I</b>	<b>N</b>										
7	Reserved Field	<b>I</b>	<b>N</b>										
8	Sender Tel. Number	<b>I</b>	<b>N</b>										
9	Characteristics of Sender	<b>I</b>	<b>N</b>										

10	Receiver ID	<b>I</b>	<b>N</b>		
11	Comments or Special Instructions	<b>I</b>	<b>N</b>		
12	Processing	<b>I</b>	<b>N</b>		
13	ASTM Version No.	<b>I</b>	<b>N</b>		
14	Date and Time	<b>I</b>	<b>A</b>	14	Date and Time of transmission: formatted as YYYYMMDDHHMMSS. For example: 3:35 PM on March 1, 1995 would be represented using the following characters: 19950301153500.
<b>Legend:</b>		<b>R</b> Required	<b>D</b> Down Load		
		<b>O</b> Optional	<b>U</b> Up Load		
		<b>I</b> Ignored	<b>N</b> Never		
		<b>A</b> Always	<b>S</b> Sometimes		

**Example Header Record Layouts (H)**

Download	
<b>Host</b>	H I \ ^ I I I HOST I I I I I I I I 19950301153500<CR>
Upload	
<b>Analyzer System</b>	H I \ ^ & I I I SHAnalyzer I I I I I I I I 19950301154000<CR>

**8.4 PATIENT RECORD (P)**

Field	Field Title	Down Load	Up Load	Max Len	Description and Valid Values
1	Record Type ID	<b>R</b>	<b>A</b>	1	This is a required field that contains a “P” identifying it as a patient record.
2	Sequence Number	<b>R</b>	<b>A</b>	3	This field starts with a “1” for the patient and is incremented by 1 for each additional patient within the transmission.
3	Practice Assigned Patient ID	<b>R</b>	<b>A</b>	15	This field can be assigned by the instrument with no corresponding download.
4	Laboratory Assigned Patient ID	<b>I</b>	<b>N</b>		
5	Patient ID No. 3	<b>I</b>	<b>N</b>		

6	Patient Name	<b>O</b>	<b>S</b>	36	This field has two components: <ul style="list-style-type: none"> <li>• Last Name (up to 20 characters)</li> <li>• First Name (up to 15 characters).</li> </ul>
7	Mother's Maiden Name	<b>I</b>	<b>N</b>		
8	Birth Date	<b>O</b>	<b>S</b>	8	Formatted as YYYYMMDD: For example, a birth date of December 1, 1980 would be represented as: 19801201
9	Patient Sex	<b>R</b>	<b>A</b>	1	The valid values are: <ul style="list-style-type: none"> <li>• M for Male</li> <li>• F for Female</li> </ul>
10	Patient Race/Ethnic Origin	<b>I</b>	<b>N</b>		The <i>Analyzer</i> System will ignore this field at launch.
11	Patient Address	<b>O</b>	<b>S</b>	60	For <i>Analyzer</i> , this is a four- component field: <ul style="list-style-type: none"> <li>• Address (25 characters)</li> <li>• City (25 characters)</li> <li>• State (2 characters e.g.: NY, IT)</li> <li>• Zip (5 characters)</li> </ul>
12	Reserved Field	<b>I</b>	<b>N</b>		
13	Patient Tel Number	<b>I</b>	<b>N</b>		
14	Attending Physician ID	<b>I</b>	<b>N</b>		
15	Special Field 1	<b>I</b>	<b>N</b>		
16	Special Field 2	<b>I</b>	<b>N</b>		
17	Patient Height	<b>I</b>	<b>N</b>		
18	Patient Weight	<b>I</b>	<b>N</b>		
19	Patient Known or Suspected Diagnosis	<b>I</b>	<b>N</b>		
20	Patient Active Medications	<b>I</b>	<b>N</b>		
21	Patient's Diet	<b>I</b>	<b>N</b>		
22	Practice Field No. 1	<b>I</b>	<b>N</b>		
23	Practice Field No. 2	<b>I</b>	<b>N</b>		
24	Admission Date and Discharge Date (if desired)	<b>O</b>	<b>S</b>	8	Admission date only. Formatted as YYYYMMDD.
25	Admission Status	<b>I</b>	<b>N</b>		
26	Location	<b>O</b>	<b>S</b>	20	

27	Nature of Alternative Diagnostic Code and classifiers	<b>I</b>	<b>N</b>		
28	Alternative Diagnostic Code and classification	<b>I</b>	<b>N</b>		
29	Patient Religion	<b>I</b>	<b>N</b>		
30	Marital Status	<b>I</b>	<b>N</b>		
31	Isolation Status	<b>I</b>	<b>N</b>		
32	Language	<b>I</b>	<b>N</b>		
33	Hospital Service	<b>I</b>	<b>N</b>		
34	Hospital Institution	<b>I</b>	<b>N</b>		
35	Dosage Category	<b>I</b>	<b>N</b>		

**Legend:** **R** Required      **D** Down Load  
**O** Optional      **U** Up Load  
**I** Ignored      **N** Never  
**A** Always      **S** Sometimes

#### Example Patient Record (P)

	Download
<b>Host</b>	P   1   B108K     MW5910^Smith    19861002   M    Park Avenue^New York^NY^10002       20020923    Hematology
<b>Analyzer System</b>	P   1   B108K     MW5910^Smith    19861002   M    Park Avenue^New York^NY^10002       20020923    Hematology

### 8.5 TEST ORDER RECORD (O)

Field	Field Title	Down Load	Up Load	Max Length	Description and Valid Values
1	Record Type ID	<b>R</b>	<b>A</b>	1	This is required field that contains an “O” identifying it as an order
2	Sequence Number	<b>R</b>	<b>A</b>	3	This field starts with “1” for the first Test Order Record and is incremented by 1 for each additional Test Order Record within the record.  This will be reset to “1” whenever another patient record

					is transmitted.
3	Specimen ID	<b>R</b>	<b>A</b>	15	Although the operator can manually edit this field at any time, the value of this field is usually assigned by the laboratory computer before down loading. The <i>Analyzer</i> uses and reports its results based on the assigned specimen ID.
4	Instrument Specimen ID	<b>I</b>	<b>N</b>		
5	Universal Test ID	<b>I</b> <b>I</b> <b>I</b> <b>R</b>	<b>N</b> <b>N</b> <b>N</b> <b>A</b>	9	This is a four-component field: <ul style="list-style-type: none"> <li>• Universal Test ID Code (not used)</li> <li>• Universal Test ID Name (not used)</li> <li>• Universal Test ID Type (not used)</li> <li>• Manufacturer's or local code (6 characters):</li> </ul> This is the code defined in the <i>Analyzer</i> .
6	Priority	<b>I</b>	<b>N</b>		
7	Request Ordered Date/Time	<b>I</b>	<b>N</b>		
8	Specimen Collected Date/Time	<b>I</b>	<b>N</b>		
9	Collection End Time	<b>I</b>	<b>N</b>		
10	Collection Volume/Units	<b>I</b>	<b>N</b>		
11	Collector ID	<b>I</b>	<b>N</b>		
12	Action Code	<b>I</b>	<b>N</b>		
13	Danger Code	<b>I</b>	<b>N</b>		

14	Relevant Clinical Info.	<b>I</b>	<b>N</b>		
15	Date/Time Specimen Received	<b>I</b>	<b>N</b>		
16	Specimen Type	<b>R</b>	<b>A</b>	1	This is a numeric field indicating the type of specimen: The Imm. System uses the following ASCII characters: <b>0</b> = Serum <b>1</b> = Urine
17	Ordering Physician	<b>I</b>	<b>N</b>		
18	Physician Tel. Number	<b>I</b>	<b>N</b>		
19	User Field No. 1	<b>I</b>	<b>N</b>		
20	User Field No. 2	<b>I</b>	<b>N</b>		
21	Lab Field No. 1	<b>I</b>	<b>N</b>		
22	Lab Field No. 2	<b>I</b>	<b>N</b>		
23	Date /Time Result Reported Last or Modified	<b>I</b>	<b>N</b>		
24	Instrument Charge	<b>I</b>	<b>N</b>		
25	Instrument Section ID	<b>I</b>	<b>N</b>		
26	Record Type	<b>I</b>	<b>A</b>	1	The field indicates the direction of the transmission: <b>O</b> - Down Loading <b>F</b> - Up Loading
27	Reserved Field	<b>I</b>	<b>N</b>		
28	Location or Ward of Specimen Collection	<b>I</b>	<b>N</b>		
29	Nosocomial Infection Flag	<b>I</b>	<b>N</b>		
30	Specimen Service	<b>I</b>	<b>N</b>		

31	Specimen Institution	<b>I</b>	<b>N</b>		
<b>Legend:</b>		<b>R</b> Required	<b>D</b> Down Load		
		<b>O</b> Optional	<b>U</b> Up Load		
		<b>I</b> Ignored	<b>N</b> Never		
		<b>A</b> Always	<b>S</b> Sometimes		

**Example Test Order Record Layouts (O)**

<b>Download</b>	
<b>Host</b>	O   1   AR102     ^^GLU                   0                 O
<b>Analyzer System</b>	O   1   AR102     ^^GLU                   0                 F

**8.6 RESULTS RECORD (R)**

Field	Field Title	Down Load	Up Load	Max Len	
1	Record Type ID		<b>A</b>	1	This is a required field that contains an “R” identifying it as a Results Record.
2	Sequence Number		<b>A</b>	3	This field starts with “1” for the first result and is incremented by 1 for each additional result within the record. This will be reset to “1” when the results from another test order record are transmitted to the laboratory computer.



3	Universal Test ID	<b>I I I R</b>	<b>N N N A</b>	9	This is a four-component field: <ul style="list-style-type: none"> <li>• Universal Test ID Code (not used)</li> <li>• Universal Test ID Name (not used)</li> <li>• Universal Test ID Type (not used)</li> <li>• Local Manufacturer's or local code (6 characters) this is the code defined in the <i>Analyzer</i>.</li> </ul>
4	Data or Measurement value		<b>A</b>	10	'Data' is a 10-character, floating point field that includes the decimal point. The number of precision point digits will vary according to the test and is configurable on the <i>Analyzer</i> .
5	Units of Measure		<b>A</b>	6	This is a field for up to 6 characters that the operator defines for analytic measurement.
6	Reference Ranges		<b>A</b>	21	This field has two components; one giving the lower limit and the other the upper limit of the range. The format for this field is N^N.
7	Result Abnormal Flags		<b>A</b>	2	This field indicates the normal status of the result. The following codes are valid values:  L - Below Low normal H - Above High normal LL - Below Panic normal HH - Above Panic normal  < - Below absolute low (under linearity) > - Above absolute high (over linearity) N - Normal A - Abnormal E – Edited
8	Nature of Abnormality Testing		<b>N</b>		
9	Result Status		<b>A</b>	1	The Imm. System currently implements only two valid values:  F - final results;

					V - operator verified/approved result.
10	Date of Change in Instrument Normative Values or Units		N		
11	Operator ID		N		
12	Date/Time Test Started		N		
13	Date/Time Test Completed		A	14	Date and Time of test completion: formatted as YYYYMMDDHHMMSS.
14	Instrument ID		N		
<b>Legend:</b>		<b>R</b> Required <b>O</b> Optional <b>I</b> Ignored <b>A</b> Always	<b>D</b> Down Load <b>U</b> Up Load <b>N</b> Never <b>S</b> Sometimes		

**Example Result Record Layouts (R)**

Upload	
<b>Analyzer System</b>	R   1   ^^GLU   70.97   UL   70^105   N     F       20020923114302

**8.7 MESSAGE TERMINATOR RECORD (L)**

Field	Field Title	Down Load	Up Load	Max Len	Description and Valid Values
1	Record Type ID	R	A	1	This is a required field that contains an “L” identifying it as an Message Terminator Record.
2	Sequence Number	R	A	1	For a message terminator, this message should always be “1”.
3	Termination Code	R	A	1	This indicates the cause of termination. The following codes are valid values for the <i>Analyzer</i> : Null or N-normal termination
<b>Legend:</b>		<b>R</b> Required <b>O</b> Optional <b>I</b> Ignored <b>A</b> Always	<b>D</b> Down Load <b>U</b> Up Load <b>N</b> Never <b>S</b> Sometimes		

**Example Message Terminator Record Layout (L)**

<b>Host</b>	<b>L I I I N</b>
<b>Analyzer System</b>	<b>L I I I N</b>

**8.8 REQUEST INFORMATION RECORD (Q):**

<b>Field</b>	<b>Field Title</b>	<b>Down Load</b>	<b>Up Load</b>	<b>Max Len</b>	<b>Description and Valid Values</b>
1	Record Type ID		<b>A</b>	<b>1</b>	This is a required field that contains a "Q" identifying it as a request.
2	Sequence Number		<b>A</b>	<b>1</b>	It is always "1".
3	Starting Range ID Number		<b>A</b>	<b>31</b>	This field can either be:  "ALL" - to mean all demographics and tests being ordered should be sent to the instrument at this time,  or can have two components: <ul style="list-style-type: none"> <li>• Computer system patient ID No. (up to 15 characters);</li> <li>• Computer system specimen ID No. (up to 15 characters).</li> </ul>
4	Ending Range ID Number		<b>N</b>		
5	Universal Test ID		<b>N</b>		
6	Nature of Request Time Limits		<b>N</b>		
7	Beginning Request Results Date and Time		<b>N</b>		
8	Ending Request Results Date and Time		<b>N</b>		
9	Requesting Physician Name		<b>N</b>		
10	Requesting Physician Telephone Number		<b>N</b>		
11	User Field No. 1		<b>N</b>		
12	User Field No. 2		<b>N</b>		
13	Request Information Status Codes		<b>A</b>	<b>1</b>	It is always "O" (requesting test orders and demographics only).

**Example Request Information Record Layouts (Q)**

Download To	
<i>Analyzer System</i>	H   \^&       SHAnalyzer               20020927100402
	Q   1   ALL               O
	L   1   N

**CHAPTER 09**

---

**- ERROR SIGNALING AND TROUBLESHOOTING -**

**INDEX**

9.1 ERROR SIGNALING.....2

9.1.1 SYSTEM ERRORS..... 2

9.1.2 RESULT FLAGS ..... 3

9.1.3 DESCRIPTION OF RESULT FLAGS ..... 4

9.2 TROUBLESHOOTING GUIDE ..... 8

## 9.1 ERROR SIGNALING

This chapter is dedicated to a description of the **error** signaling which may occur during the programming or carrying out of the various operations. Error signals can be divided into the following two groups:

- ? System Errors
- ? Result Flags

### 9.1.1 SYSTEM ERRORS

Whenever a system error is detected, it is signaled via the activation of the below-listed warning lights/buttons:

? **Warning Light:** a triangle-shape located in the upper, right-hand portion of the System Monitor screen. It lights up yellow when activated.

If this yellow triangle does light up, the operator need simply click on it to access the relative visual text warning message indicating the cause of the signaled anomaly (said window will open in that area dedicated to viewing data).

Following is a list of possible visual text “Warning!” messages:

- Liquid Alarm I (Rinse Solution)
- Liquid Alarm II (Distilled Water)
- Liquid Alarm III (Cleaning Solution)
- Temperature Out-of-Range
- Host Serial Port cannot be opened

? **Fatal Error:** an “X”-shape located in the upper, right-hand portion of the System Monitor mask. It lights up red when activated.

If this red “X” lights up, the operator need simply click on it to access the relative visual text message explaining the cause of the Fatal Error (said window will open in that area of the mask dedicated to viewing data).

Following is a list of possible Fatal Error visual text messages:

- Programming Error
- Internal Arm Error
- External Arm Error

- Filter Error
- “Z” Axis Error
- Cuvette Block Error
- Wash System Error
- Diluter Error
- Block due to Instrument Lid being open(ed) during functioning
- Macro Error
- Maximum Time Error
- Instrument either not connected or turned off
- Decode Unit cannot be launched
- Host Link Module cannot be launched
- Host Link Module launching Error
- Work Executor cannot be launched
- Error during execution of “%1” command

### 9.1.2 RESULT FLAGS

Result Flags are categorized under the following group headings, Each group is identified by a symbol, as listed here-below:

<b>X</b>	<b>PHYSICAL ERRORS</b>
<b>R</b>	<b>CONCENTRATION ERRORS</b>
<b>C</b>	<b>CALIBRATION ERRORS</b>
<b>A</b>	<b>OPTIC DENSITY ERRORS</b>
<b>E</b>	<b>RESULTS EDITED MANUALLY</b>
<b>I</b>	<b>ISE MODULE ERRORS</b>
<b>?</b>	<b>PROGRAMMING ERRORS</b>

Every Result Flag, signaling an error, is accompanied by a symbol representing the group it is part of. The operator need only click on the small red square ( ) next to the Result Flag symbol to access the visual text message explaining the cause(s) of the signaled error.

In the central column of the following tables, the user will find those symbols which signal the type of error encounter, as used in the print-out of the final report. These symbols can be modified by the operator, in the **Parameters** section, under **print options**.

**WARNING!** The use of the Error Symbols in the print-out of the results (inclusion and/or exclusion) **IS UNDER THE DIRECT AND SOLE RESPONSIBILITY OF THE USER**

### 9.1.3 DESCRIPTION OF RESULT FLAGS

#### ➤ **X Physical Errors**

<b>Temperature Error:</b>	<b>T</b>	Reaction temperature (Reaction Plate) is out-of-range.
<b>No Sample</b>	<b>S</b>	Either no sample or sample serum quantity below minimum or above maximum level for the declared container.
<b>No Reagent:</b>	<b>R</b>	Either no reagent or reagent level below minimum or above maximum level for the declared container.
<b>No Rack:</b>		* No Rack present during sampling.

#### ➤ **R Concentration Errors**

<b>Very Low and Very High:</b>	<b>L-H</b>	Flags determined by test results out-of-range as setup in the Methods.
<b>Low Alert and High Alert:</b>	<b>A</b>	Flags determined by test results out-of-range as setup in the Methods.
<b>Low and High Linearity Limit:</b>	<b>G</b>	Flags determined by test results out-of-range as setup in the Method.
<b>Calculation Error:</b>	<b>C</b>	Concentration calculation error due to foreseeable causes (Asymptote).



## ➤ C Calibration Errors

<b>RBL missing:</b>	No Reagent Blank Level.
<b>Calibration missing:</b>	*No Standard or no Calibration curve.
<b>STD Replicate insufficient:</b>	*Insufficient number of valid Standard Replicates.
<b>STD Replicate outside CV%:</b>	*Coefficient of Variation Percentage in the Standard Replicates over the set value.
<b>Invalid Calibration:</b>	*Calibration curve not valid – either because it is not monotonic or because the Fit is above the set value.

## ➤ A Optic Density Errors

<b>Inversion:</b>	<b>I</b>	Reaction direction not in line with that set-up.
<b>End Point Limit:</b>	<b>P</b>	Values over the limits setup in the Methods Parameters
<b>Depletion Limit:D</b>		Values over the limits setup in the Methods Parameters.
<b>First Limit:</b>		*Values over the limits setup in the Methods Parameters.
<b>FIT:F</b>		Values over the limits setup in the Methods Parameters.
<b>RBL out-of-range:</b>		*Reagent Blank Levels outside the range Setup.
<b>Sample outside Standard:</b>		# Sample absorbance outside the calibration curve.

## ➤ **E** **Results Edited Manually**

**Results Edited:** **E** This symbol automatically appears whenever the operator has manually modified the obtained results. This operation annuls all the symbols indicating errors which, in this case, will not be viewed.

## ➤ **I** **ISE Module Errors**

**Air in the ISE Module:** \*Air present in the ISE Module hydraulic circuit.

**Calibrate A drift:** \*Calibrate A drift in the ISE Module.

**Noise in the ISE Module:** \*Background noise present – values not reliable.

**Out-of-Range in ISE Module:** \*Linearity values out-of-range.

## ➤ **?** **Programming Errors**

All those software errors which are deemed unforeseeable are indicated using this symbol.

## ➤ **Result Asterisks**

In the event that the following situations and/or error conditions occur, the system will notify the user via a visualization of two asterisks (\*\*) in the Results field:

### PHYSICAL ERRORS

1 – very little or no sample

2 – very little or no reagent

3 – no rack

## CONCENTRATION ERRORS

- 4 – calculation error

## PROGRAMMING ERROR

- 5 – all cases

**In the following situations and/or error conditions, the results field will contain a “0”:**

- 1 – if the result is less than zero
- 2 – if the result is equal to zero

## OPTIC DENSITY ERRORS

- 3 – Flag signaling depletion limit
- 4 – Flag signaling inversion

## 9.2 TROUBLESHOOTING GUIDE

PROBLEM	? POSSIBLE CAUSE	❖ SOLUTION
Repeatability of results insufficient	<p>? Sample Probe dirty</p> <p>❖ Clean the Sample Probe as described in Chapter 07 Maintenance</p> <p>? Hydraulic leak and/or air bubble in the hydraulic circuit (sampling)</p> <p>❖ Check Sample Probe fit</p> <p>❖ Check the fit of the hydraulic tubes and their connection: if necessary, substitute the tubes and/or adjust their connections/fittings</p> <p>? Wash solution is contaminated. If the wash solution contains contaminating particles (e.g.: mold, dust, lint), these micro-particles can cause errors during the WBL running.</p> <p>❖ Change the wash solution</p> <p>❖ Clean the Wash Solution Bottle(s) and carry out the Hydraulic Circuit wash procedure as described in Chapter 07 – Maintenance</p> <p>? Deterioration of the Reagent(s)</p> <p>❖ Substitute the bad reagent(s)</p> <p>? Reaction cuvettes not dried correctly after being washed</p> <p>❖ Check the Tip used to dry the cuvettes after washing to make sure it is in good working condition. If necessary, clean the Tip or change it, following the procedure described in Chapter 07 – Maintenance</p> <p>? Light bulb not stable</p> <p>❖ Light bulb nearing the end of its 2000-hour life cycle duration, or premature deterioration. In both cases, change the light bulb following the substitution procedure described in Chapter 07 – Maintenance.</p>	
Insufficient volume/quantity of the various Rinse, Cleaning and Distilled Water solutions.	<p>? This type of problem can present itself either when the instrument is first turned on due to a lack of liquid in one or more bottles, or during test running whenever the involved liquid has been finished</p> <p>❖ Refill the required bottle. If the instrument is running tests, it will automatically <b>Pause</b>. Wait until sampling is suspended, then refill the involved bottle with the required liquid (Rinse/Cleaning/ Distilled Water). Press <b>Start</b> to continue running the rest of the programmed tests. This operation will not adversely affect the already run tests.</p>	

# CHAPTER 12

---

## – ISE MODULE -

### INDEX

12.1	ELECTROLYTE MEASUREMENT SYSTEM.....	2
12.1.1	PRODUCT DESCRIPTION.....	2
12.2	ELECTROLYTE MEASUREMENT SYSTEM FEATURES AND BENEFITS .....	3
12.3	TECHNICAL SPECIFICATION .....	4
12.4	SYSTEM OPERATION – ISE THEORY .....	5
12.5	MECHANICAL FEATURES .....	6
12.6	ELECTRODES.....	7
12.7	FLUID MANAGEMENT .....	8
12.7.1	REAGENT USED .....	8
12.7.2	REAGENTS, CALIBRATION AND SAMPLE PROCESSING.....	9
12.8	ELECTRONICS .....	11
12.8.1	GENERAL DESCRIPTION .....	11
12.9	MAINTENANCE .....	11
12.10	TROUBLE SHOOTING GUIDE .....	12

## 12.1 ELECTROLYTE MEASUREMENT SYSTEM

### 12.1.1 PRODUCT DESCRIPTION

Electrolyte Measurement System includes a small, simple, and reliable ISE module and two peristaltic pumps designed to be mounted within a existing chemistry analyzer. The module measures the concentration of Sodium, Potassium and Chloride. The module contains an integral sample entry port, positioned at the top. The compact design allows for small sample size and fast operation. The modules requires only a 70 µl sample.

The module housing contains snap-in, snap-out ISE sensor which, through simple connectors, connect directly to an electronic board in close proximity to the module. This eliminates the need for cables and maximizes noise immunity. Samples and calibrators are positioned within the module by two snap-in, snap-out peristaltic pump cassettes. The waste pump positions the sample in front of the sensor for measurement. After sample measurement, a calibrant/wash solution is pumped in front of the sensor for a single point calibration.

Provision for two-point sensor calibration and cleaning are made through use on the analyzer sample tray which serve as reservoirs for the second calibrator and cleaning solution.

The module is completely self-contained. All sample and calibrant positioning within the module is controlled by an integral microprocessor, which assure reliable electrode operation and maximum lifetime.

## 12.2 ELECTROLYTE MEASUREMENT SYSTEM FEATURES AND BENEFITS

<u>Features</u>	<u>Benefits</u>
Integral Sample Entry port	Minimal Sample Carry Over
Small Sample Size	70 ul
Electrodes Mounted Close to Electronics	Minimal electronic noise improves precision
Sample Port Adjacent to ISE Sensor	Minimal Sample Size and Carry Over
Rapid Operation (30 Sec Cycle Time)	Rapid electrolyte results
Easy to access sensors	Simple Maintenance
No membranes	Maintenance performed by lab personnel
Easily accessed pumps	Simple maintenance
Two bubbles detectors	Reliable sample processing assured
Two point calibration	High accuracy and precision
One point calibration with every sample	High accuracy and precision
Maintenance free electrodes (6 month or 10.000 sample warranty)	Low cost per test
Disposable reference electrode, no addition of filling solution required	Convenient maintenance

## 12.3 TECHNICAL SPECIFICATION

**Sample:** Whole Blood, Serum, Plasma or Urine

(Urine requires dilution 1:10)

**Sample Size:** 70 ul, (3 channels) serum; 160 ul undiluted urine

**Reproducibility:** Maximum imprecision (within run) Typical Carry Over, % (Serum)

Sodium	CV<1,5% (100 - 160 mmol/L)	<0,5%
Potassium	CV<2% (3,00 - 6,00 mmol/L)	<1,5%
Chloride	CV<2% (80,0 - 120,0 mmol/L)	<1,0%

**Analysis Time:** Serum – 30 seconds, including one point calibration

Urine – 60 seconds, including one point calibration

**Throughput:** Serum – 120 sample/hour

Urine – 60 sample/hour

**Reagents:** Calibrant “A”

Calibrant ”B”

Cleaning Solution

Urine Diluent

**Maximum Environnemental Temperature:** 38°C

(Host Analyzer requires working temperature within:  $18^{\circ}\text{C} \leq T \leq 32^{\circ}\text{C}$ )



## 12.4 SYSTEM OPERATION – ISE THEORY

Electrolyte Measurements System measure sodium, potassium and chloride in biological fluids, using ion selective electrode technology. The flow-through sodium electrode selective PVC membrane tubing, specially formulated to be sensitive to sodium ions. The potassium and chloride electrodes employ similar designs with appropriate selective membrane materials. The potential of each electrode is measured relative to a fixed, stable voltage established by the double-junction silver/silver chloride reference electrode. An ion selective electrode develops a voltage that varies with the concentration of the ion to which it responds. The relationship between the voltage developed and the concentration of the sensed ion is logarithmic, as expressed by the Nerst equation:

$$E = E^{\circ} + RT \log (e C)/nF$$

Where:

$E$  = The potential of the electrode in sample solution

$E^{\circ}$  = The potential developed under standard conditions

$RT/nF$  = A temperature dependent “constant”, termed the slope(s)

Log = Base ten logarithm function

$e$  = Activity coefficient of the measured ion in the solution

$C$  = Concentration of the measured ion in the solution

A comparative method of measurement is utilized. First, the ISE module measures the potential developed when the sample is positioned in the electrodes. Next, Calibrant A is positioned in the electrodes. The difference in the two potentials is related logarithmically to the concentration of the measured ions in the sample divided by their respective concentrations in the Calibrant solution. Since the difference in potentials and the concentration of the sodium, potassium or the other ions in the Calibrant solution are known, the computer can calculate the concentration of the ions in the sample solution, in accordance with the Nerst equation, rewritten as:

$$E - E^{\circ} = S \log (C_x/C_s) \text{ or } C_x = C_s \times 10^{[(E - E^{\circ})/S]}$$

Where:

$E$  = ISE potential developed in sample solution

$E^\circ$  = ISE potential developed in the Calibrant solution

$S$  = Electrode slope calculated during calibration

$C_x$  = Concentration of ion in the sample

$C_s$  = Concentration of ion in the Calibrant solution

“S”, the slope, is determined during calibration using Calibrants A and B, which have known levels of sodium, potassium, and chloride. When an automatic calibration is initiated, the slope is calculated between the second Calibrant A reading and the Calibrant B reading.

Excessive drift or noise reading will be flagged and the appropriate error message sent to the host analyser from the ISE Module.

## 12.5 MECHANICAL FEATURES

The electrode housing contains each of the ion-selective electrodes, as well as the reference electrode.

Two bubble detectors are also included at both the top and bottom of the electrode chain. These are used to properly position the sample for measurement. A sample port is positioned directly above the chain of electrodes on the top of the module.

An electronic signal processing board is attached to the electrode housing. This board includes high input impedance operational amplifiers to detect the ISE signals and additional digital processing circuitry serving as an A/D converter and providing an ASCII signal output to the chemistry analyser.

Each of the electrodes can be easily removed from the front of the housing.

## 12.6 ELECTRODES

The electrodes are maintenance-free and are warranted on a prorated basis for up to 10,000 samples or 6 months, whichever occurs first. Cleaning Solution, aspirated from an operator designed sample cup, is used at least once a day at the end of the day in order to minimize protein build-up in the fluid lines. A two-point calibration of the ISE module is also done at least once a day at the beginning of the first sample run. If the user is running more than 50 samples a day, both cleaning and calibrant must be performed after 8 hours by the host analyzer.

The entire double-junction reference electrode is disposable. The reference electrode is filled with sufficient KCL so that no filling solution must be added during the lifetime of the electrode.

Electrodes require calibrant sampling at 30 minutes intervals for reliable operation, but this is completely controlled by the Electrolyte Measurement System without any need for control by the host analyzer or the operator.

The electrodes require a 10 times sample dilution for measurement of urine.

The ISE module depends on the host analyzer to perform the dilution function.

It is not necessary to regulate the electrode housing temperature provided that its environmental temperature does not exceed 38° C. However, the electrode module should not be subject to changes greater than plus or minus 8° C without recalibrating.

## 12.7 FLUID MANAGEMENT

### 12.7.1 REAGENT USED

The sample is aspirated from a sample cup and dispensed into the sample port at the top of the ISE module. The sample is then positioned in front of the sensor using the bubble detector and the Waste Pump.

Four reagents are needed to operate the ISE module.

#### **Calibrant “A”:**

Used as wash solution and single-point calibrator. Calibrant A is pumped into the sample port by the Calibrant A pump and then positioned in front of the sensors

#### **Calibrant “B”:**

Used as the second point in two-point calibration. Calibrant B is aspirated from a cup on the analyzer at least once a day or every 8 hours, depending upon the laboratory schedule.

A volume of 500 ul is sufficient for one day’s requirements.

The Calibrant “B” must be placed on the system just before use to prevent a change in values from evaporation.

#### **Cleaning Solution:**

Should be run once a day to prevent protein built up or at 8 hours intervals if the ISE Module performs greater than 50 samples per day

Cleaning Solution may be aspirated from a sample cup; 500 ul is sufficient for one day’s requirements.

#### **Urine Diluent:**

This is required for urine samples. Urine samples must be diluted manually by a factor of 10 to perform urine measurement.

## 12.7.2 REAGENTS, CALIBRATION, AND SAMPLE PROCESSING

The sequence of use of Calibrant A and patient samples during processing is as follows:

1. Sample deposited into ISE Module sample port by host analyzer;
2. Sample positioned in front of electrodes of Electrolyte Measurement System by the waste pump;
3. Sample equilibration and reading occurs during 7 second period;
4. Calibrant A pumped into electrode module;
5. Calibrant A equilibration and reading occurs during 7 second period;
6. Results transmitted to the host analyzer;
7. ISE module ready for next cycle.

Pumpings of a small amount of Calibrant A are performed when the Electrolyte measurement System is in “Standby” or when it is not being used in the “Sample” mode. This significantly improves the performance of the electrodes. The Electrolyte measurement System must always be supplied with power so that “pumping” can occur. Pumping occurs automatically, without prompting by the host analyzer, beginning 30 minutes after the last sample or calibration was performed.

During sample processing, a volume of 200 µl of Calibrant A solution is used for one point calibration, sample, wash, and cleaning. A volume of 120 µl is used for each sip.

The Electrolyte Measurement System should perform a two points calibration at the beginning of the sample run. If the user is running more than 50 sample a day, both cleaning and calibration must be performed after 8 hours by the host analyzer. 140 µl of Calibrant B solution are used during two point calibration. During two points calibration, electrode calibration slopes are transmitted by the module for QC purpose and may be used by the operator to diagnose module performance. The slope is defined as:

$$\text{Slope} = (E_B - E_A) / \log (C_B / C_A)$$

Where:

$C_A$  = Calibration A concentration in mmol/L

$C_B$  = Calibration B concentration in mmol/L

$E_A$  = ISE potential developed in Cal A solution in mV

$E_B$  = ISE potential developed in Cal B solution in mV

These slopes are checked by the module's electronic processor and an error code will be transmitted if they are outside the required range.

Typical slopes are approximately 55mV/decade for Na and K and 45mV/decade for Cl.

Acceptable limit slopes are:

<u>Slope (mV/decade)</u>		<u>Range (mmol/L)</u>		
		Serum		Urine
Na	50-63	Na	20-200	20-1000
K	50-63	K	0,2-20,0	1-50
Cl	40-59	Cl	25-200	20-500

## 12.8 ELECTRONICS

### 12.8.1 GENERAL DESCRIPTION

ISE module electronics include all pre-amplifiers and perform microprocessor control of the fluid pumps, A/D conversion and RS-232C communications. The microprocessors apply proprietary mathematical algorithms to electrolyte sensor output voltage, converting them to clinical units of mmol/L.

#### Sensor Inputs:

- Na Electrode;
- K Electrode;
- Cl Electrode;
- Reference Electrode;
- Upper Bubble Detector;
- Lower Bubble Detector.

## 10.9 MAINTENANCE

The Electrolyte Measurement System has been designed to require very little operator maintenance. The only daily maintenance required is to run the Cleaning Solution after the last sample of the day.

REPLACEMENT/PART	3	6	9	12
	MO	MO	MO	MO

Pump Tube	X			
Na Electrode	X			
K Electrode	X			
Cl Electrode	X			
Reference Electrode	X			
Reagents	Refill REAGENT as required by testing needs			

Recommended Maintenance/Replacement Schedule (more than 100 samples per day)

REPLACEMENT/PART	3	6	9	12
	MO	MO	MO	MO

Pump Tube	X			
Na Electrode	10.000 samples			
K Electrode	10.000 samples			
Cl Electrode	10.000 samples			
Reference Electrode	10.000 samples			
Reagents	Refill REAGENT as required by testing needs			

## 12.10 TROUBLE SHOOTING GUIDE

SYMPTOM	PROBLEM	CORRECTION
System does not respond	1. RS232 cable is disconnected or damaged.	Reconnect or replace cable
	2. Module connector has been damaged.	Replace board
	3. Component failure on board.	Replace board.
<b>Low Slope:</b> Na or K <45mV/decade, Cl<35 mV/decade or  <b>High Slope:</b> Na or K>63mV/decade, Cl>60 mV/decade	1. Misalignment of sensors.	Remove and replace sensors to re-seat.
	2. Deterioration of calibrator solutions	Replace Cal B first and retest. If still low, replace Cal A and retest.
	3. Deterioration of sensing electrode.	Replace problem sensor and test
	4. Air bubble on Reference Electrode membrane.	Remove electrode, tap to dislodge bubble, replace, and recalibrate.
	5. Deterioration of reference electrode.	Replace reference electrode and retest.
	6. Interaction between sensing electrodes	Replace Cl electrode only and retest.
	7. Temperature of the module or liquid higher than 37° C.	Check the temperature and if necessary change location
Noise error Flag Single electrode	1. Deterioration of sensing electrode	Replace problem sensor and test
	2. Electrical noise spike from environmental source	a) Check that motor or solenoid valve near module is not activated during the read portion of the module cycle. b) Component failure on module board. Replace board.
Noise Error Flag Multiple electrodes	1. Deterioration of reference electrode	Replace reference electrode and retest.
	2. Electrical noise spike from environmental source	a) Check that motor or solenoid valve near module is not activated during the read portion of the module cycle. b) Component failure on module board. Replace board.



Drift Error Flag Single electrode	1. Deterioration of sensing electrode	Replace problem sensor and test.
	2. May occur when new sensor or new bottle of Cal A is installed on system.	Purge the Cal A and recalibrate the module. If the sensor is new it may initially drift as it rehydrates over the course of 15 minutes.
Drift Error Flag Multiple electrodes	1. Deterioration of Reference electrode	Replace reference electrode and retest
	2. Electrical spike from environmental source	a) Check that motor or solenoid valve near module is not activated during the read portion of the cycle.
		b) Component failure on module board. Replace the board.
	3. May occur when new sensor or new bottle of Cal A is installed on system.	Purge the Cal A and recalibrate the module.
Air in Sample	1. Insufficient sample pipetted into module sample port	a) Host instrument must deliver 70 ul. Increase dispensed volume b) Insufficient sample in sample cup for all tests programmed.
	2. Sample not positioned properly.	a) Pumps not connected properly. b) Pump tubing obstructed or tubing length is excessive.
Air in Sample and Cal A	1. Sample and Cal A are segmented with air.	a) Sensors are not properly compressed. Check compression plate, spring and seal. b) Ensure that all sensor and o-rings are in place.
	2. Fibrin or salt is plugging the sensor flow path	c) Use Cleaning procedure <CLEN> for module d) Disassemble module and clean or replace sensor with plugged flow path.
	3. Air Bubble Detector failure	Replace air bubble detector
	4. Waste Pump failure	Replace the waste pump.

Air in Cal B and in Cal A	1. Cal B and Cal A are segmented with air.	a) Sensors are not properly compressed. Check compression plate, spring and seal. b) Ensure that all sensor and o-rings are in place.
	2. Fibrin or salt is plugging the sensor flow path	a) Use Cleaning procedure <CLEN> for module b) Disassemble module and clean or replace sensor with plugged flow path.
	3. Waste Pump failure.	Replace the waste pump.
	4. Air Bubble Detector failure.	Replace air bubble detector.
Air in Cal B	1. Insufficient Cal B piped into module sample port	a) Host Instrument must deliver 70 µl. Increase dispensed volume b) Insufficient sample in sample cup for all tests programmed.
	2. Sample not positioned properly.	a) Pumps not connected properly. b) Pump tubing obstructed or tubing length is excessive for spring and seal.
Air in Cal A (no “Air” errors reported for Sample or Cal B)	1. Cal A bottle is empty	Replace Cal A bottle with a new one and recalibrate.
	2. Tubing is disconnected.	Reconnect or replace tubing
	3. Cal A pump is not working properly.	a) Check electrical connections. b) Replace pump cassette. c) Replace motor.
	4. Tubing is plugged, split or crimped.	Replace tubing.

## Chapter 13

---

# Spare Parts

---

This chapter contains the part lists for all the servicing aids and replaceable parts for all the ILab 300 Plus systems.

For easier identification of the parts, this section includes the following documentation:

List Sorted by IL Part Number

Pag 2

List Sorted by Description

Pag 5

<b>I.L Part Number</b>	<b>Description</b>
23010120100	Gauge Peristaltic Pump
23010122400	Joint Bottles (Rinse,Water, Cleaning)
23050040300	Micro Pump 4
23050040400	Micro Pump 2
23050040500	Connection Joint Assembly
23050041300	Interlock Assy
23050044700	Micro Pump 3
23050045901	Stat Sample Rack
23050051900	Reagent Panel
23050052000	Sample Panel
23100010800	Assy Photometer
23100020201	Electronic Controller Assy
23100028601	Peristaltic Pump Assembly
23100039900	Motors Assy for: S.Arm+ P.Pump
23135000500	Belt (arm horizontal)
23300010700	Pre-Ampl./ADC Board Main Channel
23300010703	Pre-Ampl./ADC Board from S/N 0205011 Ref Channel
23300010705	Pre-Ampl./ADC Board up to S/N 0205010 Ref Channel
23300018801	Stepper Motor Drivers Bd.
23300019201	Sampling Bd
23300019301	Reaction Tray Interface Bd.
23300019501	Analytical Control Bd.
23300034100	Motor Plate Interface Bd.
23300034200	Control Panel
23300041900	Hydraulic Interface Bd.
23300044600	Racks Identification Bd.
23350049300	Interlock Switch
23350076100	Extractor for connectors Lumberg
23350081400	Sodium Electrode
23350081500	Potassium Electrode
23350081600	Chloride Electrode
23350081700	Reference Electrode
23350081800	Ise Control Board
23350081900	Bubble Detector
23350082700	Peristaltic Pump Waste Assy
23350082000	Peristaltic Pump Cal A Assy
23350082800	Pump Motor Cal A
23350082900	Pump Motor Waste
23350083000	Head Tube Perilstatic Pump
23350083100	Ise Connection Tube 6mt
23350083200	O-Ring 4 Pcs
23350083300	Motor Cable
23350083400	Electrode Compression Plate
23350083500	Alignment Electrode Tool
23350083600	Inlet Sample Pot Assy
23350083700	Output Connector Assy
23350083800	Lower Kit Gasket
23350083900	Gasket (compression)

<b>I.L Part Number</b>	<b>Description</b>
23500034600	Cable W1
23500034700	Cable W18
23500034701	Cable W2
23500034800	Cable W3
23500034900	Cable W4
23500035000	Cable W5
23500035100	Cable W6
23500035200	Cable W7
23500035201	Cable W8
23500035202	Cable W9
23500035203	Cable W19
23500035300	Cable W11
23500035600	Cable W14
23500035700	Cable W13
23500036500	Cable W12
23500047000	Cable W20
23550120400	Kit Measurament Temperature
23650067800	Bar Codes kit ILLab 300
23900125300	Tygon Tube M.1
23901003800	Drying Pad
23901055800	Special Key
23905004200	Cannula "A0"
23905004201	Cannula "A1"
23905004300	Cannula "B"
23905004400	Cannula "C"
23905004500	Cannula "D"
23905006300	Washing Well Assembly
23905006401	Sampling Probe
23905007200	Reagent Rack 3 (from 14-23)
23905007300	Reagent Rack 4 (from 24-33)
23905007400	Reagent Rack 1 (from 1-4 + std/ctrl)
23905007500	Reagent Rack 2 (from 5-13)
23905007800	Bottle level Sensor assembly
23910000600	Washing Station Probe Assembly
23910001300	Cover Lock assembly
23910001601	Stepper Motor Assembly (motor with pulley)
23910001700	Stepper Motor (washing probes)
23910001901	Stepper Motor Assembly (Photometer)
23910002200	Temperature Sensor
23910002300	Horizontal Home Sensor Assembly (filter wheel, reaction plate, washing station)
23910002400	Home Sensor Assembly - vertical
23910002801	Micro Pump Assembly (Probe rinse)
23910003301	Front Panel Switch assembly
23910006000	Cuvette Holder Solenoid
23910006101	Reaction Plate Assy
23910006201	Complete Probe Assembly
23910006300	Level Sensor and Pre-heater Assy

<b>I.L Part Number</b>	<b>Description</b>
23915001800	Ise Module
23915002000	Arm Assembly
23915002400	Refrigerated Reagent Rack assembly
23930000400	Lamp P.W.S Board
23930001400	Label Index Board
23935001600	Halogen Lamp
23935002201	Reaction Plate Heater
23935003100	Belt (arm vertical)
23935003200	Belt for Reaction Plate
23935003401	Power Supply
23935003500	Fluidic Valve 2 ways W.P
23935003600	Fluidic Valve 3 ways W.P
23935004000	Bar Code Reader
23935004100	Bottle 2lt
23935004200	Main Power Socket
23935005501	Interconnecting Serial Cable 9 Pin 5m
23950001000	Cable
23950003000	Cable with Button (SW1) (sample&reagent Door)
23965001900	kit Ise Module
23965002701	Tubes Kit
23965002900	Interferential filters Kit (Photometer)
23965003100	Reaction Cuvette (60 pieces)
23965003200	Fuse Holder Kit
23965003500	Lamps kit (5 pieces)
23965003601	Caps Kit (yellow, green, red)
23965003700	Button Kit
23965004000	Kit Tubing peristaltic pump
23990090100	Bearing Reaction Plate
23990090600	Kit Diluter Joints
23650022702A	PC Boards Service kit ILab 300
23AS620005	Power Cord
23C10048300	Tubes Adapter (20 pcs)
23C101004201	Head Diluter Joint
23C1010122200	Joint Diluter 10 Pcs
23C101023000	Metal Rings for lev. Sens. & pre-heater fixing (10 pcs)
23C101023200	Locknut (10 pcs)
23C101023600	Level Sensor & pre-heater springs (10 pcs)
23C101057200	Level Sensor & pre-heater assy lock. Key (10 pcs)
23C135003700	Fuse 6,3x32 (10 pcs)
23C1900125400	Joint E.V (Rinse) 10 Pcs
23F35001800	Diluter 1000 microl
23F35001900	Solenoid Valve - 2 way
23F50002000	Bottle Sensor Cable

<b>I.L Part Number</b>	<b>Description</b>
23350083500	Alignment Electrode Tool
23300019501	Analytical Control Bd.
23915002000	Arm Assembly
23100010800	Assy Photometer
23935004000	Bar Code Reader
23650067800	Bar Codes kit ILab 300
23990090100	Bearing Reaction Plate
23135000500	Belt (arm horizzontal)
23935003100	Belt (arm vertical)
23935003200	Belt for Reaction Plate
23935004100	Bottle 2lt
23905007800	Bottle level Sensor assembly
23F50002000	Bottle Sensor Cable
23350081900	Bubble Detector
23965003700	Button Kit
23950001000	Cable
23500034600	Cable W1
23500035300	Cable W11
23500036500	Cable W12
23500035700	Cable W13
23500035600	Cable W14
23500034700	Cable W18
23500035203	Cable W19
23500034701	Cable W2
23500047000	Cable W20
23500034800	Cable W3
23500034900	Cable W4
23500035000	Cable W5
23500035100	Cable W6
23500035200	Cable W7
23500035201	Cable W8
23500035202	Cable W9
23950003000	Cable with Button (SW1) (sample&reagent Door)
23905004200	Cannula "A0"
23905004201	Cannula "A1"
23905004300	Cannula "B"
23905004400	Cannula "C"
23905004500	Cannula "D"
23965003601	Caps Kit (yellow, green, red)
23350081600	Chloride Electrode
23910006201	Complete Probe Assembly
23050040500	Connection Joint Assembly
23300034200	Control Panel
23910001300	Cover Lock assembly
23910006000	Cuvette Holder Solenoid
23F35001800	Diluter 1000 microl
23901003800	Drying Pad
23350083400	Electrode Compression Plate

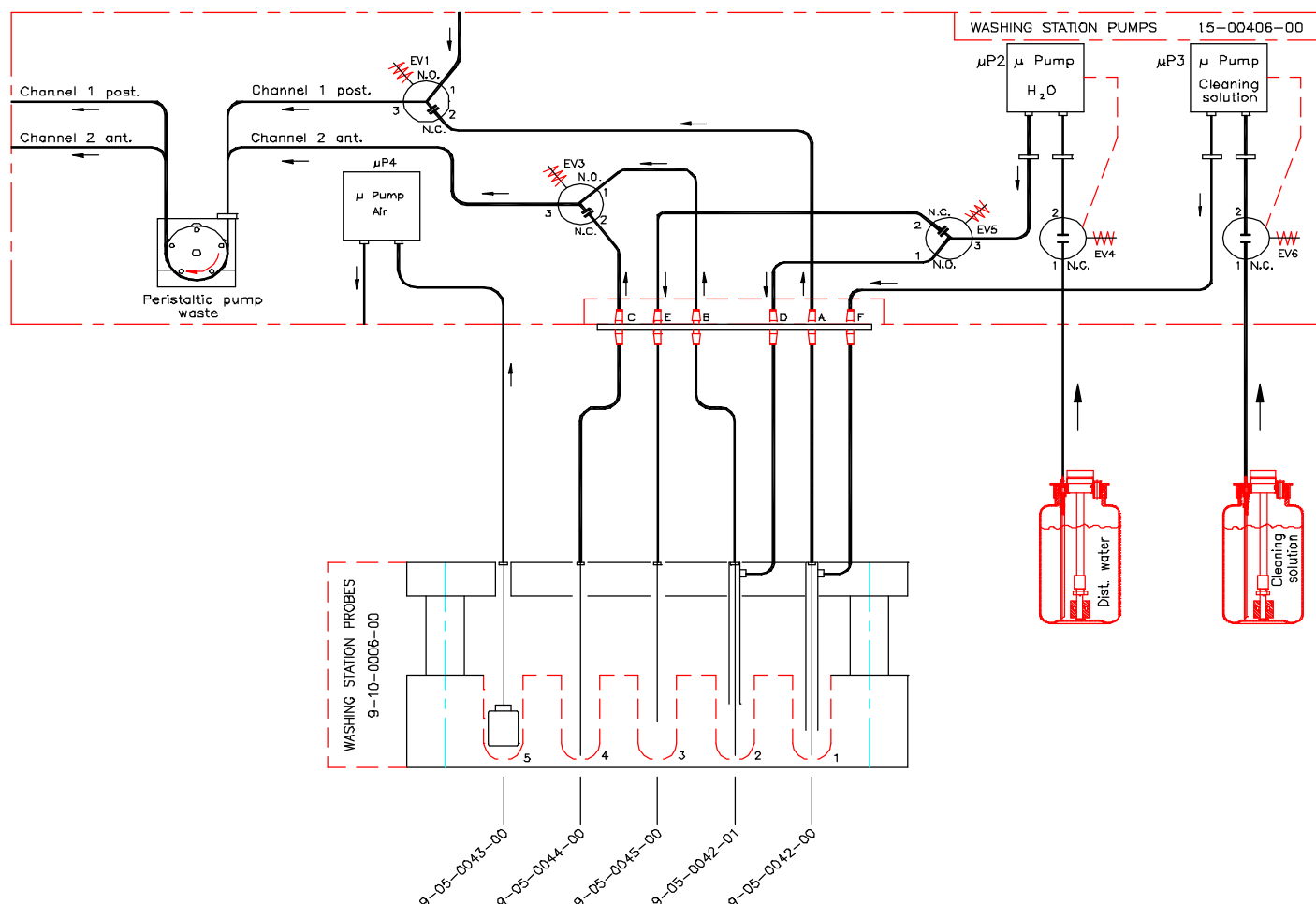
<b>I.L Part Number</b>	<b>Description</b>
23100020201	Electronic Controller Assy
23350076100	Extractor for connectors Lumberg
23935003500	Fluidic Valve 2 Ways W.P
23935003600	Fluidic Valve 3 Ways W.P
23910003301	Front Panel Switch assembly
23C135003700	Fuse 6,3x32 (10 pcs)
23965003200	Fuse Holder Kit
23350083900	Gasket (compression)
23010120100	Gauge Peristaltic Pump
23010122400	Joint Bottles (Rinse, Water, Cleaning)
23C1010122200	Joint Diluter 10 Pcs
23C1900125400	Joint E.V (Rinse) 10 Pcs
23935001600	Halogen Lamp
23C101004201	Head Diluter Joint
23350083000	Head Tube Perilstatic Pump
23910002400	Home Sensor Assembly - vertical
23910002300	Horizzontal Home Sensor Assembly (filter wheel, reaction plate, washing station)
23300041900	Hydraulic Interface Bd.
23350083600	Inlet Sample Pot Assy
23935005501	Interconnecting Serial Cable 9 Pin 5m
23965002900	Interferential filters Kit (Photometer)
23050041300	Interlock Assy
23350049300	Interlock Switch
23350083100	Ise Connection Tube 6mt
23350081800	Ise Control Board
23915001800	Ise Module
23990090600	Kit Diluter Joints
23550120400	Kit Measurament Temperature
23965001900	kit Ise Module
23965004000	Kit Tubing peristaltic pump
23930001400	Label Index Board
23930000400	Lamp P.W.S Board
23965003500	Lamps kit (5 pieces)
23C101057200	Level Sensor & pre-heater assy lock. Key (10 pcs)
23C101023600	Level Sensor & pre-heater springs (10 pcs)
23910006300	Level Sensor and Pre-heater Assy
23C101023200	Locknut (10 pcs)
23350083800	Lower Kit Gasket
23935004200	Main Power Socket
23C101023000	Metal Rings for lev. Sens. & pre-heater fixing (10 pcs)
23050040400	Micro Pump 2
23050044700	Micro Pump 3
23050040300	Micro Pump 4
23910002801	Micro Pump Assembly (Probe rinse)
23350083300	Motor Cable
23300034100	Motor Plate Interface Bd.
23100039900	Motors Assy for: S.Arm+ P.Pump
23350083200	O-Ring 4 Pcs



<b>I.L Part Number</b>	<b>Description</b>
23350083700	Output Connector Assy
23650022702A	PC Boards Service kit ILab 300
23100028601	Peristaltic Pump Assembly
23350082000	Peristaltic Pump Cal A Assy
23350082700	Peristaltic Pump Waste Assy
23350081500	Potassium Electrode
23AS620005	Power Cord
23935003401	Power Supply
23300010700	Pre-Ampl./ADC Board Main Channel
23300010703	Pre-Ampl./ADC Board from S/N 0205011 Ref. Channel
23300010705	Pre-Ampl./ADC Board up to S/N 0205010 Ref Channel
23350082800	Pump Motor Cal A
23350082900	Pump Motor Waste
23300044600	Racks Identification Bd.
23965003100	Reaction Cuvette (60 pieces)
23935002201	Reaction Plate Heather
23300019301	Reaction Tray Interface Bd.
23050051900	Reagent Panel
23910006101	Reaction Plate Assy
23905007400	Reagent Rack 1 (from 1-4 + std/ctrl)
23905007500	Reagent Rack 2 (from 5-13)
23905007200	Reagent Rack 3 (from 14-23)
23905007300	Reagent Rack 4 (from 24-33)
23350081700	Reference Electrode
23915002400	Refrigerated Reagent Rack assembly
23050052000	Sample Panel
23300019201	Sampling Bd
23905006401	Sampling Probe
23350081400	Sodium Electrode
23F35001900	Solenoid Valve - 2 way
23901055800	Special Key
23050045901	Stat Sample Rack
23910001700	Stepper Motor (washing probes)
23910001601	Stepper Motor Assembly (motor with pulley)
23910001901	Stepper Motor Assembly (Photometer)
23300018801	Stepper Motor Drivers Bd.
23910002200	Temperature Sensor
23900125300	Tygon Tube 1 m
23C10048300	Tubes Adapter (20 pcs)
23965002701	Tubes Kit
23910000600	Washing Station Probe Assembly
23905006300	Washing Well Assembly

# APPENDIX

## Washing Station's cycle :



**Rate of flow of the Peristaltic Pump = 1000ul/second**

\*\*\*\*\*

### Probe N° 1

**ON = EV1** (emptying of the channel of the Probe1 is enabled), **uP4** (ON uPump Probe N°5) , **P.M0** (ON peristaltic pump, remains ON until the end of the cycle).

**Washing Station goes down (570mS).**

**Liquid into the cuvette is aspirated.**

**Aspirated and dispenced 400uL(+50uL) of Cleaning solution.** While emptying with the Peristaltic Pump

**Wait 100mS.**

**Aspirated and dispenced 400uL(+50uL) of Cleaning solution.** . While emptying with the Peristaltic Pump

**Wait 300mS.**

**Washing Station goes up (570mS).**

**OFF = EV1** (emptying of the channel of the Probe1 is disabled), **uP4** ( OFF uPump 4 Probe5) ,  
**P.M0** (OFF Pump peristaltic).

**Cycle lasts 3790mS used 800ul(+/-100ul) of Cleaning solution.**

\*\*\*\*\*  
\*\*\*\*\*

## **Probe 2**

**ON = EV1** (emptying of the channel of the Probe1 is enabled), **uP4** (ON uPump Probe5) , **P.M0**  
ON peristaltic pump, remains ON until the end of the cycle

**Washing Station goes down** (570mS).

**Liquid into the cuvette is aspirated.**

**Aspirated and dispensed 400uL(+/-50uL) of Dist.Water.** While emptying with the Peristaltic  
Pump

**Wait 100mS.**

**Aspirated and dispensed 400uL(+/-50uL) of Dist.Water.** While emptying with the Peristaltic  
Pump

**Wait 100mS.**

**Aspirated and dispensed 400uL(+/-50uL) of Dist.Water.** While emptying with the Peristaltic  
Pump

**Wait 300mS.**

**Washing Station goes up** (570mS).

**OFF = EV1** (emptying of the channel of the Probe1 is disabled), **uP4** (OFF uPump Probe5) ,  
**P.M0** (OFF Pump peristaltic).

**Cycle lasts 3790mS used 1200ul(+/-150ul) of Dist.Water.**

\*\*\*\*\*

## **Probe 3**

**ON = EV1** (emptying of the channel of the Probe1 is enabled), **uP4** (ON uPump Probe5) , **P.M0**  
ON peristaltic pump, remains ON until the end of the cycle

**Washing Station goes down** (570mS).

**ON EV5** (filling of the channel of the Probe3 is enabled)

**Erogo 400uL(+/-50uL) of Dist.Water.**

**Washing Station goes up** (570mS).

**OFF = EV1** (emptying of the channel of the Probe1 is disabled), **uP4** (OFF uPump Probe5) ,  
**P.M0** (OFF Pump peristaltic), **EV5**(filling of the channel of the Probe3 is disabled)

**Cycle lasts 3790mS used 400ul(+/-50ul) of Dist.Water**

\*\*\*\*\*

**Before the Probe 4, the cuvette is moved in front of the photometer to be verified.**

\*\*\*\*\*

## **Probe 4**

**ON = EV1** (emptying of the channel of the Probe1 is enabled), **uP4** (ON uPump Probe5) , **P.M0**  
(ON Pump peristaltic). remains ON until the end of the cycle

**Washing Station goes down** (570mS).

**ON EV3** (emptying of the channel of the Probe4 is enabled)

**Wait 600mS**

**Washing Station goes up** (570mS).

**OFF = EV1** (emptying of the channel of the Probe1 is disabled), **uP4** (ON uPump Probe5) , **P.M0** (OFF Pump peristaltic), **EV3** (emptying of the channel of the Probe4 is disabled).

**Cycle lasts 3790mS**

\*\*\*\*\*

\*\*\*\*\*

## **Probe 5**

**ON = EV1** (emptying of the channel of the Probe1 is enabled), **uP4** (ON uPump Probe5) , **P.M0** (ON Pump peristaltic). remains ON until the end of the cycle

**Washing Station goes down** (570mS).

**Wait 2100mS**

**Washing Station goes up** (570mS).

**OFF = EV1** (emptying of the channel of the Probe1 is disabled), **uP4** (ON uPump Probe5) , **P.M0** (ON Pump peristaltic).

**Cycle lasts 3790mS**

\*\*\*\*\*